W/O: 429234	FILE NAME: 429234_MiaSmart_GBL_UG_2227647	PROJECT: Project Bouquet	BRAND: Mia Smart
COMPONENT: User Guide	LANGUAGE SET: EN/FR/DE/SP/IT/KO/SC/TC/PT/RU	DIE NO: PHX-21952743-01-PPR-2C	
REVISION HISTORY: R1 12/12/17: Created mechanical [pd]. 12/22/17: New hi-res art [pd]. R2 01/02/18: ESKO changes [pd]. 01/12/17: ESKO revisions [pd]. R3 01/24/18: ESKO revisions [pd]. R4 02/08/18: ESKO revisions [pd].			


clarisonic

MIA SMART™

 EN User Guide FR Mode d'emploi DE Benutzerhandbuch SP Guía de uso IT Guida all'utilizzo KO 사용 설명서 SC 用户指南 TC 使用者指南 PT Manual do utilizador

EN	English	3 - 4
FR	Français	5 - 7
DE	Deutsch	8 - 10
SP	Español	11 - 13
IT	Italiano	14 - 16
KO	한국어	17 - 18
SC	简体中文	19 - 22
TC	繁體中文	23 - 24
PT	Português	25 - 26

MIA SMART™

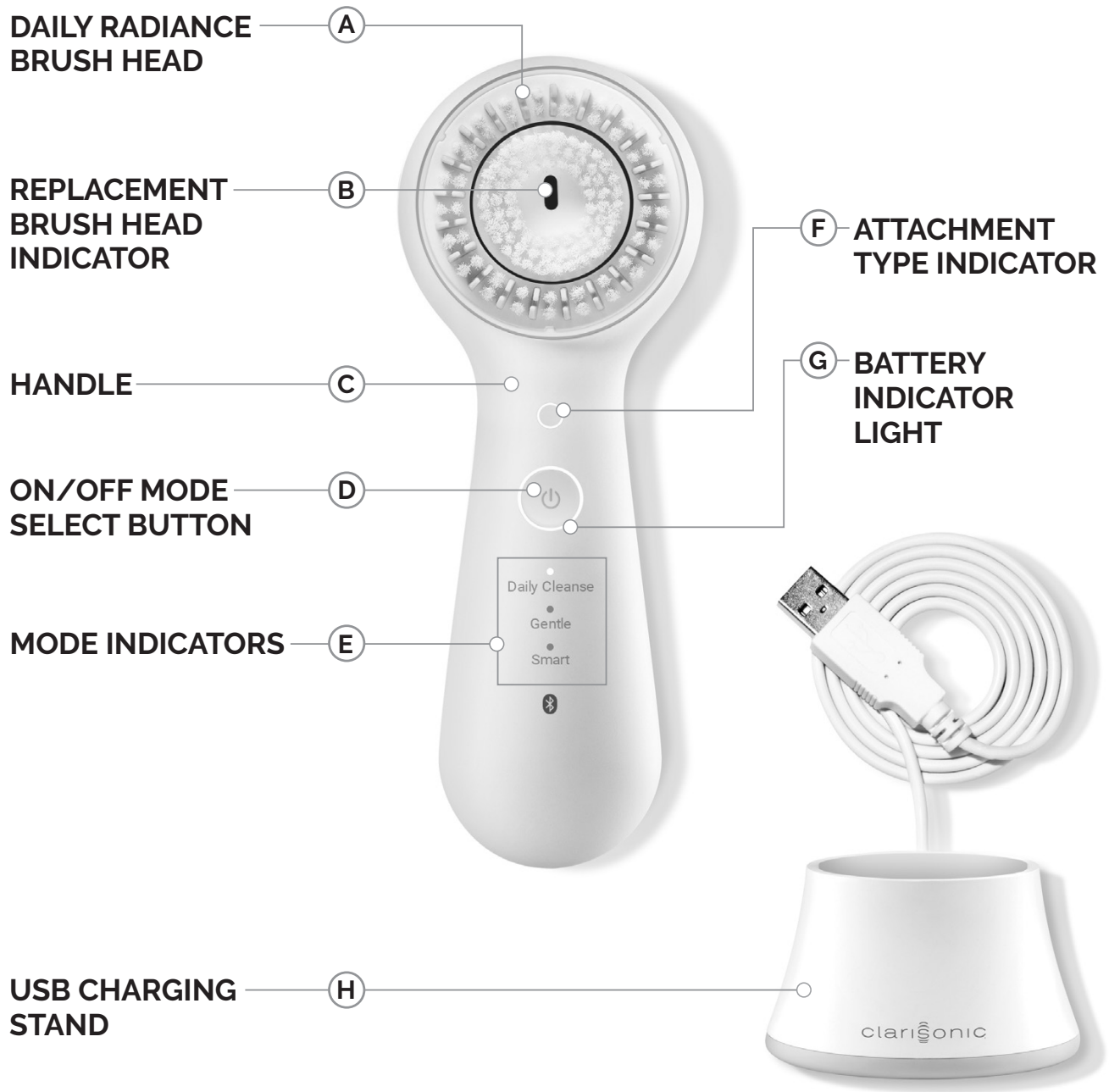- A**— EMBOUT DE BROUSSE RADIANCE QUOTIDIEN
DAILY RADIANCE BÜRSTENKOPF
CABEZAL DE CEPILLO RADIANCE DIARIO
TESTINA RADIANCE PER IL DISPOSITIVO
데일리 래디언스 브러시 헤드
每日焕彩刷头
日常使用的音波亮白刷頭
CABEÇA DE ESCOVA PARA
RADIANCE DIÁRIA
- B**— INDICATEUR DE REMPLACEMENT DE
L'EMBOUT DE BROUSSE
BÜRSTENKOPF-WECHSELANZEIGE
INDICADOR DE SUSTITUCIÓN DEL CABEZAL
DEL CEPILLO
INDICATORE SOSTITUZIONE TESTINA
브러시 헤드 교체 표시등
更換刷头指示灯
更換刷頭指示燈
INDICADOR DE SUBSTITUIÇÃO DA ESCOVA
- C**— MANCHE
GRIFF
MANGO
UNITÀ PRINCIPALE
기기
手柄
機身
APARELHO
- D**— BOUTON INTERRUPTEUR MARCHÉ/ARRÊT
EIN/AUS-KNOPF UND MODUSAUSWAHL
BOTÓN DE SELECCIÓN DE MODO
ENCENDIDO/APAGADO
PULSANTE ON/OFF PER SELEZIONARE
LA MODALITÀ
ON/OFF 버튼
开/关模式选择按钮
開關模式選擇按鈕
BOTÃO DE SELEÇÃO ON/OFF
- E**— INDICATEURS DES DIFFÉRENTS MODES
MODUSANZEIGE
INDICADORES DE MODO
INDICATORI DELLA MODALITÀ
모드 표시등
模式指示灯
模式指示
INDICADORES DE MODO
- F**— INDICATEUR DE TYPE D'EMBOUT
AUFSATZANZEIGE
INDICADOR DE TIPO DE CABEZAL
INDICATORE DEL TIPO DI ACCESSORIO
헤드 유형 표시등
附件类型指示灯
配件種類指示燈
INDICADOR DO TIPO DE CABEÇA
- G**— TÉMOIN LUMINEUX DE LA BATTERIE
AKKUANZEIGE
LUZ INDICADORA DE BATERÍA
INDICATORE LUMINOSO DELLA BATTERIA
배터리 표시등
电量指示灯
電池指示燈
LUZ INDICADORA DE PILHA
- H**— SOCLE DE RECHARGEMENT USB
USB-LADESTATION
BASE DE CARGA USB
BASE PER LA RICARICA USB
스탠드형 USB 충전기
USB 充电基座
USB 充電座
SUPORTE PARA CARREGAR VIA USB

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including but not limited to the following:

READ ALL INSTRUCTIONS BEFORE USE

DANGER – To reduce the risk of electrocution:

1. **Do not** place or store the charger or adapter where it can fall or be pulled into a bathtub or sink.
2. **Do not** place the charger in or drop the charger or adapter into water or other liquids. Charger unit is not for immersion or for use in the shower.
3. **Do not** reach for a charger or adapter that has fallen into water. Unplug immediately.

WARNING – To reduce the risk of burns, electrocution, fire, or injury to persons:

1. The USB wall adapter must be certified to UL 60950-1 and IEC 60950-1 by appropriate safety agencies in your country with an output marked Class 2 or the limited power source (LPS) requirements. Electrical Input Requirements: 5 V --- 0.5 A.
2. Close supervision is necessary when this product is used by, on, or near children or physically challenged persons.
3. This product is for use in home, spas, salons and similar locations. Use this product only for its intended use as described in this manual.
4. Keep the cord away from heated surfaces.
5. **Do not** use the charger outdoors.
6. Use only the brush head, applicators and accessories supplied by the manufacturer.
7. To charge your Clarisonic device, connect the charger into USB adapter.
8. **Do not** use attachments not recommended by the manufacturer.
9. Never operate the Clarisonic device if it is damaged. **Do not** operate the charger if it has a damaged cord or adapter, if it is not working properly, if it has been dropped or damaged, or dropped into water.

SAVE THESE INSTRUCTIONS

- The device, charger and power adapter have no serviceable parts.

IMPORTANT CAUTIONS AND WARNINGS

- If you are under a doctor's care for any skin condition or if you think you may have a skin condition, please consult your physician before using the Clarisonic device.
- We recommend that you **do not** share heads.
- **Do not** use the Clarisonic device with jagged exfoliating particles such as shells, micro-beads, and substrates.
- Using the Clarisonic device to remove makeup may stain the bristles, but will not affect performance. If discoloration occurs, clean the brush head with soap and water or makeup brush cleanser.
- This product uses magnets.
- DO NOT twist or directly pull on silicone tip when removing massage head.
- **DO NOT USE** the massaging head directly on the eye, eye lid or soft tissue under the eyes as shown in the image below.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC (Federal Communications Commission) Rules. This equipment complies with Canadian standard ICES-003.

U.S. FCC Statement: This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 and Part 18 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna. Increase the separation between the equipment and receiver. Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by Clarisonic could void the user's authority to operate the equipment.

Canada IC Statement: This device complies with industry standards of Canada license exempt RSS standards. Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Treatment of electrical and electronic equipment applicable at the end of life in the European Union and other European countries with selective collection systems: This product should not be treated as household waste but should be handled at a designated collection point for the recycling of electrical and electronic equipment. Please contact your local authority for further details of your nearest designated collection point.

RED statement for EU: This device uses Bluetooth Low Energy technology (BLE v 4.2 only and not Bluetooth Classic). Frequency Band: 2.4 GHz (2.402-2.480 GHz, Max Power: 0 dBm).

CHARGING

Allow approximately 18 hours for a full charge. A fully charged device averages 60 minutes of use depending on the type of brush head used.

CARING FOR YOUR DEVICE

To keep your Clarisonic device and accessories in optimum condition, care for it as follows:

Do not use abrasive cleansers, chemicals or exfoliators to clean any part of the Clarisonic device or accessories.

Do not put any part of the Clarisonic device or accessories in the dishwasher.

The device is sealed and fully waterproof. Never attempt to open the device for any reason. If you are not using the device for long periods of time (weeks/months), unplug the charger and store the device in a cool dry place away from direct sunlight. Charge the device every 6 months to maintain its charging capabilities.

DEVICE DISPOSAL AT END OF PRODUCT LIFE: Your Clarisonic device is sealed and contains rechargeable NiMH batteries. These batteries are non-replaceable and the sealed device is non-repairable. When your device reaches end of life: it no longer holds a charge, does not function and is out of warranty; this device must be recycled or disposed of per your local regulations. Do not dispose of it in fire.

CUSTOMER SERVICE AND WARRANTY INFORMATION: www.clarisonic.com/service

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All rights reserved/ tous droits réservés/ Alle Rechte vorbehalten/ Todos los derechos reservados/
Pacific Bioscience Laboratories, Inc.
© 2018 Pacific Bioscience Laboratories, Inc.
L'Oréal Australia Pty Ltd, 564 St Kilda Road, Melbourne, VIC 3004, Australia
L'Oréal Hong Kong Limited, 34th Floor, Tower One, Times Square,