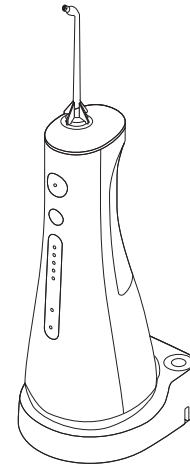
Panasonic

Operating Instructions
(Household) Rechargeable Oral Irrigator
Model No. EW1511

English	3	Português	123	Česky	245
Deutsch	23	Svenska	143	Slovensky	265
Français	43	Dansk	163	Magyar	285
Italiano	63	Norsk	183	Română	305
Español	83	Suomi	203	Türkçe	325
Nederlands	103	Polski	223		

Operating Instructions (Household) Rechargeable Oral Irrigator Model No. EW1511

Contents

Safety precautions	6	Frequently asked questions.....	17
Intended use.....	9	Troubleshooting.....	18
Parts identification	10	Battery life.....	20
Charging.....	10	Removing the built-in rechargeable battery.....	20
How to use.....	12	Specifications	21
After use	15	Changing the nozzle (Sold separately)...	21
Wall mounting.....	17		

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

Warning

English

- Do not immerse the AC adaptor and the charging stand in water or wash it with water.
Doing so may cause electric shock or fire due to a short circuit.
- Always unplug the adaptor from a household outlet when cleaning the AC adaptor and the charging stand.
Failure to do so may cause electric shock or injury.
- The supply cord cannot be replaced. If the cord is damaged, the AC adaptor should be scrapped.
Failure to do so may cause an accident or injury.
- Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.

- Oral irrigators shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions.
- Children shall not play with the appliance.
- Do not use anything other than the supplied AC adaptor for any purpose. Also, do not use any other product with the AC adaptor. (See page 10.)
Doing so may cause burn or fire due to a short circuit.
- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type reference of power supply unit is marked near the symbol.

Safety precautions

English

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

DANGER

Denotes a potential hazard that will result in serious injury or death.

WARNING

Denotes a potential hazard that could result in serious injury or death.

CAUTION

Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe types of instructions to be observed.

This symbol is used to alert users to a specific operating procedure that must not be performed.

This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

WARNING

► Preventing accidents

Do not store within the reach of children or infants. Do not let them use it.

- Doing so may cause an accident or injury due to accidental ingestion of removable parts.

► Power supply

Do not connect or disconnect the adaptor or the appliance plug to a household outlet with a wet hand.

- Doing so may cause electric shock or injury.

Do not place the AC adaptor and the charging stand over or near water filled sink or bathtub.

Do not use when the AC adaptor is damaged or when the fitting into the household outlet is loose.

Do not damage, modify, forcefully bend, pull, twist, or bundle the cord. Also, do not place anything heavy on or pinch the cord.

- Doing so may cause electric shock or fire due to a short circuit.

Do not use in a way exceeding the rating of the household outlet or the wiring.

- Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.

WARNING

Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the AC adaptor. Fully insert the adaptor.
- Failure to do so may cause fire or electric shock.

► In case of an abnormality or malfunction

Immediately stop using and remove the adaptor if there is an abnormality or malfunction.

- Failure to do so may cause fire, electric shock, or injury.

<Abnormality or malfunction cases>

- **The main unit, adaptor, charging stand or cord is deformed or abnormally hot.**
 - **The main unit, adaptor, charging stand or cord smells of burning.**
 - **There is abnormal sound during use or charging of the main unit, adaptor, charging stand or cord.**
- Immediately request inspection or repair at an authorized service centre.

► This product

This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, or leave in a high temperature environment.
- Doing so may cause overheating, ignition, or explosion.

Do not modify or repair.

- Doing so may cause fire, electric shock, or injury.
- Contact an authorized service centre for repair (battery change, etc.).

Never disassemble except when disposing of the product.

- Doing so may cause fire, electric shock, or injury.

► Cleaning

Regularly clean the power plug and the appliance plug to prevent dust from accumulating.

- Failure to do so may cause fire due to insulation failure caused by humidity. Disconnect the adaptor and wipe with a dry cloth.

CAUTION

► Preventing damage to teeth and gums

People who are unable to operate should not use this appliance.

People with weak oral senses should not use this appliance.

- Doing so may cause injury or damage to teeth and gums.

CAUTION

People who may have periodontal disease, have teeth being treated, or concerned about symptoms within their mouth should consult a dentist before use.

- Failure to do so may cause injury or damage to teeth and gums.

► This product

Never use for a purpose other than oral cleaning.

- Doing so may cause an accident or injury when the stream of water is pointed toward the eyes, nose, ears, tongue, or throat.

Do not fill the tank with water warmer than 40 °C.

- Doing so may cause burns.

Do not share nozzle with your family or other people.

- Doing so may cause infection or inflammation.

Do not place a metallic object such as a coin or a clip on top of the charging stand.

- Doing so may cause burns due to excessive heat.

► Note the following precautions

Do not allow metal objects or trash to adhere to the charging plug, the power plug or the appliance plug.

- Doing so may cause electric shock or fire due to a short circuit.

Do not drop or subject to shock.

- Doing so may cause injury.

Do not wrap the cord around the adaptor or the charging stand when storing.

- Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.

Disconnect the adaptor from the household outlet when not charging.

- Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

Disconnect the adaptor or the appliance plug by holding onto it instead of the cord.

- Failure to do so may cause electric shock or fire due to a short circuit.

After use, clean the main unit and the nozzle and dry them. (See page 15.)

- Failure to do so may cause damage to your health due to propagation of mold and bacteria.

Handling of the removed battery when disposing

DANGER

The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.

Do not charge the battery after it has been removed from the product.

Do not throw into fire or apply heat.

Do not hit, disassemble, modify, or puncture with a nail.

 Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.

Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.

Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat. Never peel off the tube.

- Doing so may cause overheating, ignition, or explosion.

WARNING

After removing the rechargeable battery, do not keep it within reach of children and infants.

 - The battery will harm the body if accidentally swallowed.

If this happens, consult a physician immediately.

 If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.

- Failure to do so may cause eye injury. Consult a physician immediately.

CAUTION

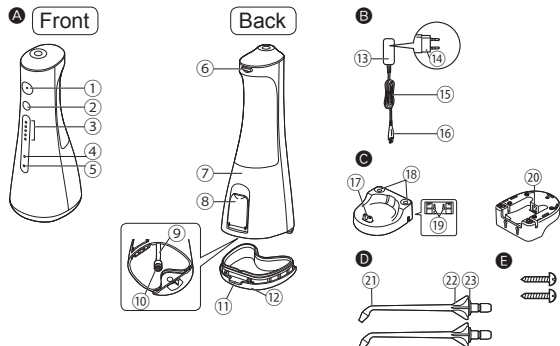
 If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.

- Failure to do so may cause inflammation.

Intended use

- When you first begin using the oral irrigator, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the oral irrigator for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the oral irrigator and consult your dentist.
- Do not put salt, mouthwash or other chemicals into the water tank. Doing so may cause malfunction.
- Do not dry the main unit with a dryer or a fan heater. Doing so may cause malfunction or deformation of the parts.
- Store the appliance in a place with low humidity after use. Leaving the appliance in a bathroom may cause malfunction.

Parts identification

A Main unit

- ① Power switch (Ⓚ)
- ② Level switch
- ③ Water pressure level LED indicator (1~5 levels)
- ④ Low battery indicator (blue)
- ⑤ Charge status indicator (red)
- ⑥ Nozzle release button
- ⑦ Water tank (Inside)
- ⑧ Water tank cap
- ⑨ Suction hose
- ⑩ Filter
- ⑪ Bottom cover release button
- ⑫ Bottom cover

B AC adaptor (RE7-87)

(The shape of the power plug differs depending on the area. One of the power supply units specified is provided with

this appliance.)

- ⑬ Adaptor
- ⑭ Power plug
- ⑮ Cord
- ⑯ Appliance plug

C Charging stand (RC8-4)

- ⑰ Charging plug
- ⑱ Nozzle stand
- ⑲ Screw holes for wall mounting
The charger can be mounted on a wall using two wood screws.

- ⑳ Stand socket

D Nozzle (x2)

- The nozzle is a consumable.
- ⑳ Nozzle
- ㉑ Nozzle Handle
- ㉒ Identification ring

Accessory

- E 2 screws

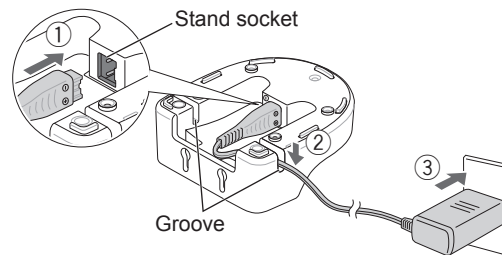
Note

There may be moisture inside the main unit (water tank and suction hose) left from the distilled water used for product testing.

Charging


1 Connect the AC adaptor.

- ① Insert the appliance plug firmly to the stand socket.
- ② Press the cord into either the left or right groove, depending on where it will be placed.
 - Place the charging stand on a flat and stable surface.
- ③ Insert the adaptor into a household outlet.

► About the charging time

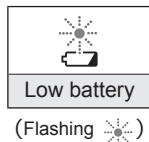- Charging is required when the low battery indicator (blue) flashes.
- Turning off the power switch “①” when the battery has little remaining charge, the low battery indicator will light off after 5 seconds.

2 Insert the main unit straight into the charging plug with the front facing outward.


► Indicator display during charging

		
Charging in progress	Charging complete	Charging abnormal (Flashing in red twice a second)
(Lights on ● Lights off ○)		(Flashing ⚡)

► How to check if charging is completed

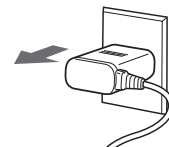- After charging is completed, if you connect and then disconnect the main unit to and from the charging stand,  (the charge status indicator) lights on and then lights off after 5 seconds.

3 For safety and energy saving, disconnect the power plug.

- Charging the battery for more than 1 hour will not cause any damage to the main unit.
- **Operating time after a full charge**
Approx. 10 minutes*
(When the water pressure level is set to 5)
- * It varies with the ambient temperature and usage conditions.
- Charging is completed after approx. 1 hour.
(The charging time may vary with the remaining battery capacity.)

Notes

- The charging time may differ when the appliance is used for the first time, or if it has not been used for more than 6 months.
In addition,  (the charge status indicator) may not light on in the first several minutes of charging, but it will light on as the charging proceeds.
- The battery will deteriorate if not used for over 6 months, so make sure to perform a full charge at least once every 6 months.
- The recommended ambient temperature for charging is 10 °C - 35 °C. It may take longer to charge, or  will flash rapidly and the battery may fail to charge properly in extremely high or low temperatures.
- To keep charging the battery will not affect battery performance.

How to use

1 Fill the water tank with water.

① Open the water supply cover.

② Add water while keeping the main unit level.

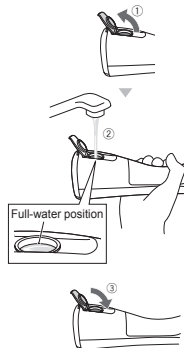
③ Close the water supply cover securely.

- Always use new tap water or lukewarm water.
- If the water feels cold, use lukewarm water with a temperature of 40 °C or less.


► If the water supply cover comes off

① Align the washer of the water supply cover with the water supply port.


② Press in until a “click” is heard.


2 Install the nozzle.

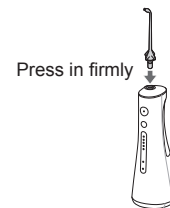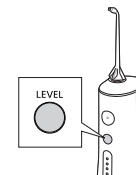3 Press the level switch, and select the water pressure level.

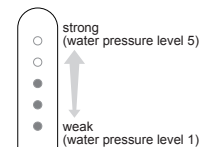- When you press the level switch, the water pressure level that was used last time is selected. Then each time you press the level switch, the water pressure level is switched as follows: 1 (weak) → 2 → 3 → 4 → 5 (strong) → 1...

- **About the water pressure level memory function. The last selected water pressure level will be remembered, from which the appliance will start next time.**

- When the appliance has not been used for a long time or when the battery is used up, the water pressure level will return to 1.

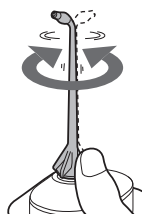
4 Clean inside your mouth with water jet.

- ① Place the nozzle in your mouth.
- ② Close your mouth gently.
- ③ Turn on the power switch “ⓘ”.

- You can switch the water pressure level by pressing the level switch, even while using the appliance.
- If the main unit is tilted too much, the water does not come out from the nozzle.
- ▶ **Operation time with full water**
Approx. 60 seconds.
(When the water pressure level is set to 5)

▶ To adjust the orientation of the nozzle

You can adjust the direction of water jet by rotating the nozzle handle.

Point 1 /
Keep your face
over the basin.

※ It can also be used in a bath.

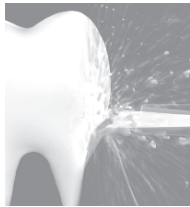
Notes

- When using the appliance for the first time, please press the power switch “ⓘ” and check the water pressure level while water is coming out. (Please start using the appliance from water pressure level.)
- The appliance may not operate in an ambient temperature approximately 5 °C or lower. The appliance may not work outside the recommended temperature range.
- Do not push the tip of the nozzle too hard against the teeth or gums.
- When you feel that the irritation is strong, turn off the power switch “ⓘ” and use the appliance from water pressure level.
- Do not turn on the power switch “ⓘ” when the water tank is empty unless you are drying the moisture inside the main unit. (Doing so may cause malfunction.)

Tips of use

► Interdental cleaning

- Direct the water jet onto the space between the teeth.
- Direct from the front and rear sides of the teeth.
- ◆ With the water jet, the space between the teeth that cannot be reached with a toothbrush is cleaned thoroughly, and the dirt trapped in periodontal pockets is removed completely. In addition, effective care is provided for the gums by comfortable intermittent stimulation.

► Periodontal pocket cleaning

- Aim the water jet at the space between teeth and gums.
- Slowly shift the water jet along the gums.
- Clean the front and back of the teeth and behind the molars as well.
- If you are concerned about sensitive gums, clean with a gentle water jet setting.

► Gums care

- Stimulate the gums with the water jet.
- If you are concerned about sensitive gums, stimulate with a gentle water jet setting.
- ◆ Since plaque adhering to the surface of teeth cannot be removed only with water jet, brushing is always necessary. It is recommended to clean inside your mouth with the water jet after brushing your teeth.

After use

Please clean the appliance after each use.

- 1 To prevent the water from splashing about, turn off the power switch "①" while keeping the nozzle in your mouth.

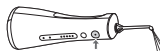
- 2 Throw away the remaining water, and wash the nozzle with water.

- ① Pour out the remaining water. Change the water every time you use the appliance.

- ② Turn on the power switch "①" to remove the remaining water. Keep it on until no water comes out.

- ③ Remove the nozzle, and wash it with water.

Pull it out while pressing.

Nozzle release button

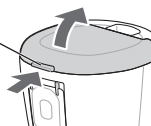
- 3 Remove the bottom cover, and wash the main unit (as a whole) and inside of the tank.

► How to remove the bottom cover


Bottom cover release button

Lift it up while pressing.

4 After drying, put the main unit upright on the charging stand.

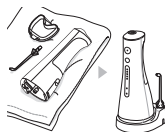- ① Wipe any water away from the nozzle and the main unit with a towel or the like, and dry them in a well-ventilated shade.
- ② Store it in the upright position
 - Attach the bottom cover, and put the main unit upright on the charging stand.

► How to attach the bottom cover

- a Insert the claws (2 places) on the front of the bottom cover into the main unit.
- b Press in the hooks (2 places) on the back.

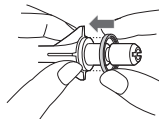
Notes

- Do not wash with lukewarm water exceeding 50 °C.
- Use neutral kitchen detergent and rinse well with water to ensure no detergent remains after cleaning.
- Never use thinners, benzene, or alcohol. (Doing so may cause malfunction, or cracking or discoloration of the parts.)

► If the identification ring comes off

Align the nozzle with the identification ring (in any direction) to install it.

- Even if the identification ring is not attached, the performance of the nozzle is not affected.


► When the water jet is poor

Remove the dirt from the filter (the tip of the suction hose).

Filter

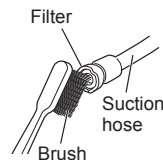Rub it gently with a soft brush.


► When dirt is a concern

Wipe with a dry cloth. (Do not rinse it with water.)

Charging stand

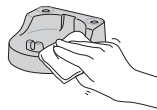AC adaptor


Wall mounting

Be sure to choose a solid wall or a crosspiece to install the charging stand.


<How to install>

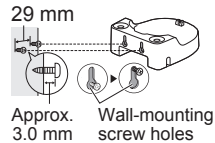- ① Tighten two commercially available round wood screws (nominal diameter: 2.7 - 3.1 mm) at intervals of 29 mm so that the screw head can float about 3.0 mm from the wall.
- ② Hook the wall-mounting screw holes on the back of the charging stand on the screws, and fix them.

<How to remove>

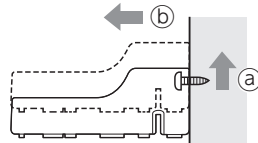- a Lift the charging stand.
- b Pull it forward.

Frequently asked questions

	Question	Answer
Battery	Will the battery be exhausted after a long period of disuse?	▶ Please charge it once every 6 months if you have not used it for over 6 months since the battery may deteriorate.
	Can the battery be charged before every use?	▶ The service life of the battery will not be affected thanks to the use of a lithium ion battery.

English

Troubleshooting

English

	Problem	Possible cause	Action
While charging	The charge status indicator does not light on.	▶ The appliance was just purchased or it has not been used for a while.	▶ Although the indicator may not light on in the first several minutes of charging, it will light on as the charging proceeds. If it does not light on after a while, please request a repair at your dealer or Service Centre.
	The charge status indicator flashes rapidly (twice a second).	▶ The main unit is not properly plugged in the charging stand.	▶ Check to see that the main unit stands upright so that its bottom comes in contact with the charging plug of the charging stand, and that the charge status indicator is on. (See page 11.)
		▶ The ambient temperature is too low or too high during charging.	▶ Please charge within the recommended temperature range of 10 °C to 35 °C. ▶ If the indicator flashes even when you charge it within the recommended temperature range, contact your dealer or Service Centre for repair.
	The appliance can be used for just several minutes even after a full charge.	▶ The main unit is not properly plugged in the charging stand.	▶ Check to see that the main unit stands upright so that its bottom comes in contact with the charging plug of the charging stand, and that the charge status indicator is on. (See page 11.)
	The low battery indicator does not flash, but the appliance stops working.	▶ The battery has reached the end of its life (approx. 3 years).	▶ Replace it with a new one. (See page 20.) Please contact your dealer or Service Centre.

	Problem	Possible cause	Action	
While using	The appliance does not operate.	▶ The battery has reached the end of its life (approx.3 years).	▶ Please remove the built-in rechargeable battery. (See page 20.) ▶ Please contact your dealer or Service Centre.	English
		▶ The appliance is used with approximately 5 °C or lower.	▶ Do not use the appliance with approximately 5 °C or lower.	
	The water jet is poor.	▶ The nozzle is clogged. The outlet of nozzle is squashed.	▶ Replace the nozzle with a new one.	
		▶ The filter is clogged.	▶ Remove the dirt from the filter. (See page 16.)	
	No water comes out.	▶ The main unit is tilted too much.	▶ Please use the appliance with the main unit upright. (See page 13.)	
		▶ It may take some time to suck up water.	▶ Please wait until water comes out.	
	Water leaks from the bottom cover.	▶ Dirt is accumulated on the rubber part of the bottom cover.	▶ Please remove the dirt. (See page 16.)	

If the problems still cannot be solved, contact the store where you purchased the unit or a service centre authorized by Panasonic for repair.

Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge.

(The operating life of the battery may differ depending on usage or storage conditions.)

Have the battery replaced by an authorized service centre.

Removing the built-in rechargeable battery

Remove the built-in rechargeable battery before disposing of the appliance.

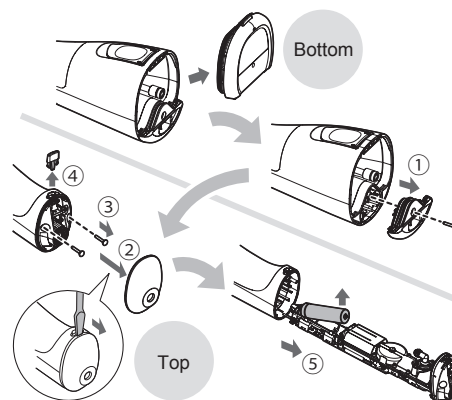
The battery must be disposed of safely.

Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it will no longer be waterproof, which may cause it to malfunction.

- Remove the appliance from the charging stand and adaptor, press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Remove the bottom cover.
- Perform steps ① to ⑤ to disassemble the main unit and remove the battery with a screwdriver.
- Insulate the ⊕ and ⊖ terminals by applying tape to them. (To prevent short circuit of terminals).

20

For environmental protection and recycling of materials

This appliance contains a Lithium-ion battery. Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

Power source	See the name plate on the product (Automatic voltage conversion)
Charging time	Approx. 1 hours
Pump frequency	Approx. 1600 pulses/minute (level 5) Approx. 960 pulses/minute (level 1)
Water pressure	LEVEL 5: Approx. 647 kPa (approx. 6.6 kgf/cm ²) LEVEL 1: Approx. 196 kPa (approx. 2.0 kgf/cm ²)
Operating time	LEVEL 5: Approx. 10 minutes* (at 20 °C when fully charged)
Tank volume	Approx. 200 ml (approx. 60 seconds*)
Dimensions	Main unit: 210 (H) x 75 (W) x 85 (D) mm (Not including nozzle) Charging stand: 40 (H) x 90 (W) x 100 (D) mm
Mass	Main unit: Approx. 270 g (Including nozzle) Charging stand: Approx. 115 g
Airborne Acoustical Noise	71 (dB (A) re 1 pW)

This product is intended for household use only.

Changing the nozzle (Sold separately)

- If the tip of the nozzle gets deformed, please replace the nozzle with a new one.
- The nozzle should be replaced once every 6 months or so due to hygiene reasons, even if it is not deformed.
- Optional nozzles with identification ring to allow several persons to share the main unit.

New nozzle (tip)

Deformed (worn) nozzle
(tip)

Part number	Replacement Nozzle for EW1511	
WEW0983	Nozzle	

Disposal of Old Equipment and Batteries
Only for European Union and countries with recycling systems

English

These symbols on the products, packaging, and/or accompanying documents mean that used electrical and electronic products and batteries must not be mixed with general household waste.

For proper treatment, recovery and recycling of old products and batteries, please take them to applicable collection points in accordance with your national legislation.

By disposing of them correctly, you will help to save valuable resources and prevent any potential negative effects on human health and the environment.

For more information about collection and recycling, please contact your local authority.

Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

Note for the battery symbol
(bottom symbol)

This symbol might be used in combination with a chemical symbol. In this case it complies with the requirement set by the Directive for the chemical involved.

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FABRİKA

Panasonic Wanbao Appliances Beauty and Living (Guangzhou) Co., Ltd.

33, Wanbao North Street, Wanbao Industry Zone, Zhongcun, Panyu District, Guangzhou City, Guangdong Province, China

YETKİLİ TEMSİLCİ

İthalatçı Firma:

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Çobançeşme Mah. Bilge Sok. No:17

34196 Yenibosna- İstanbul

Tel: 0212 454 60 70 Faks: 0212 454 60 20

www.tesan.com.tr

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Panasonic Corporation
<http://www.panasonic.com>
© Panasonic Corporation 2019

⑦ EN, GE, FR, IT, SP, DU, PT, SW, DA, NW, FI, PL, CZ, SK, HU, RO, TK

Printed in China
EW970015111503 S0319-0



Panasonic

Operating Instructions
(Household) Sonic Vibration Toothbrush
Model No. EW-DM81





(GB) English.....	5	(DK) Dansk.....	113	(CZ) Český.....	221
(D) Deutsch	23	(P) Português.....	131	(SK) Slovensky.....	239
(F) Français.....	41	(N) Norsk	149	(H) Magyar.....	257
(I) Italiano.....	59	(S) Svenska.....	167	(RO) Română	275
(NL) Nederlands.....	77	(FIN) Suomi.....	185		
(E) Español.....	95	(PL) Polski.....	203		

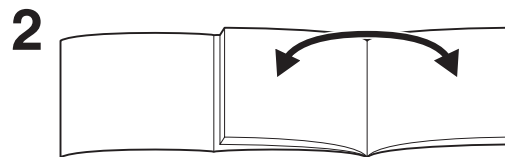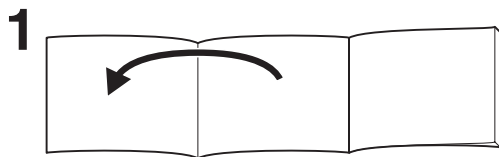Parts identification

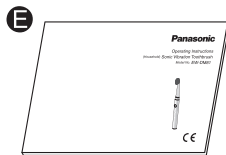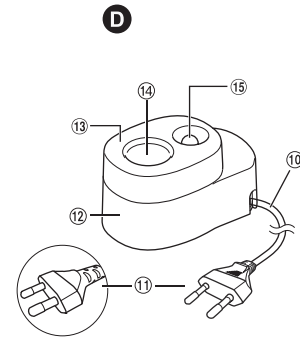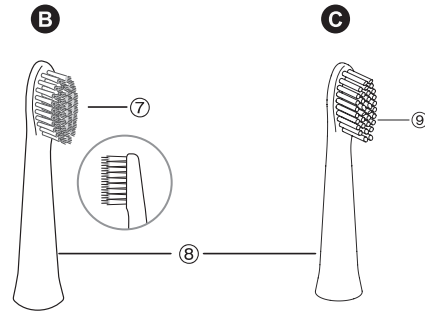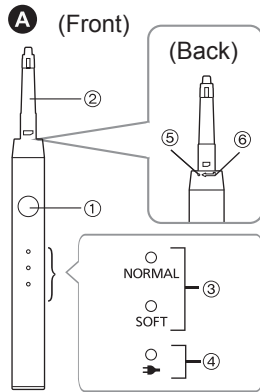
Panasonic

Operating Instructions (Household) Sonic Vibration Toothbrush

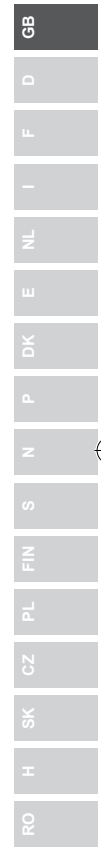Model No. EW-DM81

Contents

Safety precautions	8	Frequently Asked Questions.....	17
Intended use.....	11	Troubleshooting.....	18
Parts identification	12	Battery life.....	19
Charging the toothbrush	13	Removing the built-in rechargeable battery.....	19
Before cleaning your teeth.....	13	Specifications	20
Using the toothbrush	14	Changing the brush head (sold separately).....	20
Tips on usage	15		
After use	16		

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Warning

- Be sure to switch off the appliance and remove the power plug from a wall outlet before performing any maintenance. Failure to do so may result in electric shock.
- Do not wash the charger. Doing so may result in electric shock.
- The supply cord cannot be replaced. If the cord is damaged, the charger should be scrapped.
- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.

- Do not use anything other than the supplied charger to charge the product for any purpose. Also, do not use any other products with the charger supplied. (See page 12.)
- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type references of power supply unit and charger are marked near the symbol.

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 DANGER	Denotes a potential hazard that will result in serious injury or death.
 WARNING	Denotes a potential hazard that could result in serious injury or death.
 CAUTION	Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe the type of instructions to be observed.

	This symbol is used to alert users to a specific operating procedure that must not be performed.
	This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

WARNING

► Preventing accidents

-  **Do not store within the reach of children or infants. Do not let them use it.**
 - Doing so may cause an accident or injury due to accidental ingestion of removable parts.

► Charger

-  **Do not connect or disconnect the charger to a wall outlet with a wet hand.**
 - Doing so may cause electric shock or injury.

Do not place the charger over or near water filled sink or bathtub.

Do not immerse the charger (except for the charger cover) in water or wash it with water.

-  **Do not use when the charger is damaged or when the fitting into the wall outlet is loose. Do not damage, modify, or forcibly bend, pull, or twist the cord. Also, do not place anything heavy on or pinch the cord.**
 - Doing so may cause electric shock or fire due to a short circuit.

Do not use in a way exceeding the rating of the wall outlet or the wiring.

-  - Exceeding the rating by connecting too many power plugs to one wall outlet may cause fire due to overheating.

-  **Do not use charger near the bathtub or in the shower.**
 - Doing so may lead to electric shock and fire.

WARNING

Always ensure that the product is operated on an electric power source matched to the rated voltage indicated on the charger.

Fully insert the power plug.

- Failure to do so may cause fire or electric shock.

Always unplug the charger from a wall outlet when cleaning the charger.

- Failure to do so may cause electric shock or injury.

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- Failure to do so may cause fire due to insulation failure caused by humidity.
Disconnect the charger and wipe with a dry cloth.

► In case of an abnormality or failure

Immediately stop using and remove the charger if there is an abnormality or failure.

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<Abnormality or failure cases>

- **The main unit, charger, or cord is deformed or abnormally hot.**
- **The main unit, charger or cord smells of burning.**
- **There is abnormal sound during use or charging of the main unit, charger or cord.**
- Immediately request inspection or repair at an authorized service centre.

► The product

This product has a built-in rechargeable battery. Do not heat or throw the product into fire.

- Doing so may cause battery fluid leakage, battery overheating or explosion.

Do not modify or repair.

- Doing so may cause fire, electric shock, or injury.
→ Contact an authorized service centre for repair (battery change etc.).

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CAUTION

► The product

Never use for a purpose other than cleaning your teeth.

- Doing so may cause an accident or injury.

Do not share your brush head with your family or other people.

- Doing so may result in infection or inflammation.

► Note the following precautions

Do not drop or hit against other objects.

- Doing so may cause injury.

Do not wrap the cord tightly around the charger when storing.

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CAUTION

Disconnect the charger from the wall outlet when not charging.

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When unplugging the power plug, hold the power plug instead of the cord.

- Otherwise, the wire in the cord may break with the load, which may cause electric burns or electric shock, or may cause fire due to a short circuit.

► Protecting your teeth and gums

Do not use the brush with toothpaste that contains a high abrasive content, such as toothpaste for removing tar or for whitening.

Do not apply the brush head strongly against teeth or gums, or use it continuously against same spot.

Do not apply the plastic part of the main unit or brush head to teeth.

Do not let anyone who cannot operate the brush by themselves use it.

People with low sensitivity in their mouth should not use the brush.

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People who may have periodontal disease, have teeth being treated, or who are concerned about symptoms in their mouth should consult a dentist before use.

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Disposing of the rechargeable battery

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- Do not throw into fire or apply heat.
- Do not solder, disassemble or modify the battery.
- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewelry such as necklaces and hairpins.
- Never peel off the tube.
- Doing so may cause fluid leak, overheating, or explosion.

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- The battery will harm the body if accidentally swallowed.
- If this happens, consult a physician immediately.

If the battery fluid leaks out, take the following procedures. Do not touch the battery with your bare hands.

- The battery fluid may cause blindness if it comes in contact with your eyes.
- Do not rub your eyes. Wash immediately with clean water and consult a physician.
- The battery fluid may cause inflammation or injury if it comes in contact with the skin or clothes.
- Wash it off thoroughly with clean water and consult a physician.

Intended use

► **Precautions**

- Do not clean this product by using paint thinner, benzene, alcohol, or other chemicals. Otherwise, it may lead to product fracture or color fading. Soft cloth slightly dipped in tap water or suds can be used to clean this product's shell.
- After use, this product shall be kept away from the pool, shower or other areas exposed to water or humidity.
- This product shall be kept away from area of high temperature or direct sunlight.

GB

D

F

I

NL

E

DK

P

N

S

FIN

PL

CZ

SK

H

RO

Parts identification

A Main unit

- ① Power switch
(Power off/on/mode selection)
- ② Brush shaft
- ③ Mode Indicator (Green)
- ④ Charge status indicator (Red)
- ⑤ Fixing position mark
- ⑥ Installation position mark

Accessories

- E Operating Instructions
- F Guarantee Card

B Extra fine brush

- ⑦ Brush bristle (extra-fine bristle)
- ⑧ Toothbrush handle

C Multi-fit brush

- ⑨ Brush bristle

D Charger (RE8-55)

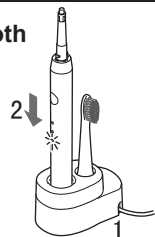
(The shape of the charger differs depending on the area.)

- ⑩ Cord
- ⑪ Power plug
- ⑫ Charging stand
- ⑬ Charger cover
- ⑭ Charge hole
- ⑮ Brush head socket


Charging the toothbrush

- 1 Place charger at a horizontal smooth location like a sink or washstand.**
- 2 Insert the main unit into the charge hole in an upright position.**
- 3 Insert the power plug into wall outlet.**
 - If the main unit is tilted, or the bottom of the main unit does not connect with charge hole, the main unit will not be charged.
 - Make sure that the charger cover is firmly connected to the charger.



Indicator in charging

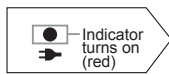Charging

Charging completed

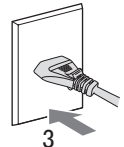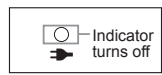If charging continuously

- It takes about 17 hours to be fully charged. When charging is completed, the charge status indicator will automatically turn off. Remaining battery capacity may shorten charging time.
- More than 17 hours of charging does not affect battery performance.
- A full charge will supply enough power for approx. 30 minutes use. (Use time may vary depending on ambient temperature and use conditions.)

4 Unplug the power plug when charging is completed.

► Notes

- If this unit is not used for more than 6 months, or is being used for the first time, the charge status indicator will not light up immediately, or the duration of use will be shorter. In this case, charge the appliance for more than 26 hours. Afterwards, the appliance will operate normally.
- Recommended temperature for charging is 0 °C-35 °C. When charging temperature is over or lower than the recommended temperature, the battery performance may be reduced.

Before cleaning your teeth

The Sonic Vibration Toothbrush is an electric toothbrush with sonic vibration* of approximately 31 000 brush strokes/minute.

* Sonic vibration in acoustic wave field

When you first begin using the toothbrush, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the toothbrush for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as pyorrhea. In this case, stop using the toothbrush and consult your dentist.

► About soft start function (For NORMAL mode only)

To prevent toothpaste from spattering and to get accustomed to the high speed vibration, the toothbrush will gradually reach the selected speed and vibration approximately 2 seconds after it starts operating.

GB
D
F
I
NL
E
DK
P
N
S
FIN
PL
CZ
SK
H
RO

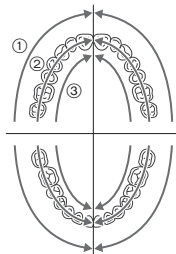
► About the 30-second brushing notification timer function

The toothbrush will pause once at 30-second intervals to indicate the brushing time. The main unit pauses twice at 2-minute intervals then starts again.

- Divide your teeth into 4 parts (upper and lower, left and right), and brush one part in order such as “① front of the teeth”, “② biting plane of the teeth”, and “③ back of the teeth” spending 30 seconds.

By brushing the above 4 parts for 30 seconds each (about 2 minutes), you can evenly brush your teeth without missing any spots in your mouth.

Using the toothbrush

- **The recommended operating temperature of electric toothbrush is 0 °C or higher. If the indoor temperature is below the recommended temperature, the electric toothbrush may stop operating.**

- Do not put the main unit in contact with hair products, hand creams, etc.
Doing so may damage the exterior.

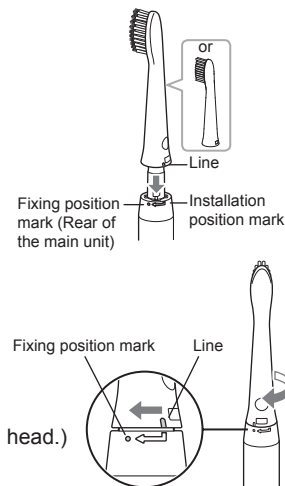
1 Firmly install the brush head onto the brush shaft.

- Insert the brush head into the brush shaft, with its line aligned to the installation mark of the main unit.

(If the line is not aligned, the brush head and main unit may malfunction.)

- Rotate the brush head to align its line to the fixing position mark of the brush shaft.

Do not rotate tightly. (Otherwise, it may cause breakage and deformation to the brush shaft and the brush head.)



2 Place toothbrush in your mouth.

- To avoid toothpaste foam from splashing, use a small amount of regular toothpaste. Before turning on the power switch, put the brush head in your mouth.

3 Press the power switch and choose a mode.

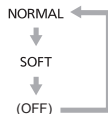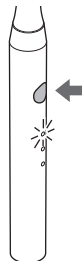- Some first-time users probably will have sense of itching and pricking because of such sound wave vibration, but after several times of use, such sense will disappear.
- Every time the power switch is turned on, the mode changes. (Each time you press the power switch, the appliance enters Normal mode.)

Mode	Description
NORMAL	Recommended for normal toothbrushing.
SOFT	When you feel Normal mode is too stimulating, selecting Soft mode is recommended.

4 Turn off the main unit when you have finished brushing.

- Press the power switch to turn off the unit.

Tips on usage

Brushing

Extra fine brush

When you want to brush mainly the periodontal pocket. (You can also brush the teeth front surface and the biting surface.)

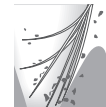
Multi-fit brush

When you want to make your teeth front surface and biting surface clean. (You cannot brush the periodontal pocket.)

- Choose the brush head depending on use and area.

Periodontal pocket

- Use the extra fine brush.
- Place the brush head at a 45° angle.
- Apply the bristles of the extra fine brush lightly into the periodontal pocket and move the toothbrush slowly.

Front surface and the biting surface of teeth

- Apply the brush head at a 90° angle.

For the back of the front teeth or misaligned teeth

- Place the brush vertically.
- Move toothbrush to match the angle of the teeth.

GB
D
F
I
NL
E
DK
P
N
S
FIN
PL
CZ
SK
H
RO

For the boundary between the teeth and gums

- Apply the brush head at a 45° angle.

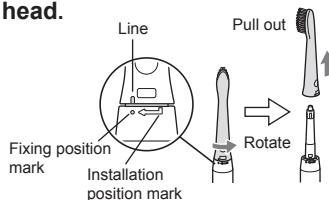
After use

Cleaning the brush head and the main unit


1 Remove the brush head.

- Rotate the brush head to align its line to the installation position mark of the main unit. (Pulling out the brush head while the brush head is locked will lead to breakage).
- Pull out the brush head vertically.


2 Clean the brush head. Wash brush head using running water.

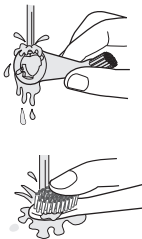- Thoroughly clean the inner and bottom of brush head.
- Do not pull the bristle on brush head. (Otherwise, it will lead to bristle falling off.)
- Do not wash it with water higher than 80 °C. (Otherwise, it may lead to brush head breakage.)

3 Clean the main unit. Wash the main unit using running water.

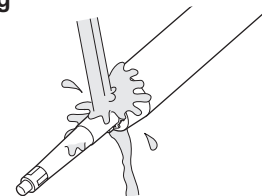- Do not immerse the main unit in water. (Otherwise, the water may seep into the main unit.)
- Do not wash it with water higher than 40 °C. (Otherwise, it may damage the main unit.)

4 Wipe off the water on brush head and main unit for storage.

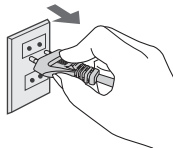
- Do not charge while it is wet. Doing so may cause soiling.


Cleaning the charger

- 1 **Unplug the power plug from the wall outlet.**

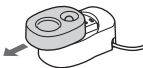
- 2 **Slide and remove the charger cover.**

- 3 **Wash the charger cover with running water before wiping remaining water off.**

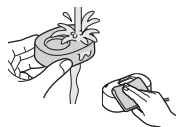- Except for charger cover, all other parts cannot be washed with water.

- 4 **Wipe off the dirt on the bottom of charger with soft cloth.**

- 5 **Wipe and clean the metal parts of the power plug with dry cloth about every 6 months.**

Frequently Asked Questions

Question	Answer
Why does it stop each 30 seconds and twice each 2 minutes?	▶ Because this product has 30-second reminder and timer, which allows toothbrush stopping at 30-second interval to remind you of brushing time. (See page 14.)
Why does each pressing of the power switch make a different operation mode?	▶ This product is in two switching modes. Each pressing will switch to another mode. (See page 15.)
Where are the replacement parts for brush head sold?	▶ Available at distributor's or an authorized service centre.

GB

D

F

I

NL

E

DK

P

N

S

FIN

PL

CZ

SK

H

RO

Troubleshooting

Problem	Diagnosis	Corrective Action
Charge status indicator does not turn on.	▶ The main unit is not properly installed on the charger.	▶ Vertically insert the main unit into charge hole, and check that the charge status indicator glows. (See page 13.)
The toothbrush can only be used for a few minutes, even if it is charged.	▶ The charger cover is not firmly connected to the charger.	▶ Connect the charger cover to the charger firmly.
The toothbrush does not operate.	▶ Rechargeable battery has reached the end of its service life (about 3 years).	▶ Replace the battery. (See page 19.) ▶ Contact the authorized service centre where purchased.
	▶ The appliance is used with approximately 0 °C or lower.	▶ The recommended operating temperature is 0 °C or higher.
Vibration weakens.	▶ Low battery.	▶ Charge for at least 17 hours. (See page 13.)
The toothbrush pauses operating with short intervals.	▶ The “30-second brushing notification timer function” is in operation. (See page 14.)	▶ This is not failure. Use it as usual.
Vibration will gradually increase after pressing power switch.	▶ The “Soft start function” is in operation. (See page 13.)	

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Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge. (The operating life of the battery may differ depending on usage or storage conditions.)

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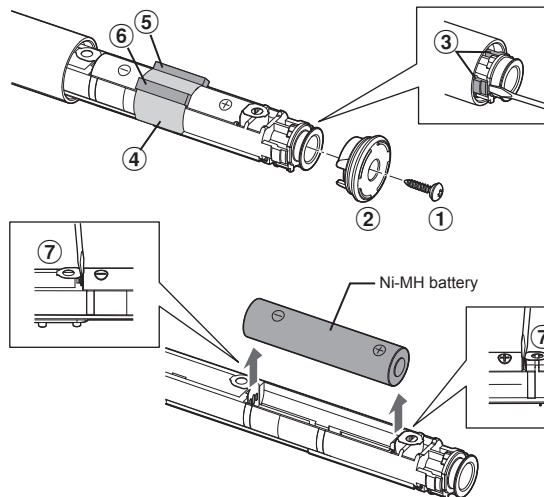
This product has a built-in rechargeable battery. Before disposing, be sure to remove the built-in battery.

If there is an officially-designated place for such battery disposal, be sure to dispose your old battery at such place.

The method as shown on below picture is only applicable to product disposal, and do not use such method to make maintenance. If you dismantle the toothbrush yourself, it will no longer be waterproof, which may cause it to malfunction.

- Press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Disassemble the toothbrush after discharging the battery completely.
- Battery must be disposed at a place as specified by the government or a recognized recycling corporation.
- Remove the battery with a flat head screw driver in the following order.
Disassemble the bottom part of the main body (①-②).
Press and hold the main body. Lever the battery clip (③).
Remove the tape from the terminals (④).
Remove the cushions (⑤-⑥).
Insert a flat head screw driver into the bottom of the welding joint, and then remove the battery (⑦).
- When disassembling, take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.

Environmental protection and resource recovery

This product employs Ni-MH battery.
If your country has relevant regulations on battery recycling, have your old battery disposed at such specified place.

Specifications

Power source	See the name plate on the product (Automatic voltage conversion)
Charging time	Approx. 17 hours (Room temperature: 0 °C- 35 °C)
Airborne Acoustical Noise	In Normal mode: 49 (dB (A) re 1 pW) In Soft mode: 44 (dB (A) re 1 pW)

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Changing the brush head (sold separately)

- Replacing the brush head each 3 months is recommended.
- Tilted or worn bristle will not produce ideal effect.
- Damaged brush head as a result of masticating will harm your teeth and gum. Replace the damaged one.
- Replacement parts are available at distributor's or an authorized service centre.
- Color identification labels are included in the brush-head replacement package. The labels can be attached on the toothbrush handle to identify users when they share the same main unit.

Part number	Replacement brush for EW-DM81	
WEW0972	Extra fine brush	
WEW0974	Multi-fit brush	

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For more information about collection and recycling, please contact your local municipality. Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

Note for the battery symbol (bottom symbol)
 This symbol might be used in combination with a chemical symbol. In this case it complies with the requirement set by the Directive for the chemical involved.

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