Oral-B®

PROFESSIONAL CARE








powered by BRAUN

500 550

CUPS EACH TOOTH
FOR A
3D
CLEAN

Type 3757

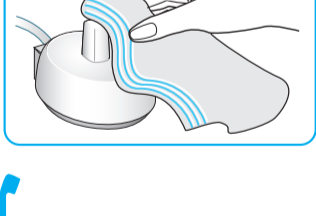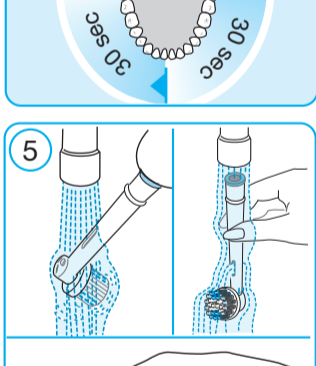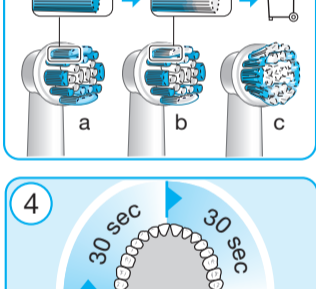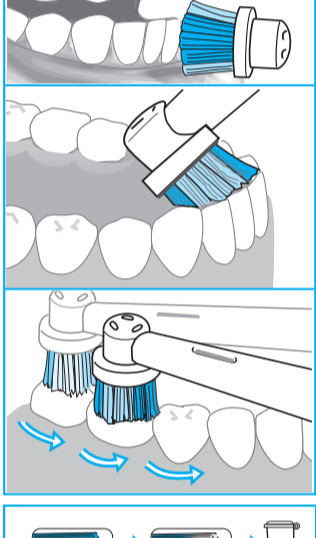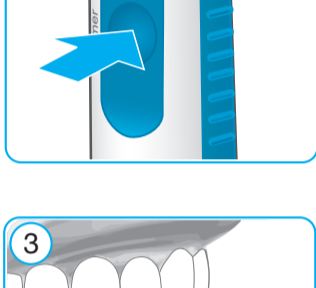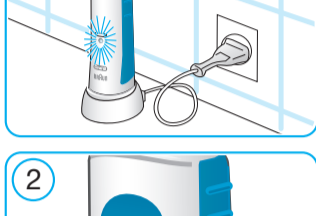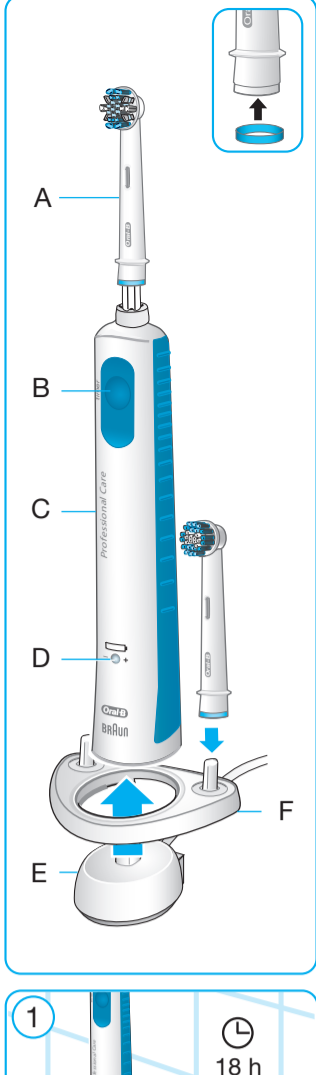- (DE) (AT) 00 800 27 28 64 63
- (CH) 08 44 - 88 40 10
- (GB) 0 800 783 7010
- (IE) 1 800 509 448
- (FR) 0 800 944 803
(appel gratuit depuis un poste fixe)
- (BE) 0 800 14 592
(9-17u/h) gratis/gratuit
- (LU) 4 00 50 51
- (ES) 901 11 61 84
- (PT) 808 20 00 33
- (IT) 800 124 600
- (NL) 0800-4 45 53 88
(9-17u) gratis
- (DK) 70 15 00 13
- (NO) 22 63 00 93
- (SE) 020 - 21 33 21
- (FI) 020 377 877
- (GR) 800 801 3457

www.oralb.com
www.braun.com
www.service.braun.com

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CH: 0844 - 88 40 10

English

Before operating this toothbrush, please read these instructions and save this manual for future reference.

IMPORTANT

- Periodically check the cord for damage. If the cord is damaged, take the charging unit to an Oral-B Braun Service Centre. A damaged or non-functioning unit must no longer be used.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.
- Children should be supervised to ensure that they do not play with the appliance.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service Centre. Do not disassemble the product except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not insert any object into any opening of the appliance/charging unit.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

A Oral-B FlossAction brush head*

B On/off button

C Handle

D Charge indicator

E Basic charger

F Brush head holder*

* depending on model

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the basic charger (E) or the assembled charging unit (E/F) into an electrical outlet and place the toothbrush handle (C) on there. For voltage specifications refer to the bottom of the charging unit.
- The green charge indicator light (D) flashes while the handle is being charged. When fully charged, the green charge indicator light will continuously shine for 5 seconds and turn off. A full charge takes at least 18 hours and enables approximately one week of regular brushing (twice a day, 2 minutes) (1).
Note: If the green light is not flashing immediately, continue charging. It will appear in 10–15 minutes.
- When the rechargeable battery is low, the red charge indicator light (D) starts flashing after turning off the toothbrush.
- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit; there is no risk of overcharging. However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Using the toothbrush

Wet brush head and apply toothpaste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on (2). Guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface. Start brushing the outsides, then the insides, finally the chewing surfaces (3). Do not press too hard or scrub.

Professional Timer

A short stuttering sound at 30-second intervals reminds you to brush all four quadrants of your mouth equally (4). A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. The elapsed brushing time is memorised when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the timer resets.

Pressure Sensor

Your toothbrush has a pressure control feature installed. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition you will hear a different sound while brushing.


Brush heads

Oral-B Precision Clean brush head

cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

Oral-B FlossAction brush head

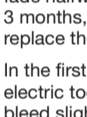features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.

Both brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, signalling it is time to replace the brush head (3a–c).


In the first few days of using any electric toothbrush, your gums may bleed slightly. In general bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

Cleaning recommendations

After brushing, rinse the brush head under running water with the handle switched on. Switch off and remove the brush head. Clean both parts separately under running water and wipe them dry. From time to time also clean the charging unit with a damp cloth (5). The brush head holder (F) is dishwasher-safe.

Subject to change without notice.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

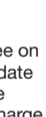Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

30-Days Money Back Guarantee

Try Oral-B Professional Care starting from the day of purchase. If you are not satisfied with this product, please return your handle, refills and charging unit in their original box with your original cash register receipt to the address below within 30 days of purchase.

Great Britain/Ireland:
«Oral-B 30-Day Challenge»
Braun Cashback (BC 024)
PO Box 504
Leicester LE94 0AE

This does not affect your statutory rights. Terms and conditions apply, for further details, please contact our customer service.

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