
déCOLITE

 the beauty specialists

ENGLISH	4
FRANÇAIS	12
DEUTSCH	20
NEDERLANDS	28
ESPAÑOL	36
PORTUGUÊS	44
ITALIANO	52

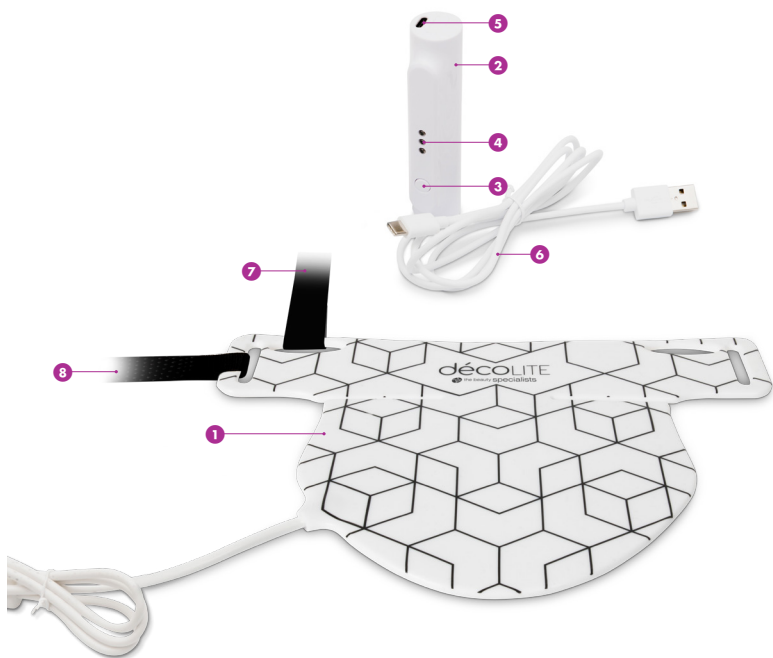WARNINGS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance. If any part of the mains adapter or cord are damaged, discontinue use. Contact Customer Care to arrange a replacement.
- Do not immerse product in water. Do NOT use solvents or abrasive cleaners.
- Do NOT use this product if it is not working properly, or if it has been dropped or damaged.
- DO NOT use décoLITE if you are taking any medication that can cause photosensitivity. If you use décoLITE and you are taking any medication that can cause photosensitivity you may cause a severe skin reaction.
- Photosensitivity is a common side effect of various medications. These can include certain antibiotics, chemotherapy drugs, and diuretics. If you are unsure about any medication you may be taking consult your healthcare provider.
- DO NOT use décoLITE to treat any other conditions apart from those listed in the indications for use.
- DO NOT use décoLITE if you are pregnant, planning on becoming pregnant or breastfeeding. décoLITE has NOT been tested on pregnant or breast-feeding women and therefore the risk to the foetus, neonate or pregnant women is unknown.

- DO NOT use décoLITE LED system if you suffer from Lupus erythematosus, photosensitive eczema or Albinism. If you use the décoLITE LED system to treat Lupus erythematosus, photosensitive eczema or Albinism you may cause a severe skin reaction.
- DO NOT use décoLITE if you suffer from any photosensitive disorder (sensitization to light).
- Before using décoLITE you must carry out a sensitization test.
- DO NOT use décoLITE in any way that is not exactly as specified in these instructions for use.
- TURN OFF and REMOVE décoLITE immediately if it feels hot or your skin feels uncomfortably warm.
- Only use décoLITE with the proprietary Hydrogel pads recommended for or supplied for use with décoLITE. Other pads are available, but have not been tested with the décoLITE system.
- DO NOT use décoLITE without first removing creams, lotions or make-up. Some ingredients in these products may cause light sensitization or may reduce the effect of the treatment.
- NEVER use décoLITE near water or where water can be spilt onto the product. DO NOT use décoLITE if you think it has come into contact with water.
- Keep décoLITE away from HOT surfaces.
- DO NOT attempt to modify décoLITE in any way. Modifying décoLITE will void your warranty. décoLITE contains no serviceable parts.
- DO NOT use décoLITE if the light pad, controller, charging cable are damaged in anyway.
- The décoLITE strap contains a mixture of polyester, nylon and lycra. If you are allergic to any of these materials, you may not be able to use the head strap provided.
- DO NOT use décoLITE if you suffer from light induced headaches.

KNOW YOUR UNIT

- | | | | |
|---|----------------------|---|--------------------|
| 1 | Light pad | 5 | USB-C socket |
| 2 | Controller | 6 | USB charging cable |
| 3 | On/Off | 7 | Long strap |
| 4 | LED indicator lights | 8 | Short strap |

PATCH TEST

Before carrying out your first treatment you must carry out a sensitivity test (patch test).

1. Place décoLITE over your inner arm.
2. If at any time your skin feels HOT or uncomfortable STOP using the device immediately.
3. After the treatment wait at least 6 hours. If the skin on your inner arm appears red and/or itchy then décoLITE is unsuitable for you.

CHARGING


To charge décoLITE insert the USB-C end of charging cable into the controller. Insert the USB connector into a USB power supply. The 3 LEDs on the control unit will start to flash to indicate charging. When fully charged all three LEDs will be illuminated.

Note: Always recharge when the last LED starts to flash. décoLITE will not operate whilst charging.


GETTING STARTED

Before each use, with the device switched off, clean the décoLITE light pad. See cleaning instructions. décoLITE comes with an adjustable strap that allows you to comfortably wear the light pad.

Use around the neck and on the chest:

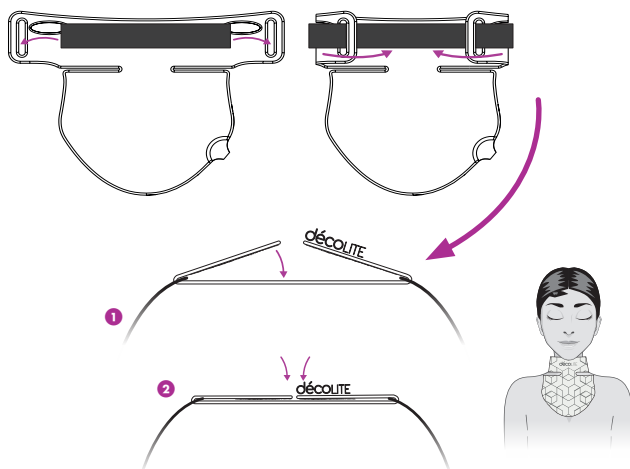Use on the chest only:

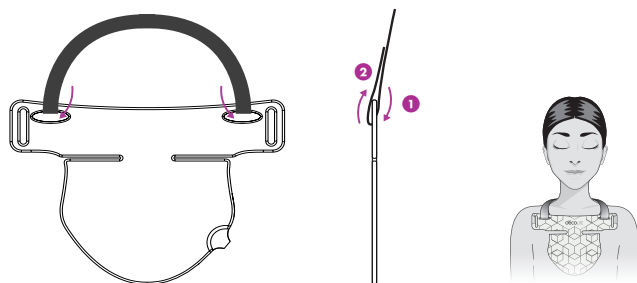TREATMENTS

For best results treat for 3-5 days a week for a minimum of 4 weeks. Maintenance treatments should be carried out as required, according to lifestyle and skin condition.

décoLITE Treatment:

1. Cleanse and exfoliate the skin to remove any traces of skincare.
2. To fit the light pad, thread the strap through the loop holes on either side of the light pad and join together using the hook and loop ends (see illustration above). Use the short strap to place around the neck or the long strap (inserting through the upper loop holes) to place over the chest.
3. Press the On/Off button to turn on. The device will automatically switch off after 10 minutes. You may stop the treatment at any time by pressing the On/Off button.
4. At the end of the treatment remove the light pad.
5. Apply your normal moisturiser and sunscreen.

Note: Immediately after treatment the skin may appear red. It should return to normal within 1-2 hours.

Hydrogel treatment (optional):

décoLITE can be combined with Hydrogel treatments for enhanced results and benefits. Please review any instructions supplied with the Hydrogel pack.

CLEANING

Before each use, with the device switched off, clean the light pad as follows:

1. Ensure that the device is NOT being charged when cleaning.
2. Clean the light pad with a water based wipe or damp clean cloth. Stubborn stains can be removed by wiping gently with 70% alcohol. DO NOT use any other cleansing products as these may damage the plastics.
3. Clean the controller with a damp clean cloth.

STORAGE

Store décoLITE at room temperature in dry conditions. Avoid storing décoLITE in dusty conditions, in direct sunlight or in damp and humid conditions.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No red or infra-red light	Not charged	Charge the controller
No red or infra-red light	Device is being charged.	Unplug charging cable from the controller

DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

SPECIFICATION

Controller Input: 5V 500mA 

Light Pad Input: 5V 280mA (CC)

-  This product complies with current EU product legislation.
-  Read these instructions carefully before use and then retain for future use.
-  This product is for indoor use only.
-  This symbol means keep the device away from direct sunlight.
-  Direct current
-  This product is NOT suitable for charging in the bathroom.
-  This symbol indicates the USB connector.
-  This symbol means that the device uses a Lithium ion Polymer battery.
-  This symbol refers to the catalogue number for the décoLITE device.
-  This symbol refers to the batch code of the décoLITE device.

Weight:

600g (controller and light pad)

LEDs:

Red: 630nm+/-10nm.

Near Infra-red: 830nm+/-17.5nm

Total Device Dimensions:

Light Pad: 300mm max (w) x 220mm max (h)

Controller: 23mm max (dia) x 104mm (l)

Programme / treatment time:

10 minutes

Operating Temperature:

0°C to 40°C

Storage and Transport Temperature:

-5°C to 55°C

Operating Humidity:

30% to 95% without condensation

Storage and Transport Humidity:

20% to 85% without condensation

CUSTOMER CARE

If you have any questions, please contact us at: online@riobeauty.com

To reorder additional accessories, please visit

www.riobeauty.com or contact your local distributor.

GUARANTEE

To register your two year guarantee online, visit:

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© The Dezac Group Ltd 2019 Manufactured by The Dezac Group Ltd PO Box 17,