PORE PERFECTION

- WIRELESS -

the beauty specialists

ENGLISH	4
FRANÇAIS	18
DEUTSCH	32
NEDERLANDS	46
ESPAÑOL	60
PORTUGUÊS	74
ITALIANO	88

WARNINGS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Keep out of the reach of children.
- Do not charge for longer than 24 hours.
- Do NOT use this appliance if it is not working properly, or if it has been dropped or damaged.
- If you suffer from a skin condition, consult your dermatologist before use.
- For hygiene reasons we recommend for single person use only.
- Do not use on cuts or broken skin.

- Do not use this appliance on loose skinned parts of your body or any areas showing signs of: varicose veins, skin diseases (for example eczema, psoriasis, dermatitis), cuts and abrasions, irritation, warts, moles, skin tags, skin rash, hypersensitivity or sunburn, phlebitis, recent scar tissue, other skin conditions. We recommend you consult your physician before treatment if you suffer from any of the above conditions.
- Do not use immediately around the eyes.
- Do not use on areas that have had cosmetic surgery.
- Do not place or store in direct sunlight. Store in a cool dry place.
- To avoid redness / bruising keep the head moving at all times and when using for the first few times, use the small Round treatment head.

KNOW YOUR UNIT

- 1 Protective pouch
- 2 Pore Perfection
- 3 Round treatment head
- 4 Oval treatment head
- 5 USB charging cable
- 6 Mode button
- 7 ON/OFF button
- 8 Charging indicator
- 9 Rubber seal
- 10 Micro USB Port

HOW IT WORKS

Pore Perfection is a non-invasive vacuum suction technology skin care system which deep cleans pores to remove dirt and blackheads, while helping to improve overall skin condition. Used with your phone the 10x magnifying camera helps identify problem areas to target and shows the actual extraction and cleansing of congested pores as it occurs.

The suction actions also helps boost micro blood circulation while the blue light targets surface bacteria which cause redness and inflammation.

TREATMENT HEADS

Oval Treatment Head

High suction, full face cleansing.


Round Treatment Head

Low suction and for difficult to reach zones (i.e. around nose area).

FITTING THE TREATMENT HEADS

1. Align the small arrow symbols and press firmly onto the unit.

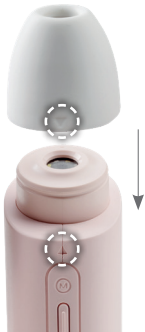HOW TO CHARGE

IMPORTANT: Before using for the first time, charge for 3.5 hours.

Charging is required if the charging indicator flashes when turning on.

1. Peel back the rubber seal on the base of the unit.
2. Plug the micro USB connector of the charging cable into the micro USB port on the unit. Then, plug the USB charging cable in to either a computer or USB (5V) power outlet.
3. The charging indicator will flash blue when charging. When fully charged, the charging indicator will remain blue.

SETTING UP THE CAMERA AND APP

1. Scan the QR code below to Install the **Pore Perfector Wireless** app or search the app store for '**Pore Perfector Wireless**'.

2. Open the **Pore Perfector** app on your phone/tablet and follow the on-screen instructions.
3. Press the ON/OFF button for one second to switch on.
4. Select the **Pore Perfector** network in your wi-fi settings to pair with the app.

NOTE: Continue to follow on-screen instructions.

WARNING:


Bruising can result if the treatment head isn't constantly moved over the skin or if too high an intensity level is used.

SETTING UP THE CAMERA AND APP

5. Place the treatment head on a patch of skin to check the camera displays on the app / phone.
6. Once the treatment is in progress, if you wish, you can use the app to take pictures of the treatment.

Pore Perfector battery health.
Touch this icon to temporarily
hide all icons. Touch
again to show icons.

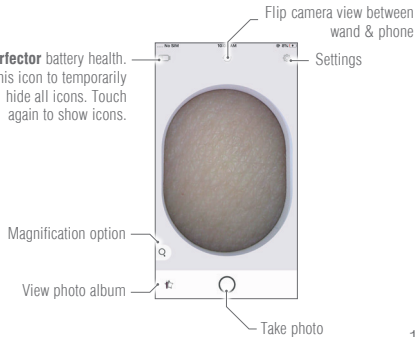PATCH TEST

IMPORTANT: Test each head on a small area on each intensity level to check for skin sensitivity. Wait for 24 hours to check for any reaction.

WARNING: Bruising can result if the treatment head isn't constantly moved over the skin or if too high an intensity level is used.

PREPARING THE SKIN

- Cleanse to remove any makeup or dirt.
- **Hold a warm cloth over the face to open the pores.**

TREATMENT

1. Select your preferred treatment head (**IMPORTANT** – always use the 'round' treatment head when using for the first time).
2. To switch on press the ON/OFF button for one second - this will automatically select suction level 1.
3. Briefly press the ON/OFF button again to select levels 2 or 3 and increase suction level.
 - Level 1 - Low suction for sensitive skin
 - Level 2 - Medium suction
 - Level 3 - Maximum suction
4. Pressing the ON/OFF button a 4th time will turn off the suction.
This is useful as the camera remains on so you can view the area of skin to be treated. Pressing again will reactivate the suction on level 1.
5. Pressing the Mode button turns the blue light on and off.

TREATMENT

NOTE: When the Mode button is pressed once, the suction treatment will stop, press the Mode button once more to continue the suction treatment at the same suction level.

6. To switch off press the ON/OFF button for 2 seconds and then release.

NOTE: During the treatment, the **Pore Perfector** will pulse on and off three times before continuing to proceed with the treatment at each 5 minute interval.

Once familiar with the functioning of the unit and having performed a patch test you are ready to treat your skin.

Always use level 1 for the first few treatments.

With the unit turned on apply the treatment head to the skin.

It is important to keep the head moving over the skin at all times to avoid bruising the skin.

Do not use over the same patch of skin for more than one second at a time.

TIP: Use the camera to guide the treatment head over the skin and assist in extracting dirt, grime or sebum plugs from the pores.

AFTER TREATMENT

- Wash skin with a mild antibacterial cleanser.
- The skin may appear red and slightly irritated after use. This should subside after 30 minutes.
- Do not use any exfoliating products before or after the treatment.
- Remove the silicon tip and thoroughly wash it with antibacterial soap or facial cleanser immediately after use.

CLEANING & MAINTENANCE

- Switch the device off and clean using a soft damp cloth.
- Do not immerse the device in water or other liquids.
- Dry thoroughly before storing.

IMPORTANT INFORMATION

There are no user serviceable parts inside the unit.

DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat

any hazardous substances and avoid the increasing landfill.

This product is classed as Electrical or Electronic Equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

CUSTOMER CARE

If you have any questions, email us at: online@riobeauty.com

To reorder additional accessories, please visit www.riobeauty.com or contact your local distributor.

GUARANTEE

To register your two year guarantee online, visit:

www.riobeauty.com/warranty_online.htm

SPECIFICATION

MODEL: DRMC2

Input: 5.0V  600mA 

Battery: Li-ion 3.7V 2000mAh

This product is for indoor use only.

Read these instructions carefully before use and then retain or future use.

This unit complies with current EU product legislation.

**UK
CA**

Direct current

Voltage: 5.0V Max

Suction: -60Kpa Max

Noise: 60dB Max

Camera Pixel: 0.3million

Camera Resolution: 640x480

Camera Zoom: 10x

White Light: 4pcs 0603 LED

Blue Light: 4pcs 0603 LED

Working Time:

Wi-fi mode On - 250 minutes

Wi-fi mode Off - 280 minutes

USB: Micro-B, 5V / 600mA input

Battery: 18650 2000mAh 3.7V with PTC

Operating Temperature / Humidity: 5~35°C / 10~90%RH

Storage Temperature / Humidity: -20~45°C / 10~90%RH

SPECIFICATION

Radio frequency output: Power: 25mW
Operating frequency band: 2.4G

Hereby, The Dezac Group declares that the radio equipment type [designation of type of radio equipment] is in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at:
<https://www.riobeauty.co.uk/facial-beauty/rejuvenation/facial-rejuvenation/pore-perfector-wireless.html>

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© The Dezac Group Ltd 2021