FACE **GLO**

1072

the beauty specialists

ENGLISH	4	ESPAÑOL	28
FRANÇAIS	10	PORTUGUÊS	34
DEUTSCH	16	ITALIANO	40
NEDERLANDS	22		

WARNINGS & PRECAUTIONS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- This appliance is only to be used with the mains adapter supplied SK01T-0500100Z. Detachable PSU IEC 60417-6181.
- If the power lead or AC/DC adapter are damaged, it must be replaced by the manufacturer or its service agent or similarly qualified person in order to avoid a hazard. Contact Customer Care to arrange a replacement.
- Do not immerse product in water. Do NOT use solvents or abrasive cleaners.
- Do NOT use this product if it is not working properly, or if it has been dropped or damaged.
- DO NOT use FaceGLO if you are taking any medication that can cause photosensitivity. If you use FaceGLO and you are taking any medication that can cause photosensitivity you may cause a severe skin reaction.
- Photosensitivity is a common side effect of various medications. These can include certain antibiotics, chemotherapy drugs, and diuretics. If you are unsure about any medication you may be taking consult your healthcare provider.
- DO NOT use FaceGLO to treat any other conditions apart from those listed in the indications for use.
- DO NOT use FaceGLO if you are pregnant, planning on becoming pregnant or breastfeeding. FaceGLO has NOT been tested on pregnant or breast-feeding women and therefore the risk to the foetus, neonate or pregnant women is unknown.
- DO NOT use FaceGLO LED system if you suffer from Lupus erythematosus, photosensitive eczema or Albinism. If you use the FaceGLO LED system to treat Lupus erythematosus, photosensitive eczema or Albinism you may cause a severe skin reaction.
- DO NOT use FaceGLO if you suffer from any photosensitive disorder (sensitization to light).
- Before using FaceGLO you must carry out a sensitization test.
- DO NOT use FaceGLO in any way that is not exactly as specified in these instructions for use.
- TURN OFF and REMOVE FaceGLO immediately if it feels hot or your skin feels uncomfortably warm.
- Only use FaceGLO with the proprietary Hydrogel masks recommended for or supplied for use with FaceGLO. Other masks are available, but have not been tested with the FaceGLO system.
- DO NOT use FaceGLO without first removing creams, lotions or make-up. Some ingredients in these products may cause light sensitization or may reduce the effect of the treatment.

- NEVER use FaceGLO near water or where water can be spilt onto the product. DO NOT use FaceGLO if you think it has come into contact with water.
- Keep FaceGLO away from HOT surfaces.
- DO NOT attempt to modify FaceGLO in any way. Modifying FaceGLO will void your warranty. FaceGLO contains no serviceable parts.
- DO NOT use FaceGLO if the mask, controller, charging cable are damaged in anyway.
- The FaceGLO strap contains a mixture of polyester, nylon and lycra. If you are allergic to any of these materials, you may not be able to use the head strap provided.
- DO NOT use the FaceGLO LED mask if you suffer from light induced headaches.
- Do not insert a plug front into a socket without the adapter connected.
- DO NOT share FaceGLO with any other person. FaceGLO is for single person use only
- DO NOT use the device to overtreat - only use per the recommend treatment protocols.

Photosensitivity is a common side effect of various medications. These can include certain antibiotics, chemotherapy drugs, and diuretics. If you are unsure about any medication, you may be taking consult your healthcare provider.

Other substances not listed above can also cause photosensitivity.

Common examples of these substances are:

St John's Wort, coal tar, deodorants, antibacterial soaps, artificial sweeteners, naphthalene (mothballs), petroleum products, brightening agents found in laundry detergent, and cadmium sulphide (a chemical injected into the skin during tattooing).

INTRODUCTION

FaceGlo is suitable for adults wishing to treat full facial wrinkles.

As we age skin becomes thinner and drier, it can develop sags and wrinkles due to the loss of collagen and elastic tissue. Aging skin is also more susceptible to the development of brown spots and pigmentation. The skin's dermis can lose from 20 %-80 % of its thickness during the aging process.

FaceGLO combines the clinically proven and trusted wavelengths in the red and near infrared region of the light spectrum (600nm to 1100 nm) that have shown in clinical studies to reduce facial wrinkles through a process called photobiomodulation.

FaceGLO has been designed to deliver light directly to the skin, ensuring that the maximum amount of light hits its cellular targets and is not wasted.

FaceGLO is an over-the-counter at home adult use device intended to emit energy in the red and near infrared region of the light spectrum for the treatment of full facial wrinkles.

KNOW YOUR UNIT

- | | | |
|---------------------|---------------------------|------------------------|
| ❶ Head strap - Hook | ❺ Adaptor with connectors | ❽ Controller |
| ❷ Head strap - Loop | ❻ Carry bag | ❾ On/Off |
| ❸ Face mask | ❼ USB-C socket | ❿ LED indicator lights |
| ❹ Eye shields | | ⓫ USB socket |

PATCH TEST

Before carrying out your first treatment you must carry out a sensitivity test (patch test).

1. Cleanse the skin on your inner forearm to remove any skincare such as sunscreens.
2. Place the FaceGLO over your inner forearm with the LEDs toward your skin.
3. Switch ON the controller and allow the FaceGLO to run for its full 10-minute treatment program.
4. If at any time your skin feels HOT or uncomfortable STOP using the device immediately.
5. After the treatment wait at least 6 hours. If the skin on your inner arm appears red and/or itchy then FaceGLO is unsuitable for you.

EYE SHIELDS

While the light used in your FaceGLO is completely safe, you may find the lights quite bright. If you do find the LED lights a little bright, the optional eye shields can be used to for enhanced levels of comfort. The eye shields are easily fitted directly into the eye holes.

* Eyewear is provided for comfort only and clips easily into the eye holes of the face mask. The use of eyewear is optional as this face mask complies to all relevant eye safety standards (IEC 62471:2006).

CHARGING

To charge FaceGLO insert the power adaptor into the USB-C socket of the controller. Plug the AC/DC adapter into a mains power supply and turn it on. The 3 LEDs on the control unit will start to flash to indicate charging. When fully charged all three LEDs will be illuminated.

Note: Always recharge when the last LED starts to flash. FaceGLO will not operate whilst charging.

GETTING STARTED


Before each use, with the device switched off, clean the FaceGLO mask. See cleaning instructions. FaceGLO comes with adjustable straps that allows you to comfortably wear the mask.

1. To fit the mask, thread the strap through the holes on either side of the mask and join together using the hook and loop ends.

2. Loop the straps around your head and secure firmly using the hook and loop ends.

3. Adjust the forehead and chin sections to achieve a close fit. For best results ensure the face mask is fitted as close to the forehead and chin as possible.

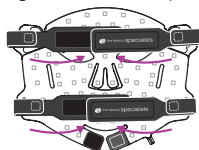TREATMENTS

For best results treat for 3-5 days a week for a minimum of 6 weeks. Maintenance treatments should be carried out as required, according to lifestyle and skin condition.

FaceGLO Treatment:

1. Cleanse and exfoliate the skin to remove any traces of skincare.
2. Place FaceGLO device over the head and onto the face, adjust the chin and forehead sections of the face mask for a close fit.
3. Press the On/Off button to turn on. The device will automatically switch off after ten minutes. You may stop the treatment at any time by pressing the On/Off button.
4. At the end of the treatment remove the mask.
5. Apply your normal moisturiser and sunscreen.

Note: Immediately after treatment the skin may appear red. It should return to normal within 1-2 hours.

After Treatment:

- Immediately after treatment your skin may appear red. The skin should return to normal within 1-2 hours.
- Some users have reported skin flaking or peeling, associated pruritis (itching), and dryness. This is transient and clears within a few days.
- Temporary headaches may occur immediately following the light treatment.

Although this device emits levels of visible and near infrared light that are considered safe for the skin, if you have an adverse reaction after using the device, you should stop using the device immediately.

Hydrogel Mask & Gel/Serum Treatments (optional):

FaceGLO can be combined with Hydrogel mask and gel/serum treatments for enhanced results and benefits. Please review any instructions supplied with the mask or treatment gel/serum.

CLEANING

Before each use, with the device switched off, clean the face mask following the cleaning instructions below.

1. Ensure that the device is NOT being charged when cleaning.
2. Clean the face mask with a water based wipe or damp clean cloth. Stubborn stains can be removed by wiping gently with 70% alcohol. DO NOT use any other cleansing products as these may damage the plastics.
3. Clean the controller with a damp clean cloth.

STORAGE

Store FaceGLO at room temperature in dry conditions. Use the carry bag provided. Avoid storing FaceGLO in dusty conditions, in direct sunlight or in damp and humid conditions.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No red or infra-red light from FaceGLO	Not charged	Plug the controller into the power adaptor making sure that the connectors are attached and leave charging until the three indicator LED's on the controller remain lit and do NOT flash.
No red or infra-red light from FaceGLO.	Device connected to the mains.	Unplug USB-C connector from controller.
LEDs do NOT turn OFF after 10 minutes.	Device is faulty.	Return to your supplier.
LEDs flicker during treatment.	Device is faulty.	Return to your supplier.

DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:

1. Check your local municipal guidelines
2. Separate the components and dispose of them correctly

SPECIFICATION

Controller Input: 5V 800mA 
Mask Input: 5V 450mA (CC)

This product complies with current EU product legislation.

This symbol indicates the USB connector.

Read these instructions carefully before use and then retain for future use.

This symbol means that the device uses a Lithium ion Polymer battery.

This product is for indoor use only.

Class III symbol (IEC 60417-5180)

This symbol means keep the device away from direct sunlight.

Class II insulation

Direct current

This unit complies with current UKCA product legislation.

Detachable PSU IEC 60417-6181

Recycle Paper

Model number: FCLT3

Program time: 10 minutes $\pm$
11 seconds

Wavelengths: • 630nm $\pm$ 10nm
• 830nm $\pm$ 10nm
• 1072nm $\pm$ 20nm

Maximum output intensity:
30mW/cm²

Operating temperature:
0°C to 30°C (32°F to 86°F)

Storage and transport temperature:
-5°C to 55°C (23°F to 131°F)

Operating humidity:
30% to 95% without condensation

Storage and transport humidity:
20% to 85% without condensation

Electrical input from power adaptor: 100-240Vac, 50/60Hz. 0.25A

Electrical output: 5V, 1A
IP22

CUSTOMER CARE

If you have any questions, please contact us at: online@riobeauty.com

To reorder additional accessories, please visit
www.riobeauty.com or contact your local distributor.

GUARANTEE

To register your two year guarantee online, visit:
www.riobeauty.com/warranty

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©The Dezac Group Ltd 2025, PO Box 17, Cheltenham Spa, England GL53 7ET
RP Regno Unito: The Dezac Group Ltd, PO Box 17,
RP UE: The Dezac Group (Europe) B.V. Azura,
UK Product RP: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, England, GL53 7ET
EU Product RP: The Dezac Group (Europe) B.V. Azura, Saturnusstraat 46-62, 2132 HB Hoofddorp, Netherlands