


EYE REFRESHING MASSAGER

 the beauty specialists

ENGLISH	4
FRANÇAIS	10
DEUTSCH	16
NEDERLANDS	22
ESPAÑOL	28
PORTUGUÊS	34
ITALIANO	40

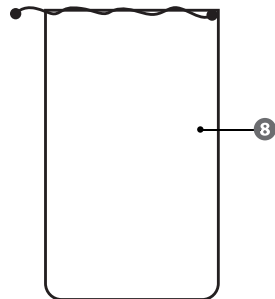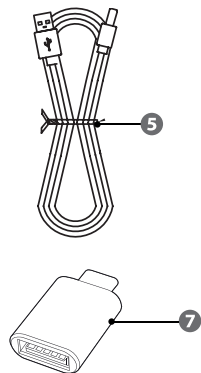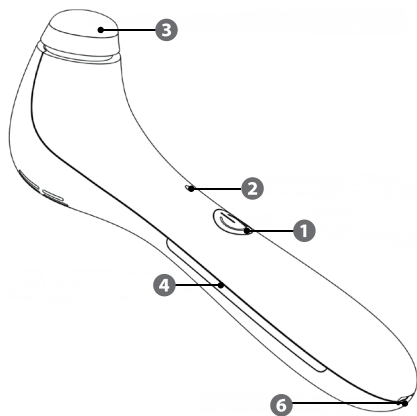WARNINGS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 12 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Do NOT use this appliance if it is not working properly, or if it has been dropped or damaged.
- This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Keep out of reach of children and pets.
- Do not use in the bath or shower.
- For external use only.
- Do not place or store in direct sunlight. Store in a cool dry place.
- If you suffer from a skin condition, consult your dermatologist before use.
- For hygiene reasons we recommend for single person use only.
- Do not use on cuts or broken skin.
- Do not use this appliance on loose skinned parts of your body or any areas showing signs of: varicose veins, skin diseases (for example eczema, psoriasis, dermatitis), cuts and abrasions, irritation, warts, moles, skin tags, skin rash, hypersensitivity or sunburn, phlebitis, recent scar tissue, other skin conditions. We recommend you consult your physician before treatment if you suffer from any of the above conditions.
- Do not use on areas that have had cosmetic surgery.

KNOW YOUR UNIT

(see page 3)

- | | |
|-------------------------|--------------------------|
| ❶ On/Off/Mode | ❸ USB-C 'In' Port |
| ❷ Indicator Light | ❹ USB-A to USB-C Adapter |
| ❺ Treatment Head | ❺ Travel Pouch |
| ❻ Silver Contact Strips | |
| ❼ USB Charging Cable | |


INTRODUCTION

Your Warming and Cooling Eye Refreshing Massager is a non-invasive at-home facial massage tool intended for cosmetic use to target puffiness, fatigue, and tension, with a synergy of heat, cooling therapy, gentle vibrations, and galvanic microcurrents.

CHARGING

Insert the USB Charging Cable into the USB C 'In' Port and connect to a USB port.

The Indicator Light will flash blue to indicate charging.

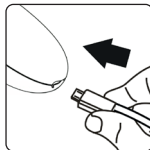When fully charged the Indicator Light will stop flashing and stay blue.

Battery Charging Time: 3 hours

Maximum Run Time: 40 minutes

NOTE: When the battery level is low the Indicator Light will flash red.

IMPORTANT: Disconnect from USB cable before using.

GETTING STARTED

Ensure the skin is clean.

Press and hold the On/Off/Mode button for 2 seconds to turn on. The Indicator Light will turn Amber, for Treatment Mode 1. Press the On/Off/Mode button again for Mode 2 (The Indicator Light will turn Red) and a third time for Mode 3 (The Indicator Light will turn Blue).

To turn off, press and hold the On/Off button for 2 seconds.

Mode 1 - Amber Light

Negatively charged Galvanic microcurrent with + 42°C warming compress.


Mode 2 - Red Light

Negatively charged Galvanic microcurrent with + 42°C warming compress and vibration massage.

Mode 3 - Blue Light

Cold compress (~ 10 °C) for 3 minutes. The cooling fan will continue to work for another minute after the treatment has finished.

Note - Each Mode will automatically stop after 4 minutes.

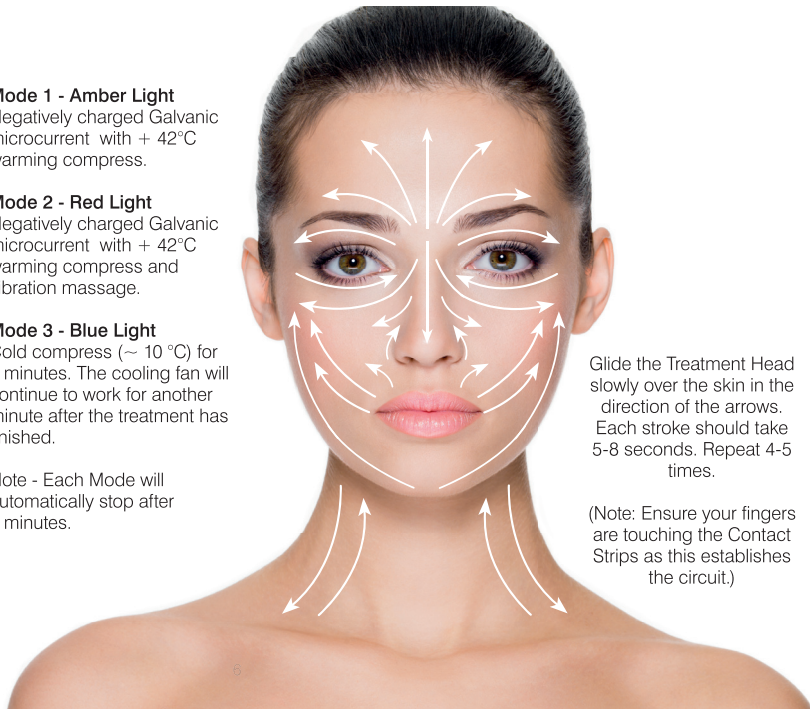Glide the Treatment Head slowly over the skin in the direction of the arrows. Each stroke should take 5-8 seconds. Repeat 4-5 times.

(Note: Ensure your fingers are touching the Contact Strips as this establishes the circuit.)

TREATMENTS EXPLAINED

Negative Potential Galvanic Microcurrent

- Helps draw out impurities
- Enhances the absorption of skincare products where the active ingredient carries a negative charge, such as hyaluronic acid.

Warming

- Increases blood circulation delivering more oxygen and nutrients to the skin.
- Opens pores making it easier to cleanse and absorb skincare products.
- Relaxes muscles to ease facial tension, which can contribute to fine lines and dullness.
- Boosts lymphatic drainage to promotes detoxification and reduces puffiness.

Cooling

- Reduces puffiness by constricting blood vessels, minimizing swelling and inflammation, especially under the eyes.
- Tightens pores, which after cleansing helps close pores, giving the skin a smoother appearance.
- Soothes irritation and calms redness which can be beneficial for sensitive or acne-prone skin.
- Enhances skin firmness, temporarily tightening then skin, creating a lifted and refreshed look.

Vibration

Gentle vibration massage can stimulate localised blood flow to the facial tissues, promote lymphatic drainage, enhance product absorption (like serums or moisturizers) and ease facial muscle tension.

AFTER USE

If you have very sensitive skin, you may notice redness after use. This should subside after 30 minutes.

CLEANING AND MAINTENANCE

Before and after use, wipe the Treatment Head with a clean damp cloth, do not immerse the unit in water.

DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:

1. Check your local municipal guidelines
2. Separate the components and dispose of them correctly

SPECIFICATION

Model: FAWC

Input: 5V  1.0A 5.0W

Battery: Li-Ion 3.7V 900mAh 3.33Wh

Charging time: 3 hours

Maximum run time: 40 minutes

Read these instructions carefully before use and then retain for future use.

This product is for indoor use only.

This product complies with current EU product legislation.

Direct current

CUSTOMER CARE

If you have any questions, please contact us at:
online@riobeauty.com

GUARANTEE

To register your two year guarantee online, visit:
www.riobeauty.com/warranty

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© The Dezac Group Ltd 2025 Manufactured by The Dezac Group Ltd PO Box 17, Cheltenham, England GL53 7ET