HAIRGROW LASER HEAD BAND

the beauty specialists

ENGLISH	4
FRANÇAIS	10
DEUTSCH	16
NEDERLANDS	22
ESPAÑOL	28
PORTUGUÊS	34
ITALIANO	40

WARNINGS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 18 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Do NOT use this appliance if it is not working properly, or if it has been dropped or damaged.
- Keep out of reach of children and pets.
- This device must only be used for its intended purpose as described in this user guide.
- If you suffer from a skin condition, consult your dermatologist before use. This may include certain types of antibiotics, cancer treatments, and diuretics. In addition to medications some common including St. John's Wort, some deodorants, antibacterial soaps, artificial sweeteners and laundry detergents, can cause photo-sensitivity.
- Do not use on cuts or broken skin.
- Do not use this appliance on areas showing signs of: varicose veins, skin diseases (for example eczema, psoriasis, dermatitis), cuts and abrasions, irritation, warts, moles, skin tags, skin rash, hypersensitivity or sunburn, phlebitis, recent scar tissue, other skin conditions. We recommend you consult your physician before treatment if you suffer from any of the above conditions.
- Do not use on areas that have had cosmetic surgery.
- Do not charge for longer than 24 hours.
- This appliance must only be supplied at safety extra low voltage (SELV) corresponding to the marking on the appliance.
- Do not use immediately around the eyes.
- Do not use on areas that have had cosmetic surgery.
- Do not leave the device unattended. Turn off before you leave the room.
- The appliance may be cleaned by wiping with a damp cloth. Do NOT use solvents or abrasive cleaners.
- If the power lead is damaged, it must be replaced by the manufacturer or its service agent or similarly qualified person in order to avoid a hazard.
- Do NOT stare directly at the light source.
- Do not use on wet or damp hair.
- Do not use where it could come into contact with water.
- Avoid direct eye exposure

KNOW YOUR UNIT

(see page 3)

- | | | | |
|----------------------------|-----------------------|----------------------|--------------|
| ❶ HairGrow Laser Head Band | ❸ Locating Band | ❺ USB Charging Cable | ❽ On/Off |
| ❷ Laser Aperture | ❹ USB-C Charging Port | ❻ Indicator Light | ❾ Travel Bag |

INTRODUCTION

There are many reasons for your hair to lose volume and lustre. The HairGrow Laser Head Band offers a natural way to counter this.


Harnessing advanced 650nm light laser diodes technology, the HairGrow Laser Head Band delivers regenerative red light deep into the skin and to the root of your hair. This gentle light energy helps increase blood circulation in the scalp and stimulates cells responsible for hair growth to produce thicker hair and encourage new hair growth.

Who Can Benefit?


Laser diode red light treatment has been demonstrated to be beneficial for androgenetic alopecia, a common cause of hair loss in men and women. For men it has been shown to promote hair growth for *Norwood-Hamilton Classifications* of II to V patterns of hair loss, and for women with *Ludwig-Savin Scale* I-1 to I-4, II-1, and II-2, or frontal patterns of hair loss.

Treatments are suitable for Fitzpatrick Skin Types I - IV (very fair skin that always burns and never tans to moderate brown skin that rarely burns and tans easily).

NORWOOD-HAMILTON CLASSIFICATION

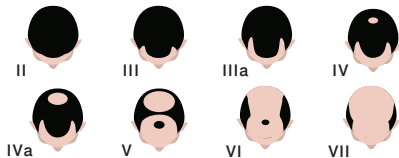LUDWIG-SAVIN SCALE

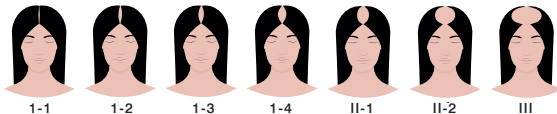CONTRADICTIONS


Consult your Doctor if:

- You are taking a medication or using an over the counter product which is known to cause photo-sensitivity. This may include certain types of antibiotics, cancer treatments, and diuretics. In addition to medications some common products including St. John's Wort, some deodorants, antibacterial soaps, artificial sweeteners and laundry detergents, can cause photo-sensitivity.
- You have any skin disorder.
- You have any eye disorder which might result in increase in photo-sensitivity.

Adverse Reactions

Skin irritation, nausea and headaches are rare but possible side effects. If experienced, stop using immediately and consult a medical professional.

CHARGING

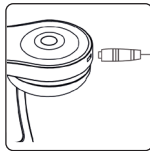Insert the USB Charging Cable into the USB-C Charging Port and connect to a USB power supply.

The Indicator Light will turn green to indicate charging and turn blue once fully charged.

Battery Charging Time: 2 hours

Maximum Run Time: 3 Hours

Note: When the battery is low, the HairGrow Laser Head Band will emit multiple distinct beeps to indicate the need for recharging.

IMPORTANT - PATCH TEST

Before using for the first time test for light sensitivity on your forearm.

Place unit over your inner forearm and run for three minutes. If you experience redness or itching which persists for more than an hour do not continue to use the product.

GETTING STARTED

Press and hold the On/Off button for 2 seconds to turn on. The unit will beep and the Indicator Light will turn blue.

Place the Laser Band over the treatment area and ensure the Locating Band is rotated clear of the Laser Band. The Laser Band should be within 3cm of the skin. The unit will emit two beeps to indicate the treatment has started. The Indicator Light around the On/Off button will flash blue.

You may feel a warm sensation during use.

If, during treatment, the laser band is too far from the skin, the unit will beep twice and the treatment laser diodes will temporarily switch off until the laser band is repositioned closer to the skin. A single beep will confirm that the laser diodes are active again.

Treatment Time:

The maximum recommended treatment time is 25 minutes. Studies have shown that longer periods are not beneficial. The unit will beep every 5 minutes and automatically turn off after 25 minutes with two beeps.

The minimum recommended time for an area is 10 minutes - if treating multiple areas you can move the band every 10 minutes

Treatment Coverage:

Depending on the area to be treated you may opt to treat different areas on alternative days.

Treatment Frequency:

Studies show that laser light therapy is effective as either a daily, or every other day treatment. Our recommendation is to use daily.

AFTER USE

If you have very sensitive skin, you may notice redness after use. This should subside after 30 minutes.

CLEANING AND MAINTENANCE

Before and after use, wipe the Treatment Head with a clean damp cloth, do not immerse the unit in water.

CUSTOMER CARE



If you have any questions, please contact us at:
online@riobeauty.com

GUARANTEE

To register your two year guarantee online, visit:
www.riobeauty.com/warranty

Treatment Positions:

The Laser Band and Location Band can be adjusted to provide full coverage of the scalp. Be careful to ensure the treatment light is not blocked by the Location Band.

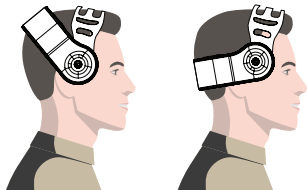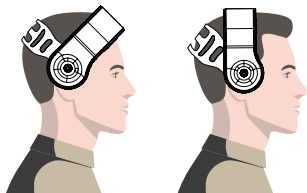DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:

1. Check your local municipal guidelines

2. Separate the components and dispose of them correctly

SPECIFICATION

Model: LAHB

Input: 5V $\equiv$ 1.5A 7.5W Max 

Battery: Li-Ion 3.7V 1500mAh 5.55Wh

Charging time: 2 hours

Maximum run time: 3 Hours

Laser Diode Wavelength: 655nm $\pm$ 10nm

Laser Power Classification: <5mW CW

Number of lasers diodes: 60

Laser Diode Life: 10,000 hours

Read these instructions carefully before use and then retain for future use.

Direct current

Recycle Paper

This product is for indoor use only.

This product complies with current UK product legislation.

Keep Dry

This product complies with current EU product legislation.

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© The Dezac Group Ltd 2025 Manufactured by The Dezac Group Ltd PO Box 17, Cheltenham, England GL53 7ET