

iWHITE
instant teeth whitening

› Sonic Toothbrush

**Charcoal-infused
toothbrush for
superior whitening**

Manual

Manufacturer:
Sylphar nv
Xavier De Cocklaan 42
B-9831 Deurle
Belgium
Tel +32 9 321 70 80
www.iwhiteinstant.com



CONGRATULATIONS ON YOUR PURCHASE!

Before using this product, read this user manual carefully to learn about the features, the do's and don'ts and how you can make best use of it.

Note: Charge the iWhite sonic toothbrush (referred to as "brush" for short in the following) fully before using it for the first time.

IMPORTANT SAFEGUARDS

WARNING: READ ALL INSTRUCTIONS BEFORE USE

1. Always check the product carefully before use. Discontinue use if the product appears damaged in any way (brush head, handle, charging pad).
2. The charging pad cord cannot be replaced. If the cord is damaged, discard the charging pad.
3. This product contains no user-serviceable parts. Do not disassemble, repair or modify the product yourself. Only certified professional technicians are qualified to fix any problems. Refer to "WARRANTY AND SUPPORT" if the product no longer works properly or needs repair.
4. Keep the product (including charging pad) out of reach of children at all times.
5. Always store the product in a cool and dry place away from direct sunlight when not using it for an extended period of time.
6. To avoid damage to the product, do not use a dishwasher to clean the brush head (including the protective cap), handle or charging pad.
7. Replace iWhite charcoal-infused brush heads every 3 months to achieve optimal results. Use only replacement iWhite charcoal-infused brush heads.
8. Do not use the brush while it is charging.
9. The product comes with a rechargeable Li-ion battery. This battery is self-contained and not consumer-replaceable. The battery must be disposed of properly when empty.
10. Always place and store the USB plug for the charging pad at a safe distance away from sources of water (i.e. sink or tub), even though the main section of the charging pad is IPX-7 waterproof.
11. Do not put the toothbrush onto a charging pad which has fallen into water. Unplug immediately.
12. Never use a charging pad with a damaged cord.
13. Do not use while bathing or showering in warm/hot water. Although the product is IPX7 waterproof, it cannot be protected effectively from hot water or steam.
14. This product is designed to clean your teeth, gums and tongue only. Use this product only for its intended use as described in this user manual or as recommended by your dental professional.
15. This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the electronic device by a person responsible for their safety.
16. This product is not recommended for use by children under the age of 12.
17. This product is a personal care device and is not intended for use on multiple patients in a dental office or institution.
18. Use a UL-listed computer, USB hub or USB duplex receptacle when charging the brush handle.
19. In order to effectively charge your brush, do not place your charging pad on a magnetic metal surface.
20. For connection to a supply not in your country, use an adapter of the proper configuration for the power outlet.
21. Consult your dental professional before using this product if you have had oral or gum surgery in the previous 2 months.
22. Contact your dental professional if excessive bleeding occurs after using this product, or if bleeding continues after 1 week of use.
23. Discontinue use of this product and contact a physician/dentist if you experience discomfort or pain.
24. This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician prior to use.
25. Consult your physician if you have other medical concerns.
26. This brush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers).
27. It is normal for your gums to feel itchy, numb or to experience slight bleeding after you use the brush for the first time. Your gums will get used to the sonic vibration in about two days and the discomfort will lessen over the course of a week. Consult your physician if discomfort continues.

BOX CONTENTS

- 1 iWhite Sonic Toothbrush
- 1 USB charger
- 1 replaceable brush head*

**please store with the transparent cap which can be used to protect the head bristles when you travel.*


Note:

There is no adaptor in the box. Any smartphone or other electronic adaptor with an output voltage of 5V can safely be used with the brush charging pad if necessary.

SPECIFICATIONS

Model No.	IWE-200-0000
Size	251X29X29mm
Frequency	31,000~40,000VPM
Brushing modes	Clean, Sensitive, Gum Care, White, Polish
Working voltage	3.3~4.2V
Battery	2,000mAh Li-ion, a single charge lasts 90 days
Charging	About 30 hours to full charge
Low battery warning	Flashing red + vibrating motor when switched on
Button	Dual Button (Power Button + Mode Button)
Indicator on charging	Charging: Solid red Fully charged: Solid green for 30s and then off
Waterproof	IPX-7
Smart memory	The brush will remember your latest brushing mode and activate it at the next use.
Brush head	Food-grade bristles with protective cap
Noise	<60dB
Weight	139g (with brush head)

PRODUCT

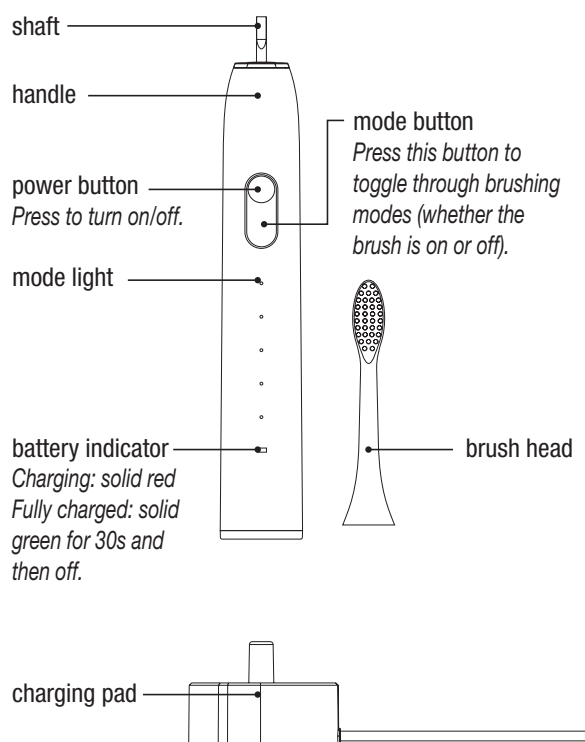


HOW TO INSTALL THE BRUSH

Align the brush head so the bristles point in the same direction as the front of the handle.
Firmly press the brush head down onto the metal shaft until it stops.


It is normal to see a slight gap between the brush head and the handle.

WEEKLY MAINTENANCE

Rinse the top of the brush handle and the button on the toothbrush head once a week.

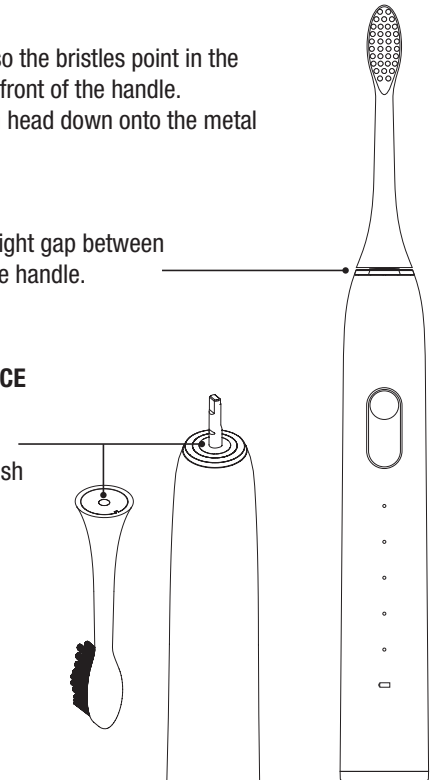

BRUSHING MODES

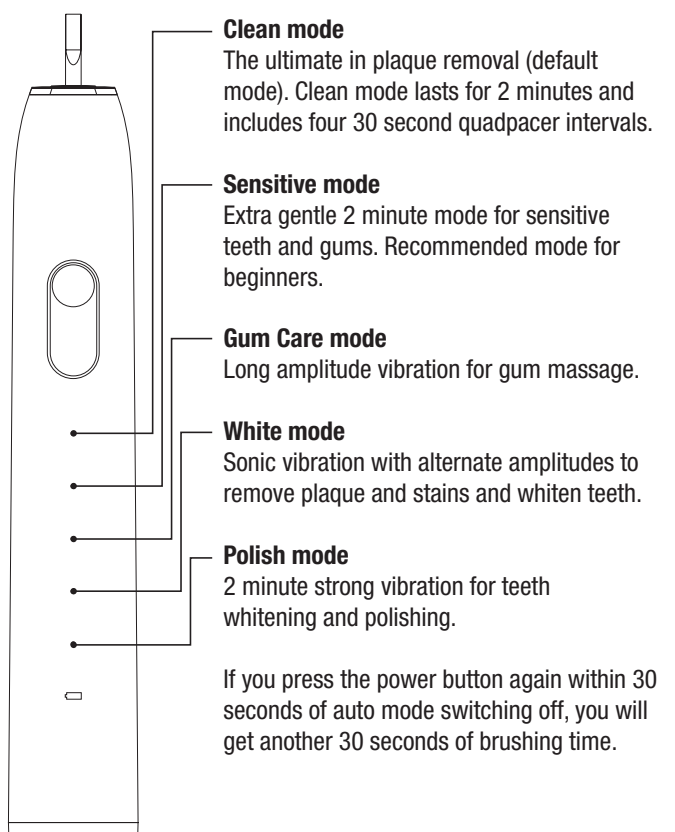


HOW TO TURN ON/OFF THE BRUSH

- ① Press the mode button to select the brushing mode.
- ② Wet the bristles and apply a small amount of toothpaste.
- ③ Place the toothbrush bristles against the teeth (do not press hard).
- ④ Press the power button to turn on.
- ⑤ When you have finished brushing, the brush will turn off automatically.
- ⑥ Hold the power button for 2 seconds to turn off.

NOTES OF BRUSHING TEETH

- ① A quadpacer interval every 30 seconds reminds you to brush each section of your mouth.


Notes: To ensure you brush your whole mouth evenly, divide your mouth into 4 sections using the quadpacer feature.

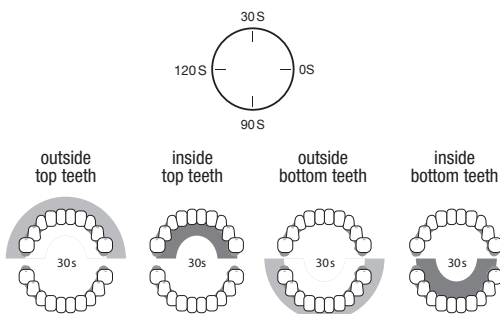

- ② For a better clean, brush your teeth at a slight angle (45 degrees).




- ③ Keep the bristles gently placed on the teeth and in the gumline. Brush your teeth using a light back and forth motion.

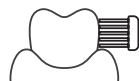

- ④ Make sure you turn off the brush before taking it out of your mouth to avoid toothpaste splashes.

- ⑤ Pull out the brush head in the direction of the arrow to replace the head.

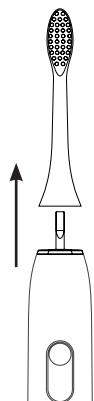

- ⑥ If you are not going to use your brush for an extended period of time, clean it and store it in a cool and dry place away from fire, heaters and direct sunlight.

- ⑦ Rinse your brush head and handle with water after each use. Wipe the entire surface of the handle with a soft tissue or cloth.

WARRANTY AND SUPPORT

We provide a 2 year warranty on the product, starting from the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at our expense, provided that convincing proof of purchase within the qualifying period is provided.

We recommend that you use genuine iWhite charcoal-infused brush heads for optimal performance. Use of other brush heads may damage your brush and limit your warranty protection.

WARRANTY EXCLUSIONS

What is not covered under warranty:

- Brush heads.
- Damage caused by use of unauthorised replacement parts or unauthorised brush heads.
- Damage caused by transport, accidental damage, misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.
- If the proof of purchase has been altered, defaced or removed.

To obtain service within the warranty period, hand in or send the complete appliance with your sales receipt to the manufacturer:

Sylphar nv
Xavier De Cocklaan 42
B-9831 Deurle
Belgium

Note: Postage costs shall not be refunded.