

GIO perfect Brush

4 in 1 LED facial brush




USER MANUAL

To avoid harm caused by misuse, please read the instructions carefully and use this product strictly in accordance with the manual.

PRODUCT INTRODUCTION

This 4-in-1 LED facial brush offers a simple and effective skincare experience. Using high-frequency vibration and heat, it promotes deep cleansing and enhances the absorption of skincare products. The replaceable silicone sleeves allow for easy maintenance and hygiene. Its soft, long-bristle brush addresses various skin concerns, leaving the skin refreshed and radiant. The ergonomic, anti-slip handle ensures comfortable use. Made from food-grade silicone, the device is safe and suitable for home use.

PARTS DESCRIPTION

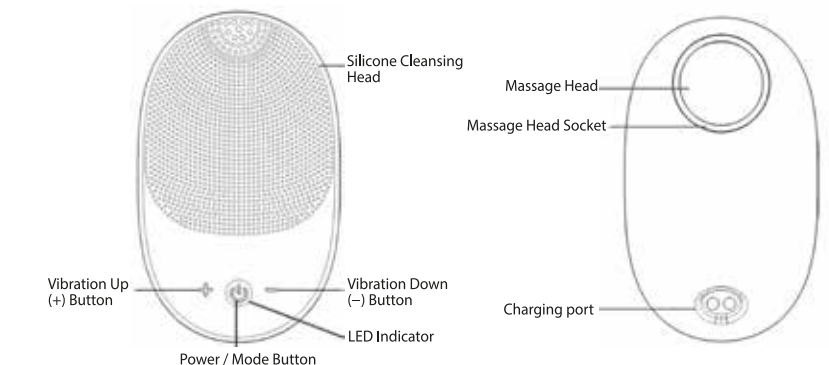

HOW TO USE

Before First Use

1. Fully charge the device using the supplied magnetic charging cable.
2. Clean the device with water and a mild facial cleanser, then rinse thoroughly with warm water. Do not use alcohol, acetone, gasoline, or corrosive cleaners.

Turning On & Selecting Modes

- Press and hold the power button for 3 seconds to turn on the device.
- The device will start in Green Light Cleansing Mode (vibration on).
- Press the function button once to switch to Blue Light Mode, and again to switch to Red Light Heating Mode.
- Modes cycle in this order: Green/ Blue /Red /Green.

Note: In Blue Light mode, the massage head may temporarily warm up after 4 minutes due to residual heat conduction.

Mode Descriptions

- Green Light – Cleansing Mode
Use with a cleanser to deeply clean pores.
- Blue Light - Whitening & Pore-Tightening Mode
Helps tighten pores and soothe the skin.
- Red Light – Heating Mode
Used with serums or essences to promote absorption and stimulate skin vitality.

Adjusting Vibration Levels

- Vibration intensity is adjustable from Level 0 (off) to Level 5.
- Default is Level 3.
- Use the "+" and "-" buttons to increase or decrease intensity.

Battery Indicator

- Blinking red light: Low battery
- Solid red light: Charging
- Light off: Fully charged

DAILY USE INSTRUCTIONS

1. Makeup Removal & Daily Cleansing



- ✓ **Step 1: Ensure the device is fully charged**

Make sure the device is fully charged. Disconnect the magnetic charging cable and ensure the brush head is clean and dry before use.

- ✓ **Step 2: Apply cleanser and activate cleansing mode**

Wet your face with warm water. Evenly apply a suitable cleansing product — such as cleansing foam, oil, or micellar water. Press and hold the power button for 3 seconds to turn on the device. It will start in Green Light Cleansing Mode.



- ✓ **Step 3: Cleanse your face using gentle motions**

Place the silicone brush flat against your skin. Gently move the brush in circular motions — from cheeks to chin, nose, and forehead. Avoid the eye area. After cleansing, rinse your face with clean water and pat dry with a soft towel.

2. Deep Cleansing (Recommended once per week)



Step 1: Check battery level

Ensure the device is sufficiently charged before starting the deep cleansing routine. Recharge if needed, then disconnect the cable.

Step 2: Wet face and apply cleanser

Rinse your face with warm water. Apply an appropriate amount of facial cleanser and gently massage it into a rich lather across your skin.



Step 3: Activate cleansing mode

Turn on the device and select Green Light Cleansing Mode or Blue Light Mode, depending on your preference. Use gentle, circular motions to deeply cleanse the skin. Focus on areas prone to clogged pores, like the T-zone.

Step 4: Rinse face and device

After cleansing, rinse your face with clean water to remove any cleanser residue. Also rinse the silicone brush head under running water.



Step 5: Apply skincare product and use heating mode

Evenly apply a serum, essence, or moisturizing product to your face. Switch to Red Light Heating Mode and gently massage the product into the skin to enhance absorption and boost skin vitality.

Step 6: Rinse the face with water

Rinse your face with clean water and pat dry with a towel.



Friendly reminder:

Use light pressure and avoid over-cleansing. Deep cleansing should not exceed 2–3 minutes in total. Always clean and dry the device after use.

CHARGING INSTRUCTIONS

1. Use the provided magnetic charging cable (no polarity).
2. Ensure the device and charging port are dry before charging.
3. Slight warming during charging is normal.
4. Do not exceed 24 hours of continuous charging.
5. When using abroad, verify that the local voltage (5V/1A) is compatible.

MAINTENANCE & CLEANING

- The device is IPX7 waterproof and can be safely rinsed with water after each use.
- After cleaning, leave it in a well-ventilated place to dry completely.
- Store in a cool, dry place, away from direct sunlight or heat sources.
- Charge the device every 3 months if unused for extended periods to preserve battery life.

CAUTIONS

1. Do not use this device if you have open wounds, skin infections, severe acne, sunburn, or have recently undergone facial surgery.
2. Do not use exfoliating scrubs, granular cleansers, or silicone-based products with this device, as they may damage the brush.
3. Do not use the device more frequently than recommended. Excessive use may irritate the skin.
4. The device is not suitable for children under 12 years of age. Minors should use it only under adult supervision.
5. Do not share the device with others for hygienic reasons.
6. Clean the brush thoroughly after each use to prevent buildup of product or bacteria.
7. Store the device in a cool, dry place away from sunlight or heat sources.

DANGER

1. Keep the device and its accessories out of reach of infants and young children.
2. Do not use or store the device near fire, heat, or in environments above 60°C.
3. Never charge the device near water (e.g., sinks, bathtubs, or in wet bathrooms).
4. Use only certified 5V/1A chargers. Do not use damaged cables or loose power sockets.
5. Do not handle the charging cable with wet hands.
6. Do not attempt to open, disassemble, puncture, or modify the device or battery.
7. Do not use the device if it is deformed, damaged, or malfunctioning.
8. Do not use the device while it is charging.

WARNING

1. If you have known skin conditions or allergies, consult a dermatologist before use.
2. Discontinue use immediately if irritation, redness, or discomfort occurs.
3. This device is intended for facial skincare only. Do not use on other body parts.
4. Keep the device dry when not in use. After cleaning, allow it to air dry fully.
5. The device contains a built-in rechargeable battery. Do not attempt to replace it.



PRODUCT SPECIFICATIONS

Product Name: GIOperfect Brush

Materials: ABS + Silicone

Product Dimensions: 100 mm × 60 mm × 40 mm

Net Weight: 86 g

Battery Capacity: 600 mAh

Charging Method: Magnetic cable

Rated Power: 3 W

Charging Voltage: 5 V / 1 A

Working Voltage: 3.7 V

Vibration Levels: 5 gears

Heating Function: 40°C – 46°C

Indicator Light: Tri-color LED

Operating modes:

- Green Light – Massage & Cleansing Mode
- Blue Light – Whitening & Pore-Tightening Mode
- Red Light – Hot Compress Mode for Enhanced Absorption

Working Time: 1–3 hours

Waterproof Rating: IPX7



Scan the QR code to download the manual



Zeskanuj kod QR, aby pobrać instrukcję



Scannen Sie den QR-Code, um das Handbuch herunterzuladen



Scannez le code QR pour télécharger le manuel



Escanee el código QR para descargar el manual



Digitalize o código QR para descarregar o manual



Scansiona il codice QR per scaricare il manuale



Scan de QR-code om de handleiding te downloaden



Скануйте QR-код, щоб завантажити інструкцію



Naskenujte QR kód a stáhněte si návod k použití



Naskenujte QR kód a stiahnite si príručku



Scanează codul QR pentru a descărca manualul



Nuskaitykite QR kodą, kad atsisiųstumėte instrukciją



Nolasiet QR kodu, lai lejupielādētu rokasgrāmatu



Lae juhend alla, skannides QR-koodi



A kézikönyv letöltéséhez olvassa be a QR-kódot.



ليزنتل لى ذى رسلا ةباجتسالال زمر حسما



innoGIO

Manufactured in China for **INNOGIO GLOBAL** sp. z o.o.
Ostródzka street 74H, 03-289 Warsaw, Poland

www.innogio.com



@innogio.official



@innogio_com