

GIO perfect Eye

Sonic LED Facial Massager with EMS Therapy



USER MANUAL

To avoid harm caused by misuse, please read the instructions carefully and use this product strictly in accordance with the manual.

IMPORTANT SAFETY INFORMATION

This product is intended for external cosmetic use only. Keep out of reach of children. Not suitable for use by individuals with medical implants or serious medical conditions without prior consultation with a healthcare professional.


Do not use if you:

- Have a pacemaker, artificial heart-lung machine, or other implanted medical electronic devices.
- Are undergoing medical treatment, have skin conditions, fever, acute illness, tuberculosis, organ disease, or telangiectasia.
- Are pregnant, menstruating, or breastfeeding (consult your doctor before use).
- Have extremely sensitive skin, metal allergies, hereditary dermatitis, or impaired temperature perception.
- Have had recent eye or facial surgery (e.g. laser vision correction, fillers, thread lifting); consult your physician.
- Are undergoing dental treatment or orthodontics, or have fillings – avoid strong stimulation near the mouth.
- Wear contact lenses when using the device around the eyes (remove lenses first).
- Have cataracts, glaucoma, retinal disorders, or any eye-related conditions – consult a physician before use.

PRODUCT PRECAUTIONS

1. Always test the device on your wrist before first use.
2. Do not immerse the product in water or expose it to moisture. It is not waterproof.
3. Clean only with a soft damp or dry cloth once the device is turned off and cooled.
4. Avoid direct sunlight and contact with corrosive materials.
5. Do not disassemble or repair the device yourself. Contact after-sales support for service.
6. Do not keep the massage head in one place on the skin for too long during use.
7. Use only by adults. Suitable for both men and women.

OPERATING INSTRUCTIONS

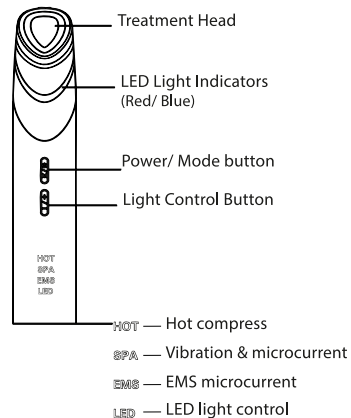

POWERING ON/ OFF

Press and hold the Power / Mode Button for 2 seconds to switch the device on. Press and hold the same button again for 2 seconds to switch the device off.

MODE SELECTION

After switching on, the device will automatically enter HOT Mode (Red LED + Heating). Each short press of the Power / Mode Button will cycle through the following modes in order: HOT Mode (Red Light + Heating) SPA Mode (Vibration + Gentle Microcurrent) EMS Mode (Blue Light + Stronger Microcurrent) The modes will repeat in the same cycle with each press.

LED Light Control

By default, the LED light is active when the device is turned on. To turn the LED lights off, press and hold the Light Control Button for 2 seconds. Press and hold again to turn the lights back on. This can be done in any mode.

Low Battery Alert

When the battery is low, the red light on the treatment head will flash several times before the device powers off automatically.

Charging Indicator

Red Light On – Device is charging
Green Light On – Battery is fully charged



HOW TO USE

1. After cleansing, apply a thin layer of gel or serum evenly around the eye area.



2. Press and hold the Power/Mode button to turn on the device. Select the desired mode and intensity level.

3. Gently glide the device from the inner to the outer corner of the eye, or use the dot pressure method: press the head to one area for 4–5 seconds, then move to the next once you feel the warming sensation.



4. After use, turn off the device, wipe the treatment head clean, and store it in a dry place.



CLEANING AND MAINTENANCE

- Always clean the treatment head immediately after use.
- Do not wash the device with water or immerse it.
- Do not use harsh or corrosive cleaning agents, thinners, or abrasive cloths.
- Store in a cool, dry, well-ventilated place, away from direct sunlight.

PRODUCT SPECIFICATION



- Product Name: GIOperfect Eye
- Product Model: GIO-720
- Product Size: 138*24*24.8mm
- Rated Power: 3W
- Operating Voltage: 3.7V
- Battery Capacity: 500mAh
- Hot Compress Function: $42 \pm 3^{\circ}\text{C}$

- 🇬🇧 Scan the QR code to download the manual
- 🇵🇱 Zeskanuj kod QR, aby pobrać instrukcję
- 🇩🇪 Scannen Sie den QR-Code, um das Handbuch herunterzuladen
- 🇫🇷 Scannez le code QR pour télécharger le manuel
- 🇪🇸 Escanee el código QR para descargar el manual
- 🇨🇿 Digitalize o código QR para descarregar o manual
- 🇮🇹 Scansiona il codice QR per scaricare il manuale
- 🇳🇱 Scan de QR-code om de handleiding te downloaden
- 🇺🇦 Скануйте QR-код, щоб завантажити інструкцію
- 🇨🇪 Naskenujte QR kód a stáhněte si návod k použití
- 🇸🇰 Naskenujte QR kód a stiahnite si príručku
- 🇷🇺 Scanează codul QR pentru a descărca manualul
- 🇵🇱 Nuskaitykite QR kodą, kad atsisiųstumėte instrukciją
- 🇵🇱 Nolasiet QR kodu, lai lejupielādētu rokasgrāmatu
- 🇵🇱 Lae juhend alla, skannides QR-koodi
- 🇮🇹 A kézikönyv letöltéséhez olvassa be a QR-kódot.
- 🇵🇱 لایلدن لایزنتک دفعہ درملا قباچتسالما زمر حسبہ



innogIO

Manufactured in China for
INNOGIO GLOBAL sp. z o.o.
Ostródzka street 74H, 03-289 Warsaw, Poland
www.innogio.com

 @innogIO.official  @innogio_com

