BEAUTIFLY

B-Hydrapeel

PRO
series

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
Specifications of B-Hydrapeel

Product colour:	white, pink gold
Weight of the instrument:	approx. 166 g
Dimensions of the instrument:	188 x 55 x 75 mm
Rated voltage:	3.7V
Supply voltage:	5V -- 1A
Rated power:	1W
Intensity levels:	3 levels
Instrument storage temperature:	0°C – 45°C
Operating environment temperature:	5°C – 35°C
Operating environment humidity:	25% – 75%

Contents of the packaging

1. Silicone cleaning head (large) x 1
2. Small cleaning head x 1
3. Large cleaning head x 1
4. Spiral cleaning head (small) x1
5. Cotton bag x1 (includes: filter cap, fleece pads, cleaning brush)
6. Facial hydrogen cleansing beauty instrument x1
7. Instructions x1
8. Charging cable x1

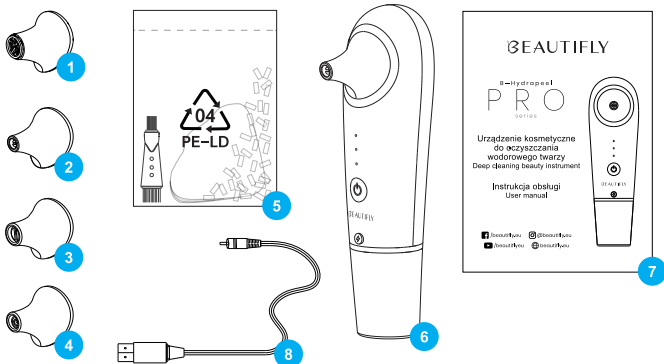
Hydrogen purification (hydrabrasion)

Hydrabrasion is a gentle cleansing, moisturising and exfoliating treatment. It is much more delicate than microdermabrasion, so there is no contraindication for it. Even people with sensitive, mature, vascular or acne-prone skin can benefit from this type of treatment. For maximum cleansing results of the treatment use special heads for particular areas of your face. They not only clean your skin even in hard-to-reach areas, but also help you get the results you want. Hydrabrasion guarantees skin as smooth as silk!

Description of the various parts of the instrument

1. Intensity level light
2. Charging socket
3. Power supply
4. Cleaning head
5. Charging light
6. Power button (On/Off switch)
7. Ventilation

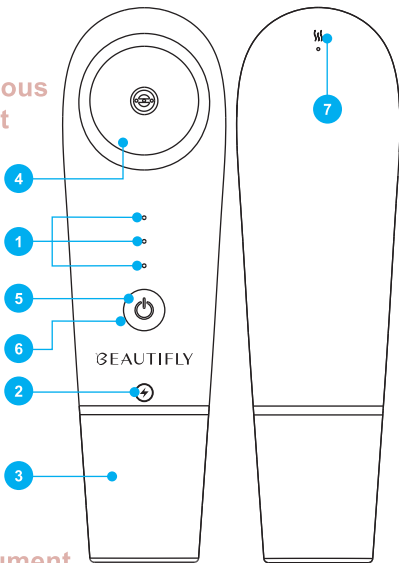Operation of the instrument

1. Press and hold the power button for 3 seconds to turn the instrument on or off.
2. Briefly press the power button to increase or decrease the intensity level.
3. The first intensity level is indicated by one green light, the second intensity level by two green lights and the third level by three green lights.
4. During charging, the charging light lights up red. When the battery is full, the red light turns off and the green light turns on.
5. It takes approximately 3.5 hours to fully charge the instrument.

Precautions during use


1. Do not place the product in an environment with a strong magnetic field.
2. Be careful not to tilt the instrument more than 45°C during use, otherwise water will spray out of the exhaust port.
3. This instrument can only use clean water. Do not use other liquids that are not suitable for this instrument.
4. Do not use the instrument while it is charging.
5. The inside of the water tank can be disinfected and sterilised.
6. This instrument should not be soaked or immersed in water for cleaning. Be careful when using the instrument and make sure that water does not penetrate inside when you use or clean the instrument.
7. After each use, drain the water that has been used from the water tank, then dry the instrument thoroughly.

How to use the instrument


1. Remove make-up from your face with a make-up remover.
2. Install the filter cap on the bottom of the instrument and fill the water tank 2/3 with clean water.
3. Install the correct cleaning head. There are different sizes of cleaning heads. Depending on your needs, you can choose the correct cleaning head for different areas.
4. After the facial cleansing treatment is complete, rinse the skin with a tonic and apply a moisturising preparation.


1

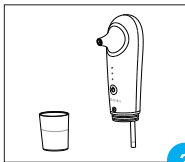2

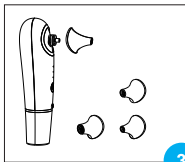
3

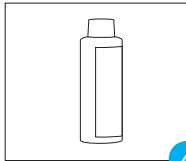4

Description of the functions of the cleaning heads

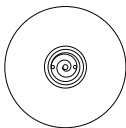
Spiral cleaning head (small)

Central suction (commonly used)

Function: dead skin, epidermis, pore cleansing (universal type)

Location: small area

Area: T zone, nasal part

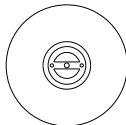
Cleaning head (large)

Large suction (commonly used)

Function: removing blackheads, dead skin, calluses, pore cleansing (universal type)

Location: large area of the face

Area: cheeks, forehead

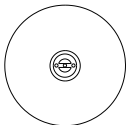Cleaning head (small)

Enhanced suction power (persistent impurities)

Function: blackheads, pimples, crusts, pore cleansing (pre-blackheads) (universal type) Location: small area Area: T-zone, nasal part

Silicone cleaning head (large)

Increased suction power (persistent impurities)

Function: acne, calluses, pore cleansing (blackheads preferred) (universal type) Location: large area of the face Area: cheeks, forehead

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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Beautifully Sp. z o.o., Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beautifully Sp. z o. o. declares that the device equipment Beautifully B-Hydrapeel is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beautifly Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.
In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@beautifly.eu

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