PL EN DE FR ES PT IT CZ SK SL EE LT LV  
NB FI HU RO AR RU EL HR SV DA NL TR BG

# BEAUTIFLY

Scrub

# PRO

series


Multifunctional skin  
scrubber

User manual

 /beautifly.eu  @beautifly.eu  beautifly.eu

 /beautiflyeu  @beautifly.eu

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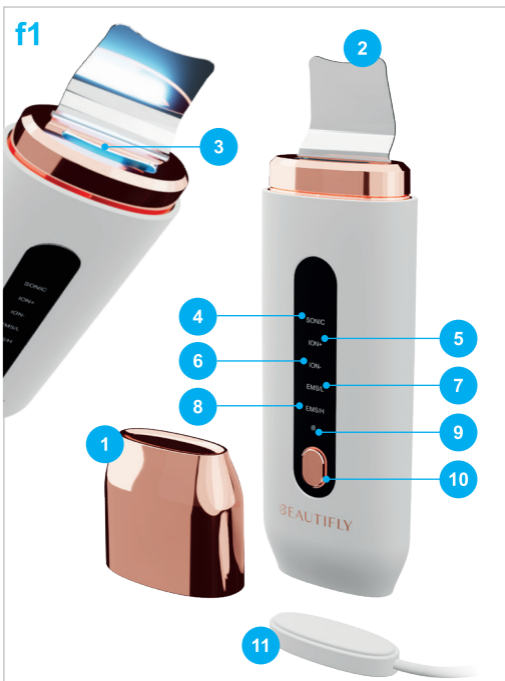rated voltage:	3.7V
vibration frequency:	29 ± 5K Hz
EMS:	low 85V / high 100V
landing:	induction base <4H
timepiece:	5 min for each mode
size:	170mm x 52mm x 26 mm
weight:	134g
material:	ABS
waterproof:	IPX5
colour:	white matt, rose gold

Part description [ see page 2 - f1 ]

1. Plug
2. Stainless steel spatula
3. LED photon light (blue / red)
4. Ultrasonic plan
5. Positive ion export plan
6. Anion import tips
7. EMS indicator (low)
8. EMS indicator (high)
9. Power indicator
10. ON / OFF / mode button
11. Induction charging base

#### Cleansing

Turn on the SONIC device with ION+, the exfoliating paddle vibrates at high frequency, the blue light is on, automatic shutdown after 5 minutes.

#### Moisturizing

Turn on the device and press the on/off/mode button on SONIC + ION- 1 time, the scrub paddle vibrates at high frequency, red light on, automatic shutdown after 5 minutes.

#### Rejuvenation

- a. Turn on the device and press the on/off/mode button 2 times under SONIC + EMS/L mode, the peeling head will vibrate in high frequency, EMS in low frequency, automatically shut down after 5 minutes.
- b. Turn on the device and press the on/off button three times under SONIC + EMS/H mode, scrub head vibrates in high frequency, EMS high intensity, automatically shuts off after 5 minutes

#### Automatic shutdown

2 beeps, light off, device off.

#### Manual shutdown

Press the on/off/mode button to turn off the device.

#### Low battery

3 beeps, button turns red and flashes.

#### Charging the device

The on/off/mode button turns red and flashes.

#### Charging the device

On/off/mode button done turns green.

#### Exfoliation of the epidermis with ultrasounds (Ultrasounds + Ionization +)

1. Before starting the treatment, thoroughly clean the skin and remove the make-up. Then moisten the face with a cotton swab or atomizer. For this purpose, you can use saline, tonic or water. The skin must be kept moist during the treatment. You cannot perform the procedure on dry skin.

2. Turn on the device and start the **Ultrasound + Ionization + mode**.
3. Gently move the spatula from bottom to top, lengthwise facial contours at a 45 degree angle. The spatula should be positioned in such a way that the photon light reaches the skin.
4. The device will turn off automatically after 5 minutes.
5. Wipe your face with toner.

Detailed, correct instructions for performing the cavitation peeling treatment are available on our YouTube channel:

#### Introducing nutrients by ultrasound (Ultrasound + Ionization-)

1. After cleansing your face, apply the product on it care (cream, serum, ampoule, gel).
2. Turn on the device and press the On/Off button once to activate the **Ultrasound + Ionization-** mode.
3. Gently move the spatula from bottom to top along the contours face. The spatula should be placed in the opposite way than in the first step.
4. The device will turn off automatically after 5 minutes.

#### EMS (Ultrasounds + EMS currents).

1. After facial skin care, turn on the device and press the on/off switch twice to activate the **Ultrasound + EMS/Weak** mode (there are two operating levels to choose from).
2. To improve skin elasticity, move the spatula in upward motions along the contours of the face at a 45-degree angle.
3. The device will turn off automatically after 5 minutes.
4. After using the device, wipe it with a damp cotton cloth.

#### Cautions

Person with electronic human organ are forbidden to use it. Person with skin surgery within 1 month or skin sunburn is forbidden to use.

Forbidden to use by person with inflammation of the skin, scar, damaged and electromagnetic sensitive is forbidden to use.

Take off watch, glasses and other metal objects before operation.

Avoid shovel head close to eyeball.

Don't operate in the same area for long time.

Do not use the device for pregnant women.

#### Cleaning

Rinse the shovel head and wipe it with a dry soft cloth after using.

Do not immerse the device into water for long time.

#### Storage

Keep it dry and clean.

Keep it away from high temperature.

Keep it out of reach of children.

More information about the company and products [page 2 - f2]

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**[EN] Rules for handling electrical and electronic equipment:** Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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**[EN] Summary declaration of conformity:** The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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## Shortened warranty terms

This Guarantee is granted by a Polish company, Beautifly Sp. z o. o. based in Warsaw, Plac Bankowy 2, 00-095 Warsaw, Poland. You can find us in the official Polish register of enterprises under the KRS number: 0000887306, NIP: 5252854218 and REGON: 388382021 We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: [www.beautifly.eu/warranty-beautifly](http://www.beautifly.eu/warranty-beautifly)