BEAUTIFLY

B-Scrub Perfume

Blush

series

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IPX5 WATERPROOF

UNIQUE DESIGN


B-Scrub Perfume from the Blush series is a device for cavitation peeling. As a result of the treatments, the skin will increase collagen production, become perfectly soft, radiant and smooth. All this, hidden in a small device, which every day is an elegant gadget for every woman. The craftsmanship was given to him by the shape of a bottle of perfume. The Beautify Blush series will make you feel extremely refined.

Technical data

rated voltage:	5VDC
rated current:	<800mA
vibration frequency:	26±5K Hz
time of surgery:	5 min.
charging time:	<4h
waterproof:	IPX5

Part description

1. Magnetic plug
2. Spatula
3. CLEAN Mode Light
4. On/Off/Mode Button
5. LIFT Mode light
6. Sensor sheet
7. Power Indicator Light
8. USB Charging port
9. Silicone Stopper
10. Type-C Cable
11. Cleaning Cloth

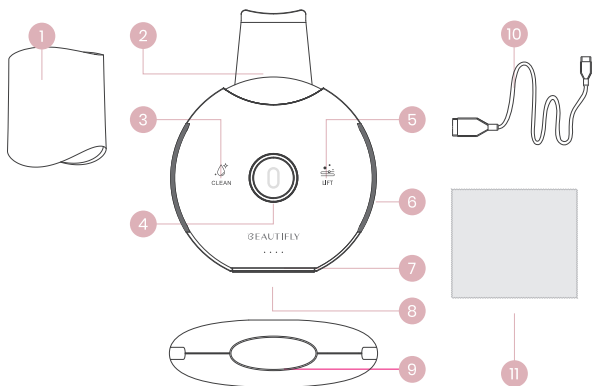This product is with built-in lithium battery and can be charged with Type-C cable (included) ONLY.

Description of the modes

	Status	Description
1.	CLEAN mode	Press the On/Off/mode Button, you can hear 1 beep and enter CLEAN Mode with the CLEAN mode light on.
2.	LIFT mode	Under CLEAN Mode, press the On/Off/Mode Button, you can hear 1 beep and enter LIFT Mode with LIFT Mode Light on.
3.	OFF	Auto-off: the device will be auto-off after 5 minutes' working along with 5 beeps. Non-auto Off: under LIFT Mode, press the On/ Off/Mode Button to turn off the device along with 1 beep.
4.	Low battery	The device would beep 3 times indicating low battery, which would repeat every 10 seconds. The device will be automatically off after 9 beeps. Note: you can check battery power when you are using the device. The fewer indicator lights are on, the lower the power is. The power is less than 25% with one indicator light on only.
5.	Charging	The Power Indicator Lights flash in turn.
6.	Full Charged	All the Power Indicator Lights are on.

CLEAN mode

1. Before starting the treatment, cleanse the skin thoroughly, wash your makeup. Then wet your face with a cotton swab or atomizer. For this purpose, you can use saline, tonic or water. The skin must be moist during the procedure. You cannot do the procedure on dry skin.
2. Turn on the device by pressing the button in the center. CLEAN mode is activated immediately.
3. Gently move the spatula from bottom to top along face contours at an angle of 45 degrees.
4. The device will turn off automatically after 5 minutes.
5. Wipe your face with tonic.

LIFT mode

1. After cleansing your face, apply a care product (cream, serum, ampoule, gel) to it.
2. Turn on the device by pressing the button twice.
3. Gently move the spatula from the bottom up, along the face contours at an angle of 45 degrees, thanks to which you will thoroughly introduce the cosmetic.

Cautions

1. A person with a pacemaker cannot perform a procedure using the B—Scrub Luxury
2. It is forbidden for people after skin surgery to use it within 1 month or for skin burns.
3. It is forbidden to use by people with dermatitis, scar, damage and electromagnetic sensitivity.
4. Avoid contact of the head with the eyeball.
5. Do not work in the same area for a long time.
6. Do not use the device by pregnant women.

Cleaning

Rinse the shovel head and wipe it with a dry soft cloth after using.
Do not immerse the device into water for long time.

Storage

1. Keep it dry and clean.
2. Keep it away from high temperature.
3. Keep it out of reach of children.

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The crossed-out wheeled bin symbol indicates that in the European Union after the end of the product's useful life it must be disposed of in a separate, dedicated collection point. Do not dispose of these products together with unsorted general household waste

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Beauty Factor Sp. z o.o.
Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beauty Factor Sp. z o. o. declares that the device equipment Beautifly B-Scrub Perfume is in compliance with Directive 2014/30/UE. The full text of the EU declaration of conformity is available at the following internet address:

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beauty Factor Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address biuro@beautifly.eu

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