BEAUTIFLY

B-Booster
PRO
series

Fale radiowe

Rejuvenation radio frequency

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User manual

[f /beautifly.eu](https://www.facebook.com/beautifly.eu) [@beautifly.eu](https://www.instagram.com/beautifly.eu) [g beautifly.eu](http://beautifly.eu)

[/beautiflyeu](https://www.youtube.com/beautiflyeu) [@beautifly.eu](https://www.tiktok.com/@beautifly.eu)

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
Technical data B—Booster

Nominal voltage:	3.7V
Operating current:	SONIC 100mA-400mA RF 200mA-2800mA EMS 8mA-30mA
Ultrasonic frequency:	3MHz±10%
Radio frequency:	1MHz
Waterproof:	IPX5
Working time:	3Min, 5Min, 6Min
Working:	<3H
Battery capacity:	1000 mAh

Note: This product has a built-in lithium battery and only supports USB-C charging cable and adapter as a charging method.

Part description

1. Transparent protective cover
2. V face roller
3. HD display
4. Switch/mode key
5. Magnetic snap ring
6. Nursing head
7. Magnetic base
8. Induction strip
9. Type-C Charging hole
10. Type-C Cable
11. Cleaning cloth

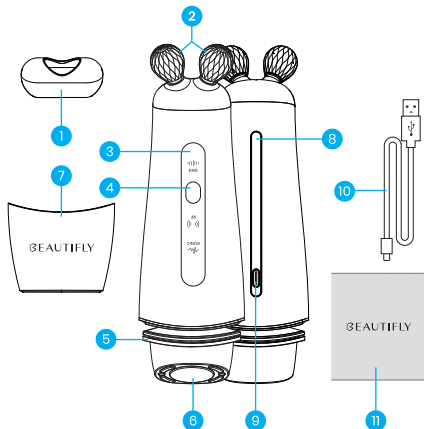Functions

Functions B-Booster Rejuvenation Radio Frequency. Strongly rejuvenating cosmetic apparatus, through ultrasound vibrations combined with positive ions, deeply exports dirt residues to the outside, preparing the skin for further treatment steps. 1MHz RF radio frequency combined with negative ions activates the production of collagen and elastin, while importing nutrients that penetrate the skin. Lifting EMS micro-currents tighten the skin, and the diamond facets of the roller deeply massage the face, revealing its perfect V-contour.

Function description

No.	Status	Indications
1.	"SONIC" mode	Press and hold the mode switching button for 1 second, the beep will sound once and the corresponding indicator will light up and turn off automatically after 3 minutes.
2.	"RF" mode	In "SONIC" mode, short press the mode switch button to switch, the buzzer will beep one beep, and the corresponding indicator will be on for a long time and automatically turn off after 5 minutes.

No.	Status	Indications
3.	"EMS" mode	In "RF" mode, short press the mode switch button to switch, the buzzer will beep one time, and the corresponding indicator will be on for a long time and automatically turn off after 6 minutes.
4.	Turning off	Manual shutdown: In any mode, press and hold the mode switch button for 1 second to turn off. Automatic shutdown will take place when the product is finished in any mode, it turns off after a beep.
5.	Not enough power	When the battery level is less than 10%, the white battery indicator will flash slowly once every 3 seconds and cyclically 3 times. Slow blinking white will change to fast blinking a third time and the buzzer will beep twice until it goes off. In radio frequency (RF) import mode, you can use the device 3 times without recharging. The low battery information is the same as the previous two modes.

No.	Status	Indications
		When the battery level is 0%, the buzzer will beep 3 times, the battery will be turned off and the white indicator will flash for 5 seconds. Then the power will be turned off. Device cannot be started.
6.	Battery charging	The white battery light blinks while charging.
7.	Battery charged	The white battery indicator is on.
8.	Memory function	After manual or scheduled shutdown, the last mode is defaulted on next start step.

How to use

Below is a demonstration process (picture and text below). You can select any mode to use alone or perform all steps sequentially as in "Function description".


SONIC mond

1. Before use, remove makeup and cleanse your face with a face cleanser.

2. Tightly clamp the cotton pad between the magnetic ring and the main head and moisten it with micellar liquid.

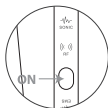3. Turn on the power and enter "SONIC mode".

4. Aim the main head at the skin of the face and move from the outside to the inside, gradually, every part of the face and neck to cleanse it thoroughly.

5. On the forehead, nose, mouth, cheeks and other parts more exposed to contamination, you can repeat the operation several times.

6. After 3 minutes, the device will turn off automatically.

7. Remove the cotton pad and clean the head with the cleaning cloth.


8. Rinse your face with clean water and apply a toner to the skin.

RF mode

1. Apply the beauty essence, serum ampoule or mask evenly on the face and neck.

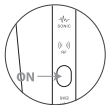2. Turn on and switch to "RF" mode.

3. Slowly move the main head along the contour of the face and apply the care evenly to the skin (do not stay in the same position on the skin for too long).

4. After 5 minutes, the device will turn off automatically.


5. Use a cloth to clean the main head.

EMS mode

1. Evenly distribute the firming care products on the face. It may be a cream that we will use to complete the care.

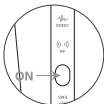2. Turn on and enter "EMS" mode.

3. Adjust the V face roller to the contour of the face, massage and lift evenly from the bottom up.

4. After 6 minutes, the product will turn off automatically.

5. Use the cleaning cloth to clean the face roller V.

Warm tips

Contraindication

- 1 Person under 18 years of age may not use the B-Booster.
- 2 It is forbidden to use the device with electronic devices inside the body. e.g. a pacemaker.
3. People who have undergone any surgical procedures or operations in the area where the procedure is to be performed during the last one. It is not recommended to use the device after exposing the skin to the sun.

4. People with damaged skin, suppuration, inflammation, fresh scars and electromagnetic sensitivity should not use the B-Booster.
5. People suffering from blood diseases, tuberculosis, celiac disease, malignant tumors, hypertension, fever, cannot use the device.
6. Do not use if you are pregnant, menstruating or breastfeeding.
7. Before proceeding with the procedure, remove watches, glasses, jewelry and other items containing metal.
8. Do not use around the eyes, mouth or throat.
9. Do not use on the same skin area for too long.
10. Do not leave this product unchecked when powered on.
11. Do not use in high temperature environment.
12. If the product is struck or dropped to the ground, check for damage before resuming use.
13. Make sure that the plug and power cord are intact to avoid the risk of electric shock.
14. Do not disassemble this device by yourself.
15. Please don't turn it on while charging.

Cleaning method

Rinse the main head and V face rollers with clean water after each use, then dry with a clean cloth. Do not use alcohol or corrosive liquids to wipe the magnetic ring.

Storage requirements

Please keep the machine in a dry and dust-free place;

Please keep away from high temperature;

Keep away from corrosive liquids;

Please keep it out of reach of children;

Do not immerse the body in water.

Reminder: If the product has not been used for a long time (6 months or more), it is recommended to use it after fully charged.

Troubleshoot

Failure phenomenon

The product does not work when turned on

Root cause analysis

- 1.The power light does not turn on (battery or motherboard failure)
2. The power light is flashing (the battery is too weak)

Solution

1. Restart after charging for 3 hours and keep turning on the light is off, please contact the after-sales service and return to factory for repair;
2. After charging for 3 hours, turn it on again.

Failure phenomenon

Poor product loading

Root cause analysis

1. The power light does not turn on (battery or motherboard failure)
2. The power light is flashing (the battery is too weak)

Solution

1. Check if the charging cable / plug is plugged in, the mode indicator flashes, it means "charging";
2. If there is no answer after filling in, please replace the product or contact the after-sales service.

Failure phenomenon

The product suddenly stops working while running

Cause analysis

1. Low battery
2. Other unknown faults may occur

Solution

1. After 3 hours of charging, restart;
2. Report to the after-sales service.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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English

Hereby, Beauty Factor Sp. z o. o. declares that the device equipment Beautifly B-Booster is in compliance with Directive 2014/30/UE.

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