PL EN DE FR ES PT IT CZ SK SL EE LT LV  
NB FI HU RO AR RU EL HR SV DA NL TR BG

# BEAUTIFLY

Scrub Slim

# Blush

series


Sonic skin scrubber

User manual

 /beautifly.global  /@beautifly.global


 /beautifly\_global  /@beautifly.global

f1

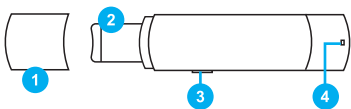f2

1

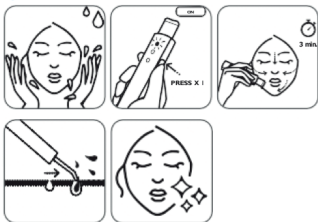2

3

f3

[facebook.com/beautifly.global](https://www.facebook.com/beautifly.global)

[youtube.com/@beautifly.global](https://www.youtube.com/@beautifly.global)

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Weight:	~43g
Size:	163*34*12mm
Rechargeable battery:	300mAh, 3.7V 6
Charging method:	Type C- USB

## Safety

Read this manual before using the product and follow the instructions in it when using the device.

## Device description [page 2 - f1]

1. Protective cover
2. Scrubber - head
3. On/Off button
4. Charge indicator

## Modes [page 2 - f2]

### 1. Peeling: Ion(+) + Ultrasonic vibration

Apply a liquid cleanser or water on your face. Switch on your device and select the mode 1 Peeling (1press). With the wave-head down, cleanse and exfoliate your skin for 3 minutes maximum. Don't press too hard, the ultrasound has the power to unclog impurities from your skin painlessly. Avoid the eye area.

### 2. Moisturizing: Ion(-) + Ultrasonic vibration

Apply your daily serum. Select the mode 2 Moisturizing (2nd press). Moisturizing your skin using the flat part of your wave-head: the ultrasonic vibrations help your cosmetic to be better absorbed. This step lasts 2 minutes maximum. Avoid the eye area.

### 3. Lifting: Ion(-) + Ultrasonic vibration

Select the mode 3 Lifting (3rd press). Slowly move the device and the flat side wave-head, from your chin up to your ears, following the contour of your face. And on the chin, from the center to the outside. For 2 minutes maximum. Avoid the eye area.

## Cautions

- Person with electronic human organ are forbidden to use it.
- Person with skin surgery within 1 month or skin sunburn is forbidden to use.
- Forbidden to use by a person with damage or inflammation of the skin, scar, or electromagnetic sensitivity.
- Take off your watch, glasses and other metal objects before use.
- Avoid shovel head close to the eyeball.
- Don't operate in the same area for a long time.
- Do not use the device for pregnant women.

## Cleaning

After using, rinse the shovel head and wipe it with a dry soft cloth.

Do not immerse the device into water for a long time.

## Storage

- Keep it dry and clean.
- Keep it away from high temperatures.
- Keep it out of reach of children.

## More information about the company and products [page 2 - f3]

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**[EN] Rules for handling electrical and electronic equipment:** Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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**[EN] Summary declaration of conformity:** The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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