
PL EN DE FR ES PT IT CZ SK SL EE LT LV
NB FI HU RO AR RU EL HR SV DA NL TR BG

BEAUTIFLY

LaserBoost

PRO

series

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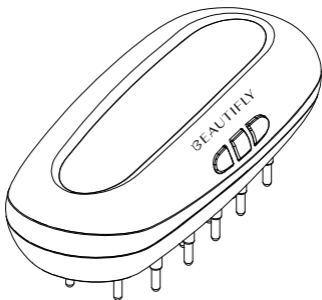Charging Type: USB-C
Supply Voltage: 3.7V
Output Power: 3.33W
Battery Capacity: 700 mAh/3.7 Wh
Charging Time: 2-3 hours
Single Mode Operating Time: 10 minutes
Standby Time: 90 days (if the device has not been used for more than 90 days, fully charge it before use)
Dimensions (W x L x H): ~136.7×58.2×44.5 mm
LED Light: 18 groups
Number of Lasers (LLLT): 6
Infrared Light Wavelength: 850±10 nm
Red Light Wavelength: 620±10 nm
Blue Light Wavelength: 465±10 nm
Laser Light Wavelength: 650±10 nm
Material: ABS + metal
Includes: device x1, USB cable x1, user manual x1, pipette x1, warranty card x1

Buttons and Indicators: [see page 2-3, f2]

1. Negative ion emission area
2. Laser light
3. EMS comb teeth
4. LED light
5. Vibration button
6. Power on/off button (level adjustment button)
7. Mode button
8. Time indicator
9. Level indicator
10. Battery capacity indicator
11. USB charging interface
12. Liquid tank (4 ml)
13. Liquid distribution teeth

Usage:

- Hold the power on/off button to enter the first level of red light mode, then press briefly to switch levels.
- Press the mode button to switch to blue light mode, then briefly press the power on/off button to adjust the level.
- Press the vibration button to turn on vibrations, and press again to turn them off.
- Hold the power on/off button to turn off the device.

Adding Hair Care Serum:

- Remove the tank from the device and lift the cover.
- Use the pipette to add oil or serum to the tank.
- Replace the cover and place the tank back in the device.

Charging:

- Make sure the device is fully charged before first use. If the battery is depleted, charge the device for at least 15-25 minutes to perform one treatment.
- Ensure the device is turned off.
- Connect the USB cable to the device and adapter (adapter not included), then insert the adapter into the outlet (the device cannot be used while charging. During charging, the battery

capacity indicator will move).

- When fully charged, the battery capacity indicator will stop moving.
- When the battery is low, the battery capacity indicator will blink three times and the device will turn off.

EN

Warnings:

1. Before using the purchased product, read the product's user manual and these safety warnings, maintenance, and storage instructions, as well as instructions for handling used electrical and electronic equipment. Follow the recommendations and use the product as intended.
2. The product is not a toy and is not intended for use by children. Only adults should use the product.
3. Do not use the product if it is not functioning, functioning improperly, or has been dropped or damaged in any way. Do not open or repair the product yourself. If the product is still under warranty, contact the seller or guarantor. If the warranty period has expired or the defect is not due to the product, have it repaired at an appropriate service center.
4. Do not use the product near heat sources, water, or moisture, and do not immerse it in water or other liquids.
5. The product is intended for use at room temperature, ranging from 10 to 40 degrees Celsius.
6. Do not expose the product to direct sunlight.
7. After use, store the product in a dry place at room temperature, ranging from 10 to 40 degrees, away from children and pets.
8. In case of dirt, disconnect the product from the power source, wipe it with a damp cloth, and let it dry before storing. Do not use abrasive or corrosive cleaning agents as they may damage the product.
9. Do not use the product if you notice a damaged power cord or plug, it emits an unpleasant smell (e.g., burning), overheats, or emits smoke. In such a situation, do not open or repair the product yourself. If the product is still under warranty, contact the seller or guarantor. If the warranty period has expired or the defect is not due to the product, have it repaired at an appropriate service center.
10. Do not leave the product unattended while it is on.
11. The power source must comply with the required product voltage.
12. After use, disconnect the product from the power source and let it cool before moving or storing it.
13. Do not twist or coil the power cord as this may damage it.
14. Accessories not from the manufacturer may cause the product to malfunction or damage it.
15. Use the device with eye protection (goggles). Do not direct the light emitted by the device towards the eyes.
16. The device is not intended for use on children.
17. Do not use the device if you are pregnant or breastfeeding.
18. The device should not be used by people with past or present:
 - skin cancer;
 - numerous moles or other skin changes;
 - collagen or vascular skin disorders, e.g., keloid scars, wound healing problems, varicose veins or vascular ectasia;

- skin sensitivity to light;
 - rash or allergic reactions on the skin where the treatment is to be applied;
 - infections, eczema, burns, folliculitis, open cuts, abrasions, herpes, wounds or injuries, or bruises where the treatment is to be applied;
 - surgical procedures where the treatment is to be applied;
 - epilepsy with sensitivity to flashing light;
 - diabetes, lupus, porphyria, or heart failure;
 - blood clotting disorders;
 - immune-lowering diseases (including HIV or AIDS infection).
19. Do not use the device in the following places:
- on the face, particularly around the eyes and near the eyebrows;
 - on nipples and areolas, labia, scrotum, or around the vagina or anus;
 - inside the nostrils and ears;
 - on or near areas containing artificial objects, such as silicone implants, pacemakers, subcutaneous injection sites (insulin dispensers), or piercings;
 - on moles, freckles, darker skin areas, scars;
 - on warts, tattoos, and places with permanent makeup;
 - where long-acting deodorants are used;
 - in sunburned areas, recently tanned areas (within the last 48 hours), and areas where self-tanners are used.
20. After the treatment, for at least 48 hours, do not expose the skin to direct sunlight.

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[EN] Used electrical and electronic equipment, as well as batteries and accumulators, must not be mixed with other waste. This is indicated by the symbol of a crossed-out bin placed on the product packaging. Electrical and electronic equipment, as well as batteries and accumulators, may contain hazardous substances, mixtures, and components harmful to the environment and human health. Improper handling or damage to them can lead to health risks or environmental pollution. Proper handling of used electrical and electronic equipment, as well as batteries and accumulators, in households positively contributes to reuse and recovery, including recycling of the used equipment. Used electrical and electronic devices, as well as batteries and accumulators, can be returned at local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and served the same functions as the equipment sold. When delivering equipment intended for households, the distributor is obliged to collect used equipment from households free of charge at the delivery location, provided that the used equipment is of the same type and served the same functions as the delivered equipment. A distributor operating a retail unit with a

sales area, as defined in Article 2, point 19 of the Act of March 27, 2003, on spatial planning and development, of at least 400 m² dedicated to the sale of equipment intended for households, is obliged to accept used equipment from households free of charge at this unit or in its immediate vicinity, provided that none of the external dimensions exceeds 25 cm, without the necessity of purchasing new equipment intended for households. The product was placed on the market after August 13, 2005.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the below product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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