BEAUTIFLY

RejuvaHair

PRO

series

Hair growth device
(LLLT)

User manual

 /beautifuly.global  /beautifuly_global

 /@beautifuly.global  /@beautifuly.global

f1


facebook.com/beautifly.global

youtube.com/@beautifly.global

f2

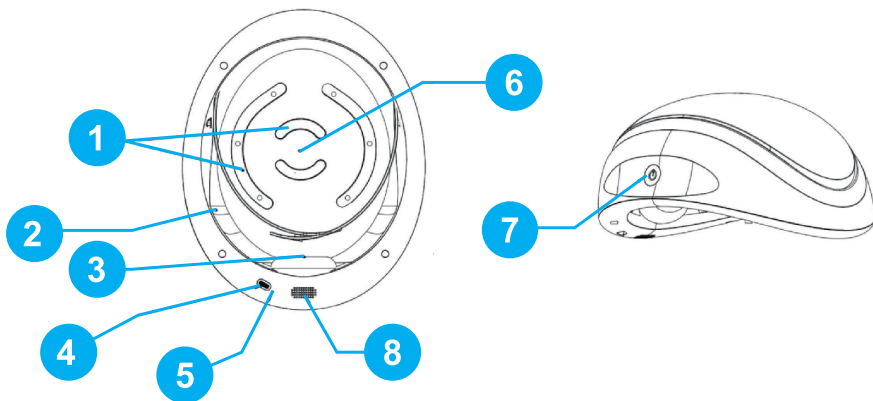f3

I - 1

I - 2

I - 3

I - 4

II - 1

II - 2

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Size: ~26cm x 21cm x 12.5cm

Net weight: ~572g

Charging time: approximately 3-4 hours (device operates for 70 minutes on a single charge)

Indication for use (see page 3, f3)

The hair growth device is indicated for aiding hair growth in women with androgenetic alopecia classified as Ludwig-Savin I-II, in men with androgenetic alopecia classified as Norwood-Hamilton IIa-V, and for both, with Fitzpatrick skin phototypes I to IV. Green areas show the type of hair loss that can be treated with the hair growth device.

1. Ludwig classifications for female pattern hair loss
2. Norwood Hamilton classifications for male pattern baldness
3. Fitzpatrick classification of skin types

Main switch light indicator

Solid red light > Error

Flashing red > Low battery level

Solid blue light > Fully charged

Flashing blue > Charging

Usage time

After turning on the device, place it on the head to start the laser session. The session will automatically end after 30 minutes. If the device is removed during the session, it will automatically stop and enter standby mode for 10 minutes.

Preparation

1. Before starting the session, clean the scalp and dry the hair.
2. For users with short hair (length should not exceed 2 cm), ensure full

exposure of the scalp in the hair loss area.

3. For users with M-type hair loss who retain hair on the top of the head, pull the hair back and secure it with pins to fully expose the scalp in the affected area.
4. For users with M-shaped hair loss and hair loss on the top of the head, pull the hair on both sides back and secure it with pins to expose the scalp in the hair loss area and prevent it from being covered again with hair combed back.
5. For users with hair loss only on the top of the head, try to move the hair away from the top of the head to expose the scalp in the hair loss area.
6. Part the hair to expose the scalp in the bald area, then put on the device.

Turning on and setting the language (see page 4, f4)

1. Turn on the device using the power button located on the side of the device.
2. Press and hold the power button for 1 second to turn on the device. You will hear a voice prompt, and the device will enter standby mode. Within 2 seconds, touch the switch five times to select the language (Chinese or English).
3. The control light will turn blue when the device is turned on.

Correct helmet wearing (see page 4, f5)

1. Switch on the right side > Correct wearing.
2. Heat dissipation holes at the back > Correct wearing.
3. Switch on the left side > Incorrect wearing.

Conducting sessions

After turning on the device and properly placing it on the head, the device will light up and emit a voice prompt. Then, the device will start the light session and automatically turn off after 30 minutes. Each laser therapy session should last 30 minutes, once every two days, for 16 weeks.

Charging

For the first use of the device, if the control light flashes red, it means the device needs to be charged before starting work. Connect the USB cable to the charging adapter, and you will hear a voice prompt. After 3-4 hours of charging, when the red light stops flashing and turns blue, the battery is fully charged. Note: The device will not function properly while charging, do not use the product during charging.

Correct connection scheme (see page 4, f6)

1. Charging port, connection with USB power cable
2. Connect the USB power cable to the charging port
3. Connect the other end of the USB power cable to the power adapter

Laser session process

1. If the device falls during the laser session process, the light will automatically stop, and it will enter pause mode. The device will emit a voice prompt, and the timer will stop.
2. After properly putting the device back on, the device will emit a message and resume operation. The timer will also continue counting.
3. If the device enters standby mode and remains inactive for 10 minutes, it will automatically turn off.
4. When the time reaches 10 minutes, the device will emit a message.
5. When the time reaches 20 minutes, the device will emit a message.
6. When the time reaches 25 minutes, the device will emit a message.
7. When the time reaches 30 minutes, it will emit a message.

Ending the session

1. After the 30-minute session, the device will emit a single beep indicating the end of the session. Then, it will emit a message and automatically turn off.
2. Press and hold the power button for 1 second. The device will emit a message and automatically turn off.

Warnings

1. Before using the product, read the user manual and these safety warnings, maintenance and storage instructions, as well as the rules for handling waste electrical and electronic equipment. Follow the recommendations contained therein and use according to its purpose.
2. The product is not a toy and is not intended for use by children. Only adults may use the product.
3. Do not use the product if it is not working or is malfunctioning, or if it has been dropped or otherwise damaged. In such a situation, do not open or repair the product on your own. If the product is still under warranty or guarantee, contact the seller or guarantor. If the warranty or guarantee period has expired, or the defect did not arise from a cause inherent in the product, take the product to an appropriate service center for repair.
4. Do not use the product near heat sources, water, or moisture, and do not immerse it in water or other liquids.
5. Do not expose the product to direct sunlight.
6. After using the product, store it in a dry place at room temperature, out of reach of children and animals.
7. Do not use the product if you notice a damaged power cord or plug, an unpleasant odor (e.g., burning), overheating or smoke emission. In such a situation, do not open or repair the product on your own. If the product is still under warranty or guarantee, contact the seller or guarantor. If the warranty or guarantee period has expired, or the defect did not arise from a cause inherent in the product, take the product to an appropriate service center for repair.
8. Do not leave the product unattended when it is turned on.
9. The power source must comply with the product's required voltage.
10. After use, disconnect the product from the power source and allow it to cool before moving or storing it.
11. Do not twist or coil the power cord, as this may damage it.
12. Accessories for the product not coming from the manufacturer may cause the product to malfunction or become damaged.
13. Use the device with eye protection (safety glasses). Do not direct the light emitted by the device towards the eyes.
14. The device is not intended for use on children.
15. Do not use the device if you are pregnant or breastfeeding.
16. The device should not be used by individuals who have or have had in the past:
 - skin cancer;
 - numerous moles or other skin changes;
 - collagen or vascular disorders of the skin, e.g., cases of keloid scarring, wound healing problems, varicose veins or telangiectasia;
 - skin sensitivity to light;
 - rashes or allergic reactions on the skin in the area to be treated;
 - infections, eczema, burns, folliculitis, open cuts, abrasions, herpes, wounds or injuries or hematomas in the area to be treated;
 - surgeries in the area to be treated;
 - epilepsy with sensitivity to flashing light;
 - diabetes, lupus erythematosus, porphyria or heart failure;
 - blood clotting disorders;
 - diseases lowering immunity (including HIV or AIDS infection).

17. Do not use the device in the following areas:

- on the face, particularly around the eyes and near the eyebrows;
 - on nipples and areolas, labia, scrotum or near the vagina or anus;
 - inside the nostrils and ears;
 - on areas or near areas containing artificial objects such as silicone implants, pacemakers, subcutaneous injection sites (insulin dispensers) or earrings;
 - on moles, freckles, darker areas of skin, scars;
 - on warts, tattoos and areas with permanent makeup;
 - in areas of long-acting deodorants;
 - on sunburned areas, recently tanned areas (within the last 48 hours) and areas with self-tanners.
18. Do not expose the skin to direct sunlight for at least 48 hours after the session.

More information can be found > see page 2, f1

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[EN] Used electrical and electronic equipment, as well as batteries and accumulators, must not be mixed with other waste. This is indicated by the symbol of a crossed-out bin placed on the product packaging. Electrical and electronic equipment, as well as batteries and accumulators, may contain hazardous substances, mixtures, and components harmful to the environment and human health. Improper handling or damage to them can lead to health risks or environmental pollution. Proper handling of used electrical and electronic equipment, as well as batteries and accumulators, in households positively contributes to reuse and recovery, including recycling of the used equipment. Used electrical and electronic devices, as well as batteries and accumulators, can be returned at local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and served the same functions as the equipment sold. When delivering equipment intended for households, the distributor is obliged to collect used equipment from households free of charge at the delivery location, provided that the used equipment is of the same type and served the same functions as the delivered equipment. A distributor operating a retail unit with a sales area, as defined in Article 2, point 19 of the Act of March 27, 2003, on spatial planning and development, of at least 400 m² dedicated to the sale of equipment intended for households, is obliged to accept used equipment from households free of charge at this unit or in its immediate vicinity, provided that none of the external dimensions exceeds 25 cm, without the necessity of purchasing new equipment intended for households. The product was placed on the market after August 13, 2005.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances. The full text of the declaration of conformity is available at the following internet address:
www.beautify.eu/deklaracja

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