PL EN DE FR ES PT IT CZ SK SL EE LT LV
NB FI HU RO AR RU EL HR SV DA NL TR BG

BEAUTIFLY

AuraSkin

PRO

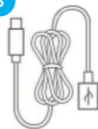
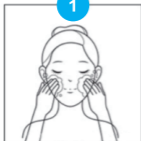
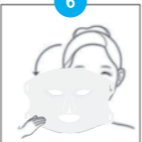
series

LED Mask
User manual

[f](#) /beautifuly.global [i](#) /beautifuly_global

[v](#) /@beautifuly.global [d](#) /@beautifuly.global

f1**1****2****3****4****5****f3****1****2****3****4****5****6****f2****1****2****3****4****f4**

[youtube.com/@beautiflyglobal](https://www.youtube.com/@beautiflyglobal)

[facebook.com/beautiflyglobal](https://www.facebook.com/beautiflyglobal)

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Moc znamionowa: 4W

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Materiał: Silikon + PP

Ilość diod LED: 70 lamp, 280 diod

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Number of LEDs: 70 lamps, 280 diodes

Luminous power: 30–50 cm²/mW

Remote control battery capacity: 1500 mAh

Operating time on a single charge: 150 minutes

Charging time: 120 minutes

Light modes: red 630nm (RL), near infrared 850nm (NIL), blue 465nm (BL), yellow 580nm (YL), green 540nm (GL), orange 610nm (OL), purple 630/465nm (PL), cyan 485nm (IL).

Safety:

Before using the product, please read this user manual and follow the recommendations contained herein during the operation of the device.

Contents of the set [page 2 - f1]:

1. Mask. / 2. Remote control. / 3. USB cable. / 4. Protective glasses. / 5. Strapping band.

Description of the remote control [page 2 - f2]:

1. Time display – shows the set operating time. / 2. Mode display – shows the currently selected light mode. / 3. Mode button – changes the light mode. / 4. Power/time setting button – turns the device on/off and sets the time.

Operating the device [page 2 - f3]:

1. Remove makeup and thoroughly cleanse your face. / 2. Put on the protective glasses. The glass marked with the letter "L" is worn on the left eye, with the letter "L" facing outward. The glass marked with the letter "R" is worn on the right eye, with the letter "R" facing outward. Attach the strapping bands. The bands do not distinguish between left and right sides. With the soft material (flannel) facing inward, thread two bands with small Velcro through the holes on the left and right sides of the mask. The Velcro is attached to the tape – the position can be adjusted left or right as needed. After putting on the mask, the large parts of the Velcro can be joined together. / 3. Connect the mask to the remote control. Put on the mask. / 4. Turn on and select the light mode. / 5/6. After use, remove and clean the mask.

Device Modes and Their Effects:

1. Red (RL, 630 nm) - stimulates collagen production, improves elasticity, and supports skin regeneration.
2. Near Infrared (NIL, 850 nm) - penetrates deeper layers of the skin, supporting metabolism and cell repair.
3. Blue (BL, 465 nm) - has antibacterial properties, soothes inflammation, and reduces the appearance of imperfections.
4. Yellow (YL, 580 nm) - brightens the complexion, evens out skin tone, and soothes sensitive skin.
5. Green (GL, 540 nm) - reduces redness and evens out skin tone.
6. Orange (OL, 610 nm) - stimulates microcirculation, restores energy and a natural radiance to the skin.
7. Purple (PL, 630/465 nm) - combines the effects of red and blue, supporting regeneration and caring for acne-prone skin.
8. Cyan (IL, 485 nm) - soothes and refreshes the complexion, improving oxygenation and the overall appearance of the skin.

How to use the device [page 2 - f2]:

1. Power/time setting button:

- Press and hold the power button to turn the device on/off.
- When turned on, the default mode is red light (RL) – the RL indicator lights up, and the default operating time is 5 minutes.
- A short press of the on/off button allows you to change the operating time to 5 or 10 minutes. The time can be set separately for each light mode.

2. Mode button:

- When the device is turned on, each subsequent press of the mode button changes the active light and activates the corresponding indicator (LED light).
3. Mode/power indicator:
- When the remote control is functioning correctly, the corresponding mode indicator lights up.
 - If the battery voltage drops below 3.4V, the indicator for that mode starts to flash.

Example:

If the red light mode (RL) is active and the voltage drops below 3.4V, the RL indicator starts to flash. When the voltage reaches 3.2V, the device automatically turns off.

Charging, maintenance, and care of the device:

The product does not support operation while charging. When the indicator is flashing, the device is in charging mode. When three brightness indicators light up simultaneously and stop flashing, charging is complete. After using the mask, wipe the silicone part with a dry or damp towel. Dry the mask in a well-ventilated area and store it in a dust-free place. Do not allow direct contact with water.

Precautions:

Keep the product away from small children. Handle the product carefully during use. It should not be used by pregnant women, those preparing for pregnancy or breastfeeding, individuals sensitive to light, or those with genetic or skin diseases (e.g., lupus erythematosus, light eczema, albinism, etc.). Do not use while charging. Do not press the device excessively to avoid deformation and malfunction. Do not repair, disassemble, or modify the device on your own. Do not drop or hit the device to avoid malfunction. Note: the remote control and cables should not come into contact with water. After each treatment, wipe the product with a damp cloth. Do not wash the device under water or immerse it. Do not use corrosive cleaning agents on the device. Store it in a cool, dry place. Regardless of the safety measures used, redness or irritation may occur on the skin due to the device's operation. In such cases, immediately stop the treatment and consult a doctor. Frequent use of the product beyond the recommended frequency may lead to irritation or inflammatory skin conditions. For at least 48 hours before using the device and for at least 48 hours after its use, do not use tanning beds or expose the skin to sunlight. The device should not be used by individuals under 18 years of age. Do not use on skin areas with an increased likelihood of skin cancer, on open wounds or rashes, on swollen, red, irritated, infected, or inflamed skin areas or skin eruptions. Use eye protection when using IPL or laser devices to avoid exposure to direct or reflected radiation. It is recommended to consult a doctor before using the device and, if necessary, to perform diagnostic tests on the area of the body undergoing treatment. The treatment should not be performed before obtaining documentation from such a consultation. Recommended use: 3–5 times a week, for 4 consecutive weeks as one cycle.

More information about the company and products [page 2 - f4]:

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[EN] Rules for handling electrical and electronic equipment: Used electrical and electronic equipment, batteries, and accumulators should not be disposed of with other waste. This is indicated by the crossed-out waste bin symbol on the product. Electrical and electronic equipment, batteries, and accumulators may contain accessories, components, and parts that are harmful to the environment and human health. Improper handling or their consequences may result in damage to the environment or the environment. Proper household handling of used electrical and electronic equipment, batteries, and accumulators is essential for the proper disposal of devices and equipment, including used electrical and electronic equipment. Used electrical and electronic equipment, batteries, and accumulators can be disposed of separately. The distributor is obligated to collect household equipment free of charge at the point of sale, provided the used equipment is of the same type and performs the same functions as commercial equipment. A distributor delivering household appliances to a purchaser is obligated to collect waste household appliances free of charge at the point of delivery, provided the waste appliances are of the same type and have served the same function as the appliances. A distributor operating a retail business with an online sales area of 400 m² of household appliances designated for sale pursuant to Article 2, Section 19 of the Act of 27 March 2003 on Spatial Planning and Development is entitled to free membership in the same entity or its immediate vicinity for discarded household appliances with external dimensions not exceeding 25 cm, without the need to purchase new household appliances. The product was registered on the market after August 13, 2005.

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