

innoGIO



GIOperfect Clean

MULTIFUNCTIONAL DEVICE
FOR MASSAGE AND SKINCARE

USER MANUAL

Introduction

Congratulations on taking the first step towards perfect, beautiful skin and buying the InnoGIO GIOperfect Clean GIO-710, a sonic facial cleansing face mask applicator which makes you feel like you're in a spa.

Before experiencing the full benefits of this advanced technology at home, please read carefully the following directions about using this device.

Read the manual before using the device and make sure you use it in line with its purpose.



WARNING: do not modify the device. Keep the device and all accessories beyond children's reach.

Storage precautions

- Proper maintenance extends the life of the device.
- Clean the device thoroughly after each use. Use soap and water, then rinse with warm water. Dry with a dust-free cloth or towel. Do not use alcohol or products which irritate the skin, such as gasoline or acetone.
- Do not expose the device to direct sunlight, high temperatures or potential sources of fire.
- Due to hygienic reasons, it is recommended that this device is used by one person only.
- Do not try to disassemble the device yourself. Store the device in appropriate conditions.

Precaution

- The device is designed for skin care. It is intended for cosmetic purposes only.
 - Do not use in case of burns, post-traumatic infections or other skin problems.
 - Use only the accessories included with the device. Using other accessories may pose a danger.
 - Do not use in the case of skin problems.
 - Discontinue use if the device is damaged.
 - Do not try to disassemble the device yourself.
- In the event the device breaks or you encounter problems using it, it is recommended to send it to a designated service centre to detect the malfunction and repair it.

Storage


If the device is not used for a long time, switch it off and store in a dry place. Avoid excessive exposure to sunlight and high temperatures.

Description of the GIO-710 Facial Cleansing Brush

A breakthrough in facial cleansing. A modern sonic brush for facial skincare provides deep cleansing while improving the condition of your skin.

The brush vibrates at a frequency of 8500 pulse per minute, which makes the silicone tips on the device vibrate. Thanks to the vibrations, the micro tips cleanse, peel and regenerate the skin with great precision. The device is made of food-grade silicone gel. It's non-sensitising, non-toxic and harmless.

The device modes allow for adjusting the frequency of vibration and heat generation. Thanks to GIOperfect Clean, your skin will be clean, radiant and regenerated. It removes dead skin cells and evens out the skin tone

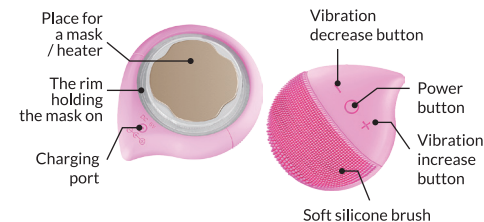

Using the GIO-710 Facial Cleansing Brush

Skin cleansing (recommended once a day):

1. Press the power button. Activate the cleansing mode and adjust the desired intensity.
2. Moisten your face with water and apply any cleansing product.
3. Massage the entire facial skin with the sonic brush by making circular movements.
4. Treatment duration: between 1 and 3 minutes.
5. Wash the face with water and then apply your preferred oil, serum or cream (do not apply skincare products which dissolve in silica, have an exfoliating or peeling effect, on the surface of the device.)
6. Done! Enjoy a fresh and cleansed face.

Soothing massage mode:

After daily facial cleaning, a soothing massage can comfort your skin. Recommended once a day. The brush comes with a heating mode (Thermotherapy).

1. After cleansing your face, press the power button, turn on the soothing massage mode and adjust the desired intensity.
2. The red multi-spectral LED light will light up for 10 seconds. Then, during the massage, the blue and red multi-spectral lights will flash alternately, making your skin more beautiful and radiant.
3. Treatment duration: between 1 and 3 minutes.
4. Done!

Face mask mode:

1. Take off the ring attached to the device. Apply the face mask on the surface of the device and fix it with the protective ring.
2. Press the power button 3 times until the blue multi-spectral light is turned on to activate the face mask mode. Adjust the desired intensity.
3. Massage the entire facial skin with the sonic brush by making circular movements.
4. The procedure takes 90 seconds. Rinse off the face mask, wash the device with clean water and dry it.

Maintenance

Clean the device thoroughly after each use. Use soap and water, then rinse with warm water. Dry with a dust-free cloth or towel. Do not use alcohol or products which irritate the skin, such as gasoline or acetone.

Charging

CAUTION: Before charging, make sure the charging port and charger are dry and there is no product residue on any part of the device. **DO NOT** charge GIOperfect Clean in places which pose a risk of the device or the charger coming in contact with water. **DO NOT** use the device while it's charging and **DO NOT** leave it connected for more than 24 hours.

- Connect the USB cable to the charging point. While the device is charging, the LED light will flash. Charging takes about 3 hours. When the battery is full, the LED light will flash continuously.
- When the battery is almost flat, the LED light will flash.

Troubleshooting

- *GIOperfect Clean does not turn on when the power button is pressed.*
The battery is flat. Charge the device.
- *GIOperfect Clean is locked.*
Press "+" and "-" at the same time to unlock the device.
- *GIOperfect Clean cannot be turned off and/or the device buttons are not responding.*
The microprocessor is temporarily malfunctioning. Plug in the charging cable to restore regular operation.
- *The LED light is not flashing when GIOperfect Clean is connected to the charger.*
The battery is full. The device is ready to use.
The battery is flat. It takes a few minutes for the LED light to flash while the device is charging.

The charging cable is not connected correctly. Check the electrical outlet/USB port and charging port.

- *The device fell into water and the charging port is wet.*
Dry the device and continue using it.

Correct disposal of the device

Disposal of waste electronic appliances (applicable in the European Union and other European countries with separate collection systems.) The crossed-out wheellie bin symbol indicates that the device must not be treated as typical household waste. It should be delivered to an appropriate collection point intended for recycling of electrical and electronic appliances.



Product specification

- Product name: GIOperfect Clean – Sonic device intended for face cleansing and application of face masks
- Model: GIO-710
- Charging voltage: 3-4.2V
- Charging time: 3h
- Usage time: 2.2Hrs
- Battery capacity: 600mAh
- Rated power: 1W
- Engine speed: 3.7V 8500rpm
- Battery voltage: 3.7V
- Input voltage: VDC5V/1.0A
- Dimensions: 9,2 x 7,7 x 3,2 cm

Disclaimer

Users of the brush take responsibility for its use. Neither the supplier nor the seller are liable for any injuries or damages, physical or otherwise, arising directly or indirectly from the use of this device.

Importer:
Innogio Sp. z o.o.
Ostródzka street 74H, 03-289 Warsaw, Poland

innogIO

Authorised Representative in Europe

www.innogio.com