





medivon.

# User manual

COLLAR MASSAGE

PL | EN | DE | FR | ES | PT | IT | CZ | SK  
SL | EE | LT | LV | HU | RO | RU | NB | FI  
EL | HR | SV | DA | NL | TR | BG | AR

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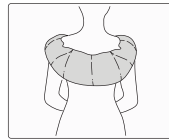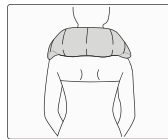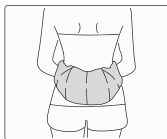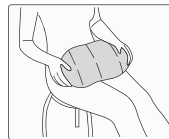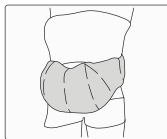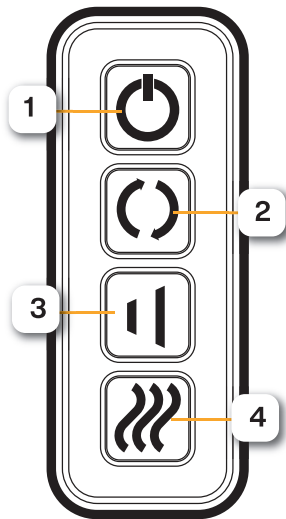### Controller support:

#### 1. Power button

Press to turn on the massager, press again to turn off the massager.

#### 2. Directional button

Press while the massager is operating to change the direction of the massage heads.

#### 3. Intensity button

Press while the massager is operating to set the intensity (high or low).

#### 4. Heating button

Press while the massager is running to turn the heating function on or off.

Note: When you turn on the massager, the default massage mode will start. The massage heads will rotate clockwise at a high intensity level along with a warm-up function.

**The device will turn off automatically after 15 minutes.**

## Safety instructions

The massager meets safety rules and technical requirements. Do not wet the device or remove the cover. The massager is not a toy, therefore children and disabled people must be under adult supervision when using it. Do not leave the massager unattended after it is plugged in. Repairs may only be performed by authorized personnel.

Improper use and unauthorized repairs may void the warranty and compromise safety. Avoid contact of the device with water, high temperatures and sunlight. Use only undamaged cables, plugs and sockets. Do not cover the ventilation opening of the device. If the plug or cord is damaged, consult the manufacturer or a specialist. Immediately disconnect the device from the network in the event of a failure. Do not use the massager if you have skin problems, wounds, swelling or inflammation.

The recommended massage time is no more than 15 minutes at a time, after this time the device will turn off automatically. Avoid massage during pregnancy in case of bruises, skin damage, inflammation and swelling. Do not use the massager while driving or sleeping.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. This device was placed on the market after 13.08.2005.

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The manufacturer hereby declares that the Collar Simple Black device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: [www.medivon.pl/deklaracje](http://www.medivon.pl/deklaracje)

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska  
Model: Collar Simple Black