
medivon.

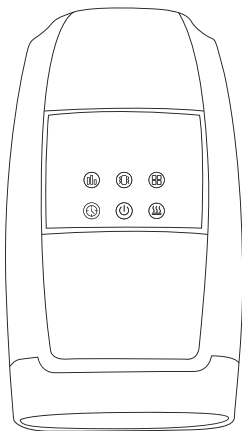User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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Device Manual

The device is very easy to use - just one press of a button to start a relaxing hand massage that will greatly improve your blood circulation and relax tense hand muscles after a hard day's work.

Put your hand into the massager and then hold the POWER button for about 2 seconds to start. The device can be turned off by pressing the POWER button again.

Changing the massage intensity - to change the massage level, just press the appropriate button on the left.

The intensity will cycle from L3 to the strongest L6 (L3 - L4 - L5 - L6 - L1) etc.

Changing the massage mode - the massager has 6 automatic massage modes defined. By pressing the MODE button, we can select the cyclically changing programs, the one that currently suits us best.

The vibration option is intended for people who appreciate even better effects of relaxation and improved circulation in the hands. using the vibration button, we can change their strength or turn it off completely.

Heating function - built-in heads act as a warming compress, this function is activated using the heating button on the right side of the control panel. After starting the warm-up function - the temperature in the device will start to rise gradually.

The time of the massage should be adjusted to your own preferences. The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems should consult a doctor before use.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. The device was placed on the market after 13.08.2005.

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The manufacturer hereby declares that the Medivon Hand Pro device complies with Directive 2014/30 / EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65 / EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the declaration of conformity is available at the internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Medivon Hand Pro