medivon®



Medivon Portia Arm

Medivon Portia Leg

Medivon Portia Center

User manual

PL | EN | DE | ES | IT | FR | CZ
SL | PT | LT | RU

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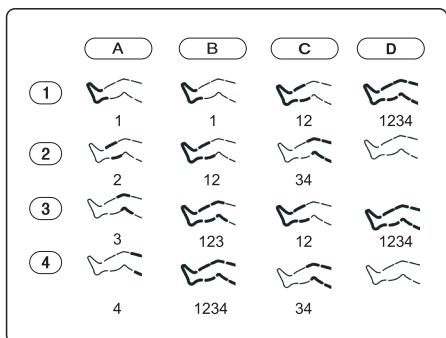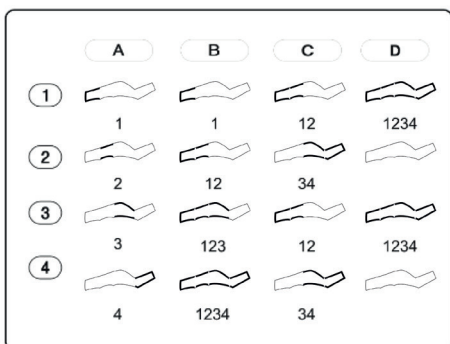Device manual

The device is ready for immediate use after purchase. Hold the ON/OFF button for about 2 seconds to start the device.

Before starting the device, connect the pump tubing to the control unit of the device, paying attention to the small recess in the tip to indicate the direction of connection. If you are using only one pumping element, use a plug to cover the other pumping port to prevent air from escaping. The power socket is located on the side of the main unit, use the power supply included in the set.

After starting the device, repeated pressing of the ON/OFF button changes the massage areas (1-2-3-4). Use the MODE button to change modes (A-B-C-D). The connection diagram of the modes and the principle of their operation are shown in the picture above. Using the TIMER button, we can freely set the sensor after which time the device will turn off automatically. The INTENSITY button changes the strength and intensity of the massage successively.

After finishing the massage, the device can be turned off by pressing the ON/OFF button again.

The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems, pacemakers, cancer or pregnant women should consult a doctor before use.

The massage area combination diagram is presented on the first page of this manual. The original Medivon pumping station is required for the correct operation of the accessory.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

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The manufacturer hereby declares that the Medivon Portia, Medivon Portia Arm, Medivon Portia Leg and Medivon Portia Center devices complies with Directive 2014/35 / EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits , Directive 2014/30 / EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65 / EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the declaration of conformity is available at the internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska

Model: Medivon Portia, Medivon Portia Arm, Medivon Portia Leg oraz Medivon Portia Center