
medivon.
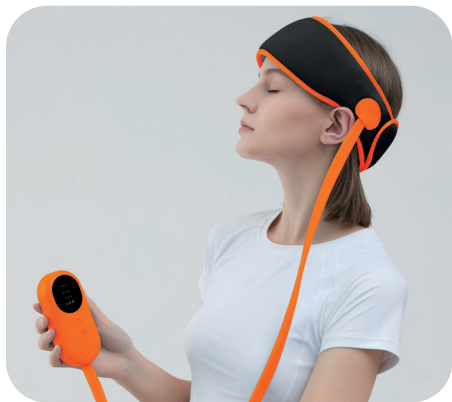
# Head Massager

Professional pressure therapy

PL | EN | DE | ES | IT | FR | CZ  
SL | PT | LT | RU

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## 1. Device operations

The device is very easy to use - just press a button to start a relaxing head massage that will allow you to relax after a hard day at work or help you fall asleep and reduce the effects of insomnia. The device is turned on using the ON/OFF button (4). Hold the button for about 3 seconds to start the massager. Before use, it is recommended to fully charge the device's battery using the USB cable. It is recommended not to fully discharge the device to extend the battery life. It is forbidden to use the massager while charging the battery. This may cause damage or fire. Before starting the massage, place the massager on your head as shown in the illustrations at the beginning of this user manual.

Improper positioning of the massager on the head during operation may cause discomfort, pain or the device may fall on the eyes, risking injury. Before starting the massage, make sure that the massager fits well to the head and does not cause discomfort. The massager has elastic Velcro that allows you to adjust the fit.

## 2. Massage and settings

After placing the massager on your head and turning it on, select the appropriate program tailored to your preferences. The device has 3 operating modes shown on the display (1). To change them, simply press the ON/OFF button (4) repeatedly.

The individual modes are: Meridian Cleanse | Reducing pain and fatigue | Acupressure

The device has a warm-up function. Just press the heating button (3) cyclically to change the heating strength on a 4-level scale, where the 4th level means turning off the heating. To turn off the device, press the ON/OFF button (4) again for a few seconds.

## 3. Dangers and warnings

Do not use spare parts that are not approved by the manufacturer. Non-professional service personnel should not disassemble and replace parts of the leg massager without permission.

If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid danger. Do not use or store the product in a place with high humidity or high temperature. After the massage is finished, unplug the power cord. Do not use the product for more than 15 minutes at a time. If any side effects occur, stop using the massager immediately. If the temperature becomes too high or the product begins to smell or emit a burning odor, stop using it immediately and unplug the power. Do not move the product by pulling the power cord to prevent the power cord from coming loose.

This product is not intended for use by persons who lack physical, sensory or mental capabilities or lack of experience and knowledge (including children).

People who are not sensitive to heat should use the heating function carefully - preferably under the supervision of a caregiver. If you have any of the following conditions, consult your doctor before use: Pregnant women, children or elderly people, People suffering from skin infections or having skin redness and burns, who use a pacemaker or similar implanted medical device; Patients with heart disease or malignant tumors etc, People who are under treatment or feel discomfort, People suffering from osteoporosis, People with high blood pressure or abnormal blood pressure.

## 4. Cleaning and care

Before cleaning, turn off the product and unplug the power cord.

Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with corrosive detergents such as benzene, thinners. Store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environments, otherwise it may cause colors to fade or damage the product.

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The device is very easy to use - just press a button to start a relaxing head massage that will allow you to relax after a hard day at work or help you fall asleep and reduce the effects of insomnia. The device is turned on using the ON/OFF button (4). Hold the button for about 3 seconds to start the massager. Before use, it is recommended to fully charge the device's battery using the USB cable. It is recommended not to fully discharge the device to extend the battery life. It is forbidden to use the massager while charging the battery. This may cause damage or fire. Before starting the massage, place the massager on your head as shown in the illustrations at the beginning of this user manual. Improper positioning of the massager on the head during operation may cause discomfort, pain or the device may fall on the eyes, risking injury. Before starting the massage, make sure that the massager fits well to the head and does not cause discomfort. The massager has elastic Velcro that allows you to adjust the fit.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

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The manufacturer declares under his sole responsibility that the HeadHug device complies with Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, as well as Directive 2011 of the European Parliament and of the Council /65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to the Directive of the European Parliament and 2011/65/EU as regards the list of restricted substances. The full text of the declaration of conformity is available at the internet address

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska  
Model: HeadHug