
medivon.

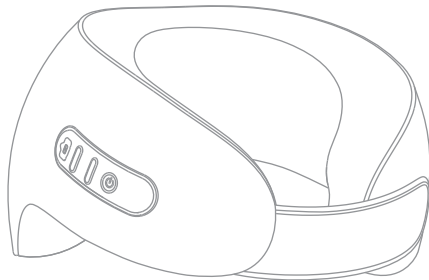Neck Massager

Deep relaxation

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1. Operation and application

Press the power button (5) on the control panel (1) for 2 seconds to start the massager and enter the operating mode. After putting on the massager, you can fasten it with elastic Velcro (2) to better fit the shape of your neck. Briefly press the mode change button (3) to adjust the massage speed, choosing between high and low speed. Briefly press the infrared heat button (4) to activate or deactivate the warm-up function. To stop using the device, long press the power button (5) for 2 seconds to turn off the product.

While charging, a white light flashes around the button, and when charging is complete, the light stays on constantly.

The massager is powered by a built-in lithium battery. If the battery runs out or is worn out, please charge the massager before using it again. Only use the included charging cable to charge, check if the power connection is correct and connect to the Type-C port. Voltage instability may occur during charging, so do not use the product while charging. The accessories included with the product are intended exclusively for this model. Accessories from other sources are not recommended. The product usage time is set to 10 minutes. After this time, the device will automatically turn off. To continue using, turn your device back on.

2. Dangers and warnings

People with medical implants, undergoing treatment, heart disease, high blood pressure, acute illness, pregnancy, osteoporosis, spine fractures, skin damage, fever above 38°C and people who are insensitive to heat should consult a doctor before use. Minors and people with intellectual disabilities may use the product only under supervision. The product is intended for use by one person, do not use it with others. After moving the product from a cold to a warm place, wait an hour before using it. Recommended use time is 10 minutes. Do not handle the product with wet hands and do not use it outdoors. Only use a power supply with the appropriate voltage. Do not smoke while using the product. Please use a qualified USB adapter for charging. Do not use in a humid environment. Do not use the product immediately after meals or simultaneously with other therapeutic devices. If there are any abnormalities during use, discontinue use immediately and consult a doctor. If the power cord is damaged, contact your distributor for replacement. When charging is complete or when not in use, unplug the plug from the outlet. Do not use the product while standing on it or damage the power cord. Do not use on an electric heating blanket or space heater. Do not disassemble the product unless you are an authorized service technician. In the event of a power outage, immediately turn off the product and unplug it from the power supply. Do not use the product while sleeping, intoxicated or feeling unwell. If there are any abnormalities (smoke, unpleasant odor, etc.), discontinue use immediately. Do not clean the product with gasoline, solvents or insecticides. In case of water spill, do not use the product. Do not use during a lightning storm. Before disposing of the product, please remove the battery. Keep these tips for safe and effective use of the product.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

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This device was placed on the market after 13.08.2005.

The manufacturer hereby declares that the Medivon Cerv device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska

Model: Medivon Cerv