


EN

GIO-711

innoGIO

# GIO perfect Care

Sonic LED facial care device with EMS therapy

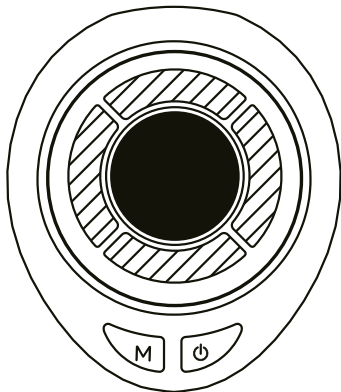

**USER MANUAL**

To avoid harm caused by misuse, please read the instructions carefully and use this product strictly in accordance with the manual.

## TABLE OF CONTENTS

Product information	1
Product structure	2
Safety Precautions	3
Operating Instructions	4
How to attach a cotton pad	4
How to use	5
Charging instructions	5
Disclaimer	5



## PRODUCT INFORMATION

### 1. Sound Wave Vibration Beauty Principle

This product uses micro-vibration and light wave therapy technologies, combined with a smart design that fits the contours of the face to deeply cleanse the skin without missing any areas. It delivers approximately 8000 high-frequency vibrations per minute to dislodge deep-seated oils and impurities from the pores, leaving the skin clean, refined, smooth, and supple—achieving a new level of effective skin cleansing.

### 2. Thermal Beauty Principle

Gentle warmth helps stimulate blood circulation in the deeper layers of the skin, accelerating its metabolic rate. This function effectively reduces dark circles and fine lines around the eyes, while enhancing skin transparency and radiance—even without makeup. When combined with micro-vibrations, the warming effect helps improve skin firmness and relieve fatigue.

#### A. Micro-vibration Effect

- Enhances skin firmness and elasticity.
- Helps prevent allergies and other skin issues.
- Shrinks pores, locks in moisture, and reduces wrinkles.
- Relieves muscle tension, skin fatigue, redness, and swelling.
- Helps prevent recurring acne.

#### B. Warming Effect

- Soothes the skin and muscles, relieves pain and soreness.
- Stimulates blood flow, opens pores, and enhances absorption of skincare products.
- Improves skin elasticity and moisture retention.
- Works with an attached cotton pad for deep cleansing, helping to extract facial oils and dirt more effectively.

### 3. Light Wave Therapy Principle

#### Red Light:

Boosts cellular activity and metabolism, stimulates collagen production, improves circulation, and enhances skin elasticity. It reduces dullness and yellow undertones, offering anti-aging and skin-repairing benefits.

#### Blue Light:

Effectively reduces inflammation during acne formation. It inhibits acne-causing bacteria without harming surrounding skin tissue and shortens the inflammatory phase.

#### Yellow Light:

Penetrates the skin's deeper layers, is absorbed by red blood cells, and promotes a photobiochemical response. This light helps break down pigmentation, promotes lymphatic drainage, and improves roughness, redness, and overall skin immunity.

### **Violet Light:**

Offers antibacterial and anti-inflammatory effects, boosts collagen production, improves elasticity and radiance, and helps reduce signs of aging and skin irritation.


### **4. Jade Therapy Effect**

Jade is believed to have natural disinfecting, detoxifying, and calming properties. It supports blood circulation, slows aging, improves sleep quality, and enhances metabolism and immunity. Jade may help inhibit a wide range of bacteria and viruses, without any side effects or harm to the human body.

### **5. EMS Microcurrent Technology**

This device uses EMS (Electrical Muscle Stimulation) microcurrents to activate skin cells, stimulating collagen production and muscle contraction. A gentle tingling sensation and a feeling of firmness may be experienced during use. EMS helps reduce fine lines and wrinkles, restoring skin's smoothness, firmness, and youthful glow—safely and painlessly.

## **PRODUCT STRUCTURE**

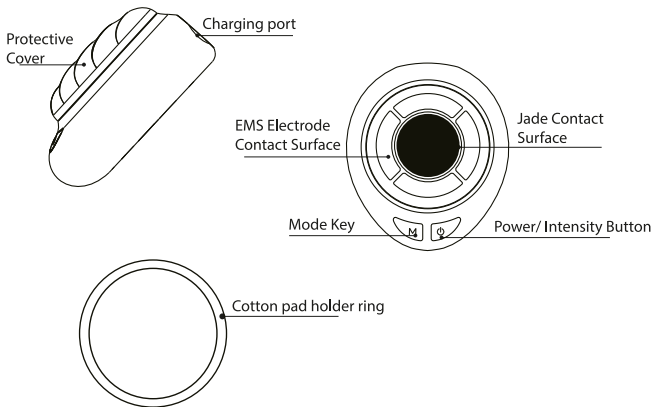

## SAFETY PRECAUTIONS


1. Not recommended for use by:
  - Individuals with implanted medical electronic devices, such as pacemakers.
  - Individuals with life-sustaining electronic medical equipment (e.g. artificial heart-lung machines) or other implanted devices such as cardiac monitors.
  - Individuals with acute illnesses or infectious diseases such as tuberculosis.
  - Individuals with high blood pressure or blood disorders.
  - Individuals with skin diseases, atopic dermatitis, allergies, or other skin conditions.
  - Individuals with pathological skin discoloration caused by physical irritation.
  - Individuals with physical, sensory, or mental impairments, or those lacking experience in using such devices.
  - Individuals with reduced sensitivity to heat or pain.
2. Do not use the device on the following areas:
  - Areas with facial implants or subcutaneous metal.
  - Inflamed, red, swollen, or irritated skin.
  - Moles that are raised or irregular, the eyelids, eyeballs, or the area directly around the eyes.
3. If you experience pain or any abnormal sensation on your skin or body, stop using the device immediately.
4. Clean the device regularly after use to prevent buildup of dust or moisture, which may reduce its lifespan.
5. Stop using the device immediately if abnormal heat is detected to prevent any risk of injury.
6. When not in use, store the device in a cool, dry, and well-ventilated place. Do not store or use the product in bathrooms or other humid environments.
7. Keep out of reach of children. Do not allow children to use or place the metal tip in their mouths, as this may lead to accidents.
8. Do not drop or strike the device. If it becomes damaged due to impact, do not use it further.
9. Do not submerge the device in water or rinse it under running water for an extended period.
10. Do not attempt to disassemble or repair the device yourself.
11. Avoid using the device on areas affected by acne, pimples, redness, or inflammation.  
For sensitive skin, it is not recommended to use the device daily—limit use to once a week.
12. Do not plug in the charging cable with wet hands. Always charge in a dry environment and ensure the device is switched off before charging.
13. If the device will not be used for more than seven days, disconnect it from the charging cable and store it with the power turned off.
14. Do not wait for the battery to fully discharge before charging. It is recommended to recharge the device as needed before each use.

## OPERATING INSTRUCTIONS


1. Press and hold [ ⏻ ] for 2 seconds to turn on the device; you will hear two beeps. Press and hold the [ ⏻ ] again for 2 seconds to turn it off; this time, you will hear a long beep.
2. Vibration intensity is set to low by default. Tap the [ ⏻ ] to adjust the level in the following order: low – medium – high. Each tap of [ ⏻ ]/[ M ] is accompanied by a short beep. Double-tap the [ ⏻ ] to turn the vibration on or off.
3. In Whitening and Rejuvenation Mode, which is activated automatically when the power is turned on, the red light, induction vibration, and heating functions are all enabled.
4. To activate Lifting Mode, tap the [ M ] this will enable the blue light, microcurrent, and induction vibration functions.
5. For Hydration and Repair Mode, tap the [ M ] again to switch to yellow light mode, which activates the induction vibration and heating functions.
6. To enter Soothing and Calming Mode, tap the [ M ] once more; this enables the violet light, induction vibration, microcurrent, and heating functions.
7. The default working light will blink to indicate low battery. During sonic vibration, the blue light will flash. The device will emit six beeps before automatically shutting down.
8. An intelligent 8-minute timer powers the device off automatically after one session.

## HOW TO ATTACH A COTTON PAD

### Assembly:

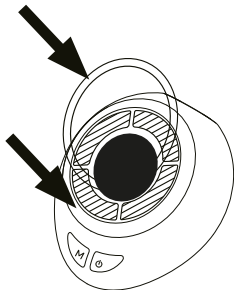

### Disassembly:

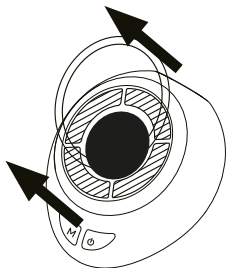

## HOW TO USE

1. Cleanse your face, remove any makeup, and dry the skin thoroughly before use. Remove all metal jewellery from your face, neck, and hands.
2. Depending on your skincare needs, apply a thin layer of conductive gel or a moderate amount of your preferred skincare product. Select the appropriate mode and intensity level.
3. Before first use, test the device on the back of your hand or behind your ear. If no skin reaction (e.g. redness, burning, irritation) occurs within 24 hours, the device can be used safely. If any abnormal reaction appears, discontinue the test and do not use the product.
4. Press and hold the power button to turn on the device. Then, short press the button to switch between modes, and short press the intensity button to adjust the level. Once the desired mode and intensity are selected and the indicator light is on, the device is ready for use.
5. Place the treatment head gently against your skin and begin to massage. Move the device slowly in upward, lifting motions, covering the entire face evenly.
6. After use, press and hold the power button until you hear a long beep, indicating the device is turned off. Wipe the treatment head with a soft, slightly damp cloth. Make sure to clean around the edges and crevices to prevent buildup that may affect future use. After cleaning, store the device in a cool, dry place.
7. Rinse your face with water and apply your usual skincare products.

### Note






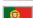











1. Always switch off the device and allow it to cool before cleaning.
2. Use non-alcoholic wipes or soft tissue to clean the treatment head and the device surface.
3. Do not use steel wool or harsh cleaning agents.
4. Keep the device dry. Do not immerse it in water or expose it to moisture for extended periods, as this may cause short-circuiting or damage.

## CHARGING INSTRUCTIONS

1. Place the device in a dry, well-ventilated area. Connect the charging cable properly to both the main unit and a power socket.
2. When the blue indicator light is blinking, the device is charging.
3. When all blue lights remain steadily on and you hear three beeps, it means the battery is fully charged. Unplug the device from the power source before switching it on and using it.
4. Do not turn on the device while it is charging.

## DISCLAIMER

This product is not intended to diagnose, treat, cure, or prevent any disease. No liability is accepted for any damage or injury resulting from improper use.

-  Scan the QR code to download the manual
-  Zeskanuj kod QR, aby pobrać instrukcję
-  Scannen Sie den QR-Code, um das Handbuch herunterzuladen
-  Scannez le code QR pour télécharger le manuel
-  Escanee el código QR para descargar el manual
-  Digitalize o código QR para descarregar o manual
-  Scansiona il codice QR per scaricare il manuale
-  Scan de QR-code om de handleiding te downloaden
-  Скануйте QR-код, щоб завантажити інструкцію
-  Naskenujte QR kód a stáhněte si návod k použití
-  Naskenujte QR kód a stiahnite si príručku
-  Scanează codul QR pentru a descărca manualul
-  Nuskaitykite QR kodą, kad atsisiųstumėte instrukciją
-  Nolasiēt QR kodu, lai lejupielādētu rokasgrāmatu
-  Lae juhend alla, skannides QR-koodi
-  A kézikönyv letöltéséhez olvassa be a QR-kódot.
-  لي ليدل ليزن تلة عيرس لة باجت سالا زم ح سا



innogIO

Manufactured in China for **INNOGIO GLOBAL** sp. z o.o.  
 Ostródzka street 74H, 03-289 Warsaw, Poland  
[www.innogio.com](http://www.innogio.com)



@innogIO.official



@innogio\_com