EN: LunaEye by GLOV – User Manual

Contents

1. General Device Information
2. Product Overview
3. Description of the Device
4. Before You Start
5. How to Use LunaEye
6. Warnings
 - Potential Benefits
 - 6.1 Tips for Use
7. Using & Caring for LunaEye
 - 7.1 General Use
 - 7.2 Cleaning
 - 7.3 Storing
 - 7.4 Troubleshooting
8. Specifications
9. Travel and International Use
10. Limited Warranty – 2 Years
11. Disposal
12. Symbols
13. Contact Information

1. General Device Information

LunaEye by GLOV is an advanced LED care device designed for the delicate under-eye area. It uses **red (630 nm)** and **infrared (830 nm)** light in **9-minute automatic sessions** to support the skin's natural appearance.

Compact, portable, and rechargeable, LunaEye fits into any daily routine. It can be used under the eyes or on other areas of the face in consecutive sessions.

⚠ Note: This is a beauty device for **personal home use only**. It is not a medical device and is not intended to diagnose, treat, or cure any medical condition.

2. Product Overview

LunaEye Case	LunaEye USB
1 × Charging Case 2 × LunaEye Patches	2 × LunaEye Patches 2 × Mini Controllers

2 × Mini Controllers 1 × USB Cable 1 × Set of Adhesive Stickers 1 × User Manual	1 × USB Cable 1 × Set of Adhesive Stickers 1 × User Manual
--	--

3. Description of the Device

- **Light Source:** Red (630 nm) + Infrared (830 nm) LEDs
- **Session Duration:** 9 minutes (auto shut-off)
- **Battery Life:** 4–6 sessions per charge
- **Design:** Ergonomic crescent-shaped patches
- **Frequency:** Safe for daily use; multiple areas can be treated consecutively

4. Before You Start

Perform a **sensitivity test** before first use:

1. Place on inner forearm.
2. Run one full session (9 min).
3. If redness or irritation occurs within 24h, discontinue use.

5. How to Use LunaEye

1. **Charge Fully** - Before first use, place the patches in the charging case or USB pod and charge for about 2 hours until fully powered.
2. **Prepare Skin**- Cleanse and dry the area you wish to treat (typically the under-eye). Ensure skin is free of oils, creams, or makeup for best contact.
3. **Apply Patches** - Attach the adhesive stickers to the back of the patches. Place gently under the eyes or on another target facial area.
4. **Start Session** - Press and hold the power button for 1.5 seconds. The device will activate and run a **9-minute automatic session**.
5. **End of Session** - The device switches off automatically after 9 minutes. Remove patches and discard adhesive stickers.
6. **After Care** - Apply your preferred serum or moisturizer to complete your routine.

△ Note: Multiple facial areas (cheeks, forehead, etc.) may be treated consecutively in the same routine.

6. Warnings

- Stop use if you experience discomfort, dizziness, or prolonged redness.
- Do not use on **wounds, ulcers, post-surgical areas, or severely irritated skin**.
- Avoid use if you have **sunburn, photosensitivity, inflamed skin**, or are on **photosensitizing medication**.
- Do not look directly into the LEDs.

- Not for children. Use only under supervision if the user has reduced physical/mental ability.
- For hygiene reasons, do not share the patches.
- Keep away from **direct sunlight, heat, or water**.
- Use only the supplied **USB cable and adapter**. Ensure charging port is dry before use.
- Do not use while charging.
- Do not disassemble the device — it contains no user-serviceable parts.
- Maximum session time is **9 minutes** (preset).
- Dispose of responsibly (see Disposal section).

7. Potential Benefits

Regular use of LunaEye may help to:

- ✨ Refresh & brighten the under-eye area
- ✨ Reduce the look of puffiness & fatigue
- ✨ Support smoother, more resilient-looking skin
- ✨ Complement daily skincare routines

7.1 Tips for Use

- Use on **clean, dry skin**.
- Apply serum or moisturizer after sessions.
- Safe for **daily use**; different areas can be treated one after another.

8. Using & Caring for LunaEye

8.1 General Use

- The device works at safe, non-invasive light levels.
- Shuts off automatically after 9 minutes.

8.2 Cleaning

- Wipe with a **damp cloth** or **70% alcohol wipe**.
- Do not use solvents or harsh cleaners.
- Never clean while charging.

8.3 Storing

- Store in the **charging case**.
- Avoid bending, crushing, or extreme temperatures.

8.4 Troubleshooting

- **Device won't turn on:** Recharge for 2 hours, then restart.

- **Device won't turn off / frozen:** Hold power button 10 seconds to reset.
- **Stops during session:** Recharge fully, clean contacts, retry.
- **Won't charge:** Use original cable, ensure port is dry, use 5V/2A adapter.

9. Specifications

LED Technology

- Wavelengths: 630 nm (red), 830 nm (infrared)
- Max. power density: $\leq 5.5 \text{ mW/cm}^2$
- Session: 9 minutes (auto shut-off)

Battery

- Type: Rechargeable Li-ion
- Capacity: 65 mAh
- Battery life: 4–6 sessions per charge
- Charging time: ~2 hours
- Power: 5V / 2A
- Standby: Up to 30 days

Physical Properties

- Materials: PC + ABS, body-safe silicone
- Colors: White/Pink
- Weight: ~24 g per patch
- Dimensions: ~70 × 40 × 12 mm

Environmental Conditions

- Storage temp: 0°C to 40°C
- Storage humidity: 30% – 85%

10. Travel and International Use

- Lightweight & compact – ideal for travel.
- Worldwide voltage (**AC100–240V, 50/60Hz**).
- Confirm airline rules before onboard use.

11. Limited Warranty – 2 Years

Phenicoptere Sp. z o.o. warrants LunaEye for **2 years** from purchase against defects in materials or workmanship under normal use.

Does **not** cover:

- Cosmetic wear and tear
- Damage from misuse, neglect, or improper storage
- Unauthorized repair or tampering

⚠ Opening or disassembling voids the warranty.

If a defect occurs, **Phenicoptere Sp. z o.o.** may repair or replace the device. Proof of purchase required. Shipping costs non-refundable.

This warranty is **in addition to your statutory consumer rights**.

Claims:

 contact@phenicoptere.com

 www.glov.co

12. Disposal

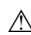
The **crossed-out bin symbol** means this product must not be disposed of with household waste. Take to an official collection point for recycling.

Battery Disposal

 **CAUTION:** Opening the device voids warranty. Remove the battery only when disposing of the device.

- Contains a lithium-ion battery — do not discard with household waste.
- Switch off and disconnect before opening housing.
- Remove and recycle the battery at an official collection point.
- Wear gloves when handling.

13. Symbols

 – Warnings / safety instructions

 – Refer to manual

 – Electronic recycling

 – Manufacturer

LOT – Batch code

 – Keep dry

EC REP – EU Representative

14. Contact Information

Manufacturer: Produced for Phenicoptere Sp. z o.o.

Address: Chocimska 14a/101, 00-791 Warsaw, Poland

Tel: +48 577 111 675

EU Representative: Phenicoptere Sp. z o.o.

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3. Description de l'appareil

Source lumineuse : LED rouge (630 nm) + LED infrarouge (830 nm)

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