
EN

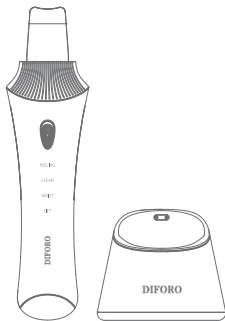DIFORO
Rives

Contents

1.	SAFETY -----	.1
2.	TECHNICAL DATA -----	.2
3.	CERTIFICATES -----	.2
4.	CLEANING AND CONSERVATION -----	.2
5.	PRODUCT DESCRIPTION -----	.3
6.	CHARGING AND USING -----	.4
7.	FACE AND BODY CARE -----	.4 i 5
8.	WARRANTY CONDITIONS -----	.6
9.	INFORMATION REGARDING DISPOSAL -----	.7

1. Safety

- People with a pacemaker should not use the device.
- During the care treatment, remember not to put them on your hands
- conductive objects, such as rings, watches, and bracelets. They can interfere with the operation of the devices.
- The spatula can damage the eyeball. Under no circumstances do not touch the eyeball.
- If your skin is irritated after sunbathing, has unhealed scars, inflamed inflammation, or any other type of damage, do not perform skin care at the site of the damaged skin.
- Pregnant women should not use the device.
- Do not carry out the procedure in one place for a long time.
- Keep the interval between treatments.
- Do not use if you have a blood clotting disorder (including HIV infection).
- Do not use with skin care products with high exfoliating properties such as liquid scrubs or granules at the same time.
- If you are unsure about any potential side effects, talk to your healthcare professional.
- The device uses a lithium battery. Do not place or store the product for a long time in a humid room or high temperature.

Technical Data

Name: Skin Scrubber

Model: Rives

Power: 2.5W

Input Voltage: 5V

Battery: 600 mAh

Charging time: 120 min.

Net weight: 150g

Waterproof Grade: IPX5

Dimensions: 197.9 * 79.9 * 45.7mm

Docking station: With a UV lamp

3. Certificates

CE / ROHS / FCC - full declarations available at diforo.pl

4. Cleaning and conservation


After each use, make sure the device is turned off, clean the ultrasonic skin care spatula and wipe the device with a paper towel or cloth. Do not immerse the device and charging base directly in water. Avoid washing with water directly so as not to cause a breakdown. Do not use steel elements, strong detergents, caustic liquids or other cleaning agents. The spatula, when placed in the docking station, will be additionally disinfected by a UV lamp.

5. Product description


Charging
connector

Spatula

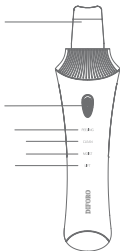Charging
connector

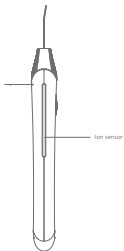6. Charging and using

The charging base and the device are connected with power pins. The charging base has a USB type C connector. It takes 2 hours to fully charge the product. When fully charged, the device works for up to 100 minutes. During charging, the LEDs blink, when fully charged, the LEDs glow continuously. When the power required is low, the current mode light flashes quickly with a beep, then the unit automatically turns off.

Basic functions after starting the device:

Power on / off: long press the switch / mode key, it will start the device with beep twice. Mode switching: press the toggle / mode key, prompt "beep" and toggle the four circulatory modes.

Skin Scrubber has a timer, after 5 minutes the device will automatically turn off.

7. Face and body care - first steps

Before use, cleanse the skin and remove any cosmetic residues.

The "PEELING" mode is the default and is used as a basic treatment.

" PEELING " ②

For use in conjunction with daily tonics etc., place an ultrasonic spatula on the face at a 30 ° angle, gently touch the skin and slowly move upward for a blackhead removal and gentle peeling effect.

" CLEAN " ③

"CLEAN" mode, for daily use with tonics. Point the spatula upwards at an angle of 30 ° gently touching the surface to be cleaned. Move slowly upwards while gently pressing for a deep cleaning effect.

" MOIST " ④

Use with lotions and creams, evenly tap cosmetics with the contoured tip of the spatula to break down the particles and penetrate the substance deeper.

" LIFT " ⑤

Use with gels and creams, moving at an angle of 30 ° from the bottom up, gently pulling the skin. Thanks to the "LIFT" mode, you will get smooth and taut skin.

After the treatment, cleanse the skin and perform natural skin care according to your preferences. You can use a serum or other care cosmetic.

Intervals between treatments and therapy.

" PEELING " Twice a week

" CLEAN " Twice a week

" MOIST " three to four times a week

" LIFT " three to four times a week

Before first use, do a skin test and wait 12 hours. This will help to eliminate any risks and adverse reaction of the skin to the treatments. Remember, the device is designed to exfoliate the epidermis, and it is still an invasive procedure. Apply a small amount of water or tonic to your forearm and run the spatula on the inside of your upper arm for about 30 seconds. If there is no strong reaction, you can safely use the product. Using it will cause the skin to appear slightly reddened or swollen. This effect will wear off after a while. This is the desired effect because this is how the peeling works. If redness does not disappear within 24 hours, consult your doctor. The skin tolerance and sensation vary from person to person and skin to complexion.

8. Warranty conditions

DIFORO provides a 2-year warranty on this device, counted from the date of the first purchase;

The warranty covers defects due to faulty workmanship or materials as well as moving parts that affect the operation of the device. If, during the warranty period, the user notices a defect and notifies the dealer, DIFORO, will replace the device free of charge or make a free repair.

The warranty does NOT cover cosmetic deterioration caused by normal wear and tear, or damage resulting from mechanical damage, misuse or neglect. Attempting to open or disassemble the device (or its accessories) will void the warranty. Claims under the warranty must be supported by reliable evidence (proof of purchase) that the claim arose within the warranty period. To validate the warranty, keep the original purchase receipt for the duration of the warranty period. The user will be provided with a Return Merchandise Authorization (RMA) number which must be handed over to the nearest DIFORO office or authorized reseller with the device. Shipping costs are not refundable. This obligation is in addition to, and in no way affects, your legal rights as a consumer.

This warranty is provided by VOBÉ sp. z o.o., ul. Rudna Mała 167A, 36-054 Mrowla. KRS: 0000800189, Registration number; NIP: 525-279-92-70, REGON: 384166183.

8. Disposal Information

Disposal of waste electrical equipment (suitable for EU countries and other European countries, with a waste recycling system).

The symbol of the crossed out wheeled bin means that the device should not be treated as municipal waste, but should be handed over to an appropriate collection point for used electrical and electronic equipment. The distributor or retailer is obligated to accept such device. If the device is used correctly, it does not pose any risk to the environment and human health, which may not be the case with inappropriate disposal of the used product. The recycling of materials also helps to conserve natural resources.

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©2022 DIFORO by VOBE sp. z o.o.. ALL RIGHTS RESERVED