
DIFORO AMIRTA
RF FACIAL BEAUTY DEVICE
INSTRUCTIONS



Contains

- Beauty device
- Charging base
- Charging cable
- Manual
- Warranty

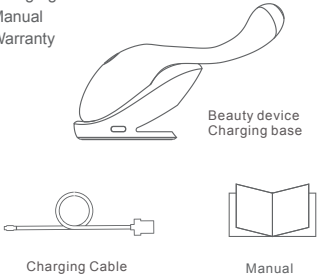Product introduction

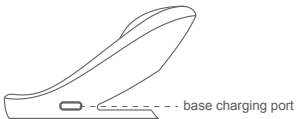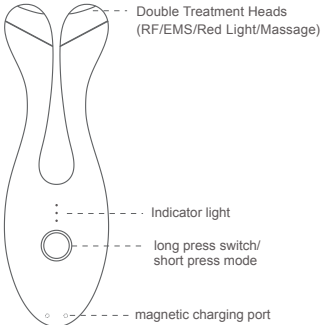Button/Mode Description

Mode Display

- Mode1: Soothing micro current+RF+Redlight
- • Mode2: Contraction microcurrent +RF+Redlight
- • • Mode3: Promote activation of microcurrent+RF+Red light

Charging Display

When the battery is low, the three charging indicators flash rapidly;
When charging, the three charging indicators flash circularly;
When full, the three charging indicators are always on.

Function Introduction

(((••))) RF(Radio Frequency)

RF(Radio Frequency) can penetrate the epidermis layer to reach the dermis layer, heat collagen cells, stimulate the growth of fibroblasts, and regenerate a large amount of collagen, thereby supporting the skin of the epidermis, so as to achieve the effect of firming the skin.

~~~~~ EMS Microcurrent

EMS microcurrent is a low-level electrical current that stimulates the face by sending gentle electrical waves through the skin and tissue to the facial muscles. Continuous tightening from the epidermis layer to the muscle layer can tighten the facial muscles. Microcurrent has been shown to stimulate the production of ATP, stimulate the vitality of facial muscles and promote the production of collagen and elastin.

~~~~~ Optical Beauty

Use 620-630nm red light. In the cells in the dermis, mitochondria can absorb a large amount of red light, accelerate the regeneration of cells, and synthesize collagen, thereby achieving the effect of plumping epidermal cells, firming the skin, diluting acne marks and melanin deposition.

~~~~~ Food Grade Silicone Body

Food-grade silica gel is an inorganic polymer colloid material made of silicic acid polycondensation. It is non-toxic and tasteless, with stable performance and soft texture. It is safe, healthy, and odor-free by selecting high-end materials. Users can use it with confidence.

~~~~~ Doubles Treatment Heads: 316L

316L is medical grade stainless steel, which is often used in medical equipment, scalpels and other medical products. It's no irritation to human epidermis that any skin can be used with confidence.

~~~~~ Memory mode, intelligent time control


The last mode you used will be remembered every time you turn it on. Automatic time control function, it will automatically shut down after 5 minutes of use, and switch the mode halfway to re-time to avoid damage to the skin caused by prolonged use.

How to Use?

It is recommended to use after cleansing the face to prevent facial oil from affecting the use experience.

① Apply the gel to the area of use. Please try it on your hands before use, and use it on your face if there is no discomfort. When using it around the eyes, there will be a flickering sensation when the eyes are opened or closed, which is the function of biological micro-electricity, which is a normal phenomenon. For better results, apply eye cream or face cream evenly. If it is used in winter, to prevent the metal from being cold when it first touches the face, you can start it on the hand for a while, and then use it on the face after it warms up.

② Double heads can be used for eye care/facial grooming/face lift, and they are perfect for different parts. The corners of the jaw, the corners of the eyes, the cheeks, the bridge of the nose and other parts can be lifted and massaged. Long-term use can effectively solve facial skin problems such as eye fine lines, nasolabial lines, neck lines, and firm and smooth skin. The use method is shown in the figure below. It is recommended to use the method from the bottom to the top to achieve the effect of lifting.

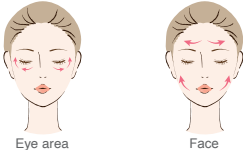Precautions

- Do not use the red light directly into the eye.
- Do not use on upper eyelids.
- This product needs to be used with gel. If there is too much gel on the guide head and the left and right guide heads are connected, there is no feeling on the face. Please wipe the excess gel between the left and right guide heads before use.
- Do not disassemble and modify this product without permission.
- The body cannot be washed.
- When using it around the eyes, there is a light flickering sensation, which is a normal biological micro-electricity effect, please do not worry.
- Due to the deviation of the resistance value of the human body, the electric shock is strong when used by individual people. Please choose whether to continue using it according to the individual tolerance level. If you feel uncomfortable, please stop using it immediately.
- After use, please wipe the guide bit with a clean cloth or wet paper towel.
- If you find any abnormality such as vibration during use, please check whether the battery is low or turn off.
- When not in use for a long time, please place the machine in a dry place to avoid direct sunlight.
- Please keep this product out of the reach of children, this product cannot be used as a toy.
- If you experience skin irritation or any discomfort during use, please stop using it immediately.

The following persons should not use

- Those equipped with medical electronic devices such as cardiac pacemakers and other implantable devices;
- When you have heart disease, cancer, or feel abnormal blood pressure;
- Suffering from diseases with fever, infectious diseases, acute diseases;
- During dental treatment;
- Patients with weakened immunity due to immunocompromised diseases such as AIDS or those who use immune preparations;
- Hemophilia, bleeding disorders. Diabetic facial neuralgia, sensory, allodynia;
- Crab foot swollen constitution, facial melanoderma, light allergy, abnormal pigment;
- Patients with microvascular dilation caused by long-term use of steroids or impaired liver skills;
- People with allergies, atopic dermatitis, skin diseases, sensitive skin or sensitive skin;
- Those with a history of epilepsy or discomfort due to flickering light;
- Those who have sequelae on the head or eyes, those who have undergone eye surgery (please consult a physician before use);
- During pregnancy, menstrual period, breast feeding;
- Young children or those who cannot express their own will, those with alcohol and hangovers;
- Those who take drowsiness-inducing drugs, those who feel extremely tired, or who are in poor physical condition;
- Those who have problems such as pigmentation and inflammation due to chemical fibers, or who have or have had dermatitis due to cosmetics;
- Those who have recently received chemical peels, and those who have had skin problems after using beauty equipment;

\* Do not use other parts except designated parts.  
\* Those who are receiving doctor's treatment should consult a doctor before using.

Product Parameters

|                  |                                |
|------------------|--------------------------------|
| Name             | Diforo RF facial Device beauty |
| Model            | AMIRTA                         |
| Size             | 142*46*44 mm                   |
| Net Wight        | 130g                           |
| Material         | 316L/Silicone/PC/ABS           |
| Battery Capacity | 400mAh                         |
| Charging time    | 120 min                        |
| Working time     | 40~60 min                      |
| Input            | 5V/0.27A                       |
| Power            | 4W                             |
| Certificate      | CE RoHS                        |

Warranty

The manufacturer guarantees its manufactured products for a period of 24 months from the date of purchase. Accessories included with the original product are warranted for a period of 12 months from the date of purchase.

For the efficient operation of products and complaints, we

- ① Study and follow the operating instructions carefully.
- ② Keep the manual, warranty card, proof of purchase.
- ③ Use the original accessories supplied with the device.
- ④ In case of problems with the device, contact the Manufacturer.
- ⑤ Do not make repairs or modifications yourself.
- ⑥ Use the product according to its intended use.

During the warranty period, all repairs are performed free of charge.

Any defects or faults will be removed within a maximum period of 14 days, calculated from the time of delivery of the device to the manufacturer's authorized service center.

Limitations

The warranty will not be recognized, in case of mechanical damage due to inattention or fortuitous events.

If the device is used contrary to its intended use and specifications.

Batteries are consumables, the warranty covers 6 months. Batteries cannot be replaced outside of an authorized service center.

Dane kontaktowe

DIFORO Brand

CUSTOMER SERVICE **+48 535 281 934**

Consultation hours : 08:00 ~ 16:00 (8 AM ~ 4 PM)  
(Monday to Friday)

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