

# medivon®
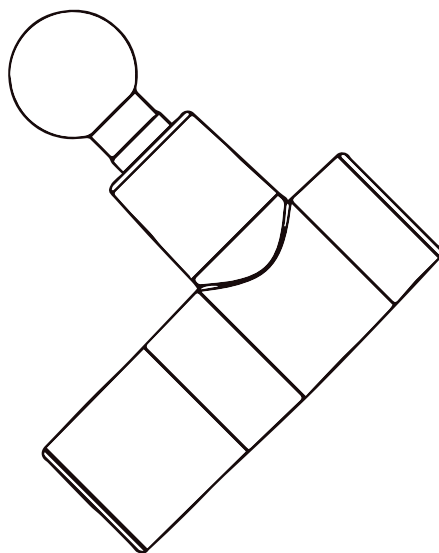
## **GUN MINI**

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# Accessory list

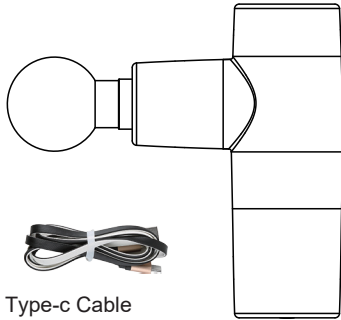Type-c Cable

Small Ball Head

For Large muscle

Six Finger  
For Hands  
Feet and Muscles

Forked Head

Suitable for Spine

Palm Heads

For All  
Kinds of Muscles



Flat Head  
Suitable for all  
parts of the body

Bullet Head

For joints  
of the body

## Illustration:

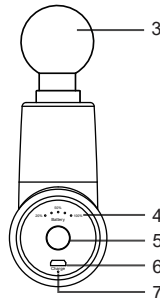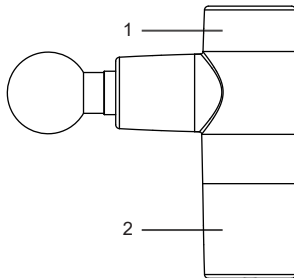
- 1.Brushless motor
- 2.Rechargeable battery
- 3.Six massage heads
- 4.Battery capacity indicator
- 5.On/Off/Mode button
- 6.Type-c charge
- 7.Charge LED

# Specification


<div>Category</div> <div>Model</div>		SPEC
Battery		Battery Type: 18650 Connection Type: Series QTY: 2 PCS Capacity: 2600 mAh Voltage: 7.4 V
Motor		Type: Brushless Maximum Speed: 3600 RPM Maximum Torque: 25 mNm Voltage: 7.4 V
Gears		5Gears 800 RPM/1200 RPM/ 1600 RPM/2000 RPM/2400 RPM
Power Supply		Common use with mobile phone charger Type-c input 5V 2A
Amplitide		8mm

# How to Use
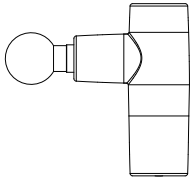
## Step 1:

Select a suitable  
massage head and install
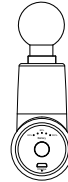
## Step 2:

Press and hold to start on, press one time  
to change one speed gear, totally 5 Gears  
Gear1-Gear2-Gear3-Gear4-Gear5  
Speed:800RPM-1200RPM-1600RPM-2000RPM  
2400RPM

# Notes

1. When you don't use the device and not switch off, 2 minutes , LED display will switch off to save power
2. When you charge the device,when the battery less than 10%,you can not use it
3. When the battery left only 5%, the device will auto stop to work to protect the battery
4. When you charge the battery, the charge led will flash, when it was full, Led will stop flashing  
And when you charge the battery, battery leds will be on ,show the capacity of battery

## Battery instructions

When the battery is empty it needs to be charged for 2.5-3 hours to be fully charged. The battery life is determined by the gear adjustment. When the battery is charging, the battery power-on indicator light (5) will light up, indicating normal charging. For your safety, the motor button cannot be activated when the battery is charging. Starting the motor during charging is prohibited.

## Notes

To reduce the risk of electric shock, fire and personal injury, this product must be used in accordance with the following instructions.

1. Only to be used by adults.
2. This equipment should not directly massage the skin surface to avoid damage to the skin. The massage part should be covered by clean, dry clothes or a clean, dry towel. It should be gently pressed and moved. Massage time on the same part should not exceed 60 seconds.
3. High speed and pressure should only be used on soft tissue parts of the body, if there is no pain or discomfort. High speed and high pressure must not be used on the head or any bones of the body.
4. Please choose the appropriate massage head, according to the massage needs of the different parts.
5. Frequent high-speed high-pressure use exerts pressure on the same body part. If abrasions, painfulness or incapacity occur, please stop using immediately.
6. Fingers, hair or any part of the body should be kept away from the telescopic rod of the device and the cooling vents of the back motor to prevent entanglement.
7. Do not put any objects into the cooling vents of the motor.

8. Do not immerse the device in water. Do not allow water to enter the device or the motor cooling vents.
9. Do not misuse the device. Dispose of the device correctly after the end of its useful life.
10. Only use the charger provided by the original manufacturer to charge.
11. Please check the equipment and battery carefully before each use.
12. Do not make any modifications to the equipment.
13. Do not run or charge the device unattended.
14. This device may not be disposed of with the normal household waste at the end of its life, but must be taken to a special collection point for the reuse of electrical and electronic devices. The materials used in the device can be recycled. By recycling used household appliances you make an important contribution to the protection of our environment. Ask your local government for the collection points.
15. The device is only intended for massaging the human body. For damage caused by the use of the device for other purposes, incorrect operation or repairs carried out by an amateur, no liability will be accepted.

## Maintenance

1. Wipe the surface of the device with a slightly wet towel and dry it with a soft piece of cloth.
2. When storing or traveling for a long time, separate the device from the battery and place it in the device's original case.

## Warning and cautions

1. Do not put the device inside water, keep the machine away from liquids and heat sources, ensure the motor cooling holes are free from dust and debris.
2. Do not remove the screws by yourself.
3. Turn off the device after one hour of operation and continue using it after resting it for 30 minutes.
4. Remove the charger from the battery after charging, and do not use the device during charging.
5. The use of this device for the injured part of the body is prohibited without consulting the doctor.

**Please read the instruction manual carefully before using this equipment, use this equipment correctly in accordance with the instructions, our company will not be liable for any damage or loss to the customer caused by the incorrect use of this equipment as explained.**


**Do not use this device without the approval of your doctor if:**

- |            |                       |                                 |
|------------|-----------------------|---------------------------------|
| -pregnancy | -pacemaker wear       | -neuropathy or retinal damage   |
| -diabetes  | -had surgery recently | -spinal inflammation            |
| -migraines | -IUD's                | -have had/have health problems. |

**You must consult your doctor first and use it with your doctor's permission.**

**We are studying the effects of the device on specific medical conditions, this may cause the above taboo symptoms to be not detailed enough, if you are unwell please stop using, follow the advice of a physician or physiotherapist.**

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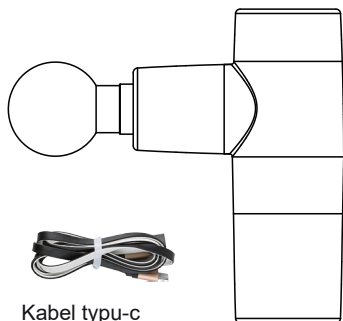
Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health damage or environmental pollution during further disposal or recycling of the equipment.

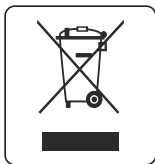Importer:  
Armare S.A.  
ul Grzybowska 87, Warszawa, Polska  
Designed by Medivon in Europe,  
produced in China

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I represent that I have read and accepted the terms and  
conditions specified in the Warranty Certificate.

.....