

medivon®

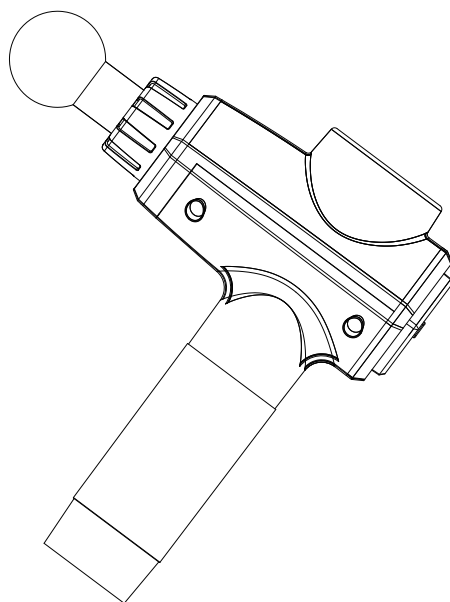GUN PRO

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Akcesoria

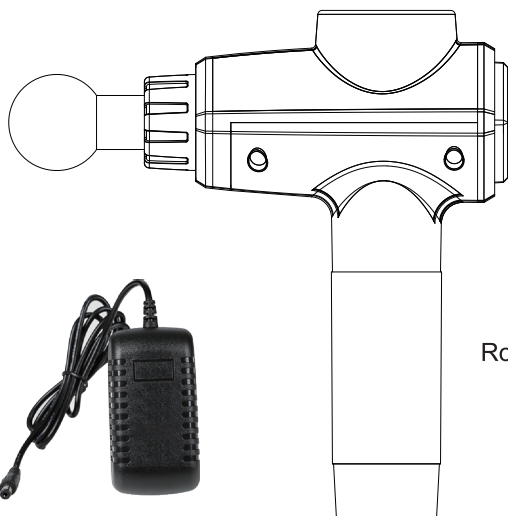Ładowarka

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
Accessory list

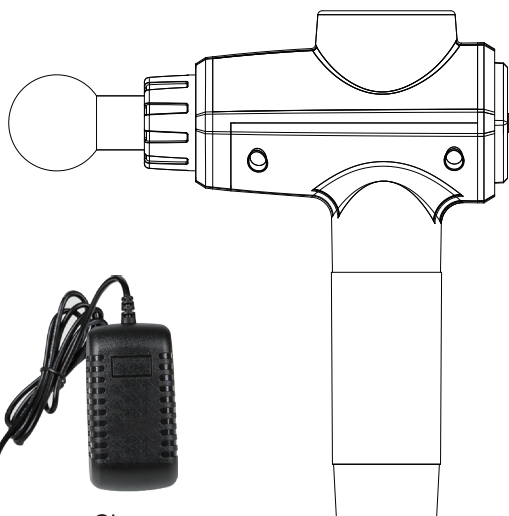Charger

Six Finger
For Hands
Feet And Muscles

Small Ball Head
For Large Muscle

Forked Head
Suitable For
Spine

Palm Heads
For All Kinds Of Muscles




Flat Head
Suitable For All Parts
Of The Body

Bullet Head
For Joints Of
The Body


Illustration:

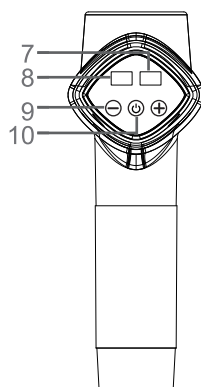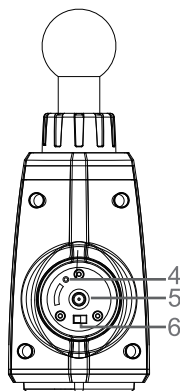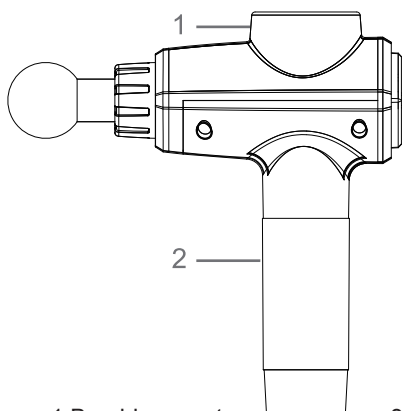- 1.Brushless motor
- 2.Rechargeable battery
- 3.Six massage heads
- 4.Charging LED
- 5.Charging jack
- 6.Power switch on/off
- 7.LED display for battery left


- 8.LED display for gears(20 Speed Version)
- 9.Touch button for speed up/down(20 Speed Version)
- 10.Touch button for stop(20 Speed Version)

How to Use

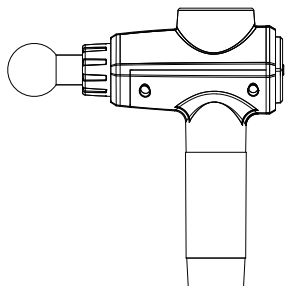
Step 1:

Select a suitable massage head and install

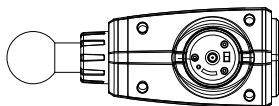Step 2:

Switch on

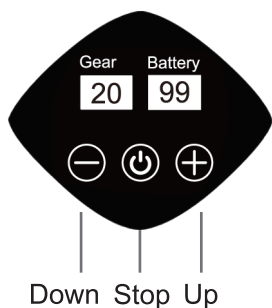Step 3(20 Gear Version):

Up/Down: one time touch only changes one gear Continuous touch to reach lower or higher gears quickly
Stop: Motor reduces to 0 speed
Up Stop Down and enters standby status

Notes

1. When you don't use the device and not switch off, 2 minutes , LED display will switch off to save power
2. Friendly definition in software
 - a. For touch version, press and hold the touch button, can fast change speed to maximum or minimum
 - b. For touch version, one button operate can reach any speed, the speed is loop out
 - c. in stand by status, touch any button reach gear 1
 - d. for button version, press and hold the button, can switch to 0 gear
3. When you charge the device, when the battery less than 10%, you can not use it
4. When the battery left only 5%, the device will auto stop to work to protect the battery
5. When you charge the battery, the charge led will flash, when it was full, Led will stop flashing
And when you charge the battery, battery leds will be on , show the capacity of battery
6. When you use the speed you like , next time when you use it will stay in the same speed

Battery instructions

When the battery is empty it needs to be charged for 2.5-3 hours to be fully charged. The battery life is determined by the gear adjustment. When the battery is charging, the battery power-on indicator light (5) will light up, indicating normal charging. For your safety, the motor button cannot be activated when the battery is charging. Starting the motor during charging is prohibited.

Safety instructions

To reduce the risk of electric shock, fire and personal injury, this product must be used in accordance with the following instructions.

1. Only to be used by adults.
2. This equipment should not directly massage the skin surface to avoid damage to the skin. The massage part should be covered by clean, dry clothes or a clean, dry towel. It should be gently pressed and moved. Massage time on the same part should not exceed 60 seconds.
3. High speed and pressure should only be used on soft tissue parts of the body, if there is no pain or discomfort. High speed and high pressure must not be used on the head or any bones of the body.
4. Please choose the appropriate massage head, according to the massage needs of the different parts.
5. Frequent high-speed high-pressure use exerts pressure on the same body part. If abrasions, painfulness or incapacity occur, please stop using immediately.

6. Fingers, hair or any part of the body should be kept away from the telescopic rod of the device and the cooling vents of the back motor to prevent entanglement.
7. Do not put any objects into the cooling vents of the motor.
8. Do not immerse the device in water. Do not allow water to enter the device or the motor cooling vents.
9. Do not misuse the device. Dispose of the device correctly after the end of its useful life.
10. Only use the charger provided by the original manufacturer to charge.
11. Please check the equipment and battery carefully before each use.
12. Do not make any modifications to the equipment.
13. Do not run or charge the device unattended.
14. This device may not be disposed of with the normal household waste at the end of its life, but must be taken to a special collection point for the reuse of electrical and electronic devices. The materials used in the device can be recycled. By recycling used household appliances you make an important contribution to the protection of our environment. Ask your local government for the collection points.
15. The device is only intended for massaging the human body. For damage caused by the use of the device for other purposes, incorrect operation or repairs carried out by an amateur, no liability will be accepted.

Maintenance

1. Wipe the surface of the device with a slightly wet towel and dry it with a soft piece of cloth.
2. When storing or traveling for a long time, separate the device from the battery and place it in the device's original case.

Warning and cautions:

1. Do not put the device inside water, keep the machine away from liquids and heat sources, ensure the motor cooling holes are free from dust and debris.
2. Do not remove the screws by yourself.
3. Turn off the device after one hour of operation and continue using it after resting it for 30 minutes.
4. Remove the charger from the battery after charging, and do not use the device during charging.
5. The use of this device for the injured part of the body is prohibited without consulting the doctor.

Please read the instruction manual carefully before using this equipment, use this equipment correctly in accordance with the instructions, our company will not be liable for any damage or loss to the customer caused by the incorrect use of this equipment as explained.

Do not use this device without the approval of your doctor if:

- | | |
|-------------------------------|---------------------------------|
| -pregnancy | -migraines |
| -diabetes | -spinal inflammation |
| -neuropathy or retinal damage | -IUD's |
| -pacemaker wear | -have had/have health problems. |
| -had surgery recently | |

You must consult your doctor first and use it with your doctor's permission.

We are studying the effects of the device on specific medical conditions, this may cause the above taboo symptoms to be not detailed enough, if you are unwell please stop using, follow the advice of a physician or physiotherapist.

We provide warranty for battery, motor, controller of one year, other wearing parts warranty of three months.

Except when:

1. Damage to this equipment caused by man-made incorrect operation or transport.
2. Unauthorized removal and repair of this equipment.
3. The instruction manual isn't being respected.
4. Overloaded use of this product design parameter.
5. Damage due to the customer's unusually poor storage or maintenance environment.
6. Failure to provide proof of purchase date.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health damage or environmental pollution during further disposal or recycling of the equipment.

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Armare S.A.