
Sönik

- EN 2-Stage Sonic Toothbrush for Babies & Toddlers
- FR Brosse à dents sonique 2 étapes pour bébés et bambins
- ES Cepillo de dientes Sonic de 2 etapas para bebés y niños
- IT Spazzolino sonico a 2 fasi per neonati e bambini piccoli
- DE 2-teilige Akustikzahnbürste für Babys und Kleinkinder
- PO Escova de dentes sônica de 2 fases para bebês e crianças pequenas

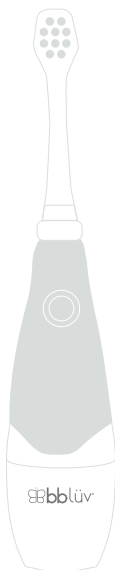SÖNIK B0116-01

Please read the instruction manual carefully before use.

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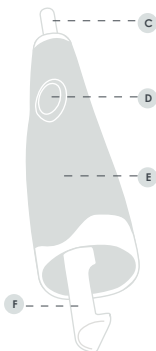
C

LED Bulb
Ampoule LED
Bombilla LED
Lampadina a LED
LED-Leuchtmittel
Lâmpada LED

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
BEFORE USE

1. Hold the toothbrush upside down. Remove the base of the brush **G** by turning it counterclockwise and pulling it upward.
2. Insert one AAA battery in the battery compartment, **F** making sure to match the indicated polarity.

3. Close the toothbrush by aligning the line markers on the base of the brush **G** and on the handle of the brush **E**.
4. Lock in the base of the brush **G** by turning it clockwise.

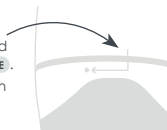
HOW TO USE YOUR SÖNIK

1. Carefully wet the bristles on the brush head **A** / **B** and apply a small amount of toothpaste.
2. Press the power button once **D** to turn on the LED light **C**.
3. Press the power button a second time **D** to activate vibration mode.
4. Use the toothbrush as recommended by your dentist or dental hygienist.
5. After brushing, press the power button **D** to exit vibration mode. Rinse the brush head **A** / **B** under running water.

NOTE: To avoid splashing, guide the brush head to the teeth before turning on vibration mode.

CHANGING THE REPLACEABLE BRUSH HEAD

1. Unlock the replaceable brush head **A** / **B** by turning it counterclockwise. Next, pull it upward to remove it.
2. Insert a new replaceable head **A** / **B** in and turn it clockwise to lock it in place.

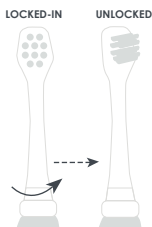NOTE: For the brush head to lock into place, the notch on the neck must align with the inner notches of the brush head.

CLEANING & MAINTENANCE

- Do not expose your Sönik directly to the sun or to high temperatures.
- Do not immerse this product in water for long periods of time.
- Clean the brush handle and the brush head separately and wipe them dry before reassembling.

WARNINGS

- The Sönik can be used by children and individuals with little experience as long as they are supervised by a person responsible for their safety.
- Children should not play with, clean, or repair the Sönik.
- Do not alter the Sönik in any way other than replacing the battery or the brush head.
- Do not try to repair the Sönik yourself.
- Insert the battery with its terminals positioned correctly.
- Do not try to recharge the non-rechargeable battery (*included*).

TECHNICAL CHARACTERISTICS

SPECIFICATIONS	
Material: DuPont Nylon/ABS/TPE	Power supply: 1 x AAA battery
Product size: 28 x 25 x 135 mm	Product standard: GB4706.59-20 08
Frequency: 16,000 movements per minute	

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