Oclean Smart Sonic Electric Toothbrush
User Manual

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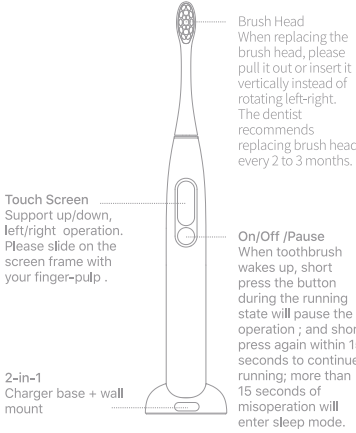
Oclean Smart Sonic Electric
Toothbrush User Manual

Instructions

- First Time Use
The toothbrush can be turned on by pressing the ON/OFF button for the first time on a brand-new toothbrush. The Bluetooth function will be activated and remain on. After installing the Oclean App, please complete the prompts to complete the connection and settings for the toothbrush.

Install App & Toothbrush Settings

- You can use the app to customize your brushing programs -automating tooth cleaning, tooth whitening and gum massaging in one brushing session. There is no need to open the Oclean App each time you brush. Simply open the app every few days to sync the brushing data.

- Direct Use without App Settings
You can also use the toothbrush without connecting to the app. It has a touchscreen display and an ON/OFF button. It offers touch control via the touch-sensitive screen.

- Brushing Plan and Performance Records
Toothbrush will start default mode with “Standard Cleaning Plan” with 2 mins if no new customized brushing plan setting in APP.

After brushing, the display will show the effective brushing time, the brushing scores and the tips of the brushing area.

- How to Charge
Remove the charging base from the wall mount and connect it to the mobile phone charger. When the bottom of the toothbrush is placed gently above the charging base, it will automatically position itself. The display will show the battery percentage when the toothbrush is charging. A full charge usually takes about two hours.

Follow the Oclean Official
Facebook Page

Oclean Facebook Page
www.facebook.com/ocleanofficial

Search "Oclean" on App Store or Google Play to download

Specifications

Product Name	Oclean Smart Sonic Electric Toothbrush
Model	Oclean X
Edition	Pro
Max Output	0dBm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V $\pm$ 1A
Battery Life	20 days (Level 1), 40 days (Level 32)
Battery Capacity	800 mAh
Charging Time	2 hours
Bluetooth	4.2 above

* Calculated by 2 minutes/time, twice/day.

(R) 201-190289

Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

Name of parts	Toxic and harmful substance or element						Environment Protection	
	Pb	Hg	Cd	Cr6+	PBB	PBDE	Term	Term
PCBA	○	○	○	○	○	○	10 Years	10 Years
Plastic parts	○	○	○	○	○	○		5 Years
Metal Hardware	○	○	○	○	○	○		
FPC	○	○	○	○	○	○		

Warranty Card

- These clauses only apply to Oclean products sold in the country or region where the purchase is made, if there is any inconsistency or omission with the national policy, the national policy shall prevail.
- For performance issues caused by the product itself, one may avail service as per below regulations with valid proof of purchase:
2.1 A customer can avail warranty service for free within 24 months from the date of purchase if there is a performance failure caused by the product itself.
- In order to avail replacement and warranty services, one should show the valid proof of purchase.
- Other product related materials (such as Packaging box, Product Manual etc.) are not covered in warranty.

Product model:	Serial Number:					
	Name:					
User Information	Tel:					
	Email:					
	Address:					
	Date of purchase:					
Problem Description:	Invoice number:					

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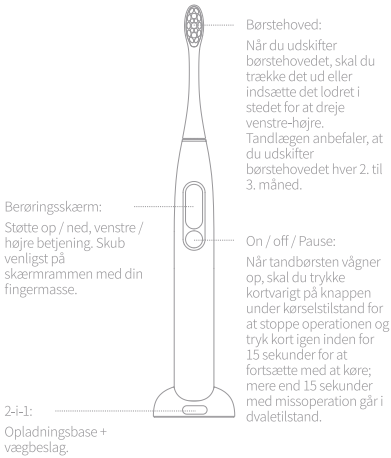
Warranty Card

- These clauses only apply to OcClean products sold in the country or region where the purchase is made, if there is any inconsistency or omission with the national policy, the national policy shall prevail.
- For performance issues caused by the product itself, one may avail service as per below regulations with valid proof of purchase:
2.1 A customer can avail warranty service for free within 24 months from the date of purchase if there is a performance failure caused by the product itself.
- In order to avail replacement and warranty services, one should show the valid proof of purchase.
- Other product related materials (such as Packaging box, Product Manual etc.) are not covered in warranty.

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Follow the Oclean Official Facebook Page

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Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- × The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- * Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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Follow the OcClean Official Facebook Page

OcClean Facebook Page
www.facebook.com/occleanofficial

GET IT ON Google Play

Download on the App Store

Search "OcClean" on App Store or Google Play to download

• Specifications	
Product Name	OcClean Smart Sonic Electric Toothbrush
Model	OcClean X
Edition	Pro
Max Output	0dBm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V ~1A
Battery Life	20 days (Level 1), 40 days (Level 32)
Battery Capacity	800 mAh
Charging Time	2 hours
Bluetooth	4.2 above

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Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- ✗ The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- * Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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