Oclean Smart Sonic Electric Toothbrush
User Manual

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
Oclean Smart Sonic Electric Toothbrush User Manual

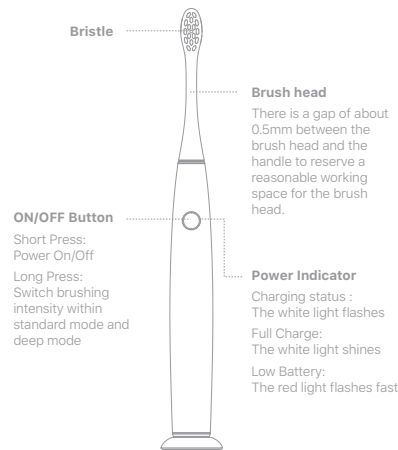Instructions

- 1.Power On**
Start the toothbrush by short pressing the power button, the toothbrush will automatically shut down after a 2-mins brushing plan is finished.
- 2.Intensity Switch**
Oclean Air2 has preset 2 types of 2mins brushing plan: Standard Cleaning(Default) and Deep Cleaning.
To switch the intensity, just long press the power button for 2s under working condition.
- 3.Battery Indication**
0%-95%, the power indicator light flashes. 95%-100%, the Power Indicator light is constantly on.

First Use : Charging time depends on the battery life demand of user.

Working Process :If battery power is less than 5%, the red indicator light will flash as a charging reminder.

Specifications	
Product Name	Oclean Air2 Sonic Electric Toothbrush
Model Name	Oclean Air2
Input Power	2.5W
Charging	Oclean USB Charger
Input	5V $\approx$ 0.5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**))
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V $\approx$ 0.5A.

Shenzhen Yunding Information Technology Co., Ltd.

32G, Tower 3, Dachong Business Building(The Third Phase), Shennan Avenue North, Tong Gu Road West, Yuehai Street Communities, Nanshan District, Shenzhen, Guangdong, China Zip Code: 518057
If there is fault occurring in the product you purchase, please contact Email: service@oclean.com

Facebook Page

www.facebook.com/ocleanofficial

[Follow the Oclean Official Facebook Page](#)

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.
If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
-- Reorient or relocate the receiving antenna.
-- Increase the separation between the equipment and receiver.
-- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
-- Consult the dealer or an experienced radio/TV technician for help.
The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

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Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- ×
- × The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.


Warranty Card

- These clauses only apply to Oclean products sold in the country or region where the purchase is made, if there is any inconsistency or omission with the national policy, the national policy shall prevail.
- For performance issues caused by the product itself, one may avail service as per below regulations with valid proof of purchase:
2.1 A customer can avail warranty service for free within 24 months from the date of purchase if there is a performance failure caused by the product itself.
- In order to avail replacement and warranty services, one should show the valid proof of purchase.
- Other product related materials (such as Packaging box, Product Manual etc.) are not covered in warranty.

Product model:	Serial Number:
User Information	Name:
	Tel:
	Email:
	Address:
	Date of purchase:
	Invoice number:
Problem Description:	

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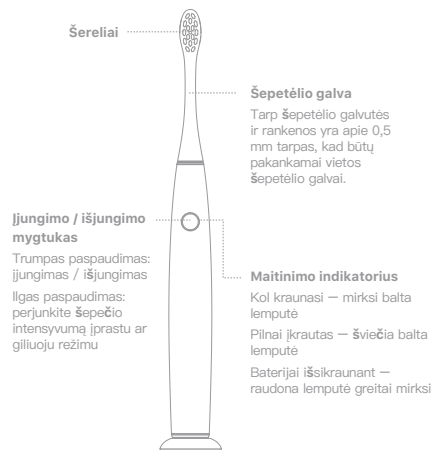Specifications

Product Name	Oclean Air2 Sonic Electric Toothbrush
Model Name	Oclean Air2
Input Power	2.5W
Charging	Oclean USB Charger
Input	5V ≈ 0.5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**)
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V ≈ 0.5A.

Shenzhen Yunding Information Technology Co., Ltd.

32G, Tower 3, Dachong Business Building(The Third Phase), Shennan Avenue North, Tong Gu Road West, Yuehai Street Communities, Nanshan District, Shenzhen, Guangdong, China. Zip Code: 518057

If there is fault occurring in the product you purchase, please contact Email: service@oclean.com

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Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
 - X The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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Facebook Page


www.facebook.com/ocleanofficial

Follow the
Oclean Official Facebook Page

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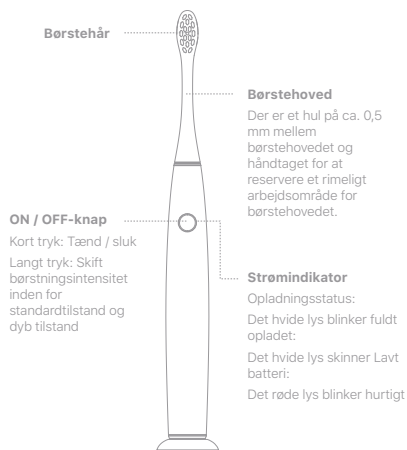Specifications

Product Name	Oclean Air2 Sonic Electric Toothbrush
Model Name	Oclean Air2
Input Power	2.5W
Charging	Oclean USB Charger
Input	5V = 0.5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**))
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V =0.5A.

Shenzhen Yunding Information Technology Co., Ltd.

32G, Tower 3, Dachong Business Building(The Third Phase), Shennan Avenue North, Tong Gu Road West, Yuehai Street Communities, Nanshan District, Shenzhen, Guangdong, China. Zip Code: 518057
If there is fault occurring in the product you purchase, please contact Email: service@oclean.com

Name of parts	Toxic and harmful substance or element						Environment Protection Term
	Pb	Hg	Cd	Cr6+	PBB	PBDE	
PCBA	x	o	o	o	o	o	10 Years
Plastic parts	o	o	o	o	o	o	
Metal Hardware	x	o	o	o	o	o	
Vibrating motor	x	o	o	o	o	o	
Accessories	x	o	o	o	o	o	5 Years
Battery	o	o	o	o	o	o	

Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
 - × The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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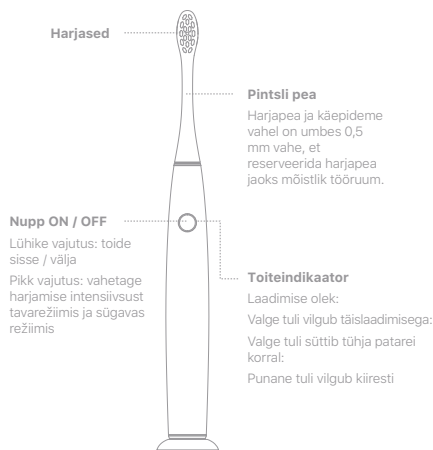
Specifications

Product Name	Occean Air2 Sonic Electric Toothbrush
Model Name	Occean Air2
Input Power	2,5W
Charging	Occean USB Charger
Input	5V = 0,5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**))
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V = 0.5A.

Shenzhen Yunding Information Technology Co., Ltd.

32G, Tower 3, Dachong Business Building(The Third Phase), Shennan Avenue North, Tong Gu Road West, Yuehai Street Communities, Nanshan District, Shenzhen, Guangdong, China. Zip Code: 518057
If there is fault occurring in the product you purchase, please contact Email: service@occean.com

Name of parts	Toxic and harmful substance or element						Environment Protection Term
	Pb	Hg	Cd	Cr6+	PBB	PBDE	
PCBA	x	o	o	o	o	o	10 Years
Plastic parts	o	o	o	o	o	o	
Metal Hardware	x	o	o	o	o	o	
Vibrating motor	x	o	o	o	o	o	
Accessories	x	o	o	o	o	o	
Battery	o	o	o	o	o	o	5 Years

Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
 - × The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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Facebook Page


www.facebook.com/occeanofficial

Follow the
Occean Official Facebook Page

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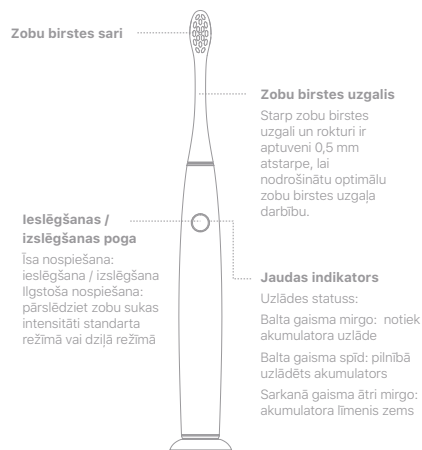Specifications

Product Name	Oclican Air2 Sonic Electric Toothbrush
Model Name	Oclican Air2
Input Power	2,5W
Charging	Oclican USB Charger
Input	5V = 0,5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**)
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V = 0.5A.

Shenzhen Yunding Information Technology Co., Ltd.

32G, Tower 3, Dachong Business Building(The Third Phase), Shennan Avenue North, Tong Gu Road West, Yuehai Street Communities, Nanshan District, Shenzhen, Guangdong, China. Zip Code: 518057
If there is fault occurring in the product you purchase, please contact Email: service@oclean.com

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Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
 - ×
 - ×
- The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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Facebook Page


www.facebook.com/ocleanofficial

Follow the
Oclican Official Facebook Page

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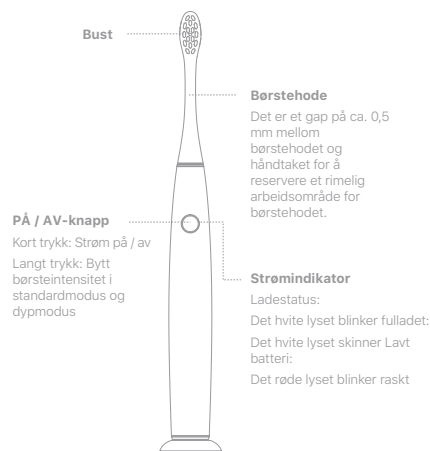Specifications

Product Name	Oclean Air2 Sonic Electric Toothbrush
Model Name	Oclean Air2
Input Power	2.5W
Charging	Oclean USB Charger
Input	5V = 0.5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**)
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V = 0.5A.

Shenzhen Yunding Information Technology Co., Ltd.

32G, Tower 3, Dachong Business Building(The Third Phase), Shennan Avenue North, Tong Gu Road West, Yuehai Street Communities, Nanshan District, Shenzhen, Guangdong, China Zip Code: 518057
If there is fault occurring in the product you purchase, please contact
Email: service@oclean.com

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Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
 - ×
 - ×
- Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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Specifications

Product Name	Oclean Air2 Sonic Electric Toothbrush
Model Name	Oclean Air2
Input Power	2.5W
Charging	Oclean USB Charger
Input	5V ≈ 0.5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**)
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V ≈ 0.5A.

Shenzhen Yunding Information Technology Co., Ltd.

32G, Tower 3, Dachong Business Building(The Third Phase), Shennan Avenue North, Tong Gu Road West, Yuehai Street Communities, Nanshan District, Shenzhen, Guangdong, China Zip Code: 518057

If there is fault occurring in the product you purchase, please contact Email: service@oclean.com

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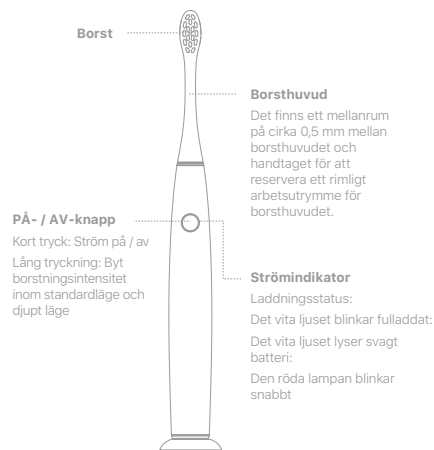
Specifications

Product Name	Oclican Air2 Sonic Electric Toothbrush
Model Name	Oclican Air2
Input Power	2,5W
Charging	Oclican USB Charger
Input	5V = 0,5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**)
Charging Time	2,5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

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*** It is suggested to use a USB charger of 5V = 0,5A.

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Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
 - ✗ The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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