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
RI 201-190289



Oclean Smart Sonic Electric Toothbrush User Manual

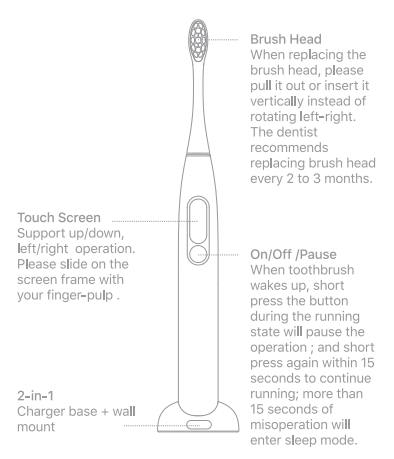First time use  
Short press on ON/OFF button to activate toothbrush and bluetooth, after being activated, the device will be constantly in stand-by status.

Brushing Plan  
Toothbrush will start default mode with “Standard Cleaning Plan” with 2 mins if no new customized brushing plan setting in APP. There will be a reminder with vibration to remind user to change the brushing area in every 30s. Toothbrush will stop operating automatically after 2 mins.

Charging  
Please wipe the toothbrush to dry before charging. Please keep an eye on the battery percentage, only charge it when battery is low, do NOT overcharge or repeat charging frequently.

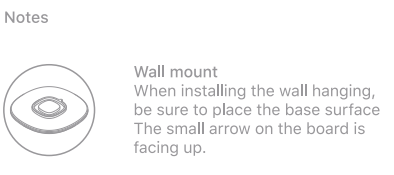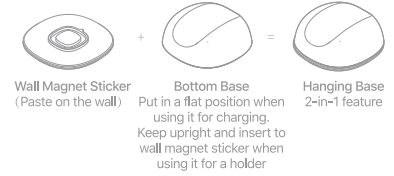Important Safety Tips

- This product is recommended for people over 12 years old.
- This product can only use OcClean brush heads, and will not bear all the consequences of non-certified third-party brush heads.
- This product is only a personal cleaning product and is not recommended for medical use or multi-user use
- This product is for indoor use only. Do not place it in direct sunlight. It is recommended to be placed in a dry and ventilated place.
- Before charging, please keep the body dry, and regularly clean the charging base with a dry cloth to prevent the deposition of stains.

Follow the OcClean Official Facebook Page

OcClean Facebook Page  
www.facebook.com/occleanofficial

GET IT ON Google Play

Download on the App Store

Search "OcClean" on App Store or Google Play to download

Welcome


Thank you for purchasing our product. We recommend you to use it with the supporting app to learn more about maintaining a healthy smile.

Not recommended for use

- Individuals who have recently undergone oral surgery. It is recommended that consumers with oral disease, such as dental cavities, periodontal (gum) disease, oral cancer, oral infectious diseases, consult with a dentist before using this toothbrush.
- Individuals with mental and/or physical disabilities should only be used under supervision.
- People who have mouth diseases, loose false teeth or gums with hemorrhage.
- People with sensitive and susceptible oral cavity, they should consult a professional dentist before using this product.
- The toothbrush is not suitable for people with reduced mental and physical capabilities. Not suitable for infant and children under the age of 12. Children who use this product should be strictly supervised by an adult.

Instructions

- First Time Use**  
The toothbrush can be turned on by pressing the ON/OFF button for the first time on a brand-new toothbrush. The bluetooth function will be activated and remain on. After installing the OcClean App, please complete the prompts to complete the connection and settings for the toothbrush.
- Install App & Toothbrush Settings**  
You can use the app to customize your brushing programs - automating tooth cleaning, tooth whitening and gum massaging in one brushing session. There is no need to open the OcClean App each time you brush. Simply open the app every few days to sync the brushing data.
- Direct Use without App Settings**  
You can also use the toothbrush without connecting to the app. It has a touchscreen display and an ON/OFF button. It offers touch control via the touch-sensitive screen.
- Brushing Performance**  
After brushing, the display will show the effective brushing time, the brushing scores and the tips of the brushing area.
- Check Brushing Records in App**  
Only after the toothbrush is paired with the app and connected to the Bluetooth, can the app get access to the brushing data stored in the toothbrush. Due to limited storage space, the brushing data can usually be stored for no more than 7 days.
- Mode Setting**  
Enter the mode-setting section by sliding up and down on the display. Each of the modes, such as Cleaning, Whitening and Massage, has 32 levels of cleaning intensity. After the new setting is completed, it will automatically register in the memory. Brushing intensity settings: Sliding up/down on the screen, the intensity will be changed by 1 level (+1 or -1); Side-sliding up/down on the screen side, the intensity will be changed by 5 levels (+5 or -5).

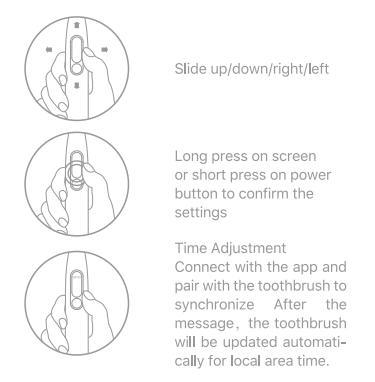- No disassembly, impact, extrusion or throwing into fire
- Please stop using the toothbrush if the battery has expanded
- Please don't put the battery in high temperature environment
- Do not use the battery if water entered the device.

Brushing Plan & App Settings

You can either choose or customize a brushing plan through the supporting OcClean App. You can also upgrade the firmware to continuously update its performance.

Low Battery Prompt

When battery percentage is lower than 20%, you'll see low battery prompt appears on your toothbrush display frequently, you will need to charge the toothbrush. When the battery percentage is lower than 1%, you will be prompted by motor vibration, and the screen will not be on.

How to Charge

Remove the charging base from the wall mount and connect it to the mobile phone charger. When the bottom of the toothbrush is placed gently above the charging base, it will automatically position itself. The display will show the battery percentage when the toothbrush is charging. A full charge usually takes about two hours.

Important Charging Tips

The toothbrush runs with a fast-charge mechanism. It is recommended that the toothbrush be placed on the charging base only when it is necessary. Please do not put the toothbrush on the charging base for a long time, as it may cause the built-in battery to decay prematurely so that the user experience is affected.

Recommended Tooth Brushing Method

Cleaning and Maintenance

Specifications	
Product Name	OcClean Smart Sonic Electric Toothbrush
Model	OcClean X
Edition	Pro
Max Output	0dBm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V=1A
Battery Life*	20 days (Level 1), 40 days (Level 32)
Battery Capacity	800 mAh
Charging Time	2 hours
Bluetooth	4.2 above

\* Calculated by 2 minutes/time, twice/day.

Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

PS E IPX7 EAC CE BC FC R 201-190289

Toxic and harmful substance or element		Environment Protection		Term	
Name of parts		Pb	Hg	Cd	Cr6+
Vibrating motor		○	○	○	○
Accessories		×	○	○	○
Battery		○	○	○	○

Toxic and harmful substance or element		Environment Protection		Term	
Name of parts		Pb	Hg	Cd	Cr6+
PCBA		○	○	○	○
Plastic parts		○	○	○	○
Metal Hardware		○	○	○	○
FPC		○	○	○	○

Guarantee Certificate of OcClean Smart Sonic Electric Toothbrush

User Information

Name: \_\_\_\_\_ Tel: \_\_\_\_\_ E-Mail: \_\_\_\_\_

Address: \_\_\_\_\_ Zip Code: \_\_\_\_\_

Product Information

Type: \_\_\_\_\_

Production Place: \_\_\_\_\_

Sales Information

Name of Retailer: \_\_\_\_\_ Tel: \_\_\_\_\_

Address: \_\_\_\_\_ Zip Code: \_\_\_\_\_

Invoice No: \_\_\_\_\_ Sales Date: \_\_\_\_\_

Fault description/Reason/Disposal		Employed equipment		Acceptance date		Serviceman	

Please ask for a maintenance report from the maintenance center and properly keep it with the warranty certificate. Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers guarantee for customers who purchase OcClean Smart Sonic Electric Toothbrush. In case of inconsistency between the service articles and the provisions of guaranteee, the provisions of guaranteee and the contents in the guarantee certificate shall govern.

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

FCC ID: 2AN5D-OCLEANX

IC WARNING

This device contains licence-exempt transmitter (s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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IC : 24982-OCLEANX

CAN ICES-3(B)/NMB-3(B)

If the product you have purchased is found to have a fault or defect, please contact us.

North America  
Email: support.us@oclean.com

Other Region  
Email: service@oclean.com

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd.

Manufacturer's address: #32G, Building 3, Dachong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China

Production address: The East of 4th Floor, Building Four, No.89, Dingshanzhong Road, Dongfu Street, Haicang District, Xiamen, China

Hereby, Shenzhen Yunding Information Technology Co., Ltd. declares that the radio equipment type Y3 is in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address: <http://www.oclean.com/en/service.aspx>

For more regulatory information and product certification/compliance logos related to OcClean X Pro, please go to Settings->About->Regulatory on your toothbrush.

SN