



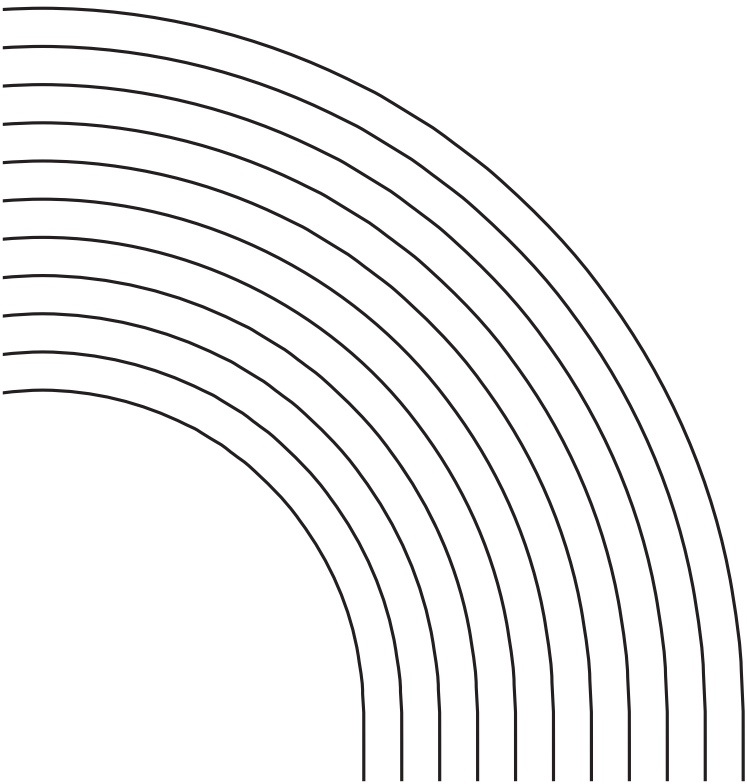


Oclean


Product User Manual

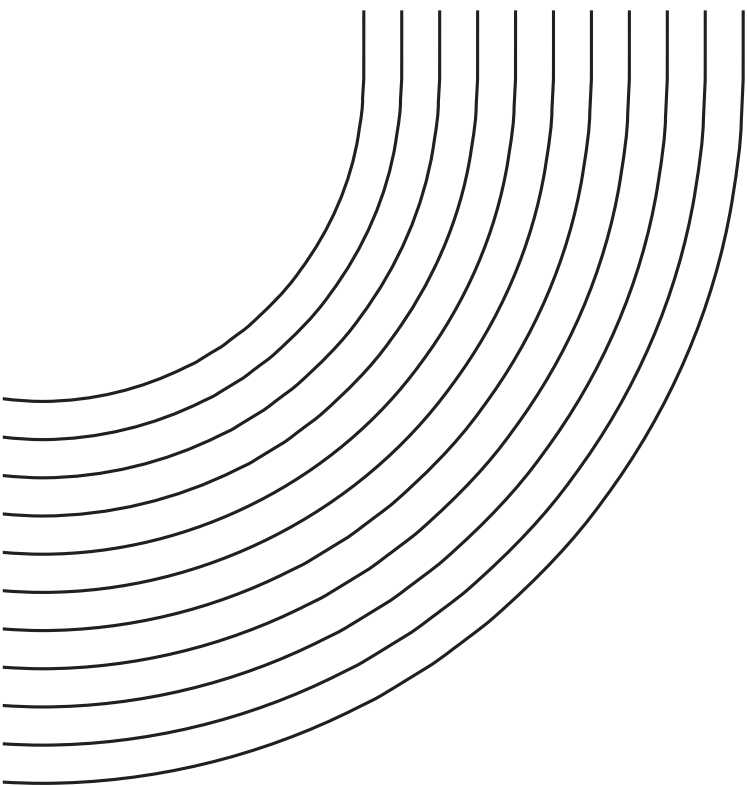


Toothbrush Body

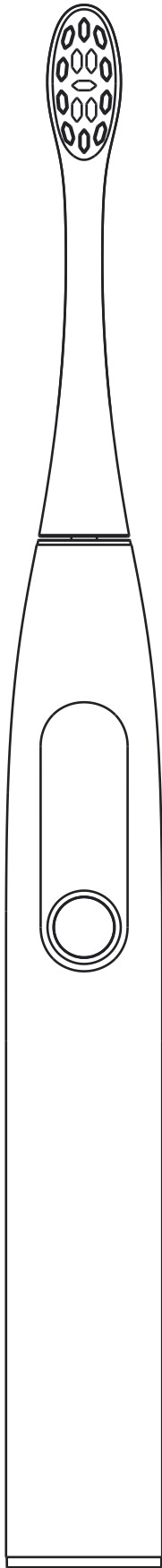


Brush Head

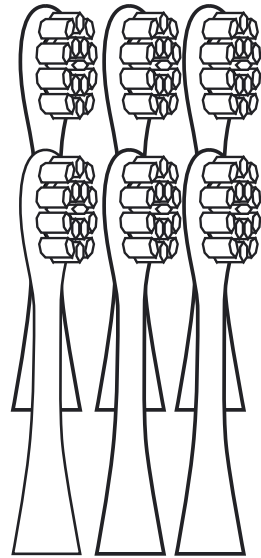


Wall-Mount

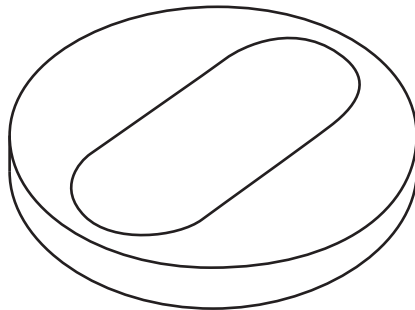

Charging Base

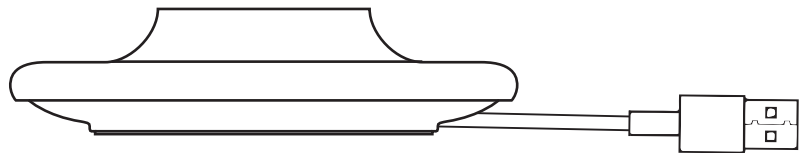


Travel Case

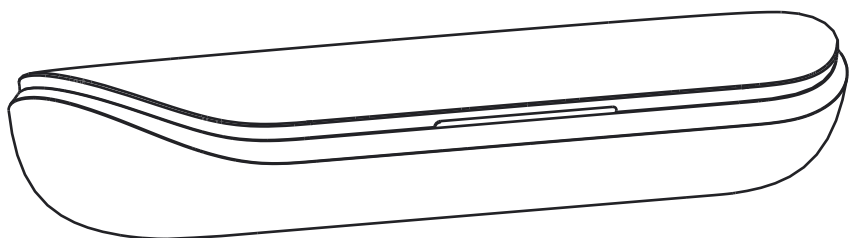

Oclean X Pro Elite Premium Set

Important Safety Tips

1. This product is recommended for individuals aged 12 and above. It only works with Oclean designated brush heads. We will not be responsible for any consequences arising from the use of non-certified third-party brush heads.
2. This product is for personal hygiene purposes only. It is not recommended for medical use or for use by multiple users. Avoid direct sunlight exposure and store it in a cool, dry, and ventilated area.
3. Ensure the device is dry before charging. Use a soft cloth to wipe it clean regularly, including the charging base, to prevent dirt buildup.
4. Individuals who have recently undergone oral surgery or are currently undergoing treatment, including those with loose dentures or significant gum bleeding, are advised against using this product.
Consult a professional dentist before use.
5. This product is not suitable for individuals with disabilities, mental illnesses, or those lacking relevant experience and knowledge, including children (0-12 years), unless under the supervision of a guardian.

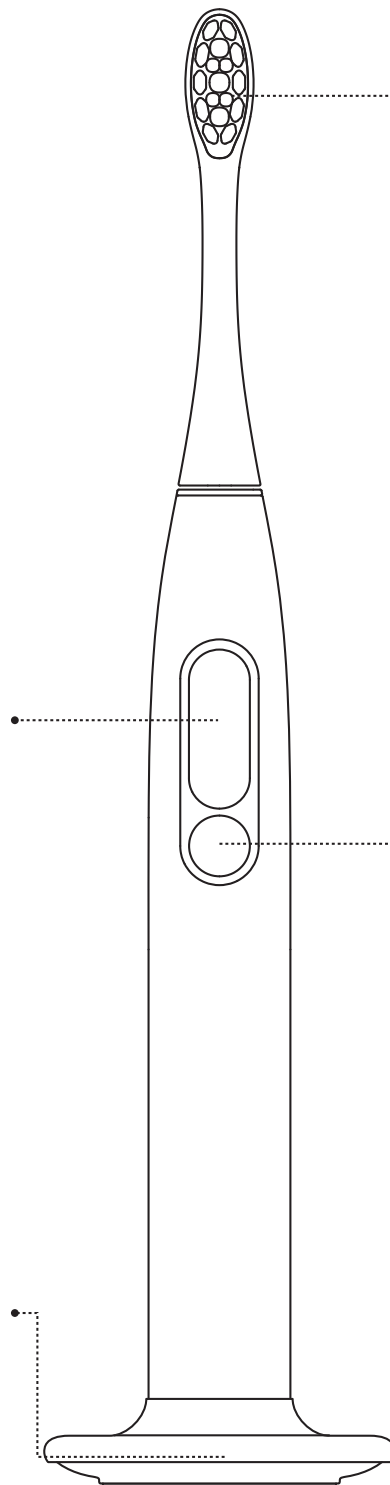

Brush Heads

When replacing the brush head, please pull it out or insert it vertically without twisting left or right. Dentists recommend replacing the brush head every 2 to 3 months.

Touch Screen

Supports vertical scrolling. Swipe up and down at the edge of the screen using your finger.

Power Button

After waking up the toothbrush, a short press during operation will pause it. Press again within 15 seconds to resume. If there's no operation for more than 15 seconds, the screen will go black and enter sleep mode.

Wireless Charging Base

The wireless charging base should be used together with the silicone holder.

User Guide

When using the toothbrush for the first time, briefly press the "Power Button" to wake it up. Bluetooth will remain on continuously. You can connect and configure the toothbrush by downloading and installing the "Oclean Care+" app. (Note: Connecting through "System Settings - Bluetooth Settings" on your phone will not work.) Once the app is installed, you can access various functions such as setting up the "Zone Change Reminder."


Additionally, you can customize your brushing plans and directly sync them to the toothbrush.

- 1 Scan the QR code or search for Oclean Care+ on your Apple Store to download.




Google Play Store & iOS APP Store


- 2 New users tap "Register account" to create an account.

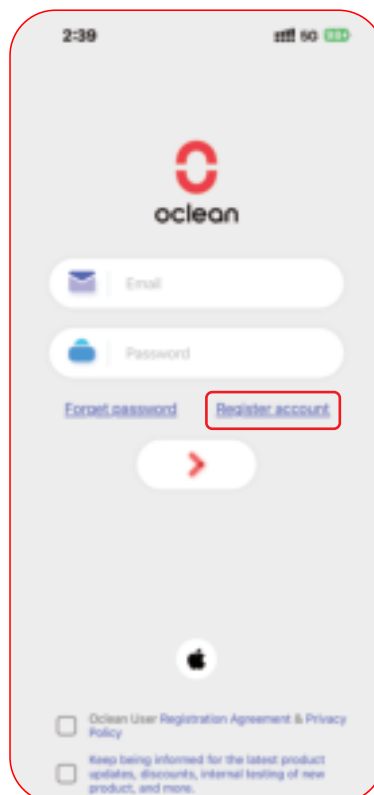

- 3 Enter your username and password, then tap ">".

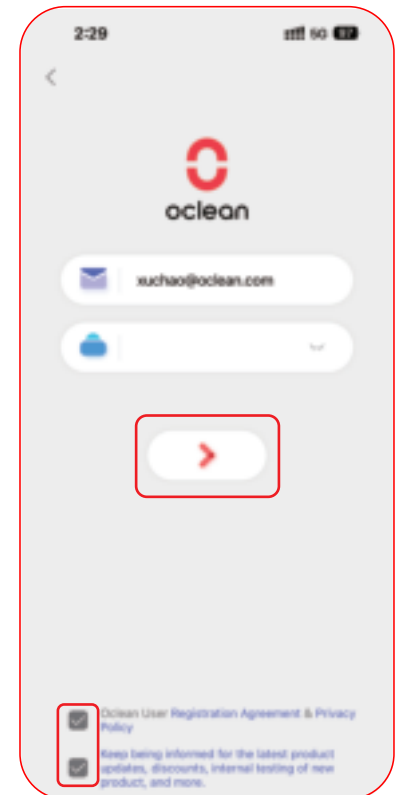


- 4 Enter the verification code, then tap "Completed".

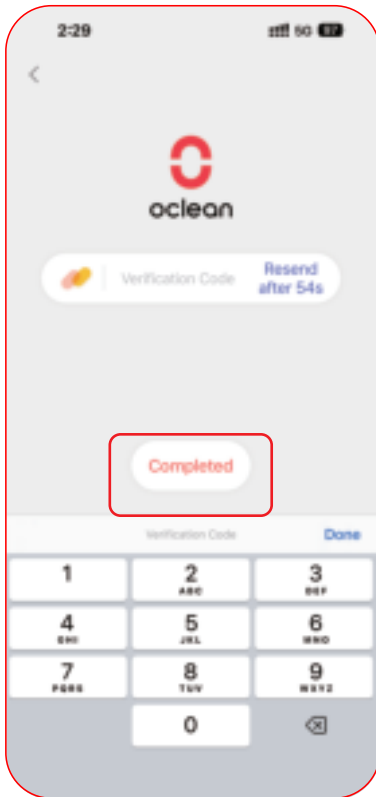

- 5 Bring the toothbrush close to your phone.




- 6 Open the app and tap "+" to go to the device connection page.

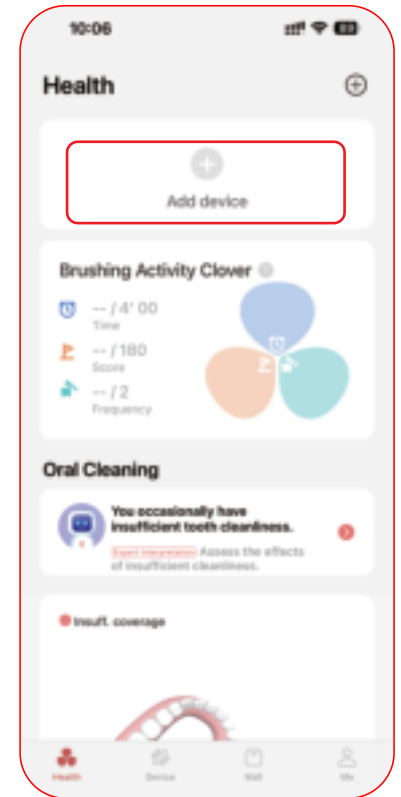


- 7 Tap "Next" to enter the search page.




- 8 Tap the first toothbrush device on the page.

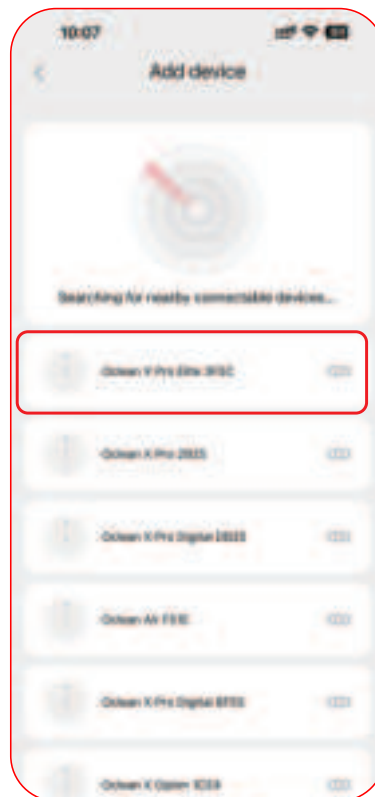


- 9 Press Set later to finish.

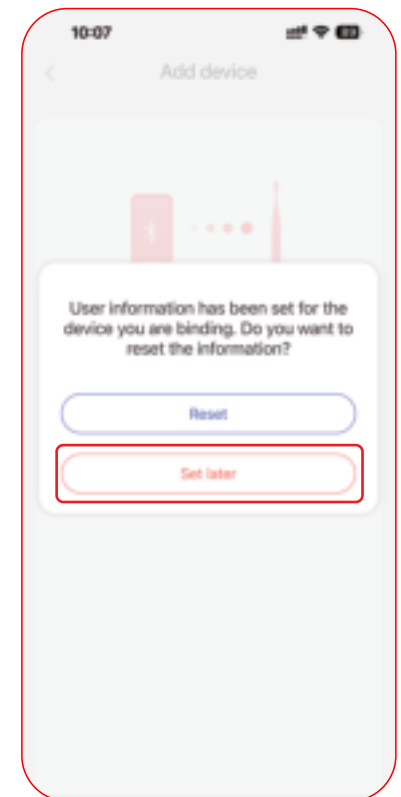

10 Tap the toothbrush device to enter the settings page.

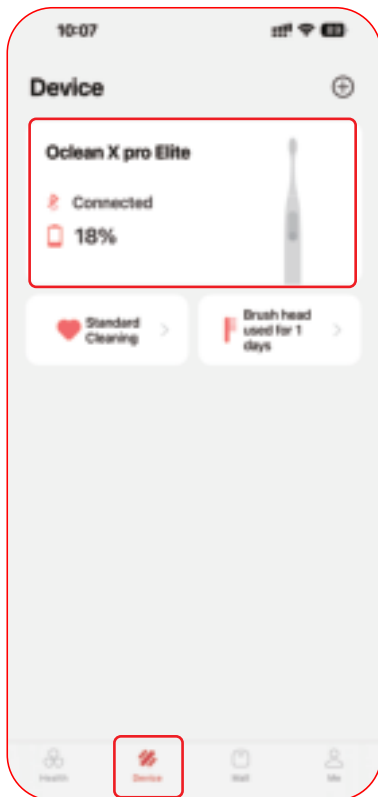


11 Enter the settings page to configure your personal preferences.

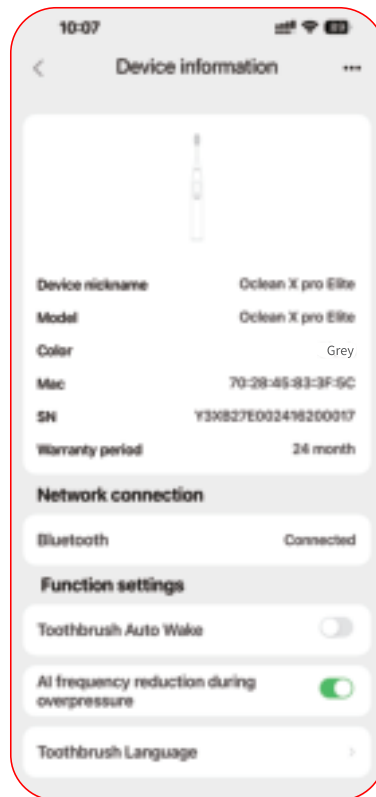

Touch Screen Toothbrush Settings

You can use the toothbrush by sliding up and down on the touch screen to select your preferred mode, duration, and intensity settings without connecting to the app. Once set, the toothbrush will automatically remember your preferences and use the new settings for future operations. You can also customize your brushing plan through the "Oclean Care+" app, where you can set your preferred mode, duration, and intensity. Once the app has sent a brushing plan to the toothbrush, it will follow these settings. If you make any changes to the mode or intensity via the touch screen after the app has sent a brushing plan, the

toothbrush will no longer follow the app's plan and will use your new settings instead.

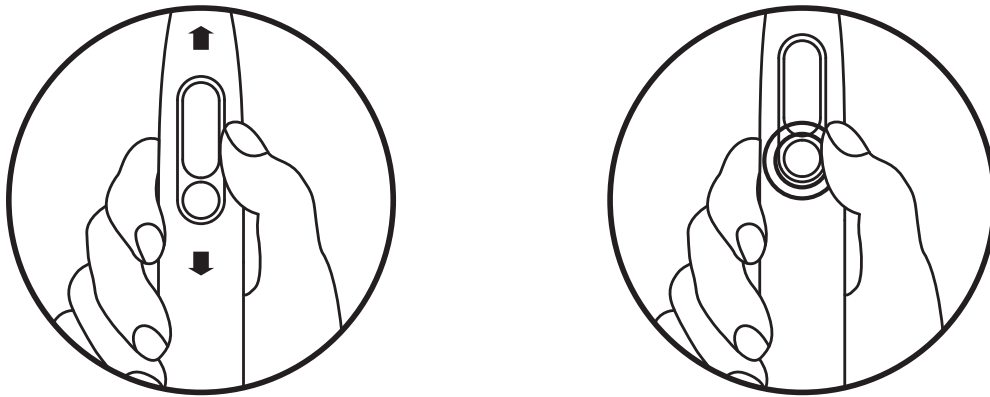

Timer and Brushing Zone Reminder

The default brushing plan has a duration of 2 minutes. Every 30 seconds, there will be a brief change in the vibration pattern to remind the user to switch to a different area of the teeth.

Brushing Report

After brushing, the product's display screen will show the duration of your brushing session, your brushing score, and any blind spots you may have missed.

Viewing Brushing Records on the App

The app can access brushing data only when the toothbrush is paired and connected to it. The toothbrush has limited storage capacity and can store brushing data for up to 7 days.


Wall Mount

The toothbrush comes with a wall mount, providing a space-saving storage solution.

Note

Wall Mount Installation:

When placing the wall mount on the wall, ensure that the arrow is pointing upwards during installation.

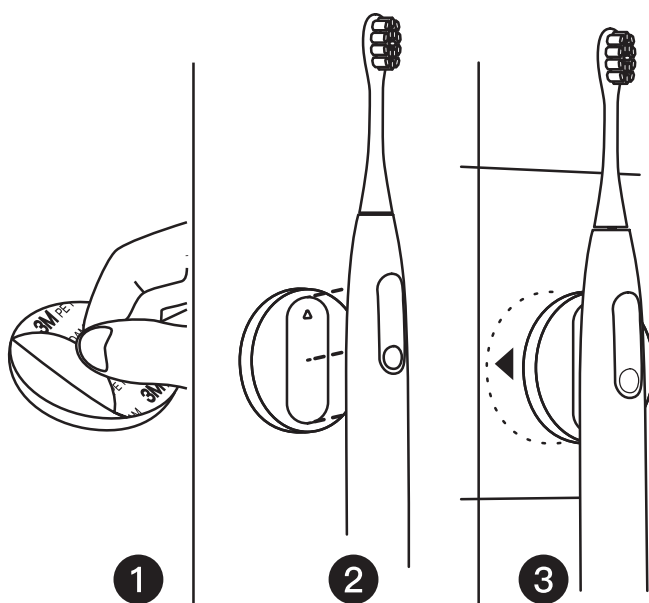

Low Battery Warning

When the toothbrush battery is below 20%, the display will frequently show a low battery warning, indicating that it needs to be charged. If the battery level falls below 1%, the toothbrush will not turn on and will alert

you with motor vibrations without lighting up the display.

Charging Instructions


1. Place the dedicated wireless charging base on a flat surface and connect the power plug. The indicator light on the base will flash when powered on.
2. Align the small silicone notch on the charging base with the bottom indicator section of the toothbrush handle and place the handle in the charging slot to begin charging.
3. Check the charging status using the indicator light on the charging base. When powered, the indicator light will flash. If charging is not occurring, the indicator light will not light up.
4. When charging is normal, the indicator light will display a steady white light. If there is an object on the wireless charger or if charging is not proceeding correctly, the indicator light will flash red rapidly.

Important Charging Notice

This product uses a wireless charging mechanism. It is recommended to place the toothbrush on the wireless charger only when it needs to be charged. When not charging, please store the toothbrush elsewhere.

Keeping the toothbrush on the powered wireless charging base for extended periods may lead to premature degradation of the internal battery, potentially affecting the user experience.

Recommended Brushing Method

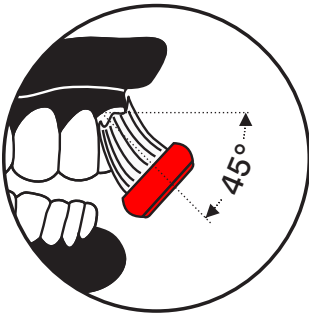

Position the brush head at a 45° angle to the gum line and move the toothbrush slowly. During this process, manual up-and-down brushing is not necessary.



Hold the toothbrush handle gently and move it in a circular motion, avoiding excessive pressure.



After brushing, rinse your mouth with clean water.

Battery Removal Instructions

1. Ensure the product is powered off before removing the battery.
2. Dispose of the battery safely after removal.
3. A professional should follow these steps to remove the battery from the product:
 - Use a heat gun to warm the charging board, then use tweezers to pry open the charging board along the edge of the bottom cover.

- Use a cross screwdriver to remove the screws securing the bottom cover and take it off.
- Use a tool to pry the clips on both sides inward and push the motor shaft forcefully until the mechanism is removed.
- Use wire cutters to cut the charging wires at the bottom of the battery, then pry open the clips and remove the battery.

Battery Safety Warnings



1. Do not disassemble, impact, crush, or expose to fire.
2. If the product shows severe swelling, do not continue to use it.
3. Do not place in high-temperature environments.
4. Do not use if the battery has been submerged in water.

Product specifications

Product Name	Oclean Smart Sonic Electric Toothbrush
Model	Y2088
IP Rating	IPX7
Power of the toothbrush body	1.5W
Charging	Wireless Charging
Input	5V $\equiv$ 0.3A
Operating time*	20 days (level 32 of intensity) - 40 days (level 1 of intensity)
Battery Capacity	800 mAh
Implementation standard	GB 4706.1-2005,GB 4706.59-2008
Bluetooth	Bluetooth 5.0 above

* Under the cleaning mode, it is calculated at 2 minutes/-time, 2 times/day working condition.