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INSTRUCTIONS FOR USE

Open the cover on the LED device
And remove the plastic safety
catch.


Attach the mouthpiece to
the LED device.

Turn the cap of the gel
before its application.

Apply ¼ of the gel on the top
surface of the mouthpiece according
to the marked line.

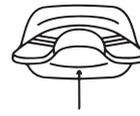Apply ¼ of the gel on the bottom
surface of the mouthpiece according
to the marked line.

Let it whiten for 10 minutes.
Rinse and repeat the next day.

Use daily for 6 consecutive days. Every gel contains amount for 2 days. For long lasting results repeat the 6-day treatment every 2 weeks.

How to do it?


①

Brush your teeth thoroughly with a
brush and dry your teeth afterwards.

②

Open the battery cover on the device,
remove the plastic safety catch
and return the batteries back.

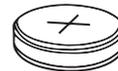
③

Mount the mouthpiece
on the device.

④

Apply ¼ of the gel on the top and
bottom rounded walls of
the mouthpiece so that the gel covers
the front side of your teeth.

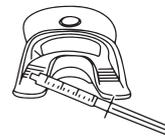⑤

Insert the mouthpiece between
your teeth and turn the device on.


⑥

The LED light will turn off
after a 10-minute cycle.

⑦

Remove the mouthpiece from the
lamp and wash it thoroughly.

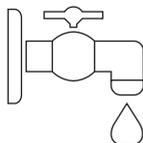
⑧

Smile!

Tips for better results

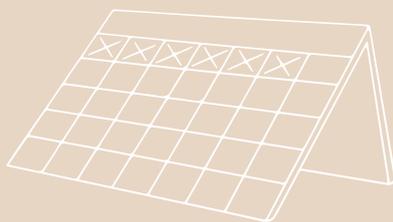Daily use is recommended
for the best results.

Reduce consumption of dark
beverages such as tea, red wine
or coffee.

During the whitening process,
keep your teeth slightly pressed together,
applying slight pressure to the mould.

EFFECTIVE SUBSTANCES

Pomegranate

Chamomile

Aloe Vera

Sodium Perborate and Hydroxyapatite

Teeth whitening gel

Ingredients: Sorbitol, Aqua, Propylene Glycol, Glycerin, Potassium Nitrate, PEG-8, Phthalimidoperoxycaproic Acid, Hydroxyapatite, Cellulose Gum, Hydroxyethylcellulose, Xanthan Gum, Menthol, Sodium Saccharin, Methylparaben, Chamomilla Recutita Extract, Aloe Barbadensis Leaf Juice, Sodium Bicarbonate, Punica Granatum Seed Extract, Propylparaben.

WARNING: Keep out of reach of children. This product is not formulated or recommended for use by children or persons under 12 years of age. All Toothy® brand products are not intended as a substitute for routine dental care. Consult your dentist for product use if you are under 18 years of age, a pregnant woman, have serious dental problems, or have had oral health related surgery. Do not swallow the gel. Avoid contact with eyes. In the case of contact, rinse eyes properly with clean water. Contact with sensitive skin may cause temporary irritation. The gel will not bleach crowns, fillings or dentures.

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