

CONTENTS Thank you for selecting this blue®m sonic+ toothbrush. Please carefully read the safety tips of the product before using it and properly keep it. Important safety tips Fitting How to use? Maintenance Charging Functional operation Instruction modes Specification Fault diagnosis / Warranty & support Hazardous substances tests Environmental protection 11

IMPORTANT SAFETY TIPS

When using electrical products, especially when children are present, basic safety precautions should alwaysbe followed,

- Unplug the charging base or charging travel case immediately if it falls into water.
- damaged cord or plug.

 Marnings To reduce the risk of burns, electrocution, fire, or
- physical injury: Do not use attachments other than those recommended by
- Voltage converters do not warrant voltage compatibility. Never force the plug into an outlet; if the plug does not easily fit into the outlet, discontinue use.
- way (brush head, handle, charging base or charging travel
- Do not use the charging base or charging travel case
- this document or as recommended by your dental
- children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning
- Children should be supervised to ensure that they do not
- patients in a dental office or institution.
- charging travel case.
- When you use a computer to recharge the toothbrush handle, the computer should be located away from sources of water (i.e., sink or tub) so as to avoid electrical shock.

- device manufacturer prior to use.

 Consult your physician/dentist if you have other medical

FITTING

including the following:

- Dangers To reduce the risk of electrocution:
 Always place and store the charging base or charging travel case at a safe distance from water.
- Never use a charging base or charging travel case with a

- Do not plug the charging base or charging travel case into an outlet with a voltage other than specified on the bottom of the charging base or charging travel case. This product is designed to operate within a range of 100 to 240 volts.

 When compared the traverse telegrape compared the little of the charging travel case compared to the charge telegrape of the charge telegrape telegrape.
- Discontinue use if the product appears damaged in any
- Do not use the charging base or charging travel case if it has dropped into water.

 Keep the cord away from heated surfaces.
- This product contains no user-serviceable parts. Refer to
- Warranty and support' (page 9) if the appliance no longer works properly or needs repair.


 Use this appliance only for its intended use as described in
- professional.

 This appliance is not intended for use by persons (including
- use of the appliance by a person responsible for their
- play with the appliance. blue m sonic+ toothbrush is a personal care device and is not intended for use on multiple
- Use UL-listed computer, USB hub, or USB duplex receptacle when you charge the toothbrush handle in the

- Medical Warnings
 Consult your dental professional before you use this product if you have had oral or gum surgery in the previous
- months.
 Contact your dental professional if excessive bleeding

- Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
 Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
 This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.



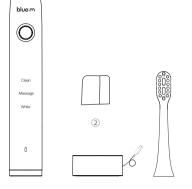

| Number | 1) | 2 | 3 | 4 |
|----------|--------------|------------------|----------------------|--------------------|
| Name | Main body | Brush head cover | Inductive Charger | Toothbrush head |
| Quantity | 1 | 2 | 1 | 2 |

Note: Please carefully check the fitting and its number. If a part is missing, contact our customer service.

HOW TO USE?

blue®m recommends BASS brushing method, brush teeth at least two times every day, each time no less than two minutes. Follow the right brushing steps and make sure to always be very careful when taking care of the hard and soft tissues in your mouth.



Put the brush head into the output axle carefully. Leave a space between the handle and the main body, which is necessary for the vibration of the brush head. toothpaste. It is relatively liquid and makes your toothbrush more effective in removing plaque. Be careful: toothpaste can splash. First put the toothbrush into your mouth, before starting.



Select the right brushing mode. Press once to clean (standard mode), press twice to massage and press three times for whitening. Place the brush head in a 45 degree angle and move the brush slowly, use minimal pressure, from teeth to teeth.

bristles reach between your teeth. Work from quadrant to quadrant, the timer will remind you to switch to each of the 4 sections in your mouth.

* Note: We strongly recommend to use the massage mode after cleaning and extend your oral care ritual with an extra minute. Press the button twice and gently indulge your gums after cleaning. Carefully, with minimal pressure, place the brush on your gums, move from one side to the other and gently stimulate your gums. The zone reminder wil give a signal to switch from upper to lower jaw or vice versa.

MAINTENANCE



TO CLEAN THE BRUSH HEAD head after using it. Let it dry naturally. • Do not pull bristles out of the brush head to avoid deformation and falling out.



TO CLEAN THE MAIN BODY

 Clean your toothbrush with soft cloth and let it air drying. Do not immerse it in water for a long time. • Do not clean with water or acid cleaning solution which is higher than 40°C.



TO CLEAN THE CHARGER



• Use a soft cloth to wipe the stains of the charging base. Before cleaning recharging bed,cut power off for safety.



TO REPLACE THE HEAD




* Regular cleaning will maintain the best condition

of the product and lengthen its service lifetime.

CHARGING

- Charge for 24 hours before your first use. Place the toothbrush in the charging base when the charging indicator is red it means your
- sonic+ is charging.
 It takes approximately 3-4 hours to recharge an empty toothbrush.
- When the indicator turns green, it means charging is completed and your sonic+ toothbrush is ready · When the toothbrush shakes continuously for three
- times it means you have low power, make sure to After a long time of non-recharging, the indicator may light up after 6-10 minutes. That is a normal

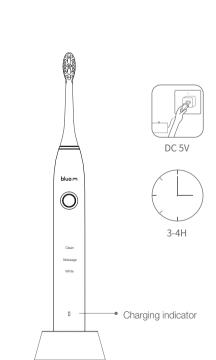

* Charging safely will also guarantee the service lifetime of the battery.


FUNCTIONAL OPERATION

Shortly press the Start button once to turn on the toothbrush. After working for 10 seconds, seconds, shortly press the button again to turn it off.

When the toothbrush is starting in a certain mode, shortly press the button once within 10 seconds to switch to

another mode.

After the toothbrush starts up, the toothbrush will automatically read the intelligent memory mode after the last brushing and start with the brushing mode you were using during your last brushing session

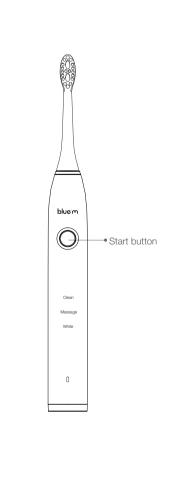

INSTRUCTION MODES

* CLEAN MODE 31000 movements/min With the high-frequency sound wave, a large of micro bubbles will be developed to generate the strong impact force, for deep cleaning the oral cavity, and easily and fully


cleaning the stubborn tartar.

* MASSAGE MODE 28000-31000 movements/min Gradually varying sonic vibration frequency, changes from weak to strong, brings comfortable massage experience to

both teeth and gums. We strongly recommend to use the massage mode after cleaning and extend your oral care ritual with an extra minute.

WHITE MODE 35000-38000 movements/min

The mutual switch of weak and strong frequency of both sound waves can easily loosen all the teeth stains on the surface of teeth to realize the polishing and cleanness under

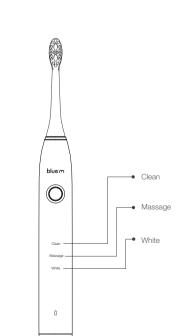

SPECIFICATION

| Name of the product: | blue°m sonic+ toothbrush |
|----------------------|--------------------------|
| Category: | Sonic Toothbrush |
| Model: | LBT-203542B |
| | |

| MAIN BODY | |
|----------------|-------|
| Rated voltage: | 3.7V= |
| Rated power: | 1.8W |
| CHARGING | |

Rated voltage:

Waterproof grade: IPX7

Execution standard: Q/LBD 1-2019 Production date: Please see the product packaging Production place: Dongguan, Guangdong, China.

* Please carefully read the parameters in the specification and use the product correctly.

FAULT DIAGNOSIS

| r dait priorioriori | Caaco anaryoro | Liiiiiiiiddioiiiiiiddiod |
|---|--|--|
| Weakening of the vibration. | ①The toothbrush is forcefully pressed on the teeth. ②The brush head is too tight to the main body. | ①Press the brush head lightly on the teeth. ②Take the brush head from the handle and then reset it on the output axle. Keep a space between the brush head and the handle. |
| The toothbrush can only be used for a short time even if it is fully charged. | The toothbrush is out of power. | ① The toothbrush isn't fully charged. ② The toothbrush has not been used for more than 6 months. ③ The battery has reached it service lifetime. |
| | | |

Fault phenomenon Cause analysis Elimination method

* Warranty & Support

There is a two-year warranty applicable on the blue®m sonic+ toothbrush. If you experience unsolvable problems, please contact your local distributor or the blue®m customer support desk via info@bluemcare.com.

HAZARDOUS

| Name an | d Co | | | h Hazard | dous Sub | stance |
|---------------------|---------------------|----|----|----------|----------|--------|
| Component | Hazardous Substance | | | | | |
| Name | Pb | Hg | Cd | Cr6+ | PBB | PBDE |
| Structural Part | 0 | 0 | 0 | 0 | 0 | 0 |
| PCBA Assemblies* | × | 0 | 0 | 0 | 0 | 0 |
| Power Adapter | × | 0 | 0 | 0 | 0 | 0 |
| Power Cord | × | 0 | 0 | 0 | 0 | 0 |
| | | | | | | |

Supporting Equipment 0 0 0 0 0 0 *PCB assemblies include the PCB and its electronic components, such as resistor, capacitor, integrated chip, etc.

This table is prepared in accordance with the provisions of SJ/T11364-2014 "Marking for the restriction of the use of hazardous substances in electrical and electronic products"

: Indicates that the content of the hazardous substance in all homogeneous materials of the component complies with the requirements of concentration limits stipulated in GB/T26572-2011 "Requirements of concentration limits for certain restricted substances in electrical

x: Indicates that the content of the hazardous substance in at least one homogeneous material of the component exceeds the requirements of concentration limits stipulated in GB/T26572-2011 "Requirements of concentration limits for certain restricted substances in electrical and electronic products".

be realized for the time being due to the global limitation of industrial technology.

10

ENVIRONMENTAL PROTECTION

This toothbrush with a built-in rechargeable battery contains substances, which may pollute the environment. Our mission is a healthier planet, therefore take proper care of disposing the blue m sonic+ toothbrush when the life cycle has ended. Please don't put dumping products and general life garbage together, you should remove the battery into the battery recycling barrel and the blue m products shall be delivered to the corresponding recycling center, which is conductive to





11