
# *beper*

IT	TAPPETINO PER AGOPRESSIONE	pag. 2
EN	ACUPRESSURE MAT	pag. 3
FR	TAPIS D'ACUPRESSION	pag. 4
DE	AKUPRESSURMATTE	pag. 5
ES	ALFOMBRA DE ACUPRESIÓN	pag. 6
GR	ΧΑΛΑΚΙ ΒΕΛΟΝΙΣΜΟΥ	pag. 7
RO	COVORAȘ DE PRESOPUNCTURĂ	pag. 8
CZ	AKUPRESURNÍ PODLOŽKA	pag. 9
NL	ACUPRESSUUR MAT	pag. 11
LV	PUNKTU MASĀŽAS PAKLĀJIŅŠ	pag. 12

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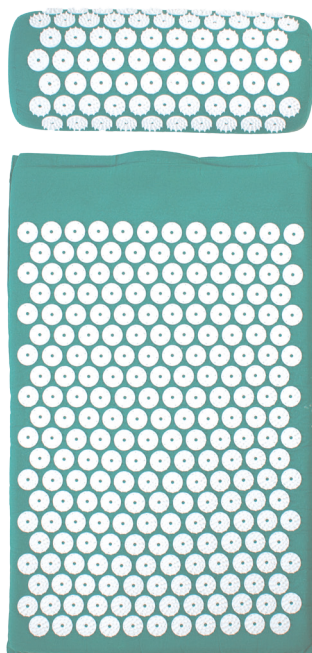---

### Tips for use

The first use on a new body area might cause a feeling of discomfort, but in time the body will adjust and the unpleasant feeling will go away.

The effectiveness of the massage depends a lot on the treated area: some body areas can require just a few minutes' treatment, while others, such as the back area, might require longer sessions to achieve the desired benefits. When treating a new body area, it is possible to regulate the intensity of the massage by varying the body weight: for example, the first time using the mat on the feet, it is possible to sit on the couch to reduce the pressure on the soles of the feet.

It is advisable to use the mat at least 3 times a week, for 15-20 minutes; frequency and duration of treatment can vary depending on the body area to be treated and on personal needs. Do not use on an inclined plane and do not move excessively or roll your body on the mat during use.

Once treatment is done, it is advisable to do some light exercise to gradually recover muscle tone and circulation.

To achieve maximum benefit, the mat must be used on bare skin; however, the sensitivity to the pain caused by the skin acupressure varies from individual to individual, and can be accentuated by certain psychophysical conditions, e.g. drugs and medicine consumption, or in case of specific medical conditions that pertain to the nervous system.

During first uses, it is thus advisable to limit the treatment with sessions of max. 10 minutes, and to wear a thin top or to lay a light cloth on the mat, to check your level of tolerance of the stimulation of the tips and the reaction of your body to the treatment, giving the body the time to adapt to the stimulation.

The use of this mat is not recommended in case of pregnancy and/or generally in case of injuries or irritations of the skin, such as sun burns, inflammation, dermatitis, eczema, psoriasis, or open wounds; chronic diseases such as rheumatoid arthritis, thrombophlebitis, or medical conditions that require use of blood thinning medicines; bleeding issues (haemophilia) or conditions that require the use of anticoagulants such as warfarin (Waran), enoxaparin and heparin (even just for the risk of skin bleeding or bruises); any heart disease or condition; low blood pressure; skin hypersensitivity of any kind.

**For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.**

### GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

**The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance**

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

### GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge. In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

Contact your distributor in your country or after sales department Beper.

E-mail [assistenza@beper.com](mailto:assistenza@beper.com) which will forward your enquiries to your distributor.

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