
CURASEPT®

SONIC

OPUSCOLO ISTRUZIONI
INSTRUCTIONS LEAFLET
NOTICE D'UTILISATION
INSTRUKCJA
NÁVOD NA POUŽITIE

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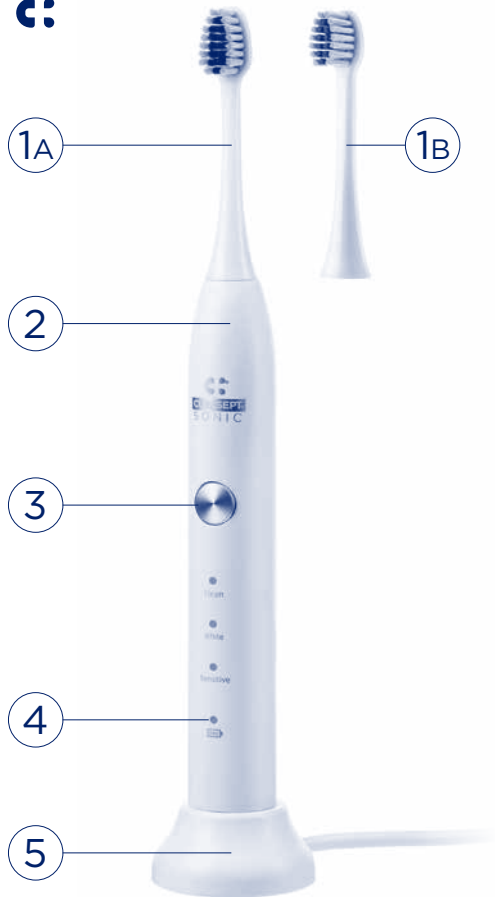Thank you for purchasing CURASEPT SONIC.

A **sonic toothbrush** is a simple and safe tool, ideal for daily oral hygiene. It operates through a fluid-dynamic and mechanical action that generates a flow of microbubbles capable of reaching all tooth surfaces.

Sonic technology toothbrushes are particularly suitable for everyone, even for people with manual coordination issues with a conventional toothbrush. They are also recommended for people who have issues with gingival recession, dentin hypersensitivity, abrasions, as - compared to conventional electric toothbrushes - it suffices to apply slight pressure on the teeth and gums.

Before using the product, please read the instructions carefully, and keep them for further reference.

The manufacturer is not liable for incorrect maintenance and use of the device.

Technical features:

Product:	Curasept Sonic
Modes:	Clean/White/Sensitive
Switch:	Single button
Battery life:	60/180 minutes
Recharge time:	9/16 hours
Voltage:	4 V
Wattage:	3 W
Recharge current:	40/70 mA
Battery:	Lithium/14500/600 mAh
Charge indicator:	Flashing light
Waterproof:	IPX7

The pack contains:

one toothbrush handle, two heads of different sizes, one charger, and one case.

Components (Image page 1):

- ①A Head Std. Soft Cleaning 1015
- ①B Head Compact Sensitive 1012
- ② Handle
- ③ Power/mode change button
- ④ Battery charge indicator
- ⑤ Stand with in-built induction charger

HEAD TYPES

The pack contains two heads; both are made with soft bristles, ideal for efficient and trauma-free cleaning.

①A **STD SOFT CLEANING 10.15:** standard-size head made of outer white ultra-soft bristles and inner blue medium-soft bristles.

①B **COMPACT SENSITIVE 10.12:** the head is smaller than the standard one and, therefore, suitable for people with small mouths or tooth sensitivity issues. In fact, it is made of outer white ultra-soft bristles and inner light blue soft bristles.

USE

- 1) Before using the toothbrush, mount the head on the upper part of the handle, making sure not to exert excessive pressure.
- 2) Before brushing your teeth, always rinse the bristles with running water and apply a small quantity of toothpaste on the brush until you have covered approximately half of the bristles' coloured part.
- 3) The toothbrush features three cleaning modes (Clean - White - Sensitive) that make it possible to regulate it based on your needs.

WHITE

Intensive action, helps to remove surface stains and prevent their formation.

48.000 Vibrations

CLEAN

Targeted and effective action for the daily hygiene of the most demanding patients.

41.000 Vibrations

SENSITIVE

Gentle action for efficient oral hygiene, even for those who experience sensitivity issues or suffer from tender gums.

31.000 Vibrations

4) The head should always be angled at 45 degrees and used with vertical movements to effectively clean all points of the dental arch. It suffices to apply very light pressure on the teeth; if the pressure is excessive, the bristles are prevented from oscillating, thus compromising the cleaning action of the toothbrush.

The correct sequence for cleaning the entire dental arch is: upper and lower outer surface, upper and lower inner surface and, finally, the biting surface.

After use, rinse with a Curasept mouthwash suitable for your needs and do not forget to rinse the toothbrush's head thoroughly, removing it from the handle.

WARNINGS FOR CLEANING THE TOOTHBRUSH

- 1) Remove the head and wash it with running water after every use.
- 2) Use a damp cloth to thoroughly clean the handle, especially where the head is inserted and at the base, in the area that connects to the battery charger.
- 3) Thoroughly dry all parts.
- 4) Before cleaning the battery charger, make sure it is disconnected from the power outlet.
- 5) Do not clean the battery charger with running water or by immersion.
- 6) Make sure that the battery charger is not wet before connecting it to the power outlet.

PROPERTIES

Saves the last mode used: The toothbrush will operate in the last cleaning mode selected.

30-second interval: The toothbrush briefly stops vibrating every 30 seconds. This makes it possible to evaluate the exact time spent brushing each surface of the dental arch.

Up to 180 minutes of autonomy: Once fully charged, the toothbrush can be used for approximately 1 month in Sensitive mode. In the other modes, autonomy is lower.

2 Minute timer: The toothbrush automatically switches itself off after 2 minutes. This helps us clean our teeth for the time recommended by the dentist.

Bristle indicator: When the bristles start losing their colour or original shape, it is recommended to replace the head. One head used correctly lasts for approximately 3 months.

Single button: The toothbrush is managed with a single button. Once switched on, pressing the button allows changing cleaning modes.

RECHARGE

- 1) Before using the toothbrush for the first time, make sure that it is fully charged.
- 2) Connect the battery charger, then insert the toothbrush's head.
- 3) The battery indicator flashes while recharging.
- 4) The indicator continues flashing until the recharge is complete.
- 5) Once charged, the device will signal that the battery is charged with one flash.
- 6) One complete battery cycle allows continuous use from 60 to 180 minutes.
- 7) When the toothbrush is charged, the battery charger may be disconnected from the power outlet (the device is in any case protected against overload).
- 8) When the battery indicator flashes rapidly, it indicates that the battery is nearly empty.

SAFETY

To prevent damage, malfunctions or injuries, it is important that you follow the basic safety rules relating to the use of electronic devices, including the following:

- 1) Before using the toothbrush for the first time, make sure you have read the instructions manual.
- 2) Before connecting the device, make sure that the voltage of the power outlet is the same as that indicated on the device.
- 3) DO NOT use the device if it is damaged in any way, if it makes unusual noises after being switched on, or if it becomes excessively warm. Do not attempt to repair the device on your own.
- 4) If the battery charger or the heads are damaged, replace them with the original spare parts to prevent any potential hazard.
- 5) When disconnecting the battery charger from the power outlet, always grab and pull the plug and not the cable.

6) Do not clean the battery charger with running water or by immersion. Make sure that the battery charger is not wet before connecting it to the power outlet.

DISPOSAL

The symbol of the crossed-through dustbin that appears on the device or its packaging indicates that, at the end of its useful life, the product must be disposed of separately from other waste so that it can be processed and recycled appropriately.

In order to dispose of the device at the end of its useful life, therefore, you must deliver it, free of charge, to the appropriate municipal centres for the separate collection of electrical and electronic waste.

Please note that the device contains rechargeable batteries that cannot be removed by the end-user: Lithium/14500/600mAh. These batteries must only be removed by qualified professionals and must not be disposed of as urban waste, as they are collected separately to prevent harm to the environment and human health, in compliance with Directive 2006/66/EC and the national Laws transposing it.

The appropriate separate collection for subsequent environmentally compatible recycling, treatment, and disposal of the discarded devices helps prevent any possible negative effects on the environment and human health and facilitates the reuse and/or recycling of the materials from which the devices are made. Illegal disposal of the product by the user entails the application of sanctions in accordance with the applicable regulatory framework on the matter.

SUPPORT

For support or information, please write to the address: info@curaseptspa.it

WARRANTY LIMITATIONS

The toothbrush has a 2-year Warranty.

The warranty does not include:

- Heads.
- Damage caused by the use of unauthorised spare parts.
- Damage due to improper use, misuse, negligence, unauthorised modifications or repairs.
- Any regular wear and tear, including chipping, scratches, abrasions, discolouration or change in colour.

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Ltd. Dongguan City, China.
Ponzini S.p.A. Lazzate (MI), Italy
Curasept S.p.A.