

# USER MANUAL


Cosmetic face device

**BeautyRelax**


**PEEL&LIFT PREMIUM**



**BR-1540**

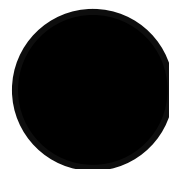

**BR-1530**

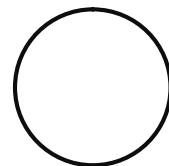



## CONTENT:

INTENDED USE	4
NOTICE	5
PACKAGE CONTENTS	5
CONTRAINDICATION	5
MAIN BENEFITS OF THE DEVICE	6
BENEFICIAL EFFECTS ON YOUR SKIN	6
TECHNICAL PARAMETERS	7
ULTRASONIC WAVES	8
GALVANIC CURRENT	9
PHOTO BLUE LIGHT	11
VIBRATING MICROMASSAGE	11
WATERPROOF PROTECTION IPX5	12
DEVICE DESCRIPTION	13
CONTROL PANEL DESCRIPTION	15
HOW TO USE	16
IMPORTANT TIPS FOR DEVICE USING	18
IMPORTANT ATTENTION	19
METHOD OF TREATMENT OF THE FACE AND NECK	20
IMPORTANT INSTRUCTIONS FOR DEVICE USING	21
MOST FREQUENT QUESTIONS AND ANSWERS	22
BATTERY CHARGING	25
CLEANING AND MAINTENANCE	27
DISPOSAL OF USED ELECTRONIC EQUIPMENT	27
CERTIFICATION	27
IMPORTER	27
DECLARATION OF CONFORMITY	29

Dear customer, thank you for your trust in the purchase of our product. Before using for the first time, please read carefully all the detailed information in the operating instructions. Pay special attention to safety warnings. Use the device exactly as instructed. All these steps will contribute to the maximum use of all functions of the device. In this guide, we will try to bring you as close as possible to the technologies used and how to make the best use of the device.

**If you encounter any ambiguities while reading the instructions, do not hesitate to contact us at the email address: [mh@santadria.cz](mailto:mh@santadria.cz) Keep the instructions for use carefully for future reference.**

## **INTENDED USE**

The BeautyRelax Peel & lift Premium cosmetic device is designed for the most demanding users. In the category of so-called "ultrasonic spatulas", it represents the absolute top. Ultrasound, vibration, galvanic current, photon light, wireless charging, water protection in IPX5 standard, luxury design with charging stand. The unique combination of functions will create a full-fledged aesthetic studio right at your home. The BeautyRelax Peel & lift Premium cosmetic device will ensure the treatment of your skin at a truly professional level!





## **NOTICE**

- Prior to initial commissioning, read these operating instructions, view all illustrations and keep these instructions for future reference.
- Do not use the appliance if your skin is affected by a serious skin condition. Then consult your doctor for the first use.
- Never operate the appliance if it has fallen, is damaged and does not work properly.
- The device is not intended for use by children or persons whose physical or mental condition does not allow the device to be used safely.
- Do not place or use the device on hot heat sources (radiator, stove, fire, etc.)
- Do not attempt to disassemble or repair the device yourself.
- Distributor shall not be liable for damages caused by improper use. He is not responsible for the consequences caused by non-compliance with the safety instructions.
- The distributor is not liable for any damage to health or property during other uses.

## **PACKAGE CONTENTS**

- Cosmetic device PEEL & LIFT PREMIUM.
- Charging stand.
- Power cable.
- Operating Instructions.

## **CONTRAINDICATION**

Under no circumstances may the device be used by persons who are affected by any of the following conditions:

- Cancer, epilepsy, viral disease
- In case of high fever
- Severe kidney or liver disease
- Heart disease, pacemaker, hypertension
- Metal implants in the body
- Vascular disorders or atherosclerosis
- Pregnancy
- Menstruation
- Acnitic keratosis, basal cell carcinoma, melanoma, spinalioma, atopic eczema, contact eczema, solar dermatitis, psoriasis
- People taking anticoagulants
- Persons who require constant supervision may only use the device after consulting a doctor

## **MAIN BENEFITS OF THE DEVICE**

- 24khz ultrasonic waves in skin cleansing and lifting modes.
- Galvanic current (positive and negative ionization).
- Photon blue light with a wavelength of 465nm.
- Vibrating micromassage.
- Stainless steel applicator for extremely long life.
- Wireless use thanks to built-in rechargeable battery.
- Wireless charging via charging station.
- Compact, small and light.
- Luxurious design.
- Smooth surface and shape thanks to smooth surface and shape.
- Water protection according to IPX5 standard (does not apply to charging station).
- High utility value.
- The highest level of certification and safety of use.

## **BENEFICIAL EFFECTS ON YOUR SKIN**

- Removal of dead skin cells and unwanted residues.
- Perfect skin peeling down to the depth of skin pores.
- Functional prevention against the formation of blackheads and clogging of skin pores.
- Effective incorporation of cosmetics, creams and serums into the skin.
- Lifting effect for fine facial wrinkles and strengthening skin elasticity.
- Non-invasive painless treatment without any discomfort.
- Local support of blood circulation to stimulate the lymphatic system of the face.



## **TECHNICAL PARAMETERS**

- Ultrasound frequency: 24khz.
- Blue light wavelength: 465nm.
- Timer: automatic switch-off after 10 minutes.
- Power: 2W.
- Operating voltage: 5V.
- Charging voltage: 15 to 500mA.
- Battery capacity: 600mAh.
- Charging time: 2 hours.
- Usage time per charge: 2 to 3 hours.
- Moisture protection standard: IPX5.
- Charging: wireless via charging station.
- The charging station is connected to the USB cable.
- Dimensions: 17.5x5.2x2cm.
- Weight: 99g.




## ULTRASONIC WAVES

The BeautyRelax Peel & lift Premium cosmetic device generates ultrasonic waves via the application head. The frequency of the ultrasound of 24khz is optimized to achieve the desired effect without injuring the skin. All dirt, makeup remnants and dead skin cells stick to the tip of the applicator. Ultrasound penetrates even the smallest folds. It adapts perfectly to the shape of the face. The other side of the applicator uses ultrasound in combination with galvanic current to help effectively incorporate cosmetics with a high degree of efficiency.



## GALVANIC CURRENT


To achieve even better results on the way to clean and radiant skin, the device is equipped with a galvanic current function. The applicator generates a galvanic current with a positive charge. Negatively charged dirt particles are attracted by the applicator like a magnet. Ultrasonic waves, in combination with a positive charge, create a highly effective tool for routine skin treatment. Regular instrumental cleansing significantly contributes to slowing down the skin aging process.

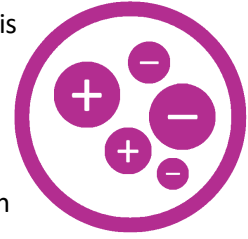

In the negatively charged application head mode, you can thoroughly incorporate the cosmetics you normally use into your skin. The active ingredients do not dry out uselessly on the surface of the skin, but penetrate effectively through the skin pores into the dermis layer. You will get the most out of commonly available cosmetics!


Differently charged ions have specific effects on the skin. The application head is used in the positive ion mode for the skin cleansing phase (**CLEAN**). In the second phase, negative ionization helps to penetrate nutrients from cosmetics many times better (**NUTRITION**).

### **Positive ionization CLEAN mode**

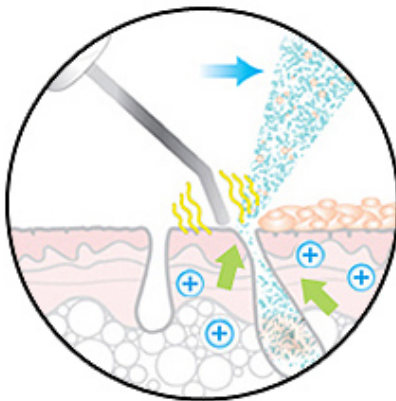



Dirt is charged with a negative charge. On the contrary, the application head creates positive ions, which attract impurities like a magnet and pull them to the surface, where it can be easily removed. Dirt is trapped not only on the surface of the skin, but deep in the pores.

### **Negative ionization NUTRITION mode**

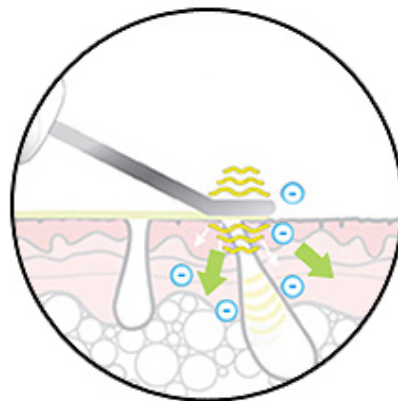

The application head generates a large amount of negative ions, which are bound to positive ions in skin cells. This creates an optimal way for the nutrients from the cream to get as deep as possible into the skin. During normal application of the cream, most of the upper layer of the skin (epidermis) absorbs. The residue dries uselessly without any further use.



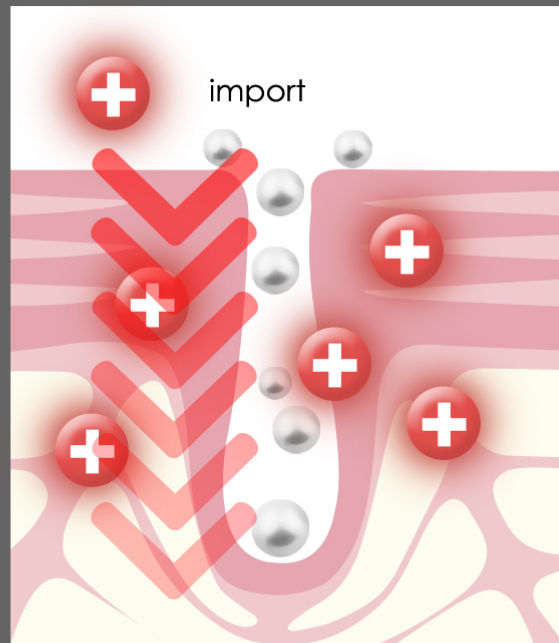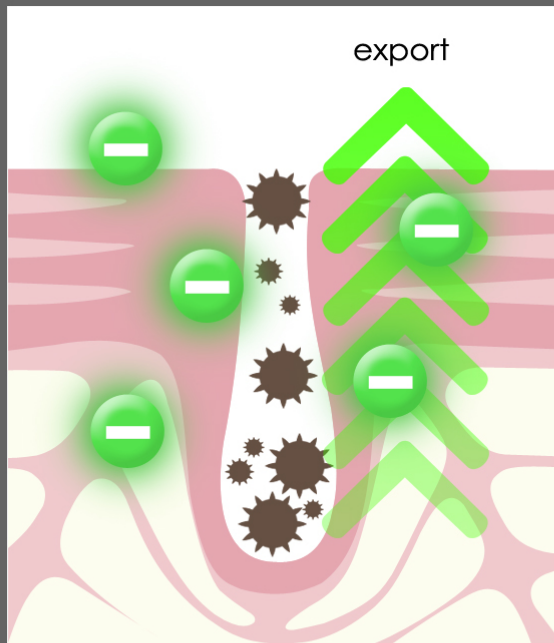

cleansing face  
cream



moisturizing skin  
serum






## **PHOTO BLUE LIGHT**

Blue photon light with a wavelength of 465nm has proven positive effects in preventing the formation of acne. Light with this specific wavelength significantly contributes to the activation of porphyrins. These simply contribute to the formation of atomic oxygen, which specifically destroys acne bacteria. It acts directly in the local area without negative effects on the surrounding tissue.



## **VIBRATING MICROMASSAGE**

The main effect is to improve blood circulation by stimulating the lymphatic system of the face. The result is smoother, fresher skin. Better tone (skin is more uniform in color). At the same time, mechanical stimulation of collagen and elastic fibers occurs.

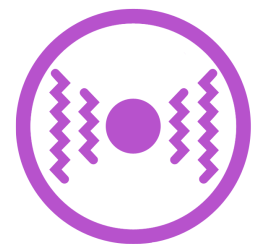



## **WATERPROOF PROTECTION IPX5**

The cosmetic device is waterproof according to the IPX5 standard. However, always make sure that it is carefully dried after use. This will significantly extend the life of the device. Waterproof does not apply to the charging stand.

# IPX5

**Protective standard**





## DEVICE DESCRIPTION



1. Application head  
2. Blue light  
3. Information panel  
**ULTRA  
CLEAN  
NUTRITION**  
BATTERY STATUS

4. **ON/OFF** button  
5. Guide side contact  
(on both sides of the device)  
6. Battery charging contacts

## **DEVICE DESCRIPTION**




- 7. Cosmetic device
- 8. Charging stand
- 9. Input for charging USB cable



- 10. Cosmetic device
- 11. Charging stand
- 12. Guide side contact (on both sides of the device)

## CONTROL PANEL DESCRIPTION

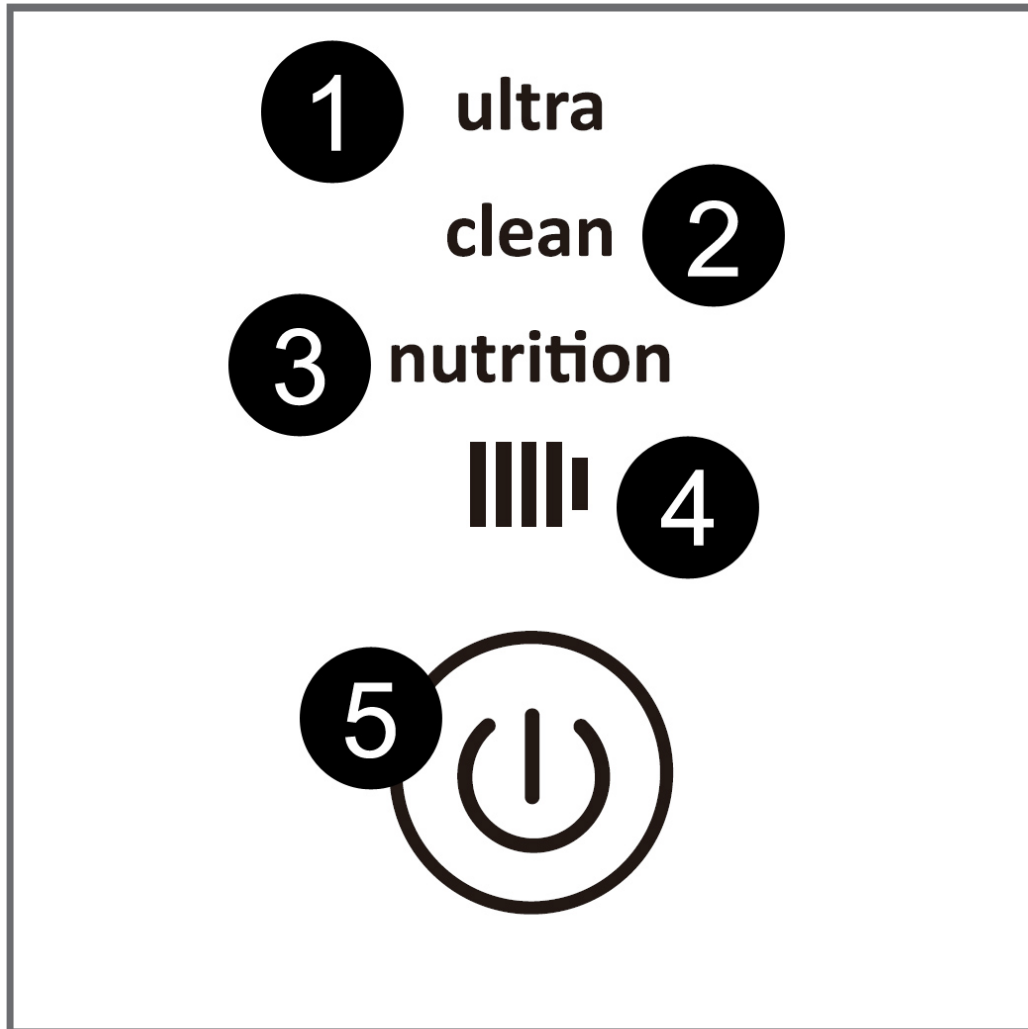

- 1. ULTRA** - ultrasonic waves work in two modes (continuous and intermittent)
- 2. CLEAN** - ultrasonic waves supplemented by positive ionization
- 3. NUTRITION** - ultrazvukové vlny doplněné o zápornou ionizaci
- 4. BATERIE** - battery status indication
- 5. ON/OFF** - on/off button, setting of individual treatment modes


## HOW TO USE

1. If you decide to treat your face with a device, remove any remnants of makeup from your body in the usual way. Rinse face with lukewarm water and wipe dry.
2. Now the skin is ideally prepared for the use of the device. Make sure that the device battery is sufficiently charged.



3. Remove the application head cover and all packing material. To switch on the device, press the **ON/OFF** button. A short beep sounds.
4. The device is now in standby mode. All three indicators are switched on (**ULTRA, CLEAN, NUTRITION**).
5. Press the **ON/OFF** button again to activate the **ULTRA** mode. The application head generates **ultrasonic waves** in continuous mode and the **blue light** turns on. The **ULTRA** indicator light is lit continuously.
6. Press the **ON/OFF** button again to switch the **ULTRA** mode to pulse mode. The application head generates ultrasonic waves intermittently. The blue light is still on. The **ULTRA** indicator light flashes.


7. Press the **ON/OFF** button again to activate the **CLEAN** mode (**ultrasonic waves** and **positive ionization**). The blue light is still on. The **CLEAN** indicator lights up continuously.
8. Press the **ON/OFF** button again to activate the **NUTRITION** mode (**continuous ultrasonic waves** and **negative ionization**). The **blue light** is still on. The **NUTRITION** indicator lights up continuously.
9. Press the **ON/OFF** button again to return to standby mode. All three indicators are lit. The ultrasonic waves turn off and the blue light turns off. You can now start setting up the functions again from point 5 above.
10. The device has an automatic timer of 10 minutes. When the time expires, an audible signal sounds and the device switches off automatically.
11. You can switch off the device at any time by long pressing the **ON/OFF** button. A long beep sounds (Beeeeeeeeeeep).

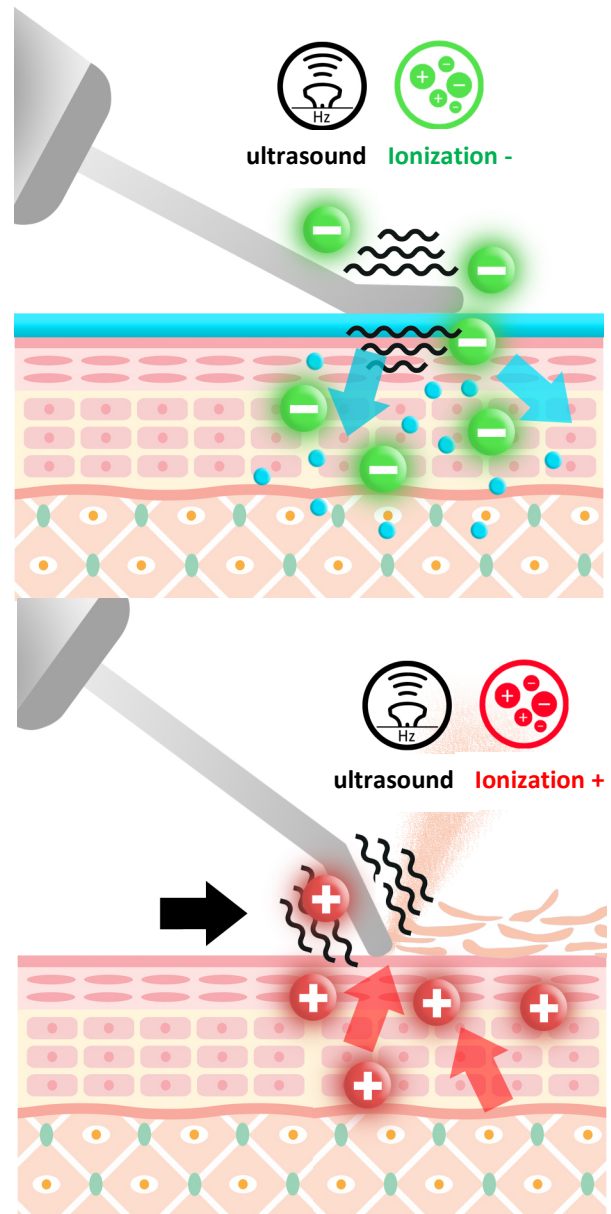

## IMPORTANT TIPS FOR DEVICE USING

1. Thanks to its special shaping, the application head is intended for skin treatment in two ways.

Two ways to treat the skin surface:

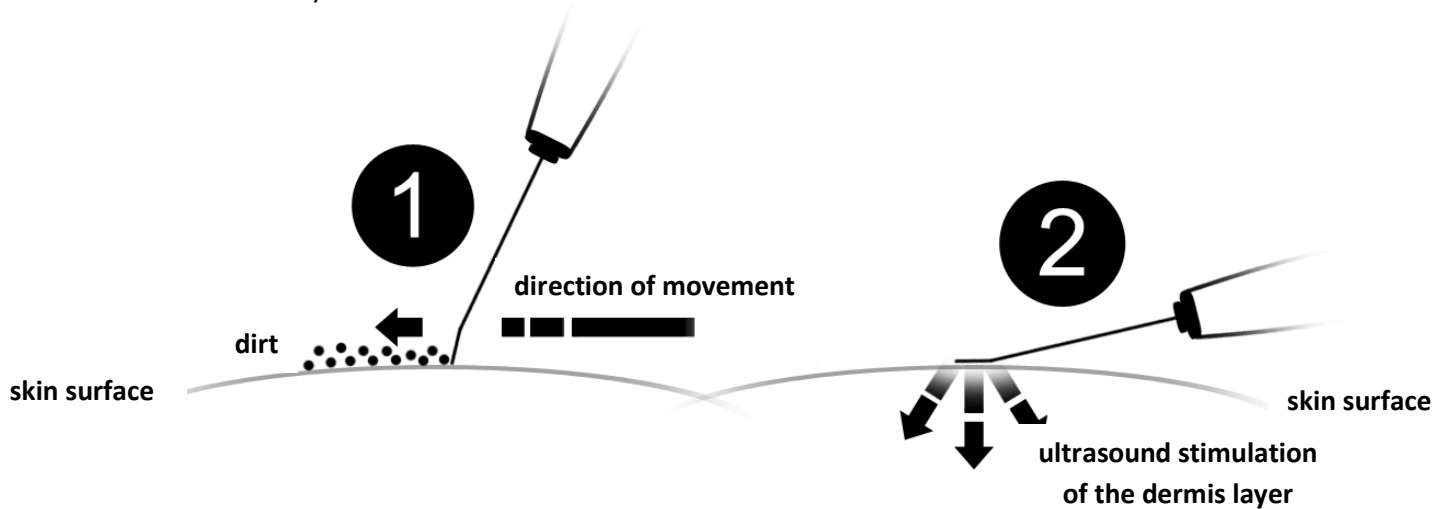

**Figure (1)** The tip of the applicator is ideal especially for intensive skin cleansing. The movement of the applicator in combination with the action of ultrasonic waves rids the skin of dead cells and residues - **CLEAN** mode.

**Figure (2)** The larger end face of the applicator is ideal for stimulating the skin with ultrasound waves in order to support local blood circulation and stimulate the formation of collagen and elastic fibers. **ULTRA** and **NUTRITION** modes are suitable for this type of treatment. However, it is always necessary to use a larger area of the applicator - see. **Figure 2**. Local blood circulation is improved and in addition you can better incorporate cosmetics. The effect of continuous ultrasonic waves is optimal.

2. If your skin is less hydrated (drier type), use sonogel, which you can buy at the pharmacy and all medical supply stores. Sonogel creates a thin layer between the application head and the skin surface. Conductivity is increased and ultrasonic waves have a better ability to penetrate through the upper layer of the skin.
3. In **NUTRITION** mode, it is ideal to apply your favorite moisturizer or serum to the skin. By the action of ultrasonic waves and negative ionization, nutrients from creams and sera will have a better ability to penetrate the layer of the dermis, where the collagen and elastic fibers are located.
4. To activate positive or negative ionization, it is essential that your fingers are in contact with the side guide contacts. Otherwise the ionization will not work properly!
5. Replace the long-life applicator guard after each use.

## **IMPORTANT ATTENTION**

To activate the galvanic current function (ionization), it is necessary for the index finger and thumb to be in contact with the side guide strips. Otherwise the function will not work!





## METHOD OF TREATMENT OF THE FACE AND NECK



**FOREHEAD:** move away from the center facing the edge, gradually cover the entire area




**NOSE:** move from the forehead around the nose to the face, from top to bottom



**FACE:** move your face in the direction bottom to top



**CHIN:** Move from the tip of the chin along the contours of the face, from bottom to top

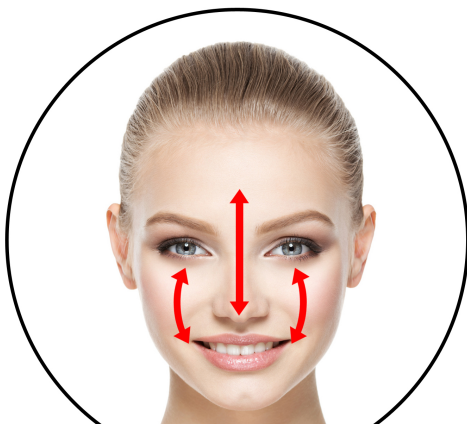

**T ZONE:** move from forehead to the tip of the nose



**NECK:** Move from the chin towards the ear and then down to the neck



## **IMPORTANT INSTRUCTIONS FOR DEVICE USING**

- Ultrasonic device PEEL & LIFT PREMIUM is a new generation from the category of cosmetic devices that are designed for non-invasive skin treatment. This means that none of the parts of the application head physically penetrate under the skin.
- Non-invasive use is a basic prerequisite for using the device at home by a normal user. Control of the device does not require the supervision of an erudite trained worker in the field of healthcare or cosmetic aesthetics.
- The PEEL & LIFT PREMIUM ultrasonic device is certified and meets all standards set for a device of a similar type.
- You can use a clear conductive sonogel during the treatment. It can be purchased at any pharmacy or health care store. It costs a few tens of crowns. Definitely enough is the cheapest. The more expensive the better. The main purpose of applying a thin layer of gel on the skin is to increase conductivity. Ultrasonic waves gain a better ability to penetrate the layer of the dermis, where the structure of collagen and elastic fibers is.
- Always clean the application head thoroughly after each treatment. It is ideal to use clear alcohol for disinfection. In the last phase, dry the applicator thoroughly. Moisture can cause oxidation of the applicator contacts. This will significantly increase the life of the application head.
- In galvanic current mode, you may feel a slight tingling sensation on your skin. Take extra care especially when treating areas of the face where the skin is very thin. In general, the thicker the layer of subcutaneous fat, the less recognizable tingling and vice versa.
- In skin cleansing mode, where the tip of the application head is in contact with the skin, do not be afraid to press lightly on the skin. The tip of the applicator has a special curvature to adhere best to the skin surface. During use, you can continuously clean the dirt from the skin with a dry cloth.
- Pay special attention to disinfecting the applicator, especially in cases where the device is used by more than one person. This will significantly prevent the mutual transmission of bacteria.

## **MOST FREQUENT QUESTIONS AND ANSWERS**

**Why is it recommended to remove makeup with a cosmetic pad before using the device? After all, the device as such is intended for cleaning the skin?**

Before using the device, it is advisable to remove makeup in the usual way, as you are used to. The main purpose is to remove the so-called "rough" impurities from the skin. With the help of the cosmetic device applicator, you can focus on really deep and thorough cleansing of the skin. For impurities that are not visible to the naked eye.



**Can I use the device on "dry" skin without using any product?**

Ultrasound waves belong to the category of non-invasive cosmetic treatment technologies. Human skin is too dry for this type of treatment. It is necessary to hydrate the skin sufficiently so that the ultrasound waves can connect to it better. When used on dry skin, the effect of ultrasound is minimal. As a last resort, you can use baby or distilled water to hydrate your skin. However, it is not very practical in practice. The water will not stay on the place you want to treat with the device. In addition, water dries quickly during use of the device. Use the device like a real professional! We recommend purchasing sonogel for the device at a pharmacy or health supply store. It is inexpensive and offers excellent properties for instrumental treatment. It is enough to apply a thin layer. The gel stays exactly where you need it. It does not run and does not dry quickly from the skin. The gooey consistency will allow smooth movement of the applicator over the skin. In this way, you can easily carry out the entire process of cleaning the skin, when you use the applicator to wipe off impurities from the skin together with the gel. After using the device, wipe off any remaining gel from the skin.

**Sonogel, cosmetic cream, serum.... So how to use or combine cosmetic products correctly? What is the best?**

It's easy. For cleaning and local lifting of the skin by ultrasound stimulation, the absolute minimum is to use at least sonogel. The basic cheap sonogel does not contain any active substances, such as aloe vera, hyaluronic acid or other folders. In the skin cleaning mode, it creates a suitable bond between the applicator and the skin. In skin massaging mode, the curved surface of the applicator can effectively stimulate fine capillary wrinkles with ultrasound. However, there are more expensive variants of sonogel on the market that contain additional ingredients for beneficial skin care. It depends on you how much you want to invest in these preparations. Whatever you decide on, be sure to monitor the proportion of these active substances. If such a product contains only a trace amount, it has de facto no effect on the skin. You can use your favorite cosmetic cream. In this case, set the NUTRITION program on the device. Ultrasonic waves hand in hand with a negative galvanic current will help the cream work better into the skin.

From the point of view of the use of cosmetic products, we would like to mention an important aspect. A cosmetic device with non-invasive skin treatment technology gets along best with sonogel or with sonogel containing other ingredients. If you use any other cream/serum, for users with very sensitive skin, the combination of ultrasound + galvanic current + cream/serum may cause an allergic reaction in the form of local skin redness. This can be caused by other substances that these preparations contain (dyes, aromas, etc.). Before facial treatment, we recommend applying the product to a small piece of skin, e.g. on your hand, using the device. If after a few minutes there is no significant redness on the skin, you can continue with the facial treatment.

**ULTRA, CLEAN, NUTRITION. What is the real meaning of these features and which one is best for me?**

You definitely do not need to use all the modes exactly one after the other when using the device. Let's analyze the features of each program in detail:

**1. ULTRA** - ultrasonic waves work in two modes (continuous and intermittent)

Ultrasonic waves are active in this mode. They work continuously or intermittently. Intermittent action is gentler. Choose if you have more time to use the device. You can use the ULTRA mode to clean the skin (with the tip of the applicator) or to stimulate specific parts of the face to achieve a lifting effect. For the ULTRA mode, the minimum is at least a thin layer of sonogel.

**2. CLEAN** - ultrasonic waves supplemented with positive ionization

CLEAN mode is very similar to ULTRA mode. In addition, the applicator switches to the positive charge of the galvanic current. It is proven that with galvanic current, cleaning is even more effective. Nevertheless, ultrasound has the biggest share in thorough peeling. Combine with a sonogel or, for example, a cosmetic cream for cleaning the skin.

**3. NUTRITION** - ultrasound waves supplemented with negative ionization

This mode is maximally adapted for the incorporation of active substances in a sonogel or cosmetic cream, provided that it contains more than a trace amount of these substances. Use the curved part of the applicator.

**How often is it recommended to use the device?**

The frequency of use should be based on the quality of the environment in which you move the most. If you spend a lot of time in a dusty environment, shorten the intervals between peeling your skin. Your skin type plays an equally important aspect. If you have skin that is more oily and tends to clog skin pores more often, shorten the intervals of use again. The skin should be smooth and soft to the touch. If it loses this feature, it is the right time to use the device.

**For which parties is the device primarily intended?**

Face, neck and décolletage.

**How to appropriately combine skin cleaning (peeling), skin lifting with the incorporation of creams into the skin?**

As mentioned above, there is no need to combine all the programs at once. The frequency of skin peeling does not have to correspond with the frequency of use of your cosmetic product. Each breed has individual needs. Try all possible combinations from the beginning. You will find out how your skin reacts to the individual possibilities of using the device. You can create your own skin care plan.

**So what are some examples of a possible skin treatment plan?**Example #1

You can only clean the skin with the device. Remove makeup residue as usual. Apply a thin layer of sonogel to the skin, choose the ULTRA or CLEAN program and clean the skin. You are done in 10 to 15 minutes.

Example #2

Remove makeup residues in the usual way (e.g. using a cosmetic pad, cleansing cream or micellar water). Focus on those parts of the face that are most quickly subject to the effects of age - the forehead, wrinkles around the mouth, hairline wrinkles around the eyes. Do not apply a layer of sonogel or other cosmetic product to these areas. If using sonogel, choose the ULTRA program - continuous or intermittent pulses. If using a high-quality cosmetic cream/serum, select the NUTRITION mode. 20 minutes and you're done.

Example #3



Remove makeup residue as usual. Apply a thin layer of sonogel to the skin, choose the ULTRA or CLEAN program and clean the skin. Remove the remains of the sonogel with a cloth. Apply a nourishing cream to the skin and work with the NUTRITION program. Split the treatment in half. 10 minutes of peeling and 10 minutes for applying the cosmetic product.

*This user manual serves only as a guide on how to work with the device. However, it also creates space for your creativity to create your own plan for using the device, which will best suit your skin. Don't be afraid to experiment a little. Engage your imagination. The technology of the device is safe for your skin and you cannot harm yourself in any way. If you don't know what to do, don't hesitate to contact us. We're here for you. The BeautyRelax team*

## **BATTERY CHARGING**

The device has a built-in rechargeable battery. It is charged using the charging stand included in the package.

When the battery is low, a flashing red battery icon appears on the display. The device switches off automatically.

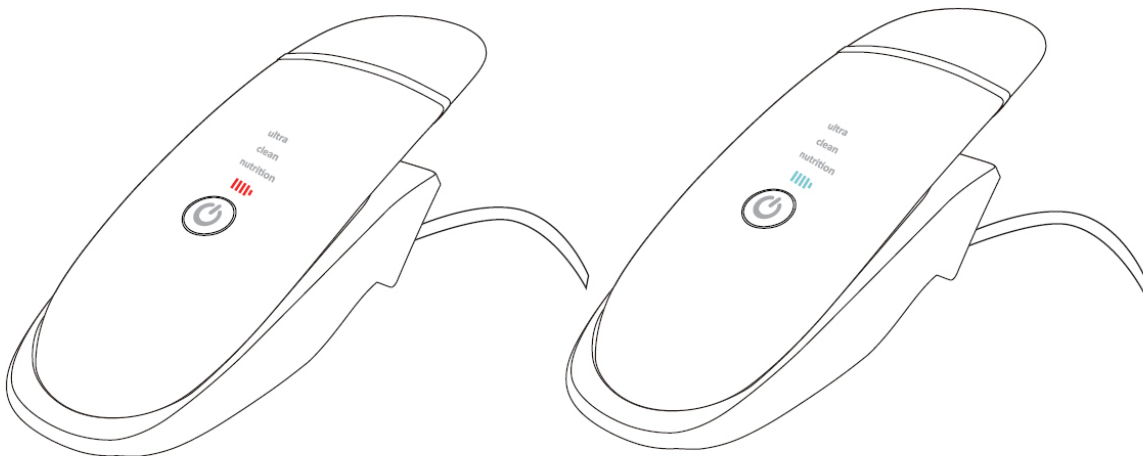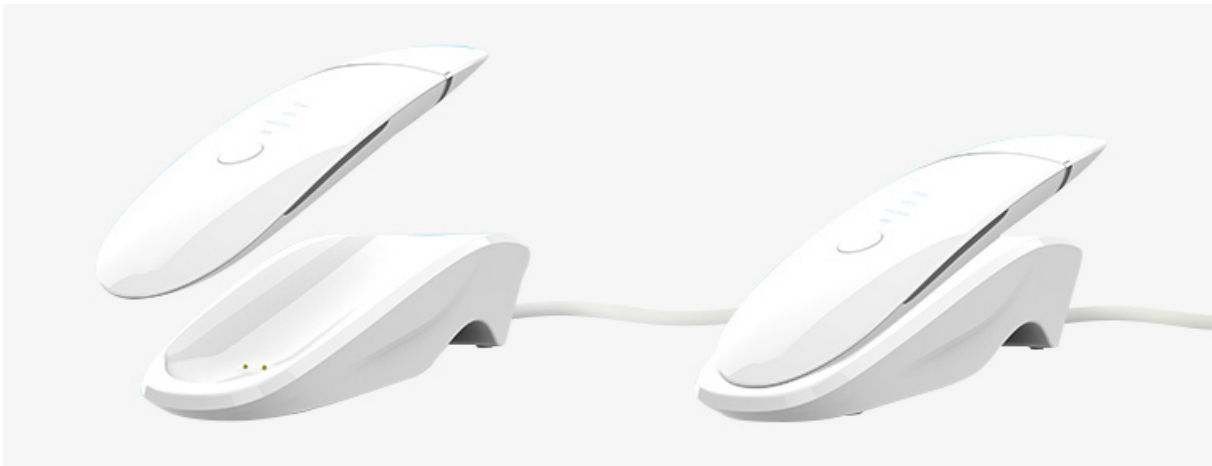


**Charging procedure:**

1. Use the included power cord to charge.
2. Connect one end of the cable to the charging stand. The second to the USB port of a computer or laptop
3. The red battery indicator flashes during charging.
4. When the battery is fully charged, the battery indicator on the device stops blinking and turns blue.



**Do not use a 220V charging adapter to charge the device. This could seriously damage the device. Connect the other end of the USB cable to a USB port on a computer or laptop.**

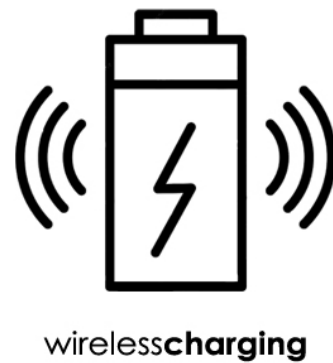

Registered



## **CLEANING AND MAINTENANCE**

Always make sure that the applicator is cleaned after the treatment. Do not use any aggressive cleaning agents. We recommend the use of special disinfectant sprays, which are commonly used, for example, in healthcare or beauty salons for cleaning tools. Classic pure alcohol will also suffice. Moisture must NEVER penetrate directly into the unit. A cotton swab or dry cloth is ideal for cleaning. Before using the device again, no gel, cosmetic cream or serum remains on the applicator. This would significantly shorten the life of the applicator. Under no circumstances use any sharp objects that could damage the surface of the applicator during maintenance.

## **DISPOSAL OF USED ELECTRONIC EQUIPMENT**



This symbol on the products or in the accompanying documents means that the used electrical and electronic products must not be added to ordinary municipal waste. For proper disposal, recovery and recycling, transfer these products to designated collection points. Alternatively, in some countries of the European Union or other European countries, you can return your products to your local retailer when buying an equivalent new product. Proper disposal of this product will help to conserve valuable natural resources and help prevent potential negative environmental and human health impacts, which could be the consequences of improper disposal of waste. Please contact your local authority or nearest collection point for further details. This symbol is valid in the European Union. To dispose of this product, ask for the proper disposal information from your local authorities or from your dealer.

## **CERTIFICATION**



This product complies with the requirements of the EU directives on electromagnetic compatibility and electrical safety.

## **IMPORTER**

**SANTADRIA s.r.o.**, Mojmírovo nám. 75/2 Brno, 612 00  
IČO: 29195454 | DIČ: CZ29195454

If you have any questions regarding the product, please do not hesitate to contact us:  
Email: [mh@santadria.cz](mailto:mh@santadria.cz) Tel: 737 196 160

**Typographical errors, properties and changes in technical parameters reserved.**

# ES DECLARATION OF CONFORMITY

Signed below, representing the company

**SANTADRIA s.r.o.**

Mojmírovo nám. 75/2, Brno 612 00, IČ: 29195454

hereby declares, as an importer in the EU, that the product with the trade name below is listed below

## **Cosmetic device BeautyRelax PEEL&LIFT PREMIUM**

placed on the market under the model name: **BR-1530/BR-1540**

EAN: **8594166116955; 8594166116962**

belonging to the product group of the brand **BeautyRelax**

**complies with the following European standards:**

**EN 55014-1: 2017**

**EN 55014-2: 2015**

**Details are given in test report number: SZAHE180319001-01**

### **Based on compliance with these guidelines, the product has been CE marked.**

The product also complies with the RoHS Directive 2002/95 / EC (RoHS). All materials are not harmful to health. They do not contain any heavy metals, such as cadmium and other toxic substances.

For safe use, it is necessary to strictly follow the instructions given in the operating instructions that are included in the package.





# VERIFIED ULTRASOUND POWER



beautyrelax **peel&lift** premium

