

USER MANUAL

Handheld Massager

BeautyRelax BR-1790B

Kineticforce Portable



CONTENT:

SAFETY NOTICE	4
PURPOSE OF USE	5
MAIN ADVANTAGES	5
4 CHANGEABLE HEADS	5
TECHNICAL PARAMETERS	5
PACKAGING CONTENT	5
DEVICE DESCRIPTION	6
USE OF INDIVIDUAL APPLICATORS	7
CONTROL PANEL DESCRIPTION	8
HOW TO USE	9
BATTERY CHARGING	9
PAY ATTENTION	9
POSSIBILITIES OF USE	10
DISPOSAL OF USED ELECTRICAL EQUIPMENT	11
CERTIFICATION	11
IMPORTER	11
ES DECLARATION OF CONFORMITY	12

Dear customer, thank you for your trust in purchasing our product. Please read all the detailed information in the instruction manual carefully before using it for the first time. Please pay particular attention to the safety warnings. Please use the device exactly as instructed. All of these steps will help to maximize the use of all functions of the device. In this manual we will try to give you the best possible idea of the technology used and how to use the device to best advantage.

If you encounter any uncertainties while reading the instructions, please do not hesitate to contact us at the following email address: mh@santadria.cz

Keep the instructions for use carefully for later reference.

SAFETY NOTICE

- Do not attempt to disassemble or repair the unit yourself in any way
- Store and use the device in a dry, closed environment
(it is forbidden to use the device in the bathroom)
- It is forbidden to use the device if you have had surgery on the cervical spine or any other part of the body for which the massage device is designed
- It is forbidden to use the device if you have serious chronic health problems with the cervical spine and if your neck area is very sensitive
- Do not use the device for too long; 5 to 20 minutes a day is sufficient for pleasant relaxation and muscle relaxation; use for longer than the specified time could cause headaches (in case of application to the cervical spine)
- If you are unsure, consult your GP or physiotherapist as a last resort; the device does not require any major mechanical maintenance
- Disconnect the unit from the power adapter before cleaning the controller; clean the surface of the controller only with a clean soft cloth
- Do not use aggressive cleaning agents; the surface of the controller could be damaged; do not leave the packaging material lying around; it could become a dangerous toy for children
- Handle the product with care; do not subject it to shocks, knocks or impacts; permissible ambient temperature: 5°C to +50°C
- Damage caused by neglect of these safety warnings and instructions for use will void any claim to the warranty of the appliance; we shall not be liable for consequential damage caused by breach of these conditions
- We do not accept any liability for personal injury caused by improper or incorrect handling of the product; in such cases, the product warranty is void

PURPOSE OF USE

The BeautyRelax Kineticforce Portable Handheld Massager is a revolutionary and innovative solution that will change the way you look at total body rehabilitation and regeneration. The device was designed by professional therapists. During the development process, emphasis was placed on top quality and maximum efficiency. There is currently no more effective device on the market of comparable quality. BeautyRelax Kineticforce Portable immediately relaxes tight muscles anywhere on the body. Treat your body to the beneficial effects of modern technology! Thanks to its compact size, the device can be with you wherever you go.

MAIN ADVANTAGES

- Immediate relaxation of stiff muscles
- Relieve stress and mental tension
- Fresh feeling during long lasting fatigue
- Local improvement of blood circulation
- Relief of swelling, especially in the limbs
- Regeneration of the mind after demanding intellectual activity

4 CHANGEABLE HEADS

- Large round applicator - suitable for large muscle groups
- Fork - area around the spine
- Flat head - suitable for all body parts
- Projectile head - knees, joints and hard to reach areas

TECHNICAL PARAMETERS

- Ergonomic handle is coated with non-slip material
- Extremely low noise operation
- 4 interchangeable application heads
- 4 massage speeds (1800 to 3200 vibrations per minute)
- Display shows battery level and massage speed
- Operating voltage: 24V
- Charging: built-in 2400mAh Li-ion battery
- Power: USB cable type C (included)
- Dimensions: 22x22x7cm
- Weight: 800g

PACKAGING CONTENT


- 1x BeautyRelax Kineticforce Portable Handheld Massager
- 4x Interchangeable applicators of different shapes
- 1x USB charging cable type C
- 1x User manual

DEVICE DESCRIPTION



1. Application head
2. Control panel
3. Handle
4. Set of changeable application heads

USE OF INDIVIDUAL APPLICATORS

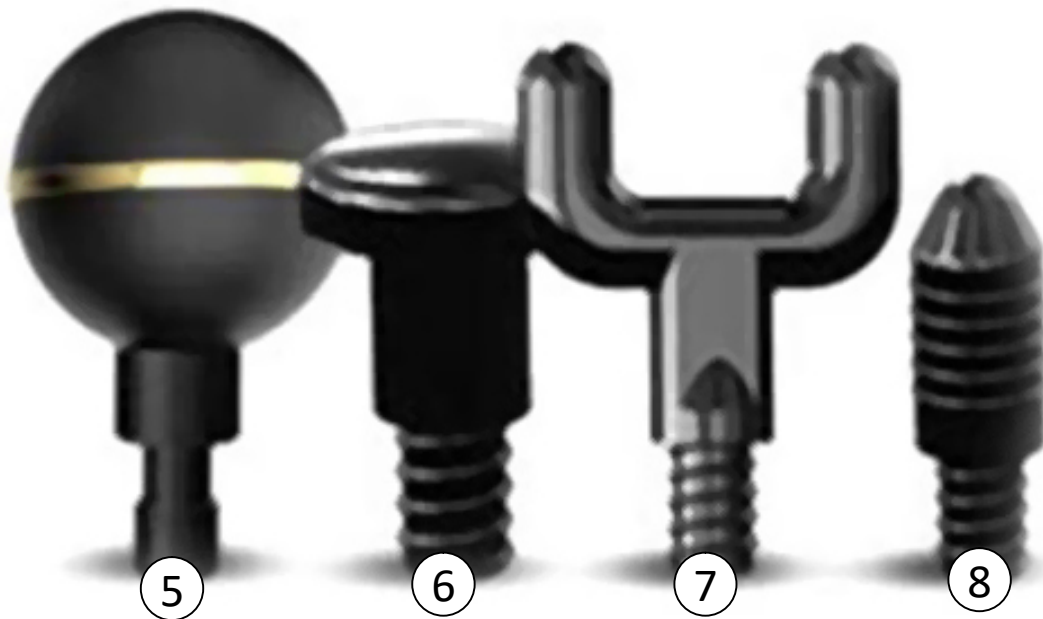

- 5. Ball** - suitable for large muscle groups
- 6. Flat head** - suitable for all body parts
- 7. Fork** - suitable for the area around the spine
- 8. Projectile** - suitable for the area around the knees, joints and hard to

CONTROL PANEL DESCRIPTION



- 1.** Power cable input
- 2.** % battery status
- 3.** Massage speed
- 4.** ON/OFF button
- 5.** Button (+) increase massage speed
- 6.** Button (-) reduce massage speed

HOW TO USE

1. Choose one of the 4 applicators depending on the purpose of use.
2. Switch the device on by pressing the **ON/OFF** button.
3. **Use the (+) and (-) buttons** to easily adjust the massage speed. Press the **(+)** button to increase the speed and the **(-)** button to decrease the massage speed. There are 20 massage levels.
4. You can keep track of the battery charge % and massage speed on the display.
5. You can switch the device off at any time using the **ON/OFF** button.


BATTERY CHARGING

1. The battery is charging with a flashing light. The charging progress is shown on the display.
2. When the battery is fully charged, the flashing stops.
3. It is not recommended to charge the battery when it is completely discharged to the lowest level.
4. The battery life per charge depends on the speed of the massage. Higher massage speeds drain the battery many times faster.
5. Do not use the device while it is charging.
6. The battery can take 3 to 6 hours to charge.

PAY ATTENTION

1. When replacing the applicator head, always make sure that it is pressed as firmly as possible into the hole provided for the applicator.
2. Do not use the device for more than 1 hour at a time. This could significantly shorten its service life. Allow at least 30 minutes between each use.
3. When the battery is charged and you decide to start using the device, always disconnect the power cord.
4. You can choose the strength of the massage yourself by the pressure you apply to each body part. Always keep in mind that the overall intensity of the massage is always determined by the force you apply to the body part and the speed you set.
5. Always start using the device at a lower massage speed initially.

POSSIBILITIES OF USE

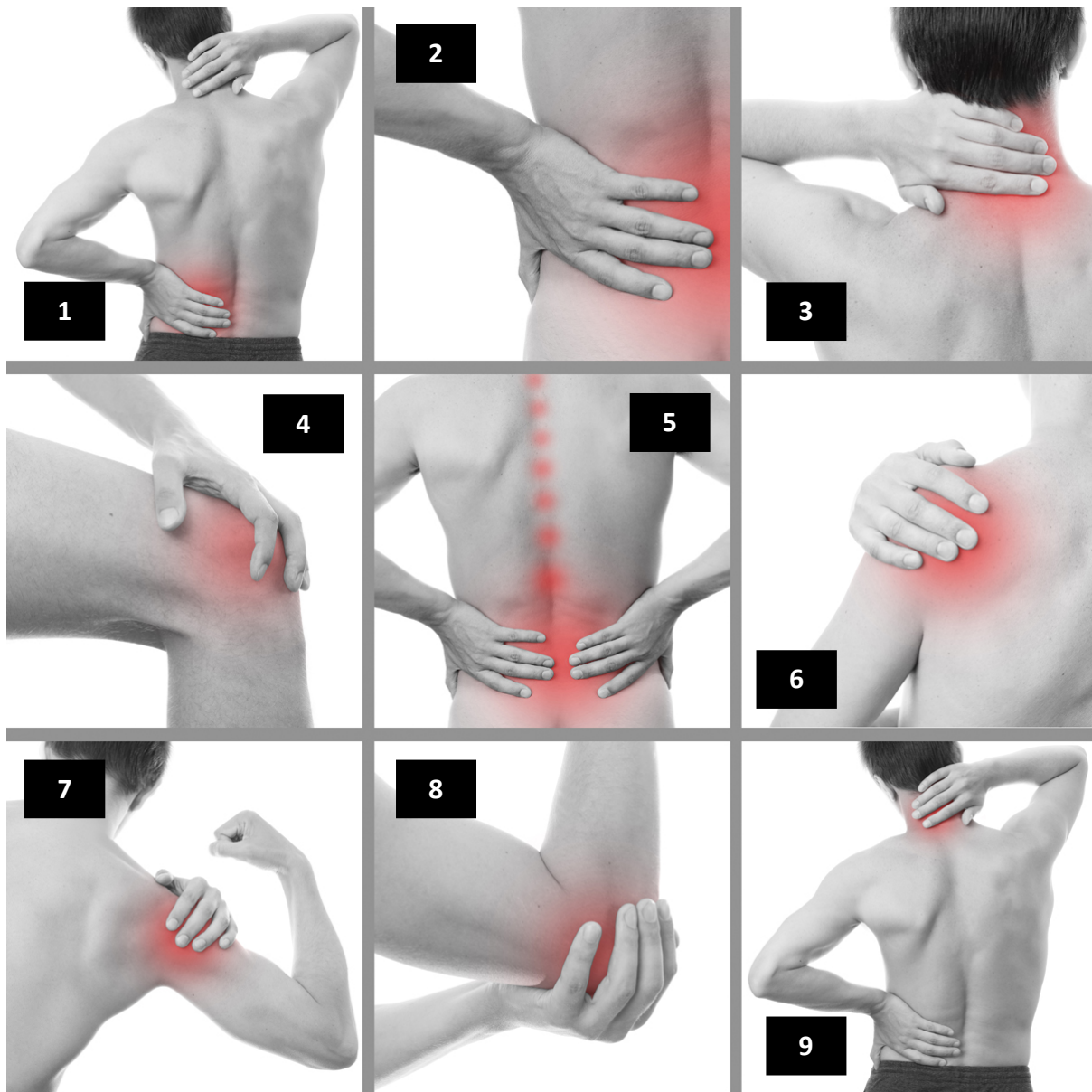

- 1. Lumbar region**
- 2. Lower back**
- 3. Cervical spine**
- 4. Knees**
- 5. The area around the spine**
- 6. Shoulder**
- 7. Shoulder**
- 8. Elbow**
- 9. Cervical spine + neck**

DISPOSAL OF USED ELECTRICAL EQUIPMENT



This symbol on products or in accompanying documents means that used electrical and electronic products must not be added to normal municipal waste. For proper disposal, refurbishment and recycling, take these products to the designated collection points. Alternatively, in some European Union or other European countries, you can return your products to your local retailer when purchasing an equivalent new product. By properly disposing of this product, you will help conserve valuable natural resources and help prevent potential negative environmental and human health impacts that could result from improper waste disposal. For further details, contact your local authority or nearest collection point. This symbol is valid in the European Union. If you wish to dispose of this product, please obtain the necessary information on the correct method of disposal from your local authorities or your dealer.

CERTIFICATION



The product is approved for use in European Union countries and bears the CE mark. The declaration of conformity and all documentation is available from the importer.

IMPORTER

SANTADRIA s.r.o. , Mojmírovo nám. 75/2 Brno, 612 00
IČO: 29195454 | DIC: CZ29195454

If you have any questions about the product, please do not hesitate to contact us:

Email: mh@santadria.cz Tel: 737196160

Printing errors, features and changes in technical parameters reserved.

ES DECLARATION OF CONFORMITY

The undersigned, representing the company

SANTADRIA s.r.o.

Mojmírovo nám. 75/2, Brno 612 00, ID: 29195454

hereby declares, as an importer into the EU, that the following product with the trade name

BeautyRelax Kineticforce Portable Massager

marketed under the model designation: **BR-1790B;BR-1790P;BR-1790R**

EAN: **8594166118539;8594166118959;8594166118966**

belonging to the **BeautyRelax** brand group of products

is in compliance with the following European standards:

EN 55014-1:2017, EN 55014-2:2015

EN 61000-3-2:2014, EN61000-3-3:2013

Based on compliance with these directives, the product has been CE marked.

The product also complies with RoHS Directive 2002/95/EC (RoHS). All materials are health-safe. They do not contain any heavy metals such as cadmium and other toxic substances.

For safe use, it is essential to strictly follow the instructions in the instruction manual included in the package.



Place of publication: Brno Day: 1.4.2020

Martin Hrbek, Managing Director