

# USER MANUAL

Cosmetic device  
**BeautyRelax BR-2080**

## EYEPEN RF





## CONTENT:

IMPORTANT SAFETY NOTICES	4
NOTICE	5
CONTRAINDICATIONS	6
USING SONOGEL	6
MAIN ADVANTAGES	7
TECHNICAL PARAMETERS	7
PACKAGING CONTENT	7
PURPOSE OF USE	8
RADIOFREQUENCY	8
VIBRACE	8
PHOTON LIGHT	8
DEVICE DESCRIPTION	9
DISPLAY DESCRIPTION	10
HOW TO USE	11
HOW TO USE DEVICE - NOTES	13
POSSIBLE AREAS OF USE	14
METHOD OF TREATMENT OF THE FACE AND NECK	15
FREQUENTLY ASKED QUESTIONS AND ANSWERS	16
BATTERY CHARGING	17
CLEANING AND MAINTENANCE	17
DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT	18
CERTIFICATION	18
IMPORTER	18

## IMPORTANT SAFETY NOTICES

- BeautyRelax EYEPEN RF is a device from the category of cosmetic devices that are designed for non-invasive skin treatment. This means that no part of the application head physically penetrates below the top layer of the skin. The device is not a medical device.
- Non-invasive use is a prerequisite for the use of the device in home conditions by the average user. Operation of the device does not require the supervision of an erudite trained medical or cosmetic aesthetics professional.
- The BeautyRelax EYEPEN RF is certified and meets all European standards for this type of device. The safety of all electronic components is ensured.
- The treatment is always associated with a clear conductive gel. It can be purchased at any pharmacy or health supply store. It costs a few tens of crowns. Definitely the cheapest is enough. The more expensive the better. The main purpose of applying a fine layer of gel to the skin is to increase conductivity. Radiofrequency waves get a better ability to penetrate the dermis layer, where the structure of collagen and elastic fibers.
- **Always proceed by holding the application head on the skin for a maximum of 1 to 2 seconds (stamping). Depending on the type of skin treatment, the technology always works when the applicator of the device is placed against the skin. Start at the lowest setting. You can gradually move up to higher levels. Skin treatments should never be unbearably painful. Try to adjust the length of contact between the applicator and the skin accordingly. Only gentle tingling to light pinching is acceptable. Each part of the face should not be treated more than twice in one use.**
- If you hold the applicator head on your skin for too long and also have the highest treatment level set, the device is powerful enough to injure your skin unpleasantly. Therefore, please approach the treatment with the utmost care and gentleness.
- The number of repetitions is very individual. It depends on the age and condition of the skin. The ideal is once or twice a week at most. The skin and especially the dermis area needs time for an autonomous regeneration process. The more frequent the use, the better the results. Neither does it apply - the greater the performance, the faster and better the result. Everyone's skin is an organ with an individual degree of sensitivity. For some, a medium level of power may already be significantly painful, while others have no problem with the highest power for longer treatments.
- Always clean the application head thoroughly after each treatment. The use of clear alcohol for disinfection is ideal. Dry the applicator thoroughly at the last stage. Moisture can cause oxidation of the applicator contacts. This will significantly increase the life of the applicator head. If you do not know what to do, do not hesitate to contact us. Phone: 737196160 or email: mh@santadria.cz
- We wish you every success in your quest for healthier and more beautiful skin.  
BeautyRelax.cz team



## **NOTICE**

- Read this manual before first use, review all illustrations and keep the manual for future reference.
- If your skin is affected by a serious skin condition, do not use the device. Consult your general practitioner for first use.
- The device is primarily intended for your personal use. Do not share it with family members or friends. If more than one person is using the device, special attention should be paid to hygiene.
- If the device is damaged, does not work properly or has fallen into the water, never put it into operation.
- If the unit is broken, do not touch the internal parts to avoid electric shock.
- If the device falls to the ground, it may be damaged. In this case, do not continue to use the device and contact the dealer's service department.
- Do not use the device near water or in places with increased humidity (bath, shower, swimming pool, etc.).
- The device is not washable. It must not be immersed in water or under a stream of running water.
- Do not leave the device unattended. Always switch off the appliance thoroughly after each use.
- If the device has been moved from a very cold to a very warm environment, wait at least 3 hours before using it.
- The appliance is not intended for use by children or persons whose physical or mental condition does not permit the safe use of the appliance.
- Do not place or use the device on hot heat sources (radiator, stove, fire, etc.)
- Do not attempt to disassemble or repair the device yourself.
- Store the device in a dry environment. Try to keep it protected from dust.
- Do not expose the device to direct sunlight.
- The distributor is not liable for damages caused by improper use contrary to the instructions in this manual. It is not liable for consequences caused by failure to observe the safety warnings.
- The distributor is not liable for any damage to health or property caused by other uses.



Dear customer, thank you for your trust in purchasing our product. Please read all the detailed information in the instruction manual carefully before using it for the first time. Please pay particular attention to the safety warnings. Please use the device exactly as instructed. All of these steps will help to maximize the use of all functions of the device. In this manual we will try to give you the best possible idea of the technology used and how to use the device to best advantage.

If you encounter any uncertainties while reading the instructions, please do not hesitate to contact us at the following email address: mh@santadria.cz

**Keep the instructions for use carefully for later reference.**

## **CONTRAINDICATIONS**

**Under no circumstances should** the device be used by persons affected by any of these conditions:

- Cancer, epilepsy, viral diseases
- In a state of high fever
- Serious kidney or liver disease
- Heart disease, pacemaker, hypertension
- Metal implants in the body
- Vascular disorders or arteriosclerosis
- Pregnancy
- Menstruation
- Acneiform keratosis, basalionomas, melanomas, spinalionomas, atopic eczema, contact eczema, solar dermatitis, psoriasis
- Persons taking anti-clotting medication
- Persons who require constant supervision may use the device only after consulting a physician
- Never move the device on the eyelids, temple area or thyroid area



## **USING SONOGEL**

Before use, a thin layer of clear conductive gel (not included) should be applied to the skin. The application gel, called SONOGEL, can be purchased at any medical supply store (can also be purchased at a pharmacy). The gel is very important for the use of the device. It forms a necessary bonding layer between the applicator and the skin. It is water-based. It cannot harm the applicator head. The use of the applicator gel is essential. This allows the radiofrequency waves to penetrate the skin more easily.



## **MAIN ADVANTAGES**

- Technology: radio frequency, red light, vibration
- Cutting-edge multipolar radio frequency technology
- Increases skin elasticity
- Very easy to achieve lifting effect
- Suppression of dark circles around the eyes
- Large LCD display for easy setting of functions
- Suitable for all skin types
- 4 different modes
- Easy operation with 2 buttons
- Built-in rechargeable battery

## **TECHNICAL PARAMETERS**

- Operating voltage: 5V
- Operating current: 1A
- Battery capacity: 500mAh
- Power: USB cable (included)
- Vibration: up to 12000 pulses per minute
- Dimensions: 15x2,8x2,2cm
- Weight: 54g

## **PACKAGING CONTENT**

- BeautyRelax BR-2080 EYEPEN RF
- USB charging cable
- Operating Instructions



## PURPOSE OF USE

BeautyRelax EyePen RF is a new generation cosmetic device that features the popular multipolar radiofrequency. Highly effective technology for lifting wrinkles in different parts of the face. Direct stimulation of collagen and elastic fibres. The skin is more elastic and looks more radiant. Radiofrequency can be combined with vibration and red photon light. Wireless use thanks to the built-in rechargeable battery.


## RADIOFREQUENCY

The application head generates energy via radio frequency waves with very specific power and frequency. The radiofrequency waves are optimized to produce a thermal gradient that significantly affects the processes in the dermis layer. The energy delivered activates the remodelling of collagen fibres and elastin. The maximum possible targeting of energy even in areas with a more pronounced fat layer and the subsequent excitation of metabolic processes has a significant impact on the very effective smoothing of sagging and deformed contours.



## VIBRACE

The proper functioning of the lymphatic system applies not only to the body's limbs, but also to our face. If our lymph and blood circulation are in order, the removal of impurities from the blood is multiply improved. Improved blood circulation = stimulation of the lymphatic system = better removal of impurities and pollutants from the body = beautiful and healthy skin. A simple equation that applies not only to the face, but to our entire body. Gentle vibration massage is technologically one of the simplest, yet very effective skin treatment techniques.

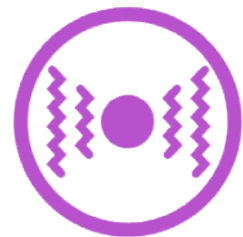

## PHOTON LIGHT

Light is simply the visible part of electromagnetic radiation of different frequencies. It reaches various wavelengths in the interval 390 to 760nm (nanometres). The human eye is capable of perceiving light with wavelengths ranging from 400 to 800nm. In photon therapy, light enters the body as a stream of energy in the form of photons. It is scientifically proven that different wavelengths of light result in specific skin care effects. Studies have verified real positive therapeutic benefits for living tissues and organisms at the cellular level. BeautyRelax EYEPEN RF offers a red light function with a wavelength of 620nm.




### ***Real positive benefits for living tissues and organisms at the cellular level***

- Photon light can enter the interior of the cell non-invasively and painlessly.
- Light activates photoreceptors in cells.
- Light therapy stimulates the healing power of cells in a natural way.
- Photon therapy does not cause any side effects or damage to the body.
- It addresses inflammatory diseases due to illness, injury or stress.



## DEVICE DESCRIPTION

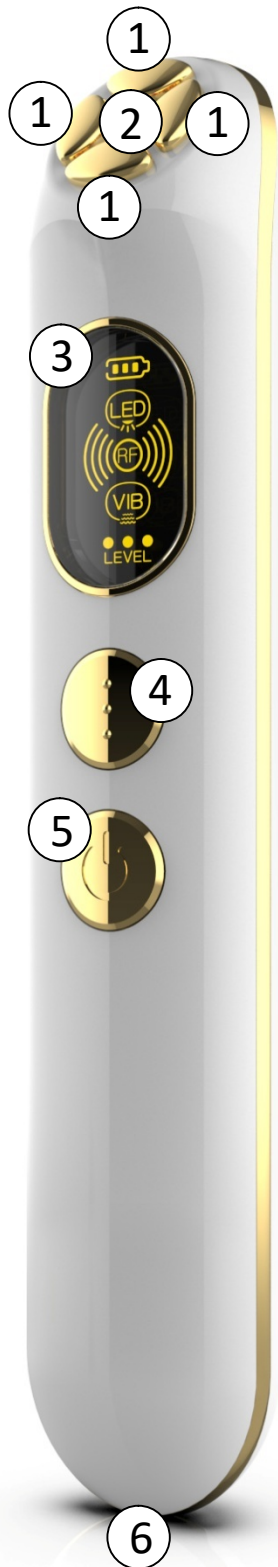

- 1. Application head (4 contact surfaces)
- 2. Red light
- 3. Display
- 4. **LEVEL** button for setting the radio frequency intensity
- 5. **ON/OFF** button
- 6. Input for USB power cable

## DISPLAY DESCRIPTION



1. Battery status indication
2. Red light indication
3. Indications for radiofrequency
4. Vibration indication
5. Intensity indication

## HOW TO USE


1. Skin should be thoroughly washed before each use. It is ideal to time your treatment after your bath. There should be no makeup, dirt or any cosmetic products on the skin. The skin should be perfectly dry.
2. Apply a fine layer of clear contact gel (not included) to the area to be treated. The contact gel, called SONOGEL, can be purchased at any medical supply store (it can also be purchased at a pharmacy). The gel is very important for the use of the device. It forms a necessary bonding layer between the application head and the skin. It is water-based. It cannot harm the application head. The use of contact gel is essential. It is not recommended to use the device without a binding agent between the skin and the applicator.
3. To switch on the device, long press the **ON/OFF** button. Device will start working immediately after switching on. The application head generates the following technologies:  
**radiofrequency + red light + vibration.**
4. With another short press of the **ON/OFF** button, the application head generates the following technologies:  
**radiofrequency + red light.**
5. With another short press of the **ON/OFF** button, the application head generates the following technologies:  
**radiofrequency + vibration.**
6. With another short press of the **ON/OFF** button, the application head generates the following technologies: **radiofrequency.**
7. **Use the LEVEL** button to adjust the intensity of the **radio frequency waves**. It is shown on the display using dots.  
**Low intensity** ●○○      **Medium intensity** ●●○      **High intensity** ●●●
8. You can switch off the device at any time by long pressing the **ON/OFF** button.

9. Instructions on how to move the applicator over the skin can be found in the passage of this manual: **METHOD OF TREATMENT OF THE FACE AND NECK.**
10. For instructions on how to properly charge the device battery, please refer to the passage in this manual: **CHARGING THE BATTERY.**
11. Instructions on how to properly care for the applicator can be found in the passage of this manual: **CLEANING AND MAINTENANCE.**

## HOW TO USE DEVICE - NOTES

1. Treat a selected part of the face one by one. Always proceed by holding the application head on the skin for a maximum of 1 to 2 seconds (stamping). Start at the lowest intensity setting. You can gradually increase the intensity according to your feeling. The skin treatment should never be unbearably painful. Try to adjust the length of contact between the applicator and the skin accordingly. Acceptable sensations are only gentle tingling to light pinching.
2. At higher power, the surface of the application head may reach higher temperatures than at lower power. If the temperature of the applicator head reaches a high temperature while the unit is in use, turn off the unit and allow the applicator head to cool down. Continue the skin treatment at a later time. Heating of the applicator head is not a malfunction. It is a side effect of the radiofrequency.
3. When using the device, always ensure that the application head (all 4 contact surfaces) is in full contact with the skin. This way the skin treatment will be as effective as possible.
4. You can alternate between different levels of radiofrequency intensity during a single skin treatment. Try to focus on the most problematic areas (fine lines on the forehead, crow's feet around the eyes, expression lines, contouring, neck lines, etc.). Use lower power on areas of thin, sensitive skin. Use higher power on areas with a higher fat layer.
5. **Please note:** Always keep the low level on when first coming into contact with skin. Apply the applicator to the part of the face where there is more subcutaneous fat. You can gradually add more to the strength. However, it always depends on your subjective feeling. If the treatment is already uncomfortable to slightly painful, adjust back to the lowest level, shorten the contact time between the applicator and the skin or stop the treatment.
6. **Tip:** The treated area should ideally be free of strong hairs. This will improve the movement of the application head over the skin. The radiofrequency waves will have a better opportunity to penetrate the subcutaneous tissue.
7. **Recommendation:** do not expose the treated skin to hot water for at least 3 days. The skin is still undergoing a recovery process after the treatment and hot water (or a stay in the sauna) can unnecessarily irritate the skin even more. The skin heats up when using the device. If it feels too warm, take a break and let your skin cool down.
8. After you have finished using the device, wipe the sonogel off your skin. If necessary, use a cream or spray to moisturize your skin. Your skin will regenerate while you sleep. Hydration will contribute significantly to faster skin recovery.

## POSSIBLE AREAS OF USE

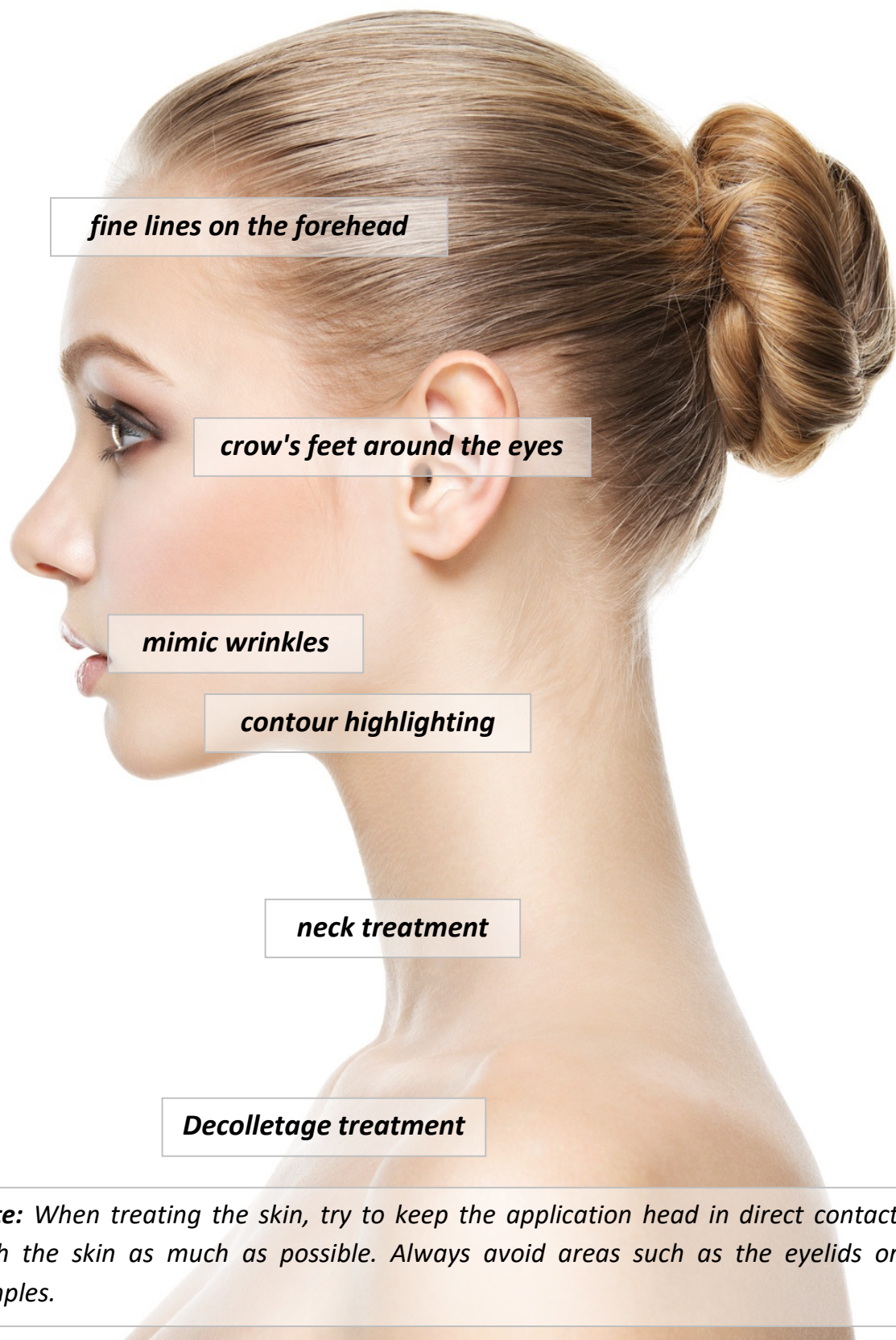

## METHOD OF TREATMENT OF THE FACE AND NECK



**FACE:** move away from the centre forehead to the edge, gradually cover the




**NOSE:** move from the forehead around the nose to the cheek, from top to bottom



**FACE:** move your face in the direction of the from the bottom up



**CHIN:** move from the tip of the chin along the contours of the face, from bottom to top

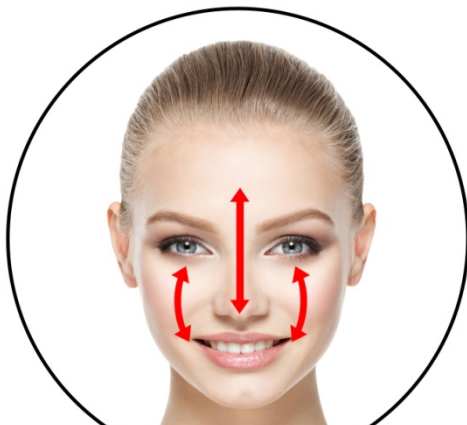

**T zone:** move from forehead to tip of nose and face, movement in both directions



**Neck:** move from the chin up to the ear and then down to the neck



# **FREQUENTLY ASKED QUESTIONS AND ANSWERS**

## **Why is there a significant price difference between devices in beauty salons and devices for home use?**

The main reason is the different frequency of use. Machines in beauty salons must be able to work continuously from morning to evening. This places increased demands on the materials used and the overall technological concept (larger power supplies, higher cooling rates, etc.) The performance of the home beauty machines is designed to meet the benefits announced to the user. Higher frequency of use = more expensive production costs.

## **What is the reason for cleaning/disinfecting the application head?**

Microscopic skin particles may adhere to the application head after the treatment. Cleaning prevents them from entering the device and ensures hygiene of operation.

## **How often can I use the device?**

It always depends on the intensity of the skin you are treating. It is important to remember that when using higher power, the skin needs some time to recover slightly. The device dries the skin slightly. It is best to use a moisturiser or ointment after each application to keep the skin hydrated.

## **How long should an individual skin treatment take?**

The recommended length of an individual treatment ranges from 10 to 15 minutes.

## **Are there any contraindications? (situations where the device cannot be used)**

The BeautyRelax EYEPEN RF is designed for all skin types without age limitation. The device should not be used if the user suffers from a serious skin condition. See the list of contraindications in the introduction of the manual.

## **Is the device a medical device?**

No, it's not. The device is not used for the treatment of serious skin diseases. The device does NOT treat: actinic keratosis, basalomas, melanomas, spinaliomas, atopic eczema, contact eczema, solar dermatitis psoriasis.

## **What should I do if the battery loses capacity after a long period of time?**

If the battery capacity decreases significantly over time, please do not hesitate to contact us at [mh@santadria.cz](mailto:mh@santadria.cz) We will send you information on how to replace the battery.

## **Does the device have any certification?**

device for standard European certification. It meets all the standards that are imposed on devices of this type.

## **Where is the device manufactured?**

We have the device manufactured in China according to our requirements and specifications. We guarantee standard warranty and post-warranty service.



**Can more than one person use the device?**

Primarily we don't recommend it. However, if family members are involved and the cleaning and disinfection process of the applicator head is followed after each treatment, multiple people can use the device.

**What cosmetic products should I combine during BeautyRelax EYEPEN RF treatment?**

During the treatment itself, it is recommended to use classic clear sonogel. When using cosmetic products or serums of a different composition, the treatment in combination with radiofrequency may cause an allergic reaction on the skin if your skin is more sensitive. The use of cosmetic products is ideal after the treatment when the skin needs to be sufficiently moisturized.

**What parts of the body can be treated with the device?**

The device is designed for the treatment of the face, neck and possibly the décolleté.

## **BATTERY CHARGING**

The device has a built-in rechargeable battery with a capacity of 500mAh.

**Note:** *Do not connect the charging cable to a 220V mains adapter. This could damage the battery and void the warranty. Use only the USB port of your computer or laptop as a power source.*

**Charging procedure:**

1. For charging, use only the power cord included in the package.
2. Connect one end of the USB cable to the cosmetic device.
3. Connect the other end to a USB port on your computer or laptop.
4. The charging progress is indicated by a flashing battery icon on the display.
5. Once the charging process is complete, the battery icon on the display will stop flashing.
6. Charging time may vary depending on the battery status.

## **CLEANING AND MAINTENANCE**

Always ensure that the application head is cleaned after the treatment is completed. Do not use any aggressive cleaning agents. We recommend the use of special disinfectant sprays, which are commonly used in e.g. healthcare or beauty salons for cleaning device. Normal plain alcohol will also suffice. Moisture must NEVER penetrate directly into the device. It is advisable to remove the application head from the device before cleaning. A cotton swab or dry cloth is ideal for cleaning. No residue of gel, cosmetic cream or serum should remain on the applicator before using the device again. Finally, always check that the applicator head is completely dry. Never use any sharp objects that could damage the surface of the applicator head. This could damage the contact surface of the head.

## **DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT**



This symbol on products or in accompanying documents means that used electrical and electronic products must not be added to normal municipal waste. For proper disposal, refurbishment and recycling, take these products to the designated collection points. Alternatively, in some European Union or other European countries, you can return your products to your local retailer when purchasing an equivalent new product. By properly disposing of this product, you will help conserve valuable natural resources and help prevent potential negative environmental and human health impacts that could result from improper waste disposal. For further details, contact your local authority or nearest collection point. This symbol is valid in the European Union. If you wish to dispose of this product, please obtain the necessary information on the correct method of disposal from your local authorities or your dealer.

## **CERTIFICATION**



The product is approved for use in European Union countries and bears the CE mark. The declaration of conformity and all documentation is available from the importer.

## **IMPORTER**

**SANTADRIA s.r.o.** , Mojmírovo nám. 75/2 Brno, 612 00  
IČO: 29195454 | DIC: CZ29195454

If you have any questions about the product, please do not hesitate to contact us:

Email: [mh@santadria.cz](mailto:mh@santadria.cz)

**Printing errors, features and changes in technical parameters reserved.**