

# USER MANUAL

Beautyrelax BR-2610 photon therapy device

## LIGHTMASK EXCLUSIVE





## CONTENT:

NOTICE	4
CONTRAINDICATIONS	5
PACKAGING CONTENT	5
TECHNICAL PARAMETERS	5
INTRODUCTION OF THE DEVICE	6
MAIN BENEFITS	6
PHOTON LIGHT THERAPY	6
BENEFICIAL EFFECTS OF PHOTON THERAPY	6
DEVICE DESCRIPTION	8
USE ACCORDING TO SKIN TYPE	9
HOW TO USE	10
BATTERY CHARGING	11
CLEANING AND MAINTENANCE	11
DISPOSAL OF USED ELECTRICAL EQUIPMENT	11
CERTIFICATION	11
SHIPPING	11
ES DECLARATION OF CONFORMITY	12

Dear customer, thank you for your trust in purchasing our product. Please read all the detailed information in the instruction manual carefully before using it for the first time. Please pay particular attention to the safety warnings. Please use the device exactly as instructed. All of these steps will help to maximize the use of all functions of the device. In this manual we will try to give you the best possible idea of the technology used and how to use the device to best advantage.

If you encounter any uncertainties while reading the instructions, please do not hesitate to contact us at the following email address: mh@santadria.cz

**Keep the instructions for use carefully for later reference.**

## **NOTICE**

- Read this manual before first use, review all illustrations and keep the manual for future reference.
- If your skin is affected by a serious skin condition, do not use the device. Consult your general practitioner for first use.
- The device is primarily intended for your personal use. Do not share it with family members or friends. If more than one person is using the device, extra attention should be paid to hygiene.
- If the device is damaged, does not work properly or has fallen into the water, never put it into operation.
- If the unit is broken, do not touch the internal parts to avoid electric shock.
- If the device falls to the ground, it may be damaged. In this case, do not continue to use the device and contact the dealer's service department.
- Do not use the device near water or in places with increased humidity (bath, shower, swimming pool, etc.).
- The device is not washable. It must not be immersed in water or under a stream of running water.
- Do not leave the device unattended. Always switch off the appliance thoroughly after each use.
- If the device has been moved from a very cold to a very warm environment, wait at least 3 hours before using it.
- The appliance is not intended for use by children or persons whose physical or mental condition does not permit the safe use of the appliance.
- Do not place or use the device on hot heat sources (radiator, stove, fire, etc.)
- Do not attempt to disassemble or repair the device yourself.
- Store the device in a dry environment. Try to keep it protected from dust.
- Do not expose the device to direct sunlight.
- The distributor is not liable for damages caused by improper use contrary to the instructions in this manual. It is not liable for consequences caused by failure to observe safety warnings.
- The distributor is not liable for any damage to health or property caused by other uses.



## **CONTRAINDICATIONS**

**Under no circumstances should** the device be used by persons affected by any of these conditions:

- Epilepsy, viral disease
- Mental illness
- States of anxiety
- Depression
- In a state of high fever
- Heart disease, pacemaker, hypertension
- Persons who require constant supervision may use the device only after consulting a physician



## **PACKAGING CONTENT**

- 1x Beautyrelax BR-2610 Facial Mask
- 1x USB charging cable
- 1x User manual

## **TECHNICAL PARAMETERS**

- Timer: yes, 28 minutes
- Operating voltage: 5V
- Operating current: 1A
- Power: 850mAh rechargeable Li-ion battery
- Charging: via USB cable (included)
- Dimensions: 19,7×19,5×13cm
- Weight: 190g

## **INTRODUCTION OF THE DEVICE**

Photon light therapy is very popular in the field of cosmetic treatment. Photon light treatment does not hurt and has very good results. The BeautyRelax Lightmask Exclusive photon therapy device provides photon therapy with full facial coverage. Your eyesight is maximally protected during use. The face mask promises a highly comfortable use compared to manual devices. You don't have to hold the device in your hand and at the same time the light covers the whole face evenly. Consistency of exposure to photon light is crucial for successful results.

## **MAIN BENEFITS**

- Even coverage of the face with light strips
- The light weight of the mask does not represent a large burden on the head
- Practical control with one button
- Easy wireless use
- Powered by a built-in rechargeable battery

## **PHOTON LIGHT THERAPY**

Light is simply the visible part of electromagnetic radiation of different frequencies. It reaches various wavelengths in the interval 390 to 760nm (nanometres). The human eye is capable of perceiving light with wavelengths in the range of 400 to 800nm. In photon therapy, light enters the body as a stream of energy in the form of photons. It is scientifically proven that different wavelengths of light result in specific skin care effects. Studies have verified real positive therapeutic benefits for living tissues and organisms at the cellular level.



## **REAL POSITIVE BENEFITS FOR LIVING TISSUES AND ORGANISMS AT THE CELLULAR LEVEL**

- Photon light can enter the interior of the cell non-invasively and painlessly.
- Light activates photoreceptors in cells.
- Light therapy stimulates the healing power of cells in a natural way.
- Photon therapy does not cause any side effects or damage to the body.
- It addresses inflammatory diseases due to illness, injury or stress.

## **BENEFICIAL EFFECTS OF PHOTON THERAPY**

Inflammatory diseases can occur as a result of illness, injury or stress. Light therapy stimulates the healing power of cells in a natural way. Light can enter the interior of the cell non-invasively and painlessly. Light activates photoreceptors in the cells. Photon therapy causes no side effects or damage to the body. Thanks to its positive effects, it has gained popularity not only in the field of dermatology, but also in the segment of cosmetic skin treatment.



### **Red light** **wavelength 660 to 665nm**

Stimulates the promotion of collagen formation. Collagen is an essential protein that serves to repair damaged and aged tissues. Improving skin quality and reducing pores. Helps to slow down skin aging.



### **Blue light** **wavelength 465 to 470nm**

The bacteria in acne contain porphyrins. Due to its high affinity with the wavelength, blue light is able to destroy porphyrins. In addition, blue light has a soothing effect that is very effective for hypersensitive skin.



### **Purple light** **wavelength 395 to 440nm**

Promotes skin cell regeneration. Targets the epidermis layer. Improves blood circulation and cell metabolism.




### **Orange light** **Wavelength 600 to 605nm**

Orange light helps lighten pigment spots. Targets the dermis layer. Promotes blood circulation and activates cells. Suppresses the formation of dark spots. Ideal for brightening the skin.

## DEVICE DESCRIPTION



- 1. ON/OFF button
- 2. Input for charging cable
- 3. Head strap
- 4. Eyeshield
- 5. Chin support
- 6. USB port for charging cable

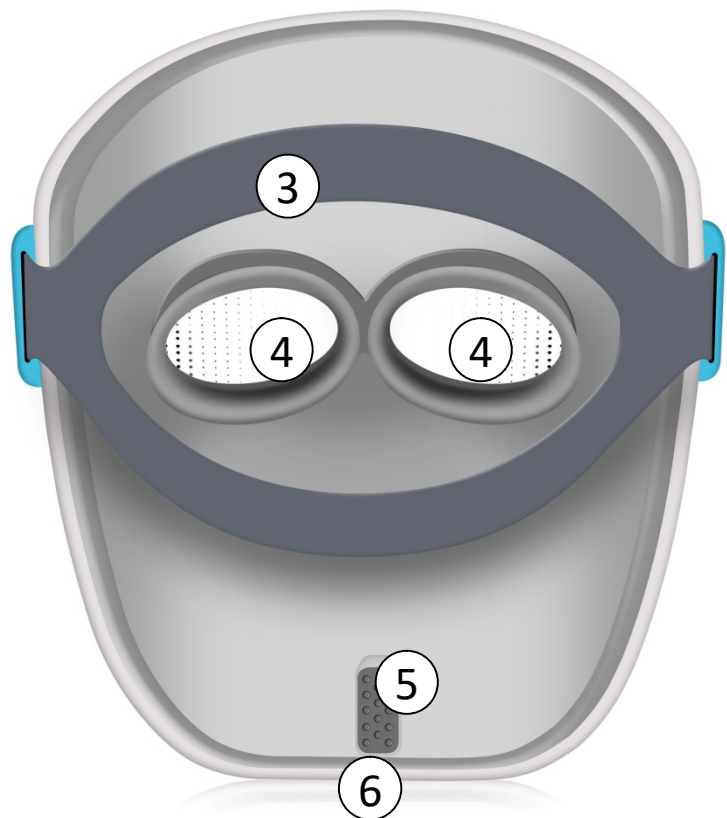



## **USE ACCORDING TO SKIN TYPE**

The gentleness of photon therapy treatments allows for the treatment of all skin types.



normal



dry



Mixed



Greasy



Sensitive



acneous

## **HOW TO USE**

1. We recommend removing makeup from the skin before ingestion. Wash and dry the skin thoroughly.
2. If necessary, you can apply your favourite cosmetic cream or serum.
3. Make sure that the battery of the device is sufficiently charged.
4. Put the device on your head in a way that is comfortable for you. The mask should not press too hard on your face. On the other hand, it should sit firmly on your face.
5. Press the **ON/OFF** button for a long time. The device will switch on. Switching on is accompanied by an acoustic signal. The red light automatically starts.
6. Each colour has a preset timer of 7 minutes. It then automatically switches to the next colour. If you want a specific colour to work, you can switch the individual colours of the light with a short press of the **ON/OFF** button.
7. To ensure comprehensive skin care with photon light, we recommend that you let the device gradually change all colors of light.
8. After 28 minutes, the device automatically switches off.
9. The ideal duration of use is 30 to 60 minutes (one or two full automatic treatments).
10. The frequency of use is optimal 2 to 3 times a week.
11. You can switch the device off at any time by long pressing the **ON/OFF** button. Switching off is accompanied by an acoustic signal.

## **BATTERY CHARGING**

The device has a built-in rechargeable battery. A sufficiently charged battery allows the device to function properly.

### **Charging procedure:**

1. For charging, use the charging adapter included in the package.
2. Once the battery is low, the device switches off and cannot be switched on.
3. Connect the other end to the charging socket.
4. The charging progress is indicated by a red light. Full charge - indication goes off.
5. The length of charge depends on the state of discharge and the overall condition of the battery.

**Please note: Only use the included adapter to charge. Using a different type of adapter may cause fatal damage to the battery and void the warranty.**

## **CLEANING AND MAINTENANCE**

Keep the surface of the device clean. Especially the area on the inside of the mask where the light therapy is applied. Do not use any aggressive cleaning agents. Moisture must NEVER penetrate directly into the device. A cotton cosmetic pad or dry cloth is ideal for cleaning. Never use any sharp objects during maintenance that could damage the surface of the device.

## **DISPOSAL OF USED ELECTRICAL EQUIPMENT**



This symbol on the products or in accompanying documents means that used electrical and electronic products must not be added to normal municipal waste. For proper disposal, refurbishment and recycling, take these products to the designated collection points. Alternatively, in some European Union or other European countries, you can return your products to your local retailer when purchasing an equivalent new product. By properly disposing of this product, you will help conserve valuable natural resources and help prevent potential negative environmental and human health impacts that could result from improper waste disposal. For further details, contact your local authority or nearest collection point. This symbol is valid in the European Union. If you wish to dispose of this product, please obtain the necessary information on the correct method of disposal from your local authorities or your dealer.

## **CERTIFICATION**



The product is approved for use in European Union countries and bears the CE mark. The declaration of conformity and all documentation is available from the importer.

## **SHIPPING**

SANTADRIA s.r.o. , Mojmírovo nám. 75/2 Brno, 612 00 IČO: 29195454 | DIC: CZ29195454

If you have any questions about the product, please do not hesitate to contact us:

Email: [mh@santadria.cz](mailto:mh@santadria.cz) Tel: 737 196 160

**Printing errors, features and changes in technical parameters reserved.**

# ES DECLARATION OF CONFORMITY

The undersigned, representing the company

**SANTADRIA s.r.o.**

Mojmírovo nám. 75/2, Brno 612 00, ID: 29195454

hereby declares, as an importer into the EU, that the following product with the trade name

**LIGHTMASK EXCLUSIVE face mask**

marketed under the model name: **BR-2610**

EAN: **8594166119482**

belonging to the **BeautyRelax** brand group of products

is in compliance with the following European standards:

**EN 55014-1:2017 EN 61000-3-2:2019**

**EN 61000-3-3:2013+A1:2019**

**EN 55014-2:2015**

**EN 61000-4-2:2009 EN 61000-4-3:2006+A1:2009+A2:2010**

**EN 61000-4-4: 2012 EN 61000-4-5: 2014+A1:2017**

**EN 61000-4-6:2014/AC:2015, EN 61000-4-11:2004+A1:2017**

**Based on compliance with these directives, the product has been CE marked.**

The product also complies with RoHS Directive 2002/95/EC (RoHS). All materials are health-safe. They do not contain any heavy metals such as cadmium and other toxic substances.

For safe use, it is essential to strictly follow the instructions in the instruction manual included in the package.



Martin Hrbek, Managing Director

Place of issue: Brno Day: 1.4.2022