

USER MANUAL

BeautyRelax BR-2620
cosmetic device for skin elasticity

EMSLIFT SMART



CONTENT:

NOTICE	4
CONTRAINDICATIONS	5
MAIN ADVANTAGES OF THE DEVICE	5
FEATURES & SPECIFICATIONS	5
PACKAGING CONTENT	6
INTRODUCTION OF THE DEVICE	6
EMS LIFTING	6
DESCRIPTION OF THE INSTRUMENT	7
HOW TO USE	8
FACIAL TREATMENT OPTIONS	9
SOLVING POTENTIAL PROBLEMS	10
FREQUENTLY ASKED QUESTIONS AND ANSWERS	11
BATTERY CHARGING	12
CLEANING AND MAINTENANCE	12
DISPOSAL OF USED ELECTRICAL EQUIPMENT	13
CERTIFICATION	13
IMPORTER	13

Dear customer, thank you for your trust in purchasing our product. Please read all the detailed information in the instruction manual carefully before using it for the first time. Please pay particular attention to the safety warnings. Please use the device exactly as instructed. All of these steps will help to maximize the use of all functions of the instrument. In this manual we will try to give you the best possible idea of the technology used and how to use the instrument to best advantage.

If you encounter any uncertainties while reading the instructions, please do not hesitate to contact us at the following email address: mh@santadria.cz

Keep the instructions for use carefully for later reference.

NOTICE

- Read this manual before first use, review all illustrations and keep the manual for future reference.
- If your skin is affected by a serious skin condition, do not use the device. Consult your general practitioner for first use.
- The device is primarily intended for your personal use. Do not share it with family members or friends. If more than one person is using the device, special attention should be paid to hygiene.
- If the device is damaged, does not work properly or has fallen into the water, never put it into operation.
- If the unit is broken, do not touch the internal parts to avoid electric shock.
- If the device falls to the ground, it may be damaged. In this case, do not continue to use the instrument and contact the dealer's service department.
- Do not use the device near water or in places with increased humidity (bath, shower, swimming pool, etc.).
- The device is not washable. It must not be immersed in water or under a stream of running water.
- Do not leave the device unattended. Always switch off the appliance thoroughly after each use.
- If the device has been moved from a very cold to a very warm environment, wait at least 3 hours before using it.
- The appliance is not intended for use by children or persons whose physical or mental condition does not permit the safe use of the appliance.
- Do not place or use the device on hot heat sources (radiator, stove, fire, etc.)
- Do not attempt to disassemble or repair the device yourself.
- Store the device in a dry environment. Try to keep it protected from dust.
- Do not expose the device to direct sunlight.
- The distributor is not liable for damages caused by improper use contrary to the instructions in this manual. It is not liable for consequences caused by failure to observe safety warnings.
- The distributor is not liable for any damage to health or property caused by other uses.
- The device is not a medical device.



CONTRAINDICATIONS

Under no circumstances should the device be used by persons affected by any of these conditions:

- Cancer, epilepsy, viral diseases
- In a state of high fever
- Serious kidney or liver disease
- Heart disease, pacemaker, hypertension
- Metal implants in the body
- Vascular disorders or arteriosclerosis
- Pregnancy
- Acneiform keratosis, basaliomas, melanomas, spinaliomas, atopic eczema, contact eczema, solar dermatitis, psoriasis
- Persons taking anti-clotting medication
- Persons who require constant supervision may use the device only after consulting a physician
- Never move the device on the eyelids, temple area or thyroid area



MAIN ADVANTAGES OF THE DEVICE

- Completely non-invasive nature of the treatment.
- Without any discomfort.
- You will only feel a slight tingling sensation when using it.
- The skin treatment is time-saving.

FEATURES & SPECIFICATIONS

- EMS skin stimulation pulses
- 4 intensity levels
- Designed for sagging and unhealthy looking skin
- Increase the elasticity of facial contours
- Intensive stimulation of collagen and elastic fibres
- Significant skin unevenness caused by age
- Largest application head compared to devices of similar type
- Completely non-invasive nature of the treatment
- Without any discomfort
- You will only feel a slight tingling sensation when using
- The skin treatment is time-saving
- Designed for face and décolleté treatment
- Compact dimensions
- Wireless (easy to use)
- Power supply: li-ion battery
- Charging: USB cable (included)
- Dimensions: 5,6×13,1cm
- Weight: 83g

PACKAGING CONTENT

- BeautyRelax BR-2620.
- USB charging cable.
- User manual.

INTRODUCTION OF THE DEVICE

BeautyRelax Emslift Smart skin elasticity device. Designed for mature and sagging skin. Maximum effect of collagen and elastic fibre stimulation. Regular use keeps the skin in great condition. Increase skin elasticity, focusing on pronounced wrinkles and uneven skin. Strengthening of facial muscles.

EMS LIFTING

The face is interwoven with a dense structure of muscles. The facial muscles are a group of muscles that result in facial expression and facial expressions. However, over the course of life, they lose their elasticity and flexibility. The power and frequency of the electrical pulses are specially adapted to enable facial treatment. The application head of the cosmetic device transmits gentle electrical pulses on contact with the skin, which penetrate through the upper part of the skin directly to the facial muscles. The muscles are contracted and thus strengthened. Collagen and elastic fibres are mechanically stimulated. In the facial area you will achieve an improvement in skin texture, skin tone and an increase in elasticity. The skin has a rejuvenated and radiant appearance. Cosmetic devices with electrical muscle stimulation function are used exclusively for facial treatments. The application head stimulates those parts of the face where the fat and muscle tissue is greatest. This improves the elasticity and firmness of the delicate facial muscles.



KEY EFFECTS OF EMS TECHNOLOGY

- support of facial blood circulation
- stimulation of the facial lymphatic system
- increase skin elasticity
- reduction of fine facial wrinkles
- remodelling of collagen fibres
- increase the elasticity of facial contours

DESCRIPTION OF THE INSTRUMENT



- 1. Application head**
- 2. button +**
- 3. Button -**
- 4. Power cable input**
- 5. Light indication - EMS pulse intensity level 1 (low)**
- 6. Light signaling - EMS pulse intensity level 2 (basic)**
- 7. Light signaling - EMS pulse intensity level 3 (high)**
- 8. Light signaling - EMS pulse intensity level 3 (highest)**

HOW TO USE

1. Before proceeding to the actual treatment, remove oil and make-up from the skin. Apply a thin layer of sonogel to the skin.
2. Press the (+) button to switch on the device. The application head generates electrical micropulses at low level 1.
3. Press the (+) button again to increase the micropulse intensity by one level to level 2.
4. Press the (+) button again to increase the micropulse intensity by one level to level 3.
5. Press the (+) button again to increase the micropulse intensity by one level to level 4.
6. Use the (-) button to reduce the intensity of the micropulses according to your needs. It is recommended to always start from the lowest value and add gradually. Each press of the button is accompanied by an acoustic signal.
7. The course of action of micropulses is accompanied by an acoustic signaling - beeping. This is an accompanying phenomenon of the device being switched on.
8. The control button operates in OFF mode - 1 - 2 - 3 - 4 - OFF.

Note 1: Always move slowly over the skin.

Note 2: Ensure that both application heads are in contact with the skin during treatment.

FACIAL TREATMENT OPTIONS



FACE: move away from the centre forehead to the edge, gradually cover the entire




NOSE: move from the forehead around the nose to the cheek, from top to bottom



FACE: move your face in the direction of the from the bottom up



CHIN: move from the tip of the chin along the contours of the face, from bottom to top

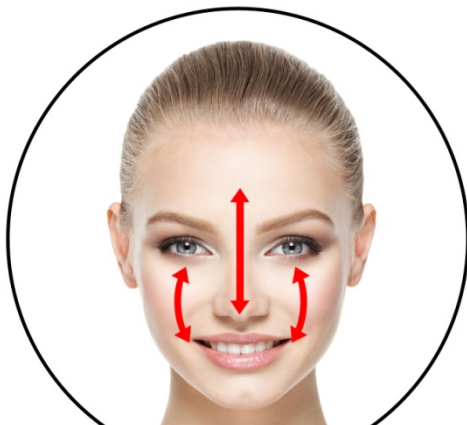

T zone: move from forehead to tip of nose and face, movement in both directions



Neck: move from the chin towards the ear and then down to the neck upwards

SOLVING POTENTIAL PROBLEMS

The Problem	Possible cause	Solution
The machine won't turn on.	The device is connected to a USB charging cable.	Disconnect the unit from the power cord before using the unit.
	The battery's dead.	Charge the instrument battery sufficiently before use.
Charging problem.	The charging cable is not sufficiently connected to the USB port on your computer or laptop.	Reconnect the end of the charging cable to the USB port on your computer or laptop.
	The charging cable is not sufficiently connected to the device.	Reconnect the end of the charging cable from the device.
	The charging cable is damaged.	Replace the charging cable with a new one.
The action of the electric micropulses is too strong.	The intensity level is too strong.	Reduce the intensity level to a lower level.
	Sonogel is on the skin in too small amounts or not at all.	Apply a larger amount of gel to the skin to increase skin conductivity.

FREQUENTLY ASKED QUESTIONS AND ANSWERS

How do I know if the EMS pulses are working properly?

As the application heads move over the skin, you may feel the presence of EMS pulses as a gentle tingling sensation. At that precise moment, all the structures of the fine facial wrinkles are stimulated. The pulses stimulate the collagen and elastic fibres of the skin.

I move the device over my skin and feel no pulses. Is that okay?

The effect of electric micropulses is very subjective and individual. It depends on several factors. The perception of electric pulses is weaker if you have a thicker fat layer on the treated area. The cause may lie in a small **amount of** cosmetic cream on the skin. You can also change the grip angle of the device during use.

How to choose the right intensity of the EMS function?

The device offers a total of 4 levels of EMS pulse intensity. Each user may have a different layer of subcutaneous fat. The fat acts as a resistance. If there is more subcutaneous fat, the electrical pulses must have more power to penetrate to the dermis layer.

Can I use the device if I have make-up on my skin?

Always remove any make-up residue before using the device. It would create a barrier between the massage rollers and the skin. The electric micropulses would not be able to penetrate the skin sufficiently.

How often and for how long to use the device?

There is no significant restriction in this area. The optimum length of treatment should be between 15 and 20 minutes. It can be done every day while resting or watching TV. In the morning before applying make-up or in the evening before going to bed. If the skin is more irritated at higher power, take a day or two off before using the device again.

Can I use the device on my body in addition to my face?

Thanks to the size of the rollers, the device can be used not only on the face, but also on the body. It is always better to use a body cream to increase the skin's conductivity and allow the electric micropulses to have the best possible effect on the skin.

Is the device high maintenance?

After each skin treatment, pay particular attention to the application heads. Always make sure there is no residual moisture on them. Use rubbing alcohol once in a while for maximum disinfection.

When can I expect some results?

The speed of achieving the beneficial effects of the device are very individual. It depends on many factors: the degree of subcutaneous fat, age, gender. Continuity of the instrumental treatment is essential to achieve results. Thanks to the compact size of the device, you can take it with you on the go. Try to make the device a regular part of your daily skincare ritual.

Is use prohibited in some cases?

There are situations where the use of the device is not recommended. Follow the passage in this manual called contraindications.

BATTERY CHARGING

The device has a built-in rechargeable battery. The input for the charging cable is located on the bottom of the device. The charging cable is included in the package. The other end plugs into a USB port on your computer or laptop.

Charging procedure:

1. For charging, use the power cord included in the package.
2. Connect one end of the cable to the charging stand. The other to a USB port on your computer or laptop.
3. The charging process is accompanied by a pulsating light signal.
4. Once the battery is charged, the light pulsing stops.
5. The length of the charge depends on the state of the battery. The charging process can take up to 2 hours.

When the battery capacity drops below 20%, the battery symbol will flash. If the battery is low, the device will automatically switch off.

Do not connect the charging cable to any type of 220V mains adapter. This could damage the device. Always use the USB port of your computer or laptop for charging.

CLEANING AND MAINTENANCE

Always ensure that the applicators are cleaned sufficiently after the treatment. Do not use any aggressive cleaning agents. We recommend the use of special disinfectant sprays, which are commonly used in e.g. healthcare or beauty salons for cleaning instruments. Normal plain rubbing alcohol will also suffice. Moisture must NEVER penetrate directly into the instrument. It is advisable to remove the application head from the instrument before cleaning. A cotton swab or dry cloth is ideal for cleaning. No residue of gel, cosmetic cream or serum should remain on the applicator before using the device again. Finally, always check that the applicator head is completely dry. Never use any sharp objects that could damage the surface of the applicator head. This could damage the contact surface of the head.

DISPOSAL OF USED ELECTRICAL EQUIPMENT



This symbol on products or in accompanying documents means that used electrical and electronic products must not be added to normal municipal waste. For proper disposal, refurbishment and recycling, take these products to the designated collection points. Alternatively, in some European Union or other European countries, you can return your products to your local retailer when purchasing an equivalent new product. By properly disposing of this product, you will help conserve valuable natural resources and help prevent potential negative environmental and human health impacts that could result from improper waste disposal. For further details, contact your local authority or nearest collection point. This symbol is valid in the European Union. If you wish to dispose of this product, please obtain the necessary information on the correct method of disposal from your local authorities or your dealer.

CERTIFICATION



The product is approved for use in European Union countries and bears the CE mark. The declaration of conformity and all documentation is available from the importer.

IMPORTER

SANTADRIA s.r.o. , Mojmírovo nám. 75/2 Brno, 612 00
IČO: 29195454 | DIC: CZ29195454

If you have any questions about the product, please do not hesitate to contact us:
Email: mh@santadria.cz Tel: 737 196 160