### EN – INSTRUCTIONS FOR USE


#### Pack Contents

- Base with UV lamp
- Brush handle
- 3 brush heads

#### Product Description


##### Brush Handle

- 1 Power button
- 2 Cleaning mode button
- 3 Mode indicator
- 4 Charging indicator
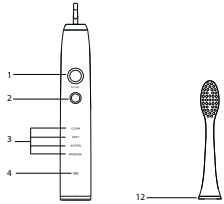
##### Base

- 5 UV lamp cover
- 6 Base
- 7 UV lamp
- 8 Sanitizer cover
- 9 Charging cradle
- 10 Power cord
- 11 Plug
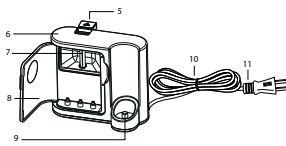
##### Brush Heads

- 12 Identification ring

#### Charging

Before charging, make sure the brush handle is turned off. Insert the plug into the socket and place the handle into the charging cradle. When charging, the charging indicator light lights red on the handle. When fully charged the indicator turns green.

1. A full charge lasts for about 24 hours. Full capacity is sufficient for 20 pre-set cycles or approximately 10 days when cleaning teeth twice a day.
2. The brush must be charged for 24 hours:
  - Before first use
  - If it hasn't been used for 3 months
  - If the battery is completely empty
3. To maintain full capacity, keep the brush in the stand when idle. The stand must be plugged into a wall socket.
4. When charging, the brush may become hot.

#### Instructions for Use

1. Attach the brush head to the brush handle. There is a gap of 0.5 mm between the head and handle.
2. Press the button to select the cleaning mode. The brush remembers the last mode used, so there is no need to reselect the last mode.
  - CLEAN – basic mode for perfectly clean teeth
  - SOFT – gentle cleaning for sensitive teeth and gums (also suitable for users not accustomed to sonic cleaning)
  - WHITEN – cleaning and polishing for whitening teeth
  - MASSAGE – gentle pulse stimulation of gums
3. Moisten the head and apply toothpaste.
4. Place the bristles on the teeth, towards the gums.
5. Turn the brush on using the ON/OFF button.
6. Move the bristles back and forth and slowly move between the teeth. Brush whilst gently applying pressure on the tooth.
7. For maximum efficiency, divide the cleaning into 4 sections. Clean each section for 30 seconds. After the 30 second cycle, the brush pauses and alerts you to move onto the next section. The brush turns itself off after 2 minutes.
8. Remove the head from the handle.
9. Rinse the head and body under running water and wipe with a cloth or towel.
10. Place the handle back into the cradle.

#### UV Lamp

- The base includes a UV lamp that disinfects the cleaning head after each use. The UV lamp can remove 99.9% of common bacteria such as E. Coli, Salmonella, Candida Albicans and the H. Simplex or H1N1 viruses.
- Only disinfect the original head with the UV lamp.
  - The cleaning cycle lasts around 10 minutes and when finished the lamp switches itself off.
  - The cycle starts automatically when the sanitizer cover is closed.

1. Clean and dry the head.
2. Open the UV lamp cover.
3. Place the head in the rack with the bristles pointing towards the UV lamp.
4. Make sure the charging base is plugged in.
5. Close the sanitizer cover and the disinfection process automatically begins.
  - The blue illumination of the UV lamp indicates that the lamp is on and the cleaning cycle is in progress. The light might not be noticeable in brightly lit spaces.
  - After the cleaning cycle is finished, the lamp automatically turns off and the light goes out.

#### Maintenance

Before cleaning, make sure the device is not plugged in.

##### Base

Wipe the dirt from the charging cradle and UV lamp with a soft cloth.

##### Brush

1. Remove the brush head.
2. Wipe the dirt off all surfaces with a soft cloth and put the brush back on the charging base. Do not immerse the handle completely in water.
3. Rinse the brush head and bristles after each use.
4. Do no wash the brush or any of its parts in the dishwasher.

#### Troubleshooting

Problem	Possible Cause	Solution
Brush does not work	Brush has never been turned on or has not been used for more than 3 months	Allow brush to charge for at least 24 hours
Weak or intermittent vibration	Brush is pressed too hard on the teeth	Reduce the amount of pressure on the teeth
Brush only works for a short amount of time or switches off immediately when disconnected from the charging base	Insufficient charging time Battery is at the end of its life cycle	Allow the brush to charge for at least 24 hours Contact customer support
UV lamp is not lit	Sanitizer cover is not closed Base not plugged in UV lamp is damaged	Close cover Plug the base into the wall socket Contact customer support

#### Safety Instructions

- Consult your dentist when using the product.
- If you experience health problems during use, stop using the product and consult your doctor.
- The UV lamp may be harmful to the eyes and skin. If it lights up after opening the cover, stop using the product and contact customer support.
- The battery capacity may diminish over time. If the brush can only last for a single cleaning cycle, contact customer support.
- Keep the lamp out of reach of children.
- If the UV lamp shatters, do not touch it and leave it for 30 minutes after which time it is safe to clear the shards and contents.

#### Responsibility of the Manufacturer

TrueLife products are subject to a warranty that covers manufacturing defects. In the event of such defects, the company elem6 s.r.o. is obliged under warranty to either repair or replace the complete product as we may choose. If any conversion, alteration or repair is undertaken by an unauthorized person, the guarantee becomes void. The warranty does not cover any damage resulting from improper use, normal wear or caused accidentally. The company elem6 s.r.o. assumes no responsibility for any accidents, injuries, death, loss or other claims relating to or arising from the use of this product. The company elem6 s.r.o. in no event shall be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Repayment or replacement of parts and/or products may be subject to shipping, handling, replacement and/or shipping charges.

#### Printing error and information changes reserved.

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