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SK	Ošetrujúca LED maska na tvár.....	5
EN	Treatment LED mask for the face.....	9
HU	Kezelés LED arcmaszk.....	13
RO	Mască LED de tratament pentru față.....	17
DE	LED-Pflegemaske für Gesicht.	21

CZ Palsar7

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TREATMENT LED MASK FOR THE FACE

PHOTONE LIGHT THERAPY

Why go for expensive chemical peels and skin cleansers, laser treatments and hyaluronic acid wrinkle fillers when you can treat your skin from the comfort of your home?

The Palsar7 treatment LED mask is an innovative way to get rid of acne, wrinkles, tired skin and other imperfections. This is thanks to photon light therapy, which is one of the most effective tools in the field of dermatology in the fight against chronic skin defects. Unlike laser treatments, photon treatment with LED lights does not cause any pain or side effects that require time to heal.

HOW PHOTONE THERAPY WORKS

Light can be simply described as a visible part of electromagnetic radiation of different frequencies. It reaches different wavelengths in the range from 390 to 760nm. In this therapy, light penetrates the body as a stream of energy in the form of photons. Different wavelengths of light have different positive effects on skin care.

LIGHT WAVELENGTHS FOR THE LED TREATMENT MASK:

Red light: **630nm**

Yellow light: **590nm**

Green light: **520nm**

Blue light: **470nm**

Using this mask is very simple, you just need to put it on, fasten it to your face with an elastic band and adjust the intensity and duration of the selected treatment with the controller. There are soft gray rubber pads around the eye holes to prevent unpleasant pressure around the eyes. The treatment LED face mask is equipped with 150 LED diodes that produce photons with different wavelengths. Each diode has a defined power and luminous flux value. Therefore, it is important to distinguish the individual colors of the diodes. One color treats acne, another treats wrinkles, scars, redness, swelling - it depends on your preference.

The LEDs on the mask do not need to be changed.

EN TREATMENT LED MASK FOR THE FACE

EFFECTS OF TREATING LED MASKS:

- Reduces wrinkles and prevents them from occurring
- Minimizes pores
- Fights acne
- Helps fight facial grease
- Improves skin elasticity
- Unifies skin tone
- Suppresses melanin production
- Improves blood circulation in subcutaneous tissue
- Increases the efficiency of absorption of active ingredients of cosmetics

MAIN FUNCTION:

- Natural light waves improve the appearance of the skin
- Each light color is effective for a specific skin treatment
- Different mode options
- Easy operation
- Easily portable

MASK DESCRIPTION:

PAD FOR FOREHEAD

ELASTIC BAND

PADS FOR EYES

LED DIODES

USB ENTRY

THE MASK IS SUITABLE FOR THOSE WHO:


- Stay up late into the night and their skin is tired
- Sit long in front of the computer and often look at his cell phone
- Who do not follow a healthy diet, which improves the condition of the skin
- Who have problematic skin: spots, acne, clogged pores, dull, yellowish, rough or sagging skin, etc.

EN TREATMENT LED MASK FOR THE FACE

INSTRUCTIONS FOR USE:

- Before using the mask, test all functions of the controller without the mask on your face. Try to avoid direct eye contact with the light coming from the mask.
- Connect the mask to the controller using the USB cable. Connect the controller to the power cable.
- The power cord must be plugged in when using the mask.
- Adjust the length of the elastic band to your face.
- Before putting on the mask, thoroughly clean the skin of the face and neck and then gently dry with a towel.
 - Press the **"power switch"** button on the controller. If the mask is on, the green **"power"** light on the top left will light up.
 - Then set the light intensity using the **"energy + and -"** buttons to your desired value (01-10). By default, the value is set to "00" = light intensity is switched off.
 - Press the **"on/off"** button and the mask will light up. If the mask doesn't light up, you don't have the light intensity set. In this case, return to Step 2, otherwise continue with Step 4.
 - Use the **"color"** button to select the color of the treatment you want to undergo. Try not to look directly into the lights of the mask when choosing a color.
 - Use the **"time"** button to set the required treatment time, the recommended time is 10-15 minutes.
 - Put the mask on your face and start the whole skin treatment process with photon therapy.

CONTROLLER DESCRIPTION:

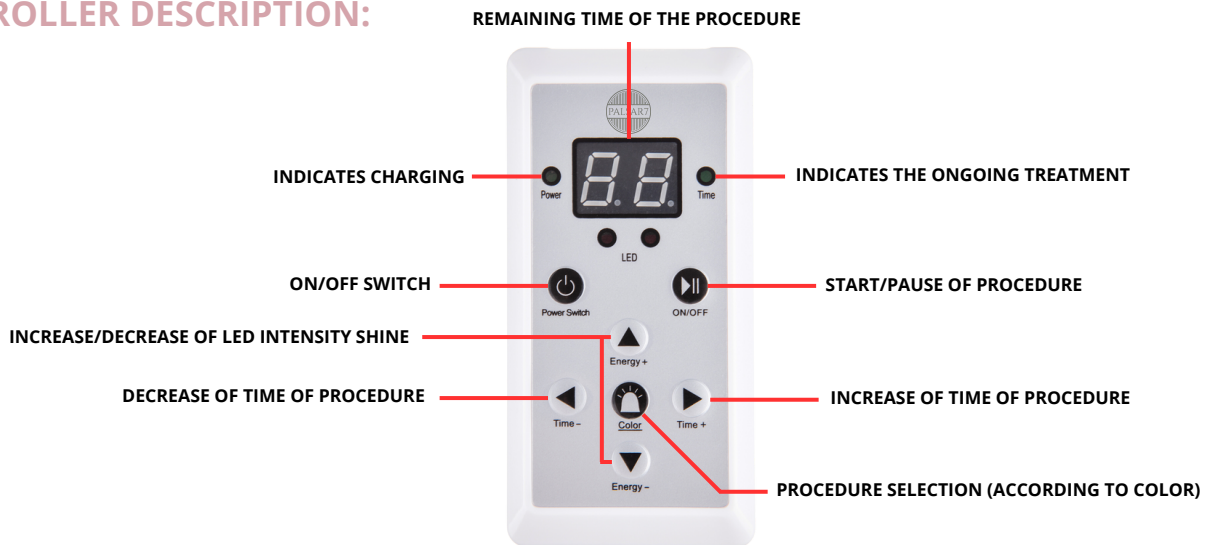The maximum time of using the mask per day is 45 minutes. If you want to try several treatments in one day, wait at least 10 minutes between them. The recommended duration of one treatment is 10-15 minutes. Ideally, use the mask 2-3 times a week. If you suffer from sensitive skin, eczema, or allergies have a cold, consult your doctor before using the mask.

IMAGES WITH DESCRIPTION ARE FOR ILLUSTRATION PURPOSES ONLY.
MADE IN CHINA.

EN TREATMENT LED MASK FOR THE FACE

EFFECTS OF INDIVIDUAL LIGHTS:

1) 630NM-LED RED LIGHT: ANTI-WRINKLE

- Bleaching effects
- For soft skin without wrinkles
- Smoothes fine wrinkles, shrinks pores
- Improves the condition of the skin
- It supports the natural production of collagen
- Increases the efficiency of absorption of active ingredients of cosmetics
- Improves blood circulation in subcutaneous tissue
- Prevents the formation of pigment

2) 520NM-LED GREEN LIGHT: ANTI-FATIGUE

- Improves the condition of the skin
- Balances skin tone
- Reduces sebum production
- Relieves mental stress
- Reduces the incidence of edema (swelling)
- It contributes to the elimination of epidermal pigment spots
- Slows down melanin production by 40 to 60%
- Has an anti-inflammatory effect (especially against the formation of acne skin)
- Suitable for all skin types

3) 470NM-LED BLUE LIGHT: AGAINST ACNE

- Effectively reduces the incidence and cures acne
- Improves the condition of the skin without leaving scars
- Activates porphyrins (start a chemical reaction that destroys harmful bacteria)
- It has antibacterial and anti-inflammatory effect

4) 590NM-LED YELLOW LIGHT: AGAINST REDNESS, TO INCREASE IMMUNITY

- Relieves skin roughness and redness
- Relieves skin diseases
- In addition, it relieves fever and increases immunity
- Supports the absorption of cosmetics
- It contributes to the improvement of cellular metabolism
- Improves detoxification properties by promoting lymphatic activity
- Particularly effective after intense mechanical peeling
- Helps slow down skin aging

5) LED PURPLE LIGHT: AGAINST ACNE SCARS

- Combination of red and blue light = combination of two treatment therapies
- Particularly effective in treating acne scars

6) LED GREEN-BLUE LIGHT: METABOLISM SUPPORT

- Gradually increases cellular energy
- Supports metabolism

7) LED WHITE LIGHT: ANTI-AGING

- Addresses age spots and sagging skin
- Reduces the appearance of fine wrinkles.

WARNING:

- Not suitable for children under 14 years.
- Prohibited to use: pregnant women, infants, children, people with colds, sore throat or flu, epileptics, allergy sufferers, patients with thyroid disease.
- It is forbidden to use by people who were ill after using the mask - please visit your doctor and consult with him when using the mask.
- When using the mask, do not watch TV, surf the Internet, read a book, etc. - there is a risk of long-term damage to the eyes.
- Avoid prolonged direct eye contact with the light source.
- If you intend to use cosmetics together with the mask, choose quality products with a natural composition.
- After opening the box, the mask may smell slightly, it is due to the fact that it is made of medical silicone, the odor will soon disappear.
- Whether you use the mask on a bed or a lounger, always keep it in a box to prevent damage.
- Do not use water to clean the mask, do not immerse the mask in water - use a damp cloth to clean it.
- Do not use any chemicals to clean the mask.
- Do not expose the mask to direct sunlight.

ACCESSORIES:

1x LED face mask, 1x remote control, 1x USB cable, 1x charging adapter, 1x user manual

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WARRANTY CARD

Product:
Model:
Order number:
Date of sale:
Warranty until:
Sale document number:
Seal of the store:

Information for Customers about Warranty Terms and Conditions:

Warranty shall apply to defects of purchased product, found during its warranty period. When purchasing a particular product, please consider carefully what properties and functions you expect your product to have. The fact that purchased product shall not satisfy your later requirements shall not give any reason for complaint. Please, prior to the first use of this product read carefully its operating instructions and follow these instructions consistently.

We recommend you to keep original product packing during the warranty period for possible transport of the product in such a way so as to prevent any product damage during its transport and handling.

For the purpose of due application of any complaint about the product we recommend you to ask Seller to fill-in and confirm this Warranty Card including identification of order number of purchased product, which shall comply with the number specified in the document proving product purchase and in the delivery note. In this connection Customer hereby takes into account that unless this Warranty Card specifies order number of purchased product, i.e., comparison of such order number with the number identified in this Warranty Card and document of product purchase/or delivery note shall not allow to confirm that particular product has come from Seller, in case of any doubts about the origin of such particular product Seller shall be entitled to automatically and without any other acts to reject the complaint of such goods.

Warranty shall not apply to:

- 1) Product defect that occurred by:
 - mechanical damage to the product or its parts
 - incorrect operating of the product
 - unprofessional intervention of a third party
 - by using the product contrary to its intended use in the enclosed operating instruction/user manual
 - improper storage or handling
 - transportation
 - natural disaster or other force measures
 - 2) Product damage suffered by unqualified installation, its modification or by static electricity
 - 3) Fault caused by external conditions (e.g. faults in the electrical network or residential installation, unsuitable charges, unsuitable operating conditions or damage by other equipment)
 - 4) A product for which the buyer (operator) detains, conceals or provides incorrect information of importance to the assessment of the defect
 - 5) External damage to the product that was clearly caused by man
 - 6) Product damage inconsistent with the operating instructions, technical standards or valid safety regulations
 - 7) Product damage resulting from normal wear and tear
 - 8) Discharged batteries, decrease in battery capacity caused by its wear and tear
 - 9) Product defects caused by the use in dusty, polluted or in any other way dissatisfactory environment pursuant to attached operating instructions
- Wear and tear of the product and its components caused by their common use (discharged batteries, loss of the original color of the product, decreased capacity of accumulator, wear and tear of mechanical or displaying components, wear and tear of supply cables, etc.) shall be subject neither to defect warranty nor quality warranty in the moment of take-over.
- Furthermore, complaints about product parameters, not included in the operating instructions or in any commercial-technical documentation of Seller, related to the product, shall not be subject to quality warranty in the moment of take-over.

Proof of product purchase, or delivery note, shall constitute an integral part of this Warranty Card.

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