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EN

Palsar7

TREATMENT LED MASK FOR THE FACE AND NECK

PHOTONE LIGHT THERAPY

Why go for expensive chemical peels and skin cleansers, laser treatments and hyaluronic acid wrinkle fillers when you can treat your skin from the comfort of your home?

The Palsar7 treatment LED mask is an innovative way to get rid of acne, wrinkles, tired skin and other imperfections. This is thanks to photon light therapy, which is one of the most effective tools in the field of dermatology in the fight against chronic skin defects. Unlike laser treatments, photon treatment with LED lights does not cause any pain or side effects that require time to heal.

HOW PHOTONE THERAPY WORKS

Light can be simply described as a visible part of electromagnetic radiation of different frequencies. It reaches different wavelengths in the range from 390 to 760nm. In this therapy, light penetrates the body as a stream of energy in the form of photons. Different wavelengths of light have different positive effects on skin care.

LIGHT WAVELENGTHS FOR THE LED TREATMENT MASK:

Red light: **630nm**

Yellow light: **590nm**

Green light: **520nm**

Blue light: **470nm**

Using this mask is very simple, you just need to put it on, fasten it to your face with an elastic band and adjust the intensity and duration of the selected treatment with the controller. There are soft gray rubber pads around the eye holes to prevent unpleasant pressure around the eyes. A second mask that treats the skin of the neck can be attached to the face mask. The treatment LED mask for the face is equipped with 138 LEDs, the treatment LED mask for the neck is equipped with 54 LEDs. These diodes produce photons of precisely different wavelengths. Each diode has a defined power and luminous flux value. Therefore, it is important to distinguish the individual colors of the diodes. One color treats acne, another treats wrinkles, scars, redness, swelling - it depends on your preference.

The treatment LED mask includes 7 galvanic electrodes (microelectronic massage knobs), which do not need to be used. If you find the galvanic electrodes uncomfortable, the mask is fully functional even without them. The massage knobs help minimize wrinkles and fine lines and promote blood circulation. **The intensity of the microelectronic treatment (so-called Galvanic mode) can be set on the controller. Before putting on the mask, we recommend reducing the intensity to a minimum. Attention! Too high an intensity can have similar effects as an electric shock. The LEDs on the mask do not need to be changed.** The controller must be connected to the usb cable and the power adapter (included in the package).

EN TREATMENT LED MASK FOR THE FACE AND NECK

EFFECTS OF TREATMENT LED MASKS:

- It reduces wrinkles and prevents their appearance
- Minimizes pores
- Fights acne
- Helps fight against facial oiliness
- Improves skin elasticity
- Evens skin tone
- It suppresses the production of melanin
- Improves blood circulation in the subcutaneous tissue
- Increases the efficiency of absorption of active ingredients of cosmetic products


MAIN FUNCTIONS:

- Natural light waves improve the appearance of the skin
- Each color of light is effective for a specific skin treatment
- Various mode options
- Easy to control
- Easily portable

NOTICE:

IN CASE OF DISCOMFORT, THE LED NECK MASK CAN ALSO BE USED WITHOUT GALVANIC ELECTRODES.

MASK DESCRIPTION:

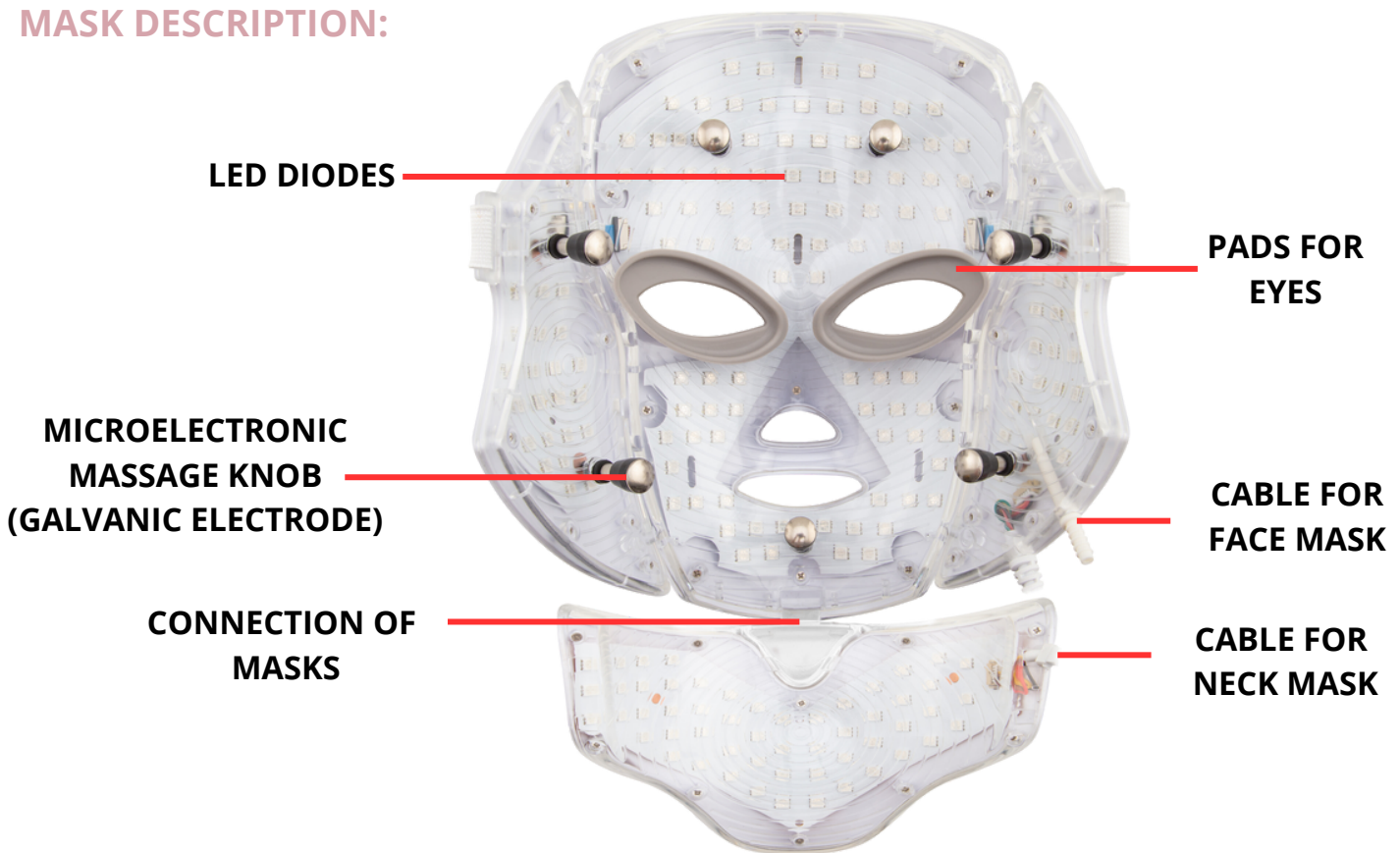THE MASK IS GREAT FOR PEOPLE WHO:


- Stay up late into the night and their skin becomes tired
- Sit in front of the computer for a long time and often look at their mobile phone
- Do not follow a healthy diet that improves skin condition
- Have problematic skin: spots, acne, clogged pores, dull, yellowish, rough or sagging skin, etc.

EN TREATMENT LED MASK FOR THE FACE AND NECK

INSTRUCTIONS FOR USE:

- Before using the mask, test all functions of the controller without the mask on your face. Try to avoid direct eye contact with the light coming from the mask.
- Before using the mask, we recommend setting all functions on the controller to minimum (see the description of the controller "default"). (too much intensity of the galvanic mode can have similar manifestations as an electric shock!). Also set a shorter treatment time (see recommended time of 10 minutes).
- Attach the neck mask to the face mask. Connect both masks to the controller using the USB cables.
- Connect the controller to the power cable. **The power cord must be plugged in when using the mask.**
- Adjust the length of the elastic band to your face and add 7 microelectronic massage tips to the mask.
- Before putting on the mask, thoroughly clean the skin of the face and neck and then gently dry with a towel.
- Start the skin treatment by setting the treatment time (using the "time" buttons) and then choose the desired treatment color (using the color select/start button).
- Put the mask on your face and start the whole skin treatment process with photon therapy.

CONTROLLER DESCRIPTION:

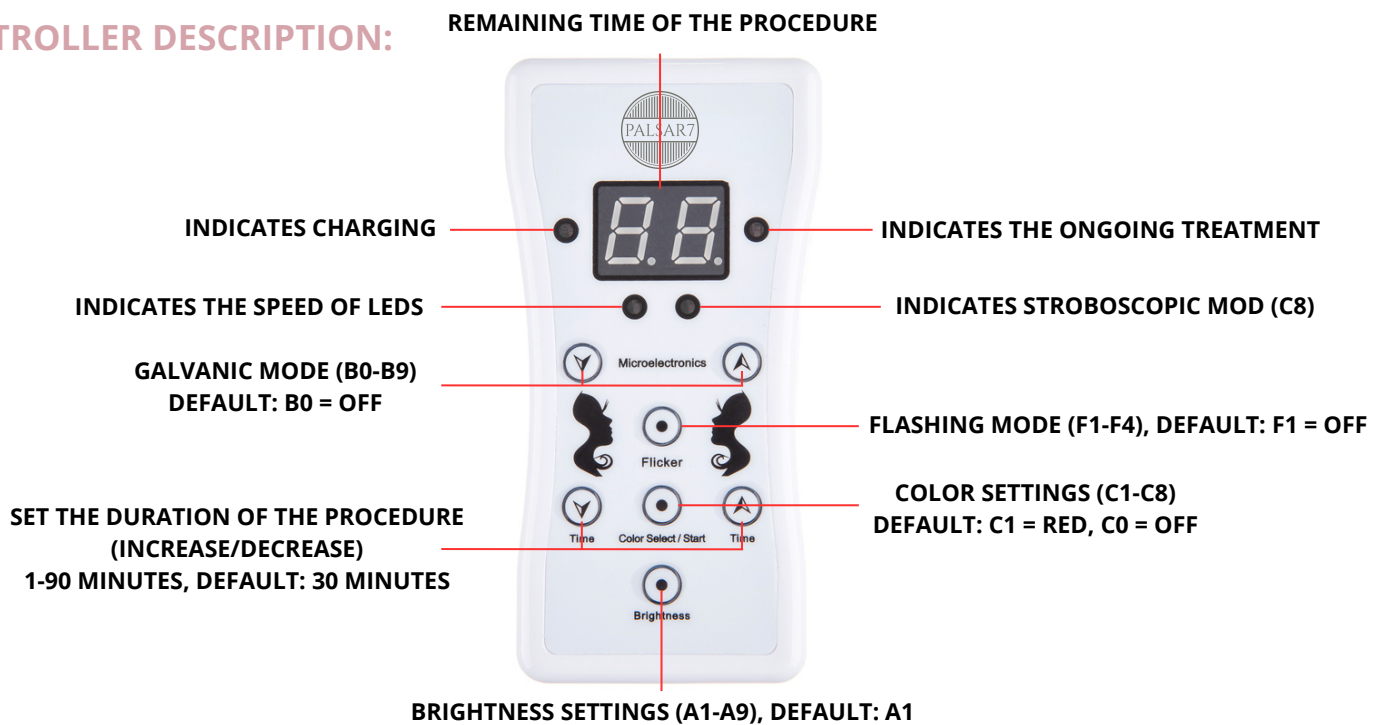The maximum time of using the mask daily is 45 minutes. If you want to try several treatments in one day, wait at least 10 minutes between them. The recommended time of one treatment is 10-15 minutes. Ideally, use the mask 2-3 times a week. If you suffer from sensitive skin, eczema, allergies or have a cold, consult your doctor before using the mask.

THE PICTURES WITH THE DESCRIPTION OF THE DRIVER AND MASK ARE FOR ILLUSTRATIVE PURPOSES ONLY.

EN TREATMENT LED MASK FOR THE FACE AND NECK

EFFECTS OF INDIVIDUAL LIGHTS:

1) 630NM-LED RED LIGHT: AGAINST WRINKLES

- Whitening effects
- For soft, wrinkle-free skin
- Smoothes fine wrinkles, shrinks pores
- Improves skin condition
- It supports the natural production of collagen
- It increases the efficiency of absorption of active substances of cosmetic preparations
- Improves blood circulation in the subcutaneous tissue
- Prevents the formation of pigment

2) 590NM-LED YELLOW LIGHT: AGAINST REDNESS, TO IMPROVE IMMUNITY

- Reduces skin roughness and redness
- Relieves skin diseases
- In addition, it relieves fever and helps to increase immunity
- Supports the absorption of cosmetics
- It contributes to the improvement of cellular metabolism
- Improves detoxification properties by supporting lymphatic activity
- Especially effective after intensive mechanical peeling
- It helps to slow down the aging of the skin

3) 520NM-LED GREEN LIGHT: AGAINST FATIGUE

- Improves skin condition
- Balances the color tone of the skin
- Reduces sebum production
- It relieves mental stress
- Reduces the occurrence of edema (swelling)
- Contributes to the elimination of epidermal pigment spots
- It slows down the production of melanin
- It has an anti-inflammatory effect (especially against the formation of acne skin)
- Suitable for all skin types

4) 470NM-LED BLUE LIGHT: AGAINST ACNE

- Effectively reduces the appearance and cures acne
- It improves the condition of the skin without leaving scars
- Activates porphyrins (they start a chemical reaction that destroys harmful bacteria)
- It has an antibacterial and anti-inflammatory effect

5) LED BRIGHT BLUE LIGHT (CALLED CYANINE): SUPPORT OF METABOLISM

- Gradually increases cellular energy
- Supports metabolism

6) LED VIOLET (PINK) LIGHT: AGAINST ACNE SCARS

- Combination of red and blue light = combination of two healing therapies
- Particularly effective in the treatment of acne scars

7) LED WHITE LIGHT: ANTI-AGING

- Treats age spots and sagging skin
- It reduces the appearance of fine wrinkles

8) STROBOSCOPE MODE - AUTOMATIC ALTERNATION OF LIGHTS

- It is not recommended for regular use

WARNING, RISKS:

- Not suitable for children under 14 years of age.
- Forbidden to use: pregnant women, infants, children, people with colds, angina or flu, epileptics, people allergic to light, patients with thyroid disease.
- **Before putting on the mask, reduce the intensity of the galvanic treatment to a minimum. WARNING! Too high an intensity can have similar effects as an electric shock.**
- Prohibited to use by people who have felt sick after using the mask - please visit your doctor and discuss the use of the mask with him
- Do not watch TV, surf the Internet, read a book, etc. while using the mask. – there is a risk of long-term damage to the eyes. Avoid prolonged direct eye contact with the light source.
- If you intend to use cosmetics at the same time as the mask, choose high-quality products with a natural composition.
- After opening the box, the mask may have a slight smell, this is because it is made of medical silicone, the smell will disappear soon.
- Use the mask in bed or on a lounger, always keep it in the box to prevent damage.
- Do not use water to clean the mask, do not immerse the mask in water - use a damp cloth to clean.
- Do not use any chemical products to clean the mask.
- Do not expose the mask to direct sunlight.

TECHNICAL PARAMETERS:

Input: 100-240VAC, 50/60HZ, Output: 12V(1A)

ACCESSORIES:

1x LED face mask, 1x LED neck mask, 7x galvanic electrode, 1x controller, 2x usb cable, 1x charging adapter, 1x user manual

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WARRANTY CARD

Product:
Model:
Order number:
Date of sale:
Warranty until:
Sale document number:
Seal of the store:

Information for Customers about Warranty Terms and Conditions:

Warranty shall apply to defects of purchased product, found during its warranty period. When purchasing a particular product, please consider carefully what properties and functions you expect your product to have. The fact that purchased product shall not satisfy your later requirements shall not give any reason for complaint. Please, prior to the first use of this product read carefully its operating instructions and follow these instructions consistently.

We recommend you to keep original product packing during the warranty period for possible transport of the product in such a way so as to prevent any product damage during its transport and handling.

For the purpose of due application of any complaint about the product we recommend you to ask Seller to fill-in and confirm this Warranty Card including identification of order number of purchased product, which shall comply with the number specified in the document proving product purchase and in the delivery note. In this connection Customer hereby takes into account that unless this Warranty Card specifies order number of purchased product, i.e., comparison of such order number with the number identified in this Warranty Card and document of product purchase/or delivery note shall not allow to confirm that particular product has come from Seller, in case of any doubts about the origin of such particular product Seller shall be entitled to automatically and without any other acts to reject the complaint of such goods.

Warranty shall not apply to:

- 1) Product defect that occurred by:
 - mechanical damage to the product or its parts
 - incorrect operating of the product
 - unprofessional intervention of a third party
 - by using the product contrary to its intended use in the enclosed operating instruction/user manual
 - improper storage or handling
 - transportation
 - natural disaster or other force measures
- 2) Product damage suffered by unqualified installation, its modification or by static electricity
- 3) Fault caused by external conditions (e.g. faults in the electrical network or residential installation, unsuitable charges, unsuitable operating conditions or damage by other equipment)
- 4) A product for which the buyer (operator) detains, conceals or provides incorrect information of importance to the assessment of the defect
- 5) External damage to the product that was clearly caused by man
- 6) Product damage inconsistent with the operating instructions, technical standards or valid safety regulations
- 7) Product damage resulting from normal wear and tear
- 8) Discharged batteries, decrease in battery capacity caused by its wear and tear
- 9) Product defects caused by the use in dusty, polluted or in any other way dissatisfactory environment pursuant to attached operating instructions

Wear and tear of the product and its components caused by their common use (discharged batteries, loss of the original color of the product, decreased capacity of accumulator, wear and tear of mechanical or displaying components, wear and tear of supply cables, etc.) shall be subject neither to defect warranty nor quality warranty in the moment of take-over.

Furthermore, complaints about product parameters, not included in the operating instructions or in any commercial-technical documentation of Seller, related to the product, shall not be subject to quality warranty in the moment of take-over.

Proof of product purchase, or delivery note, shall constitute an integral part of this Warranty Card.

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