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# EN Palsar7

**PROFESSIONAL TREATMENT LED MASK FOR THE FACE AND NECK WITH CHARGING STATION (WHITE-GOLD)**

## **PHOTON LIGHT THERAPY**

Why go for expensive chemical peels and skin cleansing, laser treatment and hyaluronic acid wrinkle filling when you can treat your skin from the comfort of your own home?

The luxurious treatment LED mask PALSAR7 is an innovative way to get rid of acne, wrinkles, tired skin and other imperfections. And this is thanks to photon light therapy, which is one of the most effective tools in the field of dermatology in the fight against chronic skin defects. Unlike laser treatments, photon treatment using LED lights does not cause any pain or side effects that need time to heal.

## **HOW DOES PHOTON LIGHT THERAPY WORK**

Light can be simply described as the visible part of electromagnetic radiation of various frequencies. It reaches different wavelengths in the interval from 415 to 635 nm. In this therapy, light penetrates the body as a stream of energy in the form of photons. Different wavelengths of light have different positive effects in skin care.

## **WAVELENGTHS OF THE LIGHTS FOR THE TREATMENT MASK:**

Red light: **635nm**

Orange light: **605nm**

Blue light: **415nm**

Using this mask is very simple, just put it on, fix it with two adjustable straps to the face and use the button on the side of the mask to set the color mode. There are soft gray rubber pads around the eye holes to prevent unpleasant pressure around the eyes. A second mask that treats the skin of the neck can be attached to the face mask. The treatment LED mask for the face and neck is equipped with a total of 92 LED diodes. These diodes produce photons of precisely different wavelengths. Each diode has a defined power and luminous flux value. Therefore, it is important to distinguish the individual colors of the diodes. Each color treats a different skin problem - it depends on your preference.

Before first use, remove the protective film from the inside of the mask.

**EN PROFESSIONAL TREATMENT LED MASK FOR THE FACE AND NECK**


## EFFECTS OF THE TREATMENT

### MASK:

- It reduces wrinkles and prevents their appearance
- Minimizes pores
- Fights acne
- Helps fight against facial oiliness
- Improves skin elasticity
- Evens skin tone
- It suppresses the production of melanin
- Improves blood circulation in the subcutaneous tissue
- Increases the efficiency of absorption of active ingredients of cosmetic products

## MAIN FUNCTIONS:

- Natural light waves improve the appearance of the skin
- Each color of light is effective for a specific skin treatment
- Various mode options
- Easy to control
- Easily portable
- Elegant and timeless design
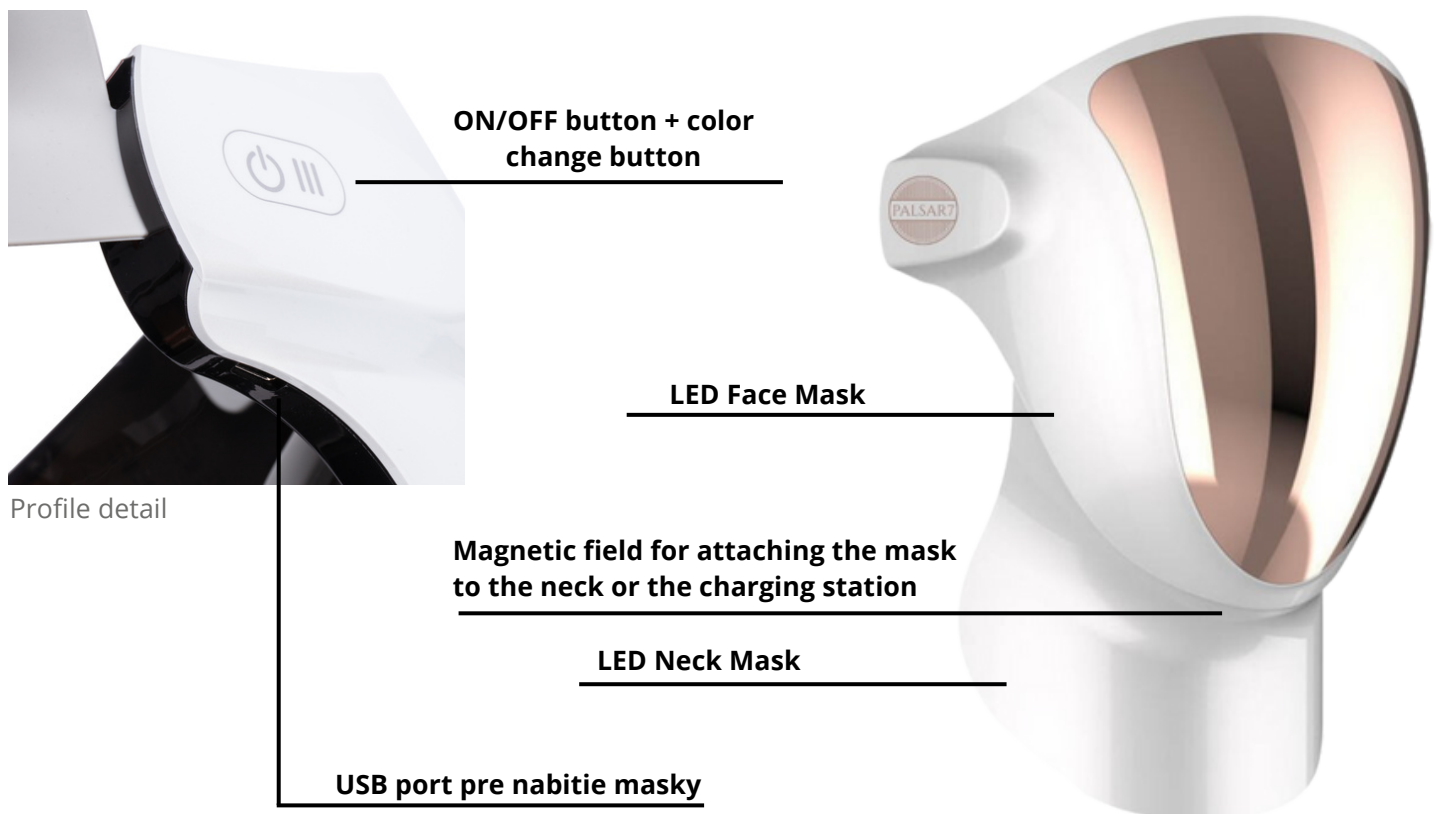
## MASK IS SUITABLE FOR PEOPLE WHO:


- Stay up so late that their skin is tired.
- Sit in front of the computer screen or mobile screen for long periods of time.
- Do not follow a healthy diet that improves skin condition
- Have problematic skin: spots, acne, clogged pores, dull, yellowish, rough or sagging skin, etc.

**EN PROFESSIONAL TREATMENT LED MASK FOR THE FACE AND NECK**

## INSTRUCTIONS FOR USE:

- **Before using the mask, remove the protective film on the inside of the mask.**
- Attach the neck mask to the face mask accordingly, the mask can also be used without the neck part.
- Before putting on the mask, thoroughly clean the skin of the face and neck and then gently dry with a towel.
- Start the skin treatment by holding the button on the side of the mask for 3 seconds. You can then switch color modes by pressing the button. (Red, Blue, Orange)
- Put the mask on your face and adjust the length of the adjustable straps to your face, so they fit you perfectly and comfortably.
- Start the whole skin treatment process with photon therapy.
- Mask color modes automatically turn off after a certain amount of time.
- Red = 5 minutes, Blue = 5 minutes, Orange = 10 minutes
- The recommended time of use is 20 minutes, so you can alternate all modes for complete care, or choose a mode according to your preferences.

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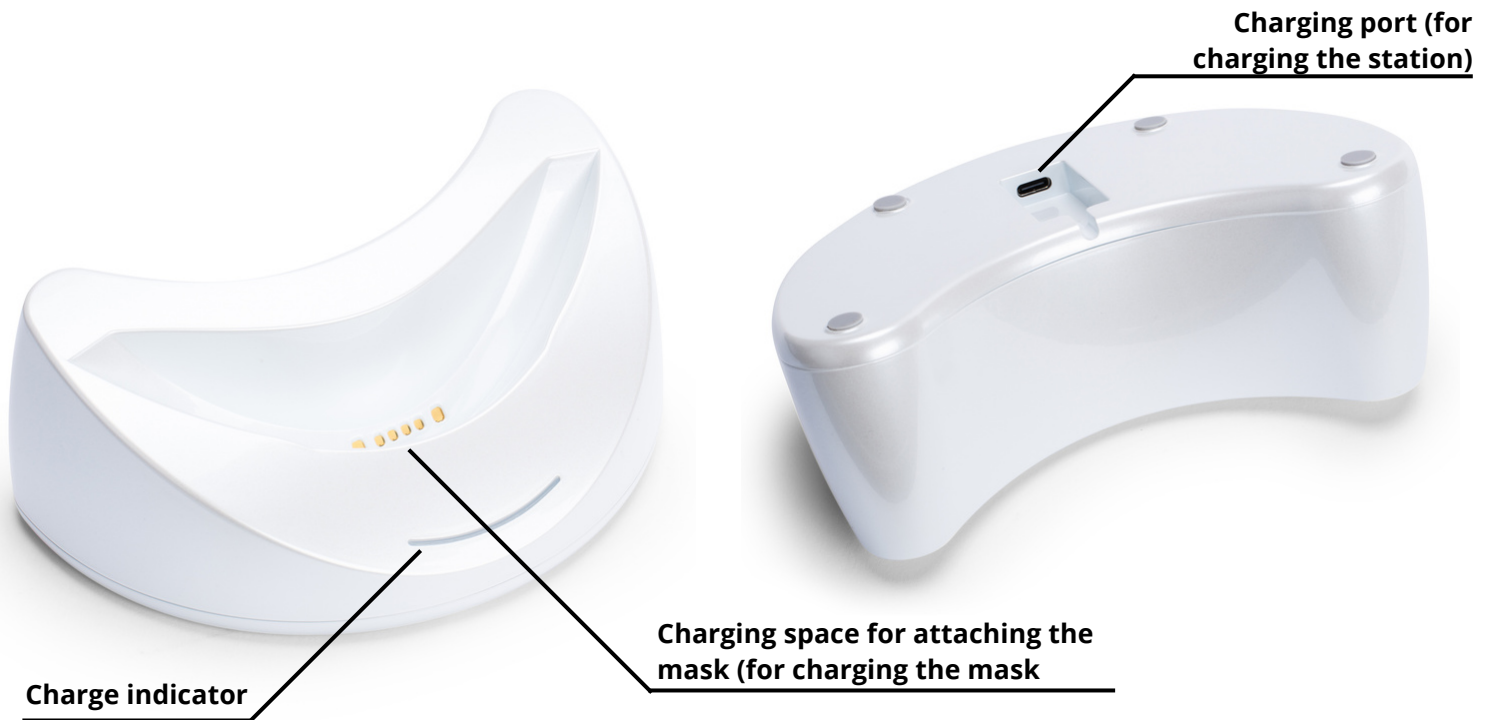**The maximum time of using the mask daily is 20 minutes. The recommended time of one treatment is 20 minutes, during one treatment you can change all color modes or choose your preferred one. Ideally, use the mask 2-3 times a week. If you suffer from sensitive skin, eczema, allergies or have a cold, consult your doctor before using the mask.**

THE PICTURES SHOWN ARE FOR ILLUSTRATION ONLY.

**EN PROFESSIONAL TREATMENT LED MASK FOR THE FACE AND NECK**

## EFFECTS OF INDIVIDUAL LIGHTS:

### 1) 635NM-LED RED LIGHT: AGAINST WRINKLES

- Whitening effects
- For soft, wrinkle-free skin
- Smooths fine wrinkles, shrinks pores
- Improves skin condition
- It supports the natural production of collagen
- Increases the efficiency of absorption of active ingredients of cosmetic products
- Improves blood circulation in the subcutaneous tissue
- Prevents the formation of pigment

### 2) 605 NM-LED ORANGE LIGHT: AGAINST FATIGUE

- Reduces skin roughness and redness
- Alleviates skin diseases
- In addition, it relieves heat and promotes an increase in immunity
- Supports the absorption of cosmetic products
- It contributes to the improvement of cellular metabolism
- Improves detoxification properties by supporting lymphatic activity
- Especially effective after intensive mechanical peeling
- Helps slow down skin aging

### 3) 415 NM-LED BLUE LIGHT: AGAINST ACNE

- Effectively reduces the appearance and cures acne
- Improves skin condition without leaving scars
- Activates porphyrins (they start a chemical reaction that destroys harmful bacteria)
- It has an antibacterial and anti-inflammatory effect

## HOW TO CHARGE THE MASK:

- The package includes a charging station that simplifies the charging process. Just detach the neck treatment part from the mask and attach the mask to the charging station on the underside.
- The charging station must first be charged using a USB-C cable. After charging the station, it can be used to charge a separate mask. The charging station will make it easier for you to handle the mask, especially when traveling.
- The mask can be charged directly, without using a charging station. Just plug the included USB-C cable into the hole near the mask launch button.
- When charging through the charging station or separately, a red light will start flashing near the charging cable. When the light turns blue, the mask is fully charged.

## TECHNICAL PAREMETERS

Input voltage: 3.7V

Battery capacity of the charging station: 1400mAh

Material: ABS, Silicone gel

Charging voltage: DC 5V - 1A

## WARNING, RISKS:

- Not suitable for children under 14 years of age.
- Forbidden to use for: pregnant women, infants, children, people with colds, angina or flu, epileptics, people allergic to light, patients with thyroid disease.
- The product is not heat-resistant, protect it from fire and caustics.
- Forbidden to use by people who felt sick after using the mask - visit your doctor and consult with him about the use of the mask.
- When using the mask, do not look into the LEDs - there is a risk of long-term eye damage. Avoid prolonged direct eye contact with the light source.
- If you intend to use cosmetic products at the same time as the mask, choose quality products with a natural composition.
- After opening the box, the mask may have a slight smell, this is because it is made of medical silicone, the smell will disappear soon.
- When the mask is not in use, keep the mask in the box or in the charging station, on a flat dry surface, away from magnetic fields. This will prevent it from being damaged.
- Do not use water to clean the mask, do not immerse the mask in water - use a damp cloth to clean.
- Do not use any chemical products to clean the mask.
- Do not expose the mask to direct sunlight.

## ACCESSORIES:

1x LED treatment mask for the face, 1x LED treatment mask for the neck, 1x usb cabel, 1x charging station, users manual

# EN PROFESSIONAL TREATMENT LED MASK FOR THE FACE AND NECK

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EU standard: CE, FCC, RoHS Certification

FCC CE RoHS