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EN PALSAR7

RECHARGEABLE SILICONE TREATMENT LED MASK FOR THE FACE (BLACK)

PHOTON LIGHT THERAPY

Why go for expensive chemical peels and skin cleansing, laser treatment and hyaluronic acid wrinkle filling when you can treat your skin from the comfort of your own home?

The luxurious treatment LED mask PALSAR7 is an innovative way to get rid of acne, wrinkles, tired skin and other imperfections. And this is thanks to photon light therapy, which is one of the most effective tools in the field of dermatology in the fight against chronic skin defects. Unlike laser treatments, photon treatment using LED lights does not cause any pain or side effects that need time to heal.

HOW DOES PHOTON LIGHT THERAPY WORK

Light can be simply described as the visible part of electromagnetic radiation of various frequencies. It reaches different wavelengths in the interval from 415 to 635 nm. In this therapy, light penetrates the body as a stream of energy in the form of photons. Different wavelengths of light have different positive effects in skin care.

WAVELENGTHS OF LIGHT FOR THE TREATMENT LED MASK:

NIR (Near infrared light): **850nm**

Red light: **630nm**

Blue light: **415nm**

Orange light: **605nm**

Using this mask is very simple, just put it on, fasten it to your face with the two adjustable straps and use the attached controller to set the desired mode. Rubber pads can be attached around the eye holes to prevent uncomfortable pressure around the eyes. The treatment LED face mask is equipped with a total of 368 LEDs. These LEDs produce photons with just different wavelengths. Each diode has a defined power and luminous flux value. It is therefore important to distinguish between the different colours of the diodes. Each colour treats a different skin problem - it depends on your preference.

EN RECHARGEABLE SILICONE TREATMENT LED FACE MASK

EFFECTS OF THE TREATMENT


MASK:

- It reduces wrinkles and prevents their appearance
- Minimizes pores
- Fights acne
- Helps fight against facial oiliness
- Improves skin elasticity
- Evens skin tone
- It suppresses the production of melanin
- Improves blood circulation in the subcutaneous tissue
- Increases the efficiency of absorption of active ingredients of cosmetic products

MAIN FUNCTIONS:

- Natural light waves improve the appearance of the skin
- Each color of light is effective for a specific skin treatment
- Various mode options
- Easy to control
- Easily portable
- Elegant and timeless design

CONTROLLER

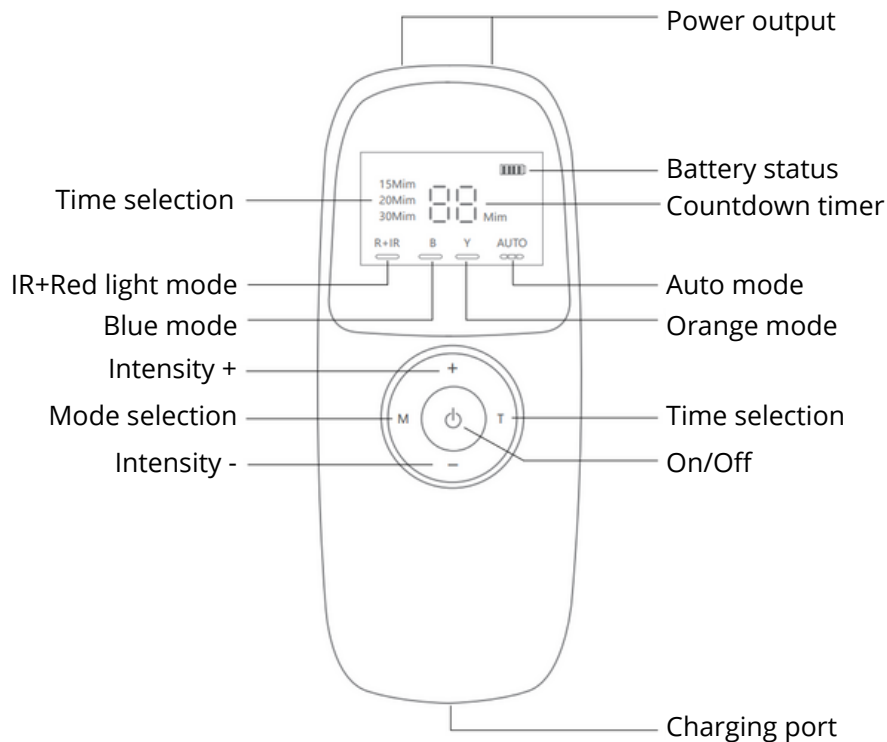MASK IS SUITABLE FOR PEOPLE WHO:

- Stay up so late that their skin is tired.
- Sit in front of the computer screen or mobile screen for long periods of time.
- Do not follow a healthy diet that improves skin condition
- Have problematic skin: spots, acne, clogged pores, dull, yellowish, rough or sagging skin, etc.

EN RECHARGEABLE SILICONE TREATMENT LED FACE MASK

INSTRUCTIONS FOR USE:

- Before putting on the mask, cleanse the face and neck thoroughly and then dry gently with a towel.
- Connect the mask to the controller using the cable. **Always connect the cable from the mask to the controller via any port on the top (both 2 ports have the same functionality). The bottom port is only for charging the controller.**
- Attach the leather straps to the mask and attach the mask to your face.
- Start the skin treatment by holding the button on the side of the mask for 3 seconds.
- Begin the entire skin treatment process with photon therapy.
- There is a basic automatic mode that alternates between all color modes. 1. Infrared and red light (5 minutes), 2. Blue light (5 minutes), 3. Orange light (10 minutes).
- You can then press the "M" button to switch color modes using the manual mode. (Red, Blue, Yellow)
- With the "T" button, the time of the selected therapy can be adjusted in manual mode.
- Use the "+" and "-" buttons to adjust the intensity of the lights.

Step 1: Remove makeup and clean your face.

Step 2: Use with a mask or a serum. Secure the device to your head using straps.

Step 3: Long press the power button for 3s to start the device. It runs automatically for 20 minutes.

Step 4: Remove the mask and apply moisturizer to the skin.

The maximum time of using the mask daily is 30 minutes. The recommended time of one treatment is 30 minutes, during one treatment you can change all color modes or choose your preferred one. Ideally, use the mask 3-5 times a week. If you suffer from sensitive skin, eczema, allergies or have a cold, consult your doctor before using the mask.

THE PICTURES SHOWN ARE FOR ILLUSTRATION ONLY.

EN RECHARGEABLE SILICONE TREATMENT LED FACE MASK

THE EFFECTS OF INDIVIDUAL LIGHTS:

1) 630NM-LED RED LIGHT + 850NM-LED INFRARED LIGHT: ANTI-WRINKLE

- Whitening effects
- For soft, wrinkle-free skin
- Smooths fine lines, shrinks pores
- Improves skin condition
- Promotes natural collagen production
- Increases the efficiency of absorption of active ingredients in cosmetic products
- Improves blood circulation in the subcutaneous tissue
- Prevents pigment formation

2) 415 NM-LED BLUE LIGHT: ANTI-ACNE

- Effectively reduces the appearance and treats acne
- Improves skin condition without leaving scars
- Activates porphyrins (starts a chemical reaction that kills harmful bacteria)
- Has antibacterial and anti-inflammatory effect

3) 605 NM-LED ORANGE LIGHT: ANTI-FATIGUE

- Relieves skin roughness and redness
- Relieves skin diseases
- In addition, it relieves hotness and promotes increased immunity
- Promotes absorption of cosmetic products
- Contributes to improved cellular metabolism
- Improves detoxifying properties by promoting lymphatic activity
- Particularly effective after intensive mechanical peeling
- Helps slow down skin ageing

WARNING, RISKS:

- Not suitable for children under 14 years of age.
- Forbidden to use for: pregnant women, infants, children, people with colds, angina or flu, epileptics, people allergic to light, patients with thyroid disease.
- The product is not heat-resistant, protect it from fire and caustics.
- Forbidden to use by people who felt sick after using the mask - visit your doctor and consult with him about the use of the mask.
- When using the mask, do not look into the LEDs - there is a risk of long-term eye damage. Avoid prolonged direct eye contact with the light source.
- If you intend to use cosmetic products at the same time as the mask, choose quality products with a natural composition.
- After opening the box, the mask may have a slight smell, this is because it is made of medical silicone, the smell will disappear soon.
- When the mask is not in use, keep the mask in the box or in the charging station, on a flat dry surface, away from magnetic fields. This will prevent it from being damaged.
- Do not use water to clean the mask, do not immerse the mask in water - use a damp cloth to clean.
- Do not use any chemical products to clean the mask.
- Do not expose the mask to direct sunlight.
- Keep all cables in a neutral position and do not twist them to avoid mechanical damage to the wires inside the protective rubber.
- Do not attempt to disconnect, rotate, or manipulate the connecting cable that is connected to the mask. The cable is firmly connected to the mask and cannot be separated from the mask.

CHARGING THE DEVICE:

- When the battery indicator starts flashing, the controller must be charged.
- Connect the charging cable to the controller and allow it to fully charge.

TECHNICAL PARAMETERS:

Input voltage: 3.7V

Rated power: 4.5W

Controller battery capacity: 2000mAh

Material: food grade silicone, plastic, synthetic leather

Charging voltage: DC 5V - 1A

ACCESSORIES:

1x LED face mask, 2x fastening leatherette strips, 1x controller, 1x cable for connecting the controller and mask, 1x USB cable, 2x eye protector, user manual

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